


FLUID MANAGEMENT

COMPONENTS

v.029

 **Nordson**
MEDICAL

Tube Size Reference

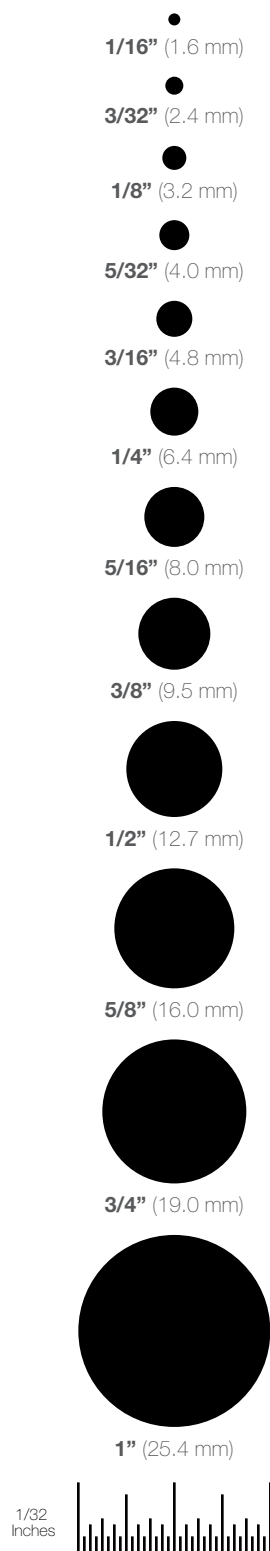

Barb Reference

Pg 138 Barb Series' Characteristics
Pg 138 Barb Dimensions
Pg 140 Barb Performance

Common Material Reference

-001 (Various, see part description)
-06 Orange Nylon
-07 Yellow Nylon
-0 Natural Nylon
-1 White Nylon
-2 Black Nylon
-3 Red Nylon
-4 Green Nylon
-5 Blue Nylon
-7 Glass Filled Natural Nylon
-8 Gray Nylon
-8B Gray Nylon
-9 Clear Polycarbonate
-10 Acetal
-40 Polysulfone
-50 Acrylic
-81 Lace White ABS
-91 White Polycarbonate
-1006 Natural Acetal
-1007 Black Acetal
-6005 Animal-free Polypropylene
-6006 White Polypropylene
-6012 Black Polypropylene
-8003 Natural ABS
-8012 White ABS
-9002 Rad. Stable Polycarbonate
-9010 Rad. Stable Polycarbonate
-9024 Polycarbonate
-B Buna-N O-ring
-CM006 Polyethylene
-CM018 Copolyester
-CM030 Dark Blue Nylon
-E PTEE/EPDM
-EP EPDM O-ring
-J1A Natural Kynar PVDF
-J1A2 Natural Kynar PVDF
-PX Platinum-cured Silicone
-U71 White Class VI Glass Filled Nylon
-U75 Blue Class VI Glass Filled Nylon
-VP1 Animal-free Polyethylene
-X2 Buna-N Rubber
-XO Stainless Steel

Resin Information

Pg 146 Chemical Resistance
Pg 147 Regulatory Status
Pg 148 Sterilization Stability
Pg 149 Typical Properties

INTRODUCTION.....	4	Swabable Break-Off Valves.....	89
COMPONENTS & TECHNOLOGIES.....	6	4-Way Swabable Stopcock.....	89
LOCATIONS.....	7	Swabable "Y" Valves.....	89
WHAT'S NEW.....	8	Swabable Straight Valves.....	89
GLOBAL PARTNERS.....	9	Bondable Valves.....	89
DESIGN CENTER.....	10	Swabable T-Port Valves.....	89
LUERS.....	12	DIVERT™ STOPCOCKS.....	90
Male.....	13	High Pressure Stopcocks.....	90
Female.....	15	Medium Pressure Stopcocks.....	91
Adapters and Couplers.....	17	Low Pressure Stopcocks.....	91
Tees.....	17	High Flow Stopcocks.....	93
Elbows.....	17	BIOPROCESS FITTINGS.....	94
Plugs and Caps.....	18	Bag Ports.....	94
Luer to Thread.....	19	Opposing Barb Ports.....	98
Panel Mounts.....	20	BioClamp®.....	99
Panel Mount Accessories.....	21	BioValve®.....	100
TUBE FITTINGS.....	22	PharmaLok™ Clamps.....	101
Straight Through.....	23	PharmaLok™ Tube Clamps.....	104
Straight Through Reducer.....	24	Spaulding™ Aseptic Disconnects.....	105
Elbow.....	27	CYLINDRAFlow™ Manifolds.....	107
Elbow Reducer.....	28	Sanitary Fittings - Mini Flange.....	108
Tee.....	29	Sanitary Fittings - Maxi Flange.....	109
Tee Reducer.....	30	Sanitary Fittings - TC Reducers.....	110
Y.....	32	Instrument Tees.....	111
Y Reducer.....	34	Gaskets.....	112
Multiport.....	34	BLOOD PRESSURE.....	114
Multiport Reducer.....	35	BPF Series.....	114
Press-In Plugs.....	36	Bayonet Male Cuff Connectors.....	115
Ultrasonically Weldable Barbs.....	37	Bayonet Panel Mount Connectors.....	115
Panel Mounts.....	37	Dinamap® Monitor Adapters.....	116
Panel Mount Accessories.....	37	Dinamap® Cuff Connectors.....	116
ISO 80369.....	38	1/4 Turn Connectors.....	117
ISO 80369-3 (Enteral Feeding) Compliant Connectors.....	39	THREADED.....	118
ISO 80369-5 (Blood Pressure) Compliant Connectors.....	39	10-32 UNF, Special Tapered Thread.....	119
NEURAXIAL CONNECTORS.....	40	1/4-28 UNF.....	122
Male Connectors.....	41	1/4-28 UNF Bottom Sealing.....	123
Plugs.....	41	Metric M6x1.....	124
Accessories.....	41	Metric M5x.8.....	124
QUICK CONNECT COUPLINGS.....	42	1/16-27 NPT.....	125
10 Series - 1/16" Flow Size.....	42	1/8-27 NPT.....	126
20 Series - 1/8" Flow Size.....	46	1/4-18 NPT Reducers.....	127
30 Series - 1/8" Flow Size.....	54	BUILD-A-PART®.....	128
33 Series - 1/8" Flow Size.....	55	Junction Blocks.....	129
40 Series - 1/4" Flow Size.....	56	Accessories.....	130
42 Series - 3/32" Flow Size.....	63	Barbs.....	130
50 Series - 1/4" Flow Size.....	63	Threads.....	131
55 Series - 1/4" Flow Size.....	68	TUBING.....	132
60 Series - 3/8" Flow Size.....	68	Linear Low Density Polyethylene.....	132
65 Series - 1/2" Flow Size.....	75	Thermoplastic Elastomer (TPE).....	132
70 Series - 1/4" Flow Size.....	75	Platinum-Cured Silicone.....	133
100 Series - 1/16" Flow Size.....	76	Pump Grade Platinum-Cured Silicone.....	133
MQC Series.....	78	Polyurethane.....	133
RQ Series.....	79	Anti-Microbial Polyvinyl Chloride (PVC).....	133
SBL Series.....	80	OTHER COMPONENTS.....	134
XQ Series.....	81	Tapered Seal.....	134
Custom Solutions.....	83	Clamps.....	137
CHECK VALVES.....	84	Injection Sites.....	137
Single Check Valves - Anti-Siphon.....	84	Spikes.....	137
High Flow Check Valves.....	85	TECHNICAL INFORMATION.....	138
Single Check Valves - One-Way.....	85	Barb Dimensions.....	138
Single Check Valves - Barbed.....	86	Barb Performance.....	140
Double Check Valves.....	87	Quick Connect Performance Ratings.....	141
NEEDLEFREE SWABABLE VALVES.....	88	BPF Series Properties and Performance Data.....	145
Vial Caps.....	88	Panel Mount Instructions.....	145
Tube End Valves.....	89	Resin Information.....	146
Swabable Bag Valves.....	89	Trademarks.....	149
Valved Male Luers.....	89	Index.....	150

Welcome to Nordson MEDICAL!

Thank you for your interest in our products. If you're a current customer, we appreciate your business and look forward to exceeding your expectations in the future. If you're seeing us for the first time, we welcome the opportunity to earn your trust.

Nordson MEDICAL is backed by the extensive resources of Nordson Corporation (Nasdaq: NDSN), the world leader in precision dispensing, fluid management and related technologies.

We are your global integrated solutions partner in the design, development, and manufacturing of complex medical devices and component technologies. We serve interventional, surgical, and specialized markets as the engine of innovations that save or enhance patients' lives.

We work with companies at any point in the product lifecycle to bring their innovative ideas to life. As a single-source partner, we enable our customers to save costs and speed time to market.

Our online store is the largest in the industry, giving customers fast access to 7,000+ components and technologies in stock and ready to ship, most within 24 hours.

Our three areas of focus provide customers with a complete solution:

Design & Development

Our 200+ technical team brings broad medical device expertise and a deep understanding of clinical applications. Services include:

- Concept generation
- Product design & engineering
- Rapid prototyping
- Project management
- New product introduction
- Testing & validation
- Quality & regulatory services

Components & Technologies

Our extensive portfolio of stock and custom components gives customers one-stop access to the critical building blocks for their complex, innovative medical devices.

- Fluid management components
- Catheters & cannulae
- Medical balloons
- Medical tubing
- Biomaterial delivery devices
- Nitinol components

Finished Device Solutions

We can leverage our in-house component and technology expertise to design, develop, and manufacture finished medical devices and subassemblies. This vertical integration helps teams reduce costs, mitigate risk, streamline supply chain, and ensure consistent quality. Our services include:

- Full-service manufacturing
- Finished devices & subassemblies
- Global manufacturing sites
- Dedicated transfer team
- Supply chain management
- Robust continuous improvement

What We Do

Our fluid management portfolio includes over 9,600 components designed specifically for flexible tubing in different style and material combinations. The Value Plastics product line includes luer fittings, check valves, bioprocessing fittings, tube fittings, threaded fittings, and blood pressure monitoring components with tubing IDs ranging from 1/16" to 1". Our LinkTech line of quick connect couplings allows you to quickly and easily connect or disconnect fluid transfer lines and is available in a wide range of flow sizes. We also offer custom component manufacturing, tubing, and select medical OEM components. Our application expertise extends to bioprocessing, blood pressure, cardiovascular, ophthalmic, and more. Our products find global application in demanding healthcare OEM, research and specialty industrial applications.

Our Quality Commitment to You

At Nordson MEDICAL, quality is designed in from the outset; it's not an add-on at the end of the process. We are an FDA registered manufacturing facility and our Quality Management System is designed as a system of interrelated processes meeting the requirements of ISO 13485. We welcome our customers to conduct an on-site supplier audit with the knowledge that our systems have satisfied the most demanding requirements of medical device manufacturers worldwide. We are dedicated to exceeding your expectations.

Clean Room Manufacturing

Most Nordson MEDICAL products are manufactured and packaged in an environment that meets the requirements of ISO 14644-1 Class 8 (Class 100,000) standards. We utilize both manual and automated assembly processes tailored to manufacture the highest quality components.

Assembly and Finishing Capabilities

Nordson MEDICAL's quality products are sophisticated and reliable. We have a solid reputation in today's medical device and biotech markets with pressure and leak testing. We also offer manual, semi-automatic, and fully automatic finished good assembly.

Design, Engineering and R&D

Nordson MEDICAL's staff of product development professionals uses the most modern equipment to design the best products available anywhere. Whether your need is for new products or custom solutions, our focus on healthcare OEMs means we understand your needs and don't need to be coached to your industry requirements. Designs incorporating human factors, finite element analysis, DFMA, rapid prototyping and concurrent engineering are just some of the tools used to bring you the components you require, when you want them.

Manufacturing Consistency

Nordson MEDICAL is equipped with multiple injection molding machines including plastic injection molding and over molding including 2-shot and insert molding. Our all-electric injection molding presses are chosen for shot-to-shot consistency and outfitted with the most up-to-date process monitoring software available. Each operator is responsible for product quality at the press and has every tool at their disposal. All molding, secondary operations and packaging takes place in a controlled environment to assure component readiness for your production floor.

Rapid Prototyping

Nordson MEDICAL's engineering department is a leader in the industry. Using rapid prototyping and an in-house, high performance 3D printer, we can utilize multiple resins to provide certain material properties to components. We can execute highly representative prototypes and build a mold in two days, allowing us to deliver solutions with testable parts in a short time.

We use a system to verify potential product performance to expedite quality product design and development. The use of several tooling technologies allows us to incorporate innovative materials and processes to provide plastic components in a fraction of the time it may take to utilize metals.

We also execute detailed processes for evaluating application requirements to ensure the design focuses on ease of use, cost efficiency, design manufacturability, and top notch performance.

Regulatory Management

Nordson MEDICAL strives to provide our customers with the most up-to-date technical information, directly from the resin manufacturer, on the resins we use. We maintain a comprehensive list of all material and regulatory information, enabling you to select the most appropriate products for your particular application. Information we provide includes:

- Regulatory data
- Material data
- MSDS
- Typical resin properties
- Sterilization stability information
- Chemical resistance information

Visit www.nordsonmedical.com for technical information.

Customer Satisfaction

As with all of our medical parts, each bag comes with a product certification and traceable lot number. Our parts are manufactured and bagged in a clean room environment to accelerate your own usage.

At Nordson MEDICAL, we are dedicated to exceeding your expectations. As such, it is important to us that when you place an order you know exactly what to expect, when to expect it and that we will do everything we can to ensure that you receive components made to the required specifications.

Nordson MEDICAL implements and regulates an internal customer support system. Our goal is to provide you exceptional service and quality solutions. We track all support issues in order to meet ISO requirements, enhance customer satisfaction, and to continuously improve our initiatives.

When a support issue concerns one of our products, we begin an immediate investigation of our available stock and quarantine any suspected parts. We then complete all required formal corrective actions and inform the customer when the issue is resolved.


To contact customer service, please call us at (888) 404-5837 or (970) 267-5200.

Contact Us

Visit us at www.nordsonmedical.com for all your low pressure fluid management needs.

Fluid Management Components

With the largest selection in the industry, we offer cost-effective solutions that connect fluid management systems. Our components join flexible tubing assemblies, allowing you to easily regulate the flow of valuable media. Process management components, like quick connect fittings, luers, and tube-to-tube fittings, effortlessly and securely connect different elements of fluid administration devices. We also offer devices like stopcocks and check valves that control the transfer of fluids by regulating directional flow.

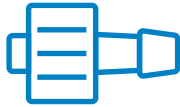


Catheters and Cannulae

Nordson MEDICAL can partner with you in all phases of catheter development, from design and development through manufacturing and assembly. In-house capabilities include custom balloon extrusions, balloon blowing, bonding, welding, and folding, coating and masking, printing, and more. Our expertise in engineered shafts includes braid- and coil-reinforced shafts. We can also utilize polymer solution casting technologies to build multiple performance characteristics and varying wall thicknesses and diameters into a single-piece construction, while maintaining tight tolerances.




Medical Balloons

We have the widest selection of medical balloons in the industry and can create any geometry with any material and virtually no limits on design. We offer the most complex medical balloons with the tightest tolerances, highest burst pressures, thinnest walls, and multiple unique design options, including coatings, textured surfaces, braided, and multilumen balloons.

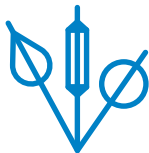


Medical Tubing

From simple medical tubing to complex, tight-tolerance extruded tubing in a wide range of materials, we can meet your needs. We're known for producing the world's thinnest, smallest, and strongest heat shrink tubing, and our tubing line has expanded to include specialized items like polyimide tubing, PTFE liners, and a full range of single-lumen, multilumen, and multilayer extruded tubing.

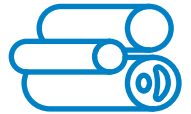

Biomaterial Delivery Devices

Nordson MEDICAL provides proven solutions for accurate and precise delivery of tissue sealants, hemostatics, adhesion barriers, and other autologous or synthetic biomaterials for surgical applications. We develop biomaterial delivery systems that accurately deliver a single biomaterial—or precisely mix and deliver a dual-component therapy—to the surgical site. Our OsteoXpress™ and OsteoPrecision™ Graft Delivery Devices deliver allograft, autograft, or synthetic bone graft materials to an orthopedic surgical site, and we provide stock, custom, and private label products to a wide range of customers.

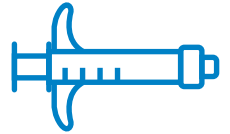

Nitinol Components

We are experts in complex manual nitinol wire braiding and forming and can achieve unique and complex geometries not possible by other methods. Symmetric and asymmetric structures are possible with spiral, wavy, or straight tip attachment options.



UNITED STATES:**Loveland, CO**

805 West 71st Street
 Loveland, CO 80538, USA
 Phone: 970.267.5200
 Toll Free: 888.404.5837
 FMF-sales@nordsonmedical.com

MEXICO:**Guaymas, Mexico**

Carretera Int No. 15 Km129 Salida Norte
 Parque Industrial Roca Fuerte, Edificio 21
 Guaymas, Sonora C.P. 85400, Mexico
 FDS-sales@nordsonmedical.com

OTHER NORDSON MEDICAL LOCATIONS**IRELAND:****Boyle, Ireland**

Abbeytown
 Boyle, County Roscommon, Ireland
 Phone: 53.71.966.3038
 EXT-sales@nordsonmedical.com

Galway, Ireland

Ballybrit Upper Industrial Estate
 Monivea Road
 H91 XP49, Galway, Ireland
 DND-sales@nordsonmedical.com

ISRAEL:**Katzrin, Israel**

20 Street
 Katzrin 1290000, Israel
 Phone: 972.4.6738600
 NWF-sales@nordsonmedical.com

Rosh Ha'Ayin, Israel

22 Hamelacha Street, Building 4-B
 Rosh Ha'Ayin 48091, Israel
 Phone: 972.76.5399300
 DND-sales@nordsonmedical.com

UNITED STATES:**Chattanooga, TN**

4701 Wilson Road, Suite A
 Chattanooga, TN 37410, USA
 Phone: 423.648.7838
 FCM-sales@nordsonmedical.com

Huntington Beach, CA

7612 Woodwind Drive
 Huntington Beach, CA 92647, USA
 Phone: 657.215.4200
 PSC-sales@nordsonmedical.com

Louisville, CO

1772 Prairie Way, Suite K
 Louisville, CO 80027, USA
 Phone: 303.440.0884
 DND-sales@nordsonmedical.com

Marlborough, MA

261 Cedar Hill Street
 Marlborough, MA 01752, USA
 Phone: 508.481.6233
 DND-sales@nordsonmedical.com

Minneapolis, MN

7125 Northland Terrace N, Suite 200
 Minneapolis, MN 55428, USA
 Phone: 763.291.4302
 DND-sales@nordsonmedical.com

Salem, NH

29 Northwestern Drive
 Salem, NH 03079, USA
 Phone: 603.327.0600
 BAL-sales@nordsonmedical.com

San Jose, CA

105 East Tasman Drive
 San Jose, CA 95134, USA
 Phone: 408.730.0771
 DND-sales@nordsonmedical.com

St. Paul, MN

1270 Eagan Industrial Road
 St. Paul, MN 55121, USA
 Phone: 651.452.1977
 BMD-sales@nordsonmedical.com

Sunnyvale, CA

610 Palomar Avenue
 Sunnyvale, CA 94085, USA
 Phone: 408.730.0771
 DND-sales@nordsonmedical.com

MRI-Compatible Valved Quick Connect Coupling

Nordson MEDICAL is excited to introduce the first MRI-compatible valved quick connect in a compact footprint for medical applications. The 33 Series push-to-connect option utilizes a metal-free valving solution. The connector is cleanroom manufactured and ideally suited for applications such as biopsy sampling, wound therapy, dialysis, hydration, etc. The built-in duck bill check valve opens to allow a maximized unidirectional flow with minimal turbulence under pressure or vacuum. These quick connects are suitable for applications needing a differentiated connection other than the standard luer. **PG-55**



4" Bag Port

Nordson MEDICAL's high tech series of bag ports offers a variety of benefits to the biopharmaceutical industry through a unique design. Value Pharma™ resin holds credentials that many manufacturers deem fundamental, such as animal derivative free qualification, USP Class VI certification, and the guarantee of heat sealing and conformance with polyethylene bags. The resin then goes even further to provide a strong, 24-hour resistance to alcohol stress cracking. Their unique design sets these ports apart from existing solutions. The most notable of the bag ports' specifications are the parabolic lead and the alignment ribs. These two features, respectively, reduce the amount of time it takes to drain a bag by up to 24 percent and improve flow dynamics through the port. **PG-96**



Spaulding™ Aseptic Disconnects

Providing a leak-free connection, the Spaulding™ Series Aseptic Disconnects eliminate the need for pinch clamps and tube welders. Smooth disconnects are achievable with the easy-to-use design, and valuable fluids are protected by valves that automatically shut off flow when the fittings are disconnected. Available in USP Class VI, ADCF Polysulfone, validation and test results confirm their suitability in pre-sterilized systems and are also gamma radiation stable. These connectors are interchangeable with similar plastic connectors. **PG-105**



Looking for Something Different?

Most of our best product ideas have come from you, our customers. While we pride ourselves on our extensive product offering, we recognize that your needs are ever-changing and that our products, too, must adapt to your needs.

If there's a new product you think Nordson MEDICAL should make, please send us an Email at designcenter@nordsonmedical.com describing your product concept in as much detail as you desire. We will contact you with any questions, and who knows - the next catalog edition could include products YOU told us to build.

See page 10 for more information on the Design Center!



Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

EUROPE:**Jean Gastaldi**

Sales Manager Europe
Phone: +33 687 753 458
jean.gastaldi@nordsonmedical.com

AUSTRIA	FRANCE	ITALY
MONACO	PORTUGAL	SPAIN
	SWITZERLAND	

Promepla

Contact: Francoise Audouy or Petra Gerner
Phone: +377 97 98 42 32
sales@promepla.com
www.promepla.com

IRELAND	UNITED KINGDOM
----------------	-----------------------

The West Group

Contact: Stuart Barnett
Phone: +44 23 9226 6031
sales@westgroup.co.uk
www.westgroup.co.uk

BELGIUM	CZECH REPUBLIC
DENMARK	FINLAND
GERMANY	HUNGARY
LUXEMBOURG	NETHERLANDS
NORWAY	POLAND
SWEDEN	TURKEY

MedNet

Contact: Matthias Heinz
Phone: +49 251 32266-0
info@medneteurope.com
Israel Contact: David Thaler
david.thaler@medneteurope.com
Israeli Email: info@mednetil.com
www.medneteurope.com

SOUTH AMERICA:

ARGENTINA	BOLIVIA	COLUMBIA
CHILE	ECUADOR	PERU
URUGUAY	VENEZUELA	

R.D.V. International

Contact: Galo H. Gutierrez
Phone: +593-3-285-4369
rdvusarep@comcast.net
galohgutierrez@comcast.net

BRAZIL:**Trusty Comércio e Representações Ltda. (TCR)**

Phone: 55 12 3132-3418
rui@tcr-brazil.com.br
www.tcr-brazil.com.br

AUSTRALIA/NEW ZEALAND**John Morris Scientific**

Phone: +61 2 9496 4200
Enquiry Phone: 1800 251 799
info@johnmorris.com.au
www.johnmorrisgroup.com

ASIA:**Steven Sun**

Account Manager, East China & Taiwan
Phone: +86 3866 9166-725
steven.sun@nordsonmedical.com

Stanley Xu

Account Manager, South & North China
Phone: +86 136 8683 8521
stanley.xu@nordsonmedical.com

CHINA:**LINK (Shanghai) Fluid Technology Co., Ltd**

Phone: +86 21 64958516
sales@linkfluidics.com
www.linkfluidics.com

Shenzhen Precise Technology Co., Ltd

Phone: +86 755 83903104
sales@precise-pipeline.com
www.precise-pipeline.com

Shanghai Lepure Biotech Co., Ltd (BioProcessing)

Phone: +86 21 57791703 667
marketing@lepurechina.com
www.lepurechina.com

Austar Lifesciences Limited (BioProcessing)

Phone: +86 311 85115868
jennaxie@austar.com.cn
www.austar.com.cn

INDIA:**Cole-Parmer India Pvt. Ltd.**

Phone: 022-6139-4444
support@coleparmer.in
www.coleparmer.in

INDONESIA:**P.T. Air Dynamics**

Phone: 62-21-29451011
sales-id@airdynamics.com.sg
www.airdynamics.com.sg

P.T. Kian Soon Hardware

Phone: +62 778 466868
sales@kiansoon.com.sg
www.kiansoon.com.sg

JAPAN:**International Scientific Instruments Supply Co. LTD (ISIS)**

Phone: +81-6-6325-1406
mail@isis-ltd.co.jp
www.isis-ltd.co.jp

Air Precision Co., Ltd

Phone: +81-45-228-2230
info@airprecision.co.jp
www.airprecision.co.jp

Bio-Pharma Support Co., Ltd. (BioProcessing)

Phone: +81-3-3263-7257
info@bp-s.com
www.bp-s.com

KOREA:**Hanmi Rubber & Plastics**

Phone: +82-2-2275-5885
order@hmrp.co.kr
www.hmrp.co.kr

Noanix Corporation

Phone: +82-43-905-0972
sales@noanix.com
www.noanix.com

BT Resources (BTR) (BioProcessing)

Phone: +82 31 8003 5885
ga@btresources.co.kr
www.btresources.co.kr

MALAYSIA:**PSI Pneumatic Control Sdn Bhd**

Phone: 60-4-6442627
sales-psi@airdynamics.com.sg
www.airdynamics.com.sg

Kian Soon Mechanical Components Sdn Bhd

Phone: +(604) 6460 318
sales@kiansoon.com.sg
www.kiansoon.com.sg

PHILIPPINES:**Air Dynamics, Inc.**

Phone: 63-2-7532881
inquiry@airdynamics.com.sg
www.airdynamics.com.sg

SINGAPORE:**Air Dynamics Pte Ltd**

Phone: 65-67466884
sales@airdynamics.com.sg
www.airdynamics.com.sg

Kian Soon

Phone: +(65) 62911177
sales@kiansoon.com.sg
www.kiansoon.com.sg

PureSys (BioProcessing)

Phone: +(65) 6291 1177
sales@puresys.com.sg
www.puresys.com.sg

TAIWAN:**Cintrade Enterprise, Inc (BioProcessing)**

Phone: +886-2-8941-3600
info@cintrade.com.tw
www.cintrade.com.tw

THAILAND:**Kian Soon CTC Mechanical Components Co., Ltd**

Phone: +(66) 2 682-9119
sales@kiansoon.com.sg
www.kiansoon.com.sg

DESIGN

Center



“It would have taken me two to three times as long and two to three times as much money without Nordson MEDICAL. This is who you want to work with. They have become part of our team.”



Fitting The World's Most Fluid Designs

How it Works:



Define Your Product Requirements:

Describe your application and requirements or provide us with a drawing or simple sketch. From here we can determine what details are needed and can help define the most ideal and cost-effective connector concept for your application.



Call the Design Center to Explore Your Connector Options:

There may be a suitable off-the-shelf solution or perhaps we can adapt an existing connector to your application. Often the modification of an existing mold can be significantly faster and less expensive.



Collaborate:

Working with your requirements, we can help walk you through the most cost effective solutions to achieve your desired performance. That may mean simply taking 2 or 3 connectors and identifying a low-cost means to roll them into a single connector.



Go-to-Market:

As we arrive at a solution together, we closely involve your team in the design and development process to ensure that as we approach production you have everything you need for qualification to get you to market both quickly and safely. The rest is up to you!



What happens when “off-the-shelf” just doesn’t work? The Design Center builds a new shelf!

Your Solution. Faster. For Less.

With project requirements and timelines becoming increasingly demanding, it's critical to partner with a team who understands fluid management solutions and the rigors of today's medical device and biotech markets. Team up with the Nordson MEDICAL Design Center for your custom connector needs, whose market leading service, quality and delivery you've already come to know.

Here's What You Can Expect From Our Design Center and the “Specialists in Fluid Connections”:

Save Time – Talk to someone who understands your application requirements. The Design Center works with you to determine the most suitable material and design options to deliver the best performing product as cost-effectively as possible. We also utilize state of the art rapid prototyping techniques and in-house testing facilities to keep your project moving quickly and economically.

Save Money – Start with our existing product offering of over 4,000 fluid connectors, or utilize some of our existing mold configurations toward perfectly customizing a part to your application. Don't forget, our technical and design assistance is FREE!

Save Trouble – Ensure that everything around your project will be done right the first time. Our dedicated and experienced staff will make sure that the finest details are covered from the time that we review your product requirements to the final product qualification.

The Design Center can also assist you in developing fluid connectors that comply to emerging standards. Nordson MEDICAL is an active member of the ISO 80369 standards group defining new standard connectors for small bore, medical device connections.

Custom Assembly – Our knowledgeable team of engineers offers capabilities to the biotechnical market through custom-made assemblies created in our design center. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

To get started today contact your Nordson MEDICAL Design Services Manager:
Toll Free - +1 (888) 404-5837
Email - designcenter@nordsonmedical.com

LUER

Fittings

Male Luers
PG-13



Female Luers
PG-15



Adapters
PG-17



Couplers
PG-17



Tees
PG-17



Elbows
PG-17



Plugs and Caps
PG-18



Luer to Thread
PG-19



Panel Mounts & Accessories
PG-20



Our Luer Tapers

Nordson MEDICAL's luer tapers are made to the ISO 594-1 and ISO 594-2 standards, and may be connected to other luers meeting these ISO standards. Please test your application to verify the performance of the interconnection.

Nordson MEDICAL's luer fittings come in over 100 configurations and several materials for the device designer, joining tubing from 1/16" (1.6 mm) ID to 1/4" (6.4 mm) ID.

Features & Benefits:

Published ISO Design Standard

- Well understood worldwide
- Often assures interconnections without need to retrofit

Rotating Lock Permits Tube Positioning Prior To Tightening

- Minimizes twisting and kinking of tubing

Tapered Cone Seal

- Reliable
- Inexpensive
- Uncomplicated

Available In Many Materials

- Facilitates optimal material selection for your application

Inexpensive

- Facilitates single-use

Easy To Understand

- Ease of use in field

Many Configuration Options

- Adaptable to wide range of applications

Security Options (Fixed Lock, Rotating Lock, No Lock)

- Tailorable to your specific application

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

MTLL210 (-1 -2 -3 -4 -5 -06 -6005 -6012 -07 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



MTLL410 (-1 -6005 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 400 Series Barb, 1/16" (1.6 mm) ID Tubing



MTLL004 (-1 -2 -3 -4 -5 -06 -6005 -07 -9 -9002)

Male Luer Integral Lock Ring to 500 Series Barb, 1/16" (1.6 mm) ID Tubing



MTLL220 (-1 -2 -3 -4 -5 -06 -6005 -6012 -07 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



MTLL420 (-1 -6005 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 400 Series Barb, 3/32" (2.4 mm) ID Tubing



MTLL007 (-1 -2 -3 -4 -5 -6005 -06 -07 -9 -9002)

Male Luer Integral Lock Ring to 500 Series Barb, 3/32" (2.4 mm) ID Tubing



MTLL230 (-1 -2 -3 -4 -5 -06 -6005 -6012 -07 -9 -9002 -J1A -CM550)

Male Luer Integral Lock Ring to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



MTLL430 (-1 -3 -6005 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 400 Series Barb, 1/8" (3.2 mm) ID Tubing



MTLL013 (-1 -2 -3 -4 -5 -06 -6005 -07 -9 -9002)

Male Luer Integral Lock Ring to 500 Series Barb, 1/8" (3.2 mm) ID Tubing



MTLL240 (-1 -2 -3 -4 -5 -6005 -6012 -06 -07 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



MTLL025 (-1 -2 -3 -4 -5 -6005 -06 -07 -9 -9002)

Male Luer Integral Lock Ring to 500 Series Barb, 5/32" (4.0 mm) ID Tubing



MTLL445 (-1 -6005 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 400 Series Barb, .170" (4.3 mm) ID Tubing



MTLL250 (-1 -2 -3 -4 -5 -06 -6005 -07 -8003 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 200 Series Barb, 3/16" (4.8 mm) ID Tubing



MTLL450 (-1 -6005 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 400 Series Barb, 3/16" (4.8 mm) ID Tubing



MTLL035 (-1 -2 -3 -4 -5 -06 -6005 -07 -9 -9002)

Male Luer Integral Lock Ring to 500 Series Barb, 3/16" (4.8 mm) ID Tubing



MTLL360 (-1 -2 -6005 -6012 -9 -9002 -J1A)

Male Luer Integral Lock Ring to 300 Series Barb, 1/4" (6.4 mm) ID Tubing



MTLL055 (-1 -2 -3 -4 -06 -6005 -07 -9 -9002)

Male Luer Integral Lock Ring to 500 Series Barb, 1/4" (6.4 mm) ID Tubing



Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

MALE LUER:

Slip, or for use with Lock Ring

FSLLR (-1 -2 -3 -4 -5 -06 -6005 -6012 -07 -9 -9002)

Finger Snap Luer Lock Ring

**MLRL004 (-1 -6005 -9)**

Male Luer to 500 Series Barb, 1/16" (1.6 mm) ID Tubing (May be used with separate rotating lock ring; FSLLR)

**MLSL004 (-1 -6005 -9)**

Male Luer to 500 Series Barb, 1/16" (1.6 mm) ID Tubing (May be used with separate stationary lock ring; FSLLR)

**MLRL007 (-1 -6005 -6012 -9)**

Male Luer to 500 Series Barb, 3/32" (2.4 mm) ID Tubing (May be used with separate rotating lock ring; FSLLR)

**MLSL007 (-1 -6005 -9)**

Male Luer to 500 Series Barb, 3/32" (2.4 mm) ID Tubing (May be used with separate stationary lock ring; FSLLR)

**MLRL013 (-1 -6005 -6012 -9)**

Male Luer to 500 Series Barb, 1/8" (3.2 mm) ID Tubing (May be used with separate rotating lock ring; FSLLR)

**MLSL013 (-1 -6005 -9)**

Male Luer to 500 Series Barb, 1/8" (3.2 mm) ID Tubing (May be used with separate stationary lock ring; FSLLR)

**MLRL025 (-1 -6005 -6012 -9)**

Male Luer to 500 Series Barb, 5/32" (4.0 mm) ID Tubing (May be used with separate rotating lock ring; FSLLR)

**MLSL025 (-1 -6005 -9)**

Male Luer to 500 Series Barb, 5/32" (4.0 mm) ID Tubing (May be used with separate stationary lock ring; FSLLR)

**MLRL035 (-1 -6005 -6012 -9 -9002)**

Male Luer to 500 Series Barb, 3/16" (4.8 mm) ID Tubing (May be used with separate rotating lock ring; FSLLR)

**MLSL035 (-1 -6005 -9)**

Male Luer to 500 Series Barb, 3/16" (4.8 mm) ID Tubing (May be used with separate stationary lock ring; FSLLR)

**MLRL055 (-1 -6005 -9)**

Male Luer to 500 Series Barb, 1/4" (6.4 mm) ID Tubing (May be used with separate rotating lock ring; FSLLR)

**MLSL055 (-1 -6005 -9)**

Male Luer to 500 Series Barb, 1/4" (6.4 mm) ID Tubing (May be used with separate stationary lock ring; FSLLR)

**It's a Lock with Nordson MEDICAL!**

Our MLRL and MLSL Series of fittings can be used both as a slip luer or as a locking luer with the addition of our FSLLR luer lock ring, which allows color coding without flow path contact. Our MLSL/FSLLR combination results in a stationary lock ring. Use a combination of MLRL and FSLLR to create a color coded rotating lock ring.



MALE LUER:

Slip

MTLS210 (-1 -6005 -9 -J1A)

Male Luer Slip to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



MSL004 (-1 -6005 -9 -9002)

Male Luer Slip to 500 Series Barb, 1/16" (1.6 mm) ID Tubing



MTLS220 (-1 -3 -5 -6005 -6012 -81 -9 -J1A)

Male Luer Slip to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



MSL007 (-1 -6005 -9 -9002)

Male Luer Slip to 500 Series Barb, 3/32" (2.4 mm) ID Tubing



MTLS230 (-1 -3 -5 -6005 -81 -9 -9002 -J1A)

Male Luer Slip to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



MSL013 (-1 -6005 -9 -9002)

Male Luer Slip to 500 Series Barb, 1/8" (3.2 mm) ID Tubing



MTLS240 (-1 -2 -5 -6005 -81 -9 -J1A)

Male Luer Slip to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



MSL025 (-1 -6005 -9 -9002)

Male Luer Slip to 500 Series Barb, 5/32" (4.0 mm) ID Tubing



MSL035 (-1 -6005 -9 -9002 -J1A)

Male Luer Slip to 500 Series Barb, 3/16" (4.8 mm) ID Tubing



MSL055 (-1 -6005 -9 -9002)

Male Luer Slip to 500 Series Barb, 1/4" (6.4 mm) ID Tubing



FEMALE LUER:

Thread Style

FTLL210 (-1 -2 -3 -4 -5 -6005 -6012 -9 -9002 -J1A)

Female Luer Thread Style to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



FTLL410 (-1 -6005 -9 -9002 -J1A)

Female Luer Thread Style to 400 Series Barb, 1/16" (1.6 mm) ID Tubing



FTLL004 (-1 -2 -3 -4 -5 -06 -6005 -07 -9 -9002 -J1A)

Female Luer Thread Style to 500 Series Barb, 1/16" (1.6 mm) ID Tubing



FTLL220 (-1 -2 -3 -5 -6005 -6012 -9 -9002 -J1A)

Female Luer Thread Style to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



FTLL420 (-1 -6005 -9 -9002 -J1A)

Female Luer Thread Style to 400 Series Barb, 3/32" (2.4 mm) ID Tubing



FTLL007 (-1 -2 -3 -4 -5 -06 -6005 -07 -9 -9002)

Female Luer Thread Style to 500 Series Barb, 3/32" (2.4 mm) ID Tubing



FTLL230 (-1 -2 -3 -5 -6005 -6012 -9 -9002 -J1A -CM550)

Female Luer Thread Style to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



FTLL430 (-1 -6005 -9 -9002 -J1A)

Female Luer Thread Style to 400 Series Barb, 1/8" (3.2 mm) ID Tubing



FTLL013 (-1 -2 -3 -4 -5 -06 -6005 -07 -9 -9002)

Female Luer Thread Style to 500 Series Barb, 1/8" (3.2 mm) ID Tubing



FEMALE LUER: Thread Style

FTLL240 (-1 -2 -3 -5 -6005 -6012 -9 -9002 -J1A)

Female Luer Thread Style to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



FTLL025 (-1 -2 -3 -4 -5 -06 -6005 -07 -9 -9002)

Female Luer Thread Style to 500 Series Barb, 5/32" (4.0 mm) ID Tubing



FTLL445 (-1 -6005 -9 -9002 -J1A)

Female Luer Thread Style to 400 Series Barb, .170" (4.3 mm) ID Tubing



FTLL250 (-1 -3 -5 -6005 -9 -9002 -J1A)

Female Luer Thread Style to 200 Series Barb, 3/16" (4.8 mm) ID Tubing



FTLL450 (-1 -6005 -9 -9002 -J1A)

Female Luer Thread Style to 400 Series Barb, 3/16" (4.8 mm) ID Tubing



FTLL035 (-1 -2 -3 -4 -5 -06 -6005 -07 -9 -9002 -J1A)

Female Luer Thread Style to 500 Series Barb, 3/16" (4.8 mm) ID Tubing



FTLL360 (-1 -3 -5 -6005 -6012 -9 -9002 -J1A)

Female Luer Thread Style to 300 Series Barb, 1/4" (6.4 mm) ID Tubing



FTLL055 (-1 -2 -3 -4 -5 -06 -6005 -07 -9 -9002)

Female Luer Thread Style to 500 Series Barb, 1/4" (6.4 mm) ID Tubing



FEMALE LUER:

Smooth Bore Thread Style

FTLLSB004 (-1 -6005 -9)

Female Luer Thread Style Smooth Bore to 500 Series Barb, 1/16" (1.6 mm) ID Tubing



FTLLSB007 (-1 -4 -6005 -9)

Female Luer Thread Style Smooth Bore to 500 Series Barb, 3/32" (2.4 mm) ID Tubing



FTLLSB013 (-1 -6005 -9)

Female Luer Thread Style Smooth Bore to 500 Series Barb, 1/8" (3.2 mm) ID Tubing



Why Choose a Smooth Bore?

Our Smooth Bore Luers have a smooth inner transition providing for easy insertion of guide wires, probes and other devices.

FEMALE LUER:

Lug Style

FTL10 (-1 -6005 -9 -J1A)

Female Luer Lug Style to Classic Series Barb, 1/16" (1.6 mm) ID Tubing



FTL210 (-1 -6005 -9 -J1A)

Female Luer Lug Style to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



FTL20 (-1 -6005 -07 -81 -9 -J1A)

Female Luer Lug Style to Classic Series Barb, 3/32" (2.4 mm) ID Tubing



FTL220 (-1 -3 -5 -6005 -6012 -9 -J1A)

Female Luer Lug Style to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



FTL30 (-1 -6005 -6012 -9 -J1A)

Female Luer Lug Style to Classic Series Barb, 1/8" (3.2 mm) ID Tubing



FEMALE LUER:

Lug Style

FTL230 (-1 -3 -5 -6005 -9 -J1A)

Female Luer Lug Style to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



FTL240 (-1 -6005 -9 -J1A)

Female Luer Lug Style to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



LUER:

Adapters, Couplers, and Lock Rings

FSLLR (-1 -2 -3 -4 -5 -06 -6005 -6012 -07 -9 -9002)

Finger Snap Luer Lock Ring



MLRLC (-1 -6005 -9)

Male Luer Coupler (May be used with separate rotating lock rings; FSLLR)



MTLCS (-1 -6005 -9 -J1A)

Male Luer Slip Coupler



FTLLC (-1 -6005 -9 -9002 -J1A)

Female Luer Thread Style Coupler




FTLC (-1 -6005 -9 -9002 -J1A)

Female Luer Lug Style Coupler



LC34 (-1 -6005 -9 -9002 -J1A)

Male Luer Integral Lock Ring to Female Luer Thread Style Coupler

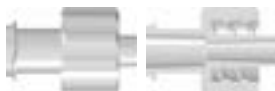

LC78 (-1 -6005 -9 -9002 -J1A)

Male Luer to Female Luer Thread Style Coupler (May be used with separate rotating lock ring; FSLLR)



LT878 (-1 -6005 -9 -9002)

Female Luer Thread Style (2) to Male Luer Tee (May be used with separate rotating lock ring; FSLLR)



LT787 (-1 -6005 -8012 -9 -9002)

Male Luer (2) to Female Luer Thread Style Tee (May be used with separate rotating lock ring; FSLLR)



FTLT (-1 -6005 -9 -J1A)

Female Luer Lug Style Tee



FTLLE (-1 -6005 -9 -J1A)

Female Luer Thread Style Elbow



FTLE (-1 -6005 -9 -J1A)

Female Luer Lug Style Elbow



LE87 (-1 -6005 -8003 -9 -9002)

Male Luer to Female Luer Thread Style Elbow (May be used with separate rotating lock ring; FSLLR)



LE7220 (-1 -2 -9 -6012)

Male Luer to 200 Series Barb Elbow, 3/32" (2.4 mm) ID Tubing (May be used with separate rotating lock ring; FSLLR)



LUER: Adapters, Couplers, and Lock Rings

LE7230 (-1 -2 -6005 -6012 -9 -9010)

Male Luer to 200 Series Barb Elbow 1/8" (3.2 mm) ID Tubing (May be used with separate rotating lock ring; FSLLR)



CA091 (-001) Gamma Stable Polycarbonate

Male Luer Elbow to 3/16" Classic Series Barb, (CF361-9002) with Rotating Lock Ring (FSLLR-9002)



LUER: Plugs and Caps

MTLLP (-1 -2 -3 -6005 -6012 -07 -9 -J1A)

Male Luer Integral Lock Ring Plug, Closed at Luer Tip



LP4 (-1 -2 -3 -4 -5 -6005 -9 -9002 -J1A)

Male Luer Integral Lock Ring Plug, Closed at Grip



VLP4 (-1 -3 -6005 -81 -J1A)

Male Luer Integral Lock Ring Vented Plug, Closed at Grip



MLRLP (-1 -6005 -9)

Male Luer Plug (May be used with separate rotating lock ring; FSLLR)



MLSLP (-1 -6005 -9)

Male Luer Plug (May be used with separate stationary lock ring; FSLLR)



MTLP (-1 -6005 -9 -J1A)

Male Luer Plug, Closed at Tether Loop (May be used with separate stationary lock ring; MLLR)



MTLSP (-1 -6005 -81 -9 -J1A)

Male Luer Slip Plug, Closed at Grip



FTLP (-1 -6005 -9 -J1A)

Female Luer Lug Style Cap with Tether Loop



FTLLP (-1 -2 -3 -5 -6005 -6012 -07 -81 -9 -9002 -J1A)

Female Luer Thread Style Cap



LP34 (-1 -3 -5 -6005 -9 -9002 -J1A)

Male Luer Integral Lock Ring to Female Luer Thread Style Plug



Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!

www.nordsonmedical.com

ABSML-6005 UNDER (-1 -2 -3 -4 -5 -06 -07 -6012)

Male Luer with Rotating Lock Ring and 5/16" Hex to 1/4-28 UNF Bottom Sealing Thread



ABSML-J1A UNDER (-1 -2 -3 -4 -5 -06 -07)

Male Luer with Rotating Lock Ring and 5/16" Hex to 1/4-28 UNF Bottom Sealing Thread



XMTLL (-1 -6005 -9 -J1A)

Male Luer Integral Lock Ring to 10-32 Special Tapered Thread



SMTLL (-1 -6005 -6012 -9 -J1A)

Male Luer Integral Lock Ring to 1/4-28 UNF Thread



18MTLL (-1 -2 -6005 -6012 -J1A)

Male Luer Integral Lock Ring with 7/16" Hex to 1/8-27 NPT Thread



14MTLL (-1 -2 -6005 -6012)

Male Luer Integral Lock Ring with 9/16" Hex to 1/4-18 NPT Thread



XMLRL (-1 -6005 -9 -9002)

Male Luer with 1/4" Hex to 10-32 Special Tapered Thread (May be used with separate rotating lock ring; FSLLR)



SMLRL (-1 -6005 -9 -J1A)

Male Luer with 5/16" Hex to 1/4-28 UNF Thread (May be used with separate rotating lock ring; FSLLR)



BSMLRL (-6005 -J1A)

Male Luer with 5/16" Hex to 1/4-28 Bottom-Sealing Thread (May be used with separate rotating lock ring; FSLLR)



Looking for Something Different?

Most of our best product ideas have come from you, our customers. While we pride ourselves on our extensive product offering, we recognize that your needs are ever-changing and that our products, too, must adapt to your needs.

If there's a new product you think Nordson MEDICAL should make, please send us an Email at designcenter@nordsonmedical.com describing your product concept in as much detail as you desire. We will contact you with any questions, and who knows - the next catalog edition could include products YOU told us to build.

See page 10 for more information on the Design Center!



Need More Details?

Part drawings, minimum order quantities, and CAD requests can be found online. Visit www.nordsonmedical.com or contact Sales at FMF-sales@nordsonmedical.com.

LUER TO THREAD:
Female

XFTLL (-1 -6005 -9 -J1A)

Female Luer Thread Style with 1/4" Hex to 10-32 Special Tapered Thread



KFTLL (-1 -6005 -9 -J1A)

Female Luer Thread Style with 1/4" Hex to 10-32 UNF Thread



SFTLL (-1 -2 -4 -06 -6005 -6012 -9 -J1A)

Female Luer Thread Style with 5/16" Hex to 1/4-28 UNF Thread



BSFTLL (-1 -6005 -J1A)

Female Luer Thread Style with 1/4" Hex to 1/4-28 UNF Bottom Sealing Thread



XFTL (-1 -6005 -9 -J1A)

Female Luer Lug Style with 1/4" Hex to 10-32 Special Tapered Thread



KFTL (-1 -6005 -9 -J1A)

Female Luer Lug Style with 1/4" Hex to 10-32 UNF Thread



SFTL (-1 -6005 -9)

Female Luer Lug Style with 5/16" Hex to 1/4-28 UNF Thread



18FTLL (-1 -6005 -6012 -J1A)

Female Luer Thread Style with 7/16" Hex to 1/8-27 NPT Thread



18LFTLL (-1 -4 -5 -6005 -6012 -J1A)

Female Luer Thread Style Elbow with 7/16" Hex to 1/8-27 NPT Thread



PANEL MOUNTS:
Female Thread Style

FTLLB210 (-1 -2 -6005 -9 -J1A)

Female Luer Thread Style Panel Mount 1/4-28 UNF to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



FTLLB220 (-1 -2 -6005 -9 -J1A)

Female Luer Thread Style Panel Mount 1/4-28 UNF to 200 Series Barb, 3/32" (2.4 mm) ID Tubing




FTLLB230 (-1 -2 -6005 -9 -J1A)

Female Luer Thread Style Panel Mount 1/4-28 UNF to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



FTLLBMLRL (-1 -6005)

Female Luer Thread Style Panel Mount 1/4-28 UNF to Male Luer (May be used with separate rotating lock ring; FSLLR)

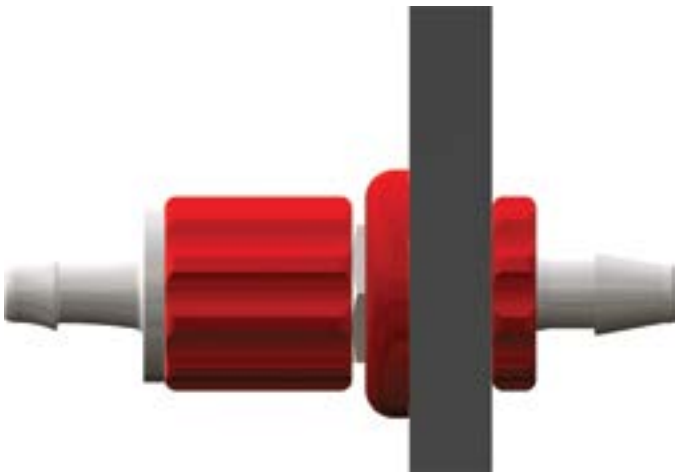

Fully Assembled Panel-Mounted Luer

Need More Details?

Additional threaded solutions are shown on **PG-118-127**. Check out all of our threaded style parts at www.nordsonmedical.com or contact Sales at FMF-sales@nordsonmedical.com.

Did You Know?

Nordson MEDICAL has an interactive PDF available online at www.nordsonmedical.com. Access part drawings, reference charts, data sheets, and more from just one location!

PANEL MOUNTS:

Female Lug Style and Male Luer

FTLB210 (-1 -6005 -9 -J1A)

Female Luer Lug Style Panel Mount 1/4-28 UNF to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



FTLB220 (-1 -6005 -9 -J1A)

Female Luer Lug Style Panel Mount 1/4-28 UNF to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



FTLB230 (-1 -6005 -9 -J1A)

Female Luer Lug Style Panel Mount 1/4-28 UNF to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



MLRLB210 (-1 -6005)

Male Luer Panel Mount 10-32 UNF to 200 Series Barb, 1/16" (1.6 mm) ID Tubing, (May be used with separate rotating lock ring; FSLLR)



MLRLB220 (-1 -6005)

Male Luer Panel Mount 10-32 UNF to 200 Series Barb, 3/32" (2.4 mm) ID Tubing, (May be used with separate rotating lock ring; FSLLR)



MLRLB230 (-1 -6005)

Male Luer Panel Mount 1/4-28 UNF to 200 Series Barb, 1/8" (3.2 mm) ID Tubing, (May be used with separate rotating lock ring; FSLLR)



PANEL MOUNT ACCESSORIES:

Lock Nuts and Color Coding Rings

CCR (-1 -2 -3 -4 -5 -06 -07)

Color Coded Ring (For use with PMS, PMK, FTLLB or FTLB panel mount fittings)



CCLR (-1 -2 -3 -4 -5 -06 -07)

Color Coded Lock Ring (For use with PMS, PMK, FTLLB or FTLB panel mount fittings)



PMKN (-XO)

Panel Mount Nut 10-32 UNF with 3/8" Hex (For use with PMK panel mount fittings)



PMSN (-XO)

Panel Mount Nut 1/4-28 UNF with 7/16" Hex (For use with PMS, FTLLB, and FTLLB panel mount fittings)





LNS (-1 -2 -3 -4 -5 -06 -07 -6005 -81 -J1A)

1/4-28 UNF Panel Mount Lock Nut (For use with PMS, FTLLB or FTLB panel mount fittings)



Select the Correct Lock Nut for Your Panel Mount:

PART	Dimension A	Use Lock Nut
FTLLB Series	.255"	PMSN or LNS
FTLB Series	.255"	PMSN or LNS
FTLLBMLRL	.255"	PMSN or LNS
MLRLB210	.193"	PMKN
MLRLB220	.193"	PMKN
MLRLB230	.255"	PMSN or LNS

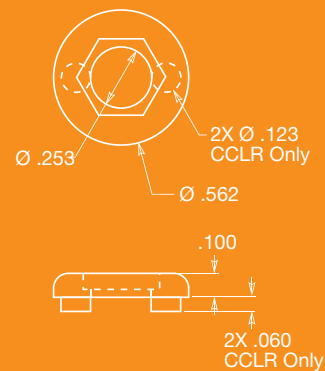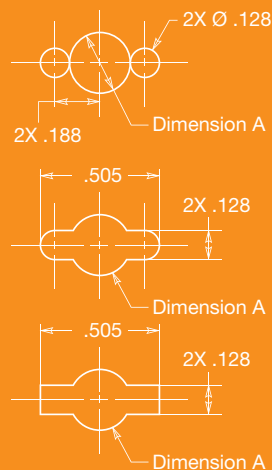

TUBE

Fittings

Straight Throughs
PG-23



Straight Through Reducers
PG-24



Elbows
PG-27



Elbow Reducers
PG-28



Tees
PG-29



Tee Reducers
PG-30



Ys
PG-32



Y Reducers
PG-34



Multiports
PG-34



Multiport Reducers
PG-35



Press-In Plugs
PG-36



Ultrasonically Weldable Barbs
PG-37



Panel Mounts & Accessories
PG-37



Nordson MEDICAL offers barbed tube fittings and tubing connectors in sizes ranging from 1/16" (1.6 mm) to 1" (25.4 mm) tubing ID in over 100 separate configurations in nylon, polypropylene, polycarbonate and PVDF, with other engineering resins available on special order. Available barbed fitting styles include straight tube connectors, elbow fittings, tee fittings, crosses, Y (wye) connectors, multi-port connectors, and multiple straight, elbow and tee reducer combinations.

Several barb styles are offered to provide designers the widest range of assembly forces, pull-off resistance and blow-off resistance – without the need for clamps. Review the descriptions for each of the styles to determine what's best for your application.

Small Bore Barb Styles

Tube ID	Classic	200 Series	400 Series	500 Series
1/16" (1.6 mm)	X	X	X	X
3/32" (2.4 mm)	X	X	X	X
1/8" (3.2 mm)	X	X	X	X
5/32" (4.0 mm)	X	X	X	X
3/16" (4.8 mm)	X	X	X	X
1/4" (6.4 mm)	X	X	X	X
5/16" (8.0 mm)				X

Large Bore Barb Styles

Tube ID	500 Series	600 Series	900 Series
3/8" (9.5 mm)	X	X	X
1/2" (12.7 mm)	X	X	X
5/8" (16.0 mm)	X	X	
3/4" (19.0 mm)	X	X	
1" (25.4 mm)	X	X	

AA (-1 -2 -6005 -9 -J1A)

Straight Through Connector with Classic Series Barbs, 1/16" (1.6 mm) ID Tubing



N210 (-1 -2 -6005 -9 -J1A)

Straight Through Connector with 200 Series Barbs, 1/16" (1.6 mm) ID Tubing



N410 (-1 -2 -6005 -J1A)

Straight Through Connector with 400 Series Barbs, 1/16" (1.6 mm) ID Tubing



N004 (-1 -2 -6005 -9 -J1A)

Straight Through Connector with 500 Series Barbs, 1/16" (1.6 mm) ID Tubing



DD (-1 -2 -6005 -9 -J1A)

Straight Through Connector with Classic Series Barbs, 3/32" (2.4 mm) ID Tubing



N220 (-1 -2 -6005 -9 -9002 -J1A)

Straight Through Connector with 200 Series Barbs, 3/32" (2.4 mm) ID Tubing



N420 (-1 -2 -6005 -J1A)

Straight Through Connector with 400 Series Barbs, 3/32" (2.4 mm) ID Tubing



N007 (-1 -2 -6005 -9 -J1A)

Straight Through Connector with 500 Series Barbs, 3/32" (2.4 mm) ID Tubing



CC (-1 -2 -6005 -9 -J1A)

Straight Through Connector with Classic Series Barbs, 1/8" (3.2 mm) ID Tubing



N230 (-1 -2 -3 -5 -6005 -9 -J1A)

Straight Through Connector with 200 Series Barbs, 1/8" (3.2 mm) ID Tubing



N430 (-1 -2 -6005 -J1A)

Straight Through Connector with 400 Series Barbs, 1/8" (3.2 mm) ID Tubing



N013 (-1 -2 -6005 -9 -J1A -9002)

Straight Through Connector with 500 Series Barbs, 1/8" (3.2 mm) ID Tubing



4040 (-1 -2 -6005 -9 -J1A)

Straight Through Connector with Classic Series Barbs, 5/32" (4.0 mm) ID Tubing



N440 (-1 -2 -5 -6005 -9 -J1A)

Straight Through Connector with 400 Series Barbs, 5/32" (4.0 mm) ID Tubing



N025 (-1 -2 -6005 -9 -J1A)

Straight Through Connector with 500 Series Barbs, 5/32" (4.0 mm) ID Tubing



5050 (-1 -2 -6005 -9 -9002 -J1A)

Straight Through Connector with Classic Series Barbs, 3/16" (4.8 mm) ID Tubing



N035 (-1 -2 -6005 -9 -J1A)

Straight Through Connector with 500 Series Barbs, 3/16" (4.8 mm) ID Tubing



6060 (-1 -2 -6005 -9 -9002 -J1A)

Straight Through Connector with Classic Series Barbs, 1/4" (6.4 mm) ID Tubing



N055 (-1 -2 -6005 -9 -J1A)

Straight Through Connector with 500 Series Barbs, 1/4" (6.4 mm) ID Tubing



N065 (-1 -2 -6005 -9 -J1A)

Straight Through Connector with 500 Series Barbs, 5/16" (8 mm) ID Tubing



LARGE TUBE FITTINGS: Straight Through



Reference part shown: N670

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

N070 (-0 -6005 -9)

Straight Through Connector with 500 Series Barbs, 3/8" (9.5 mm) ID Tubing

N670 (-6005 -9 -J1A)

Straight Through Connector with 600 Series Barbs, 3/8" (9.5 mm) ID Tubing

N970 (-9)

Straight Through Connector with 900 Series Barbs, 3/8" (9.5 mm) ID Tubing

N080 (-0 -6005 -9)

Straight Through Connector with 500 Series Barbs, 1/2" (12.7 mm) ID Tubing

N680 (-6005 -9 -J1A)

Straight Through Connector with 600 Series Barbs, 1/2" (12.7 mm) ID Tubing

N975 (-9)

Straight Through Connector with 900 Series Barbs, 1/2" (12.7 mm) ID Tubing

N090 (-0 -6005 -9)

Straight Through Connector with 500 Series Barbs, 5/8" (16.0 mm) ID Tubing

N690 (-0 -6005 -9)

Straight Through Connector with 600 Series Barbs, 5/8" (16.0 mm) ID Tubing

N100 (-0 -6005)

Straight Through Connector with 500 Series Barbs, 3/4" (19.0 mm) ID Tubing

N6100 (-6005)

Straight Through Connector with 600 Series Barbs, 3/4" (19.0 mm) ID Tubing

N110 (-0 -6005)

Straight Through Connector with 500 Series Barbs, 1" (25.4 mm) ID Tubing

N6110 (-6005)

Straight Through Connector with 600 Series Barbs, 1" (25.4 mm) ID Tubing

TUBE FITTINGS: Straight Through Reducer

AD (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



N220/210 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 200 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



N420/410 (-1 -2 -6005 -10 -J1A)

Straight Through Reduction Connector with 400 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



N007/004 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



AC (-0 -1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing



N230/210 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 200 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing



N430/410 (-1 -2 -6005 -J1A)

Straight Through Reduction Connector with 400 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing



N013/004 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing



DC (-1 -2 -6005 -9)

Straight Through Reduction Connector with Classic Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing



N230/220 (-1 -2 -6005 -07 -9 -J1A)

Straight Through Reduction Connector with 200 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing



N430/420 (-1 -2 -6005 -6012 -J1A)

Straight Through Reduction Connector with 400 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing



N013/007 (-1 -2 -6005 -6012 -9 -9002 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing



1040 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 5/32" (4.0 mm) and 1/16" (1.6 mm) ID Tubing



N025/004 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 5/32" (4.0 mm) and 1/16" (1.6 mm) ID Tubing



2040 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 5/32" (4.0 mm) and 3/32" (2.4 mm) ID Tubing



N025/007 (-1 -2 -6005 -6012 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 5/32" (4.0 mm) and 3/32" (2.4 mm) ID Tubing



3040 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 5/32" (4.0 mm) and 1/8" (3.2 mm) ID Tubing



N025/013 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 5/32" (4.0 mm) and 1/8" (3.2 mm) ID Tubing



1050 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 3/16" (4.8 mm) and 1/16" (1.6 mm) ID Tubing



2050 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 3/16" (4.8 mm) and 3/32" (2.4 mm) ID Tubing



N035/007 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 3/16" (4.8 mm) and 3/32" (2.4 mm) ID Tubing



3050 (-1 -2 -5 -6005 -81 -9 -9002 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 3/16" (4.8 mm) and 1/8" (3.2 mm) ID Tubing



N035/013 (-1 -2 -6005 -81 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 3/16" (4.8 mm) and 1/8" (3.2 mm) ID Tubing



4050 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 3/16" (4.8 mm) and 5/32" (4.0 mm) ID Tubing



N035/025 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 3/16" (4.8 mm) and 5/32" (4.0 mm) ID Tubing



3060 (-1 -2 -6005 -9 -J1A -J1A2)

Straight Through Reduction Connector with Classic Series Barbs, 1/4" (6.4 mm) and 1/8" (3.2 mm) ID Tubing



N055/013 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 1/4" (6.4 mm) and 1/8" (3.2 mm) ID Tubing



4060 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with Classic Series Barbs, 1/4" (6.4 mm) and 5/32" (4.0 mm) ID Tubing



N055/025 (-1 -2 -6005 -9 -J1A)

Straight Through Reduction Connector with 500 Series Barbs, 1/4" (6.4 mm) and 5/32" (4.0 mm) ID Tubing



TUBE FITTINGS:
Straight Through Reducer

5060 (-1 -2 -6005 -9 -9002 -J1A)
Straight Through Reduction Connector with Classic Series Barbs,
1/4" (6.4 mm) and 3/16" (4.8 mm) ID Tubing



N055/035 (-1 -2 -6005 -9 -J1A)
Straight Through Reduction Connector with 500 Series Barbs,
1/4" (6.4 mm) and 3/16" (4.8 mm) ID Tubing



N065/025 (-1 -2 -6005 -9 -J1A)
Straight Through Reduction Connector with 500 Series Barbs,
5/16" (8.0 mm) and 5/32" (4.0 mm) ID Tubing



N065/035 (-1 -2 -6005 -9 -J1A)
Straight Through Reduction Connector with 500 Series Barbs,
5/16" (8.0 mm) and 3/16" (4.8 mm) ID Tubing




N065/055 (-1 -2 -6005 -9 -J1A)
Straight Through Reduction Connector with 500 Series Barbs,
5/16" (8.0 mm) and 1/4" (6.4 mm) ID Tubing



Did You Know?

Nordson MEDICAL has an interactive PDF available online
at www.nordsonmedical.com. Access part drawings,
reference charts, data sheets, and more from just one location!

LARGE TUBE FITTINGS:
Straight Through Reducer

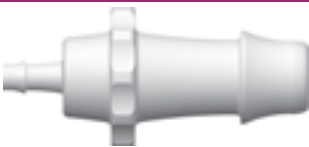

Reference part shown: N070/013

Due to size we are not able to display all large tube fittings.
Please see www.nordsonmedical.com for part reference.
Minimum package quantities vary by tube IDs, please verify
minimum package sizes when placing orders.

N070/013 (-0 -6005 -9 -9002 -J1A)
Straight Through Reduction Connector with 500 Series Barbs, 3/8" (9.5 mm)
and 1/8" (3.2 mm) ID Tubing

N070/055 (-0 -6005 -9)
Straight Through Reduction Connector with 500 Series Barbs, 3/8" (9.5 mm)
and 1/4" (6.4 mm) ID Tubing

N670/655 (-6005 -9 -J1A)
Straight Through Reduction Connector with 600 Series Barbs, 3/8" (9.5 mm)
and 1/4" (6.4 mm) ID Tubing

N080/055 (-0 -6005 -9)
Straight Through Reduction Connector with 500 Series Barbs, 1/2" (12.7 mm)
and 1/4" (6.4 mm) ID Tubing

N680/655 (-6005 -9 -J1A)
Straight Through Reduction Connector with 600 Series Barbs, 1/2" (12.7 mm)
and 1/4" (6.4 mm) ID Tubing

N080/070 (-0 -6005 -9)
Straight Through Reduction Connector with 500 Series Barbs, 1/2" (12.7 mm)
and 3/8" (9.5 mm) ID Tubing

N680/670 (-6005 -9)
Straight Through Reduction Connector with 600 Series Barbs, 1/2" (12.7 mm)
and 3/8" (9.5 mm) ID Tubing

N090/055 (-0 -6005)
Straight Through Reduction Connector with 500 Series Barbs, 5/8" (16.0 mm)
and 1/4" (6.4 mm) ID Tubing

N690/655 (-6005 -9)
Straight Through Reduction Connector with 600 Series Barbs, 5/8" (16.0 mm)
and 1/4" (6.4 mm) ID Tubing

N090/070 (-0 -6005)
Straight Through Reduction Connector with 500 Series Barbs, 5/8" (16.0 mm)
and 3/8" (9.5 mm) ID Tubing

N690/670 (-6005 -9)
Straight Through Reduction Connector with 600 Series Barbs, 5/8" (16.0 mm)
and 3/8" (9.5 mm) ID Tubing

N6100/670 (-6005)
Straight Through Reduction Connector with 600 Series Barbs, 3/4" (19 mm)
and 3/8" (9.5 mm) ID Tubing

N100/080 (-0 -6005)
Straight Through Reduction Connector with 500 Series Barbs, 3/4" (19.0 mm)
and 1/2" (12.7 mm) ID Tubing

N6100/680 (-6005)
Straight Through Reduction Connector with 600 Series Barbs, 3/4" (19 mm)
and 1/2" (12.7 mm) ID Tubing

N6110/680 (-6005)
Straight Through Reduction Connector with 600 Series Barbs, 1" (25.4 mm)
and 1/2" (12.7 mm) ID Tubing

N110/100 (-0 -6005)
Straight Through Reduction Connector with 500 Series Barbs, 1" (25.4 mm)
and 3/4" (19.0 mm) ID Tubing

N6110/6100 (-6005 -9)
Straight Through Reduction Connector with 600 Series Barbs, 1" (25.4 mm)
and 3/4" (19.0 mm) ID Tubing

L10 (-1 -2 -6005 -J1A)

Elbow Connector with Classic Series Barbs, 1/16" (1.6 mm) ID Tubing



L210 (-1 -2 -6005 -9 -J1A)

Elbow Connector with 200 Series Barbs, 1/16" (1.6 mm) ID Tubing



L410 (-1 -2 -6005 -J1A)

Elbow Connector with 400 Series Barbs, 1/16" (1.6 mm) ID Tubing



L20 (-1 -2 -6005 -J1A)

Elbow Connector with Classic Series Barbs, 3/32" (2.4 mm) ID Tubing



L220 (-1 -2 -6005 -6012 -9 -9002 -J1A)

Elbow Connector with 200 Series Barbs, 3/32" (2.4 mm) ID Tubing



L420 (-1 -2 -6005 -J1A)

Elbow Connector with 400 Series Barbs, 3/32" (2.4 mm) ID Tubing



L30 (-1 -2 -6005 -J1A)

Elbow Connector with Classic Series Barbs, 1/8" (3.2 mm) ID Tubing



L230 (-1 -2 -6005 -6012 -9 -J1A)

Elbow Connector with 200 Series Barbs, 1/8" (3.2 mm) ID Tubing



L430 (-1 -2 -6005 -J1A)

Elbow Connector with 400 Series Barbs, 1/8" (3.2 mm) ID Tubing



L40 (-1 -2 -6005 -J1A)

Elbow Connector with Classic Series Barbs, 5/32" (4.0 mm) ID Tubing



L440 (-2)

Elbow Connector with 400 Series Barbs, 5/32" (4.0 mm) ID Tubing



L50 (-1 -2 -6005 -J1A)

Elbow Connector with Classic Series Barbs, 3/16" (4.8 mm) ID Tubing



L60 (-1 -2 -6005 -J1A)

Elbow Connector with Classic Series Barbs, 1/4" (6.4 mm) ID Tubing



Need More Details?

Part drawings, minimum order quantities, and CAD requests can be found online. Visit www.nordsonmedical.com or contact Sales at FMF-sales@nordsonmedical.com.

LARGE TUBE FITTINGS:
Elbow



Reference part shown: L070

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

L070 (-0 -6005)
Elbow Connector with 500 Series Barbs, 3/8" (9.5 mm) ID Tubing

L970 (-9)
Elbow Connector with 900 Series Barbs, 3/8" (9.5 mm) ID Tubing

L080 (-0 -6005)
Elbow Connector with 500 Series Barbs, 1/2" (12.7 mm) ID Tubing

L975 (-9)
Elbow Connector with 900 Series Barbs, 1/2" (12.7 mm) ID Tubing

L090 (-0 -6005)
Elbow Connector with 500 Series Barbs, 5/8" (16.0 mm) ID Tubing

L100 (-0 -6005)
Elbow Connector with 500 Series Barbs, 3/4" (19.0 mm) ID Tubing

L110 (-0 -6005)
Elbow Connector with 500 Series Barbs, 1" (25.4 mm) ID Tubing

L6110 (-6005)
Elbow Connector with 600 Series Barbs, 1" (25.4 mm) ID Tubing

CF494 (-6005)
Low-profile Elbow Connector with Custom 600 Series Barbs, 1" (25.4 mm) ID Tubing

TUBE FITTINGS:
Elbow Reducer

L220/210 (-1 -2 -6005 -9 -J1A)
Elbow Reduction Connector with 200 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



L420/410 (-1 -2 -6005 -J1A)
Elbow Reduction Connector with 400 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



L230/210 (-1 -2 -6005 -9 -J1A)
Elbow Reduction Connector with 200 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing



L430/410 (-1 -2 -6005 -J1A)
Elbow Reduction Connector with 400 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing



L230/220 (-1 -2 -6005 -9 -J1A)
Elbow Reduction Connector with 200 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing



L430/420 (-1 -2 -6005 -J1A)
Elbow Reduction Connector with 400 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing



L430/220 (-9002)
Elbow Reduction Connector with a 400 Series Barb, 1/8" (3.2 mm) and a 200 Series Barb, 3/32" (2.4 mm) ID Tubing



LARGE TUBE FITTINGS:

Elbow Reducer



Reference part shown: **L070/055**

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

L070/055 (-0 -6005)

Elbow Reduction Connector with 500 Series Barbs, 3/8" (9.5 mm) and 1/4" (6.4 mm) ID Tubing

L080/070 (-0 -6005)

Elbow Reduction Connector with 500 Series Barbs, 1/2" (12.7 mm) and 3/8" (9.5 mm) ID Tubing

L100/080 (-0 -6005)

Elbow Reduction Connector with 500 Series Barbs, 3/4" (19.0 mm) and 1/2" (12.7 mm) ID Tubing

L110/100 (-0 -6005)

Elbow Reduction Connector with 500 Series Barbs, 1" (25.4 mm) and 3/4" (19.0 mm) ID Tubing

TUBE FITTINGS:

Tee

T10 (-1 -2 -6005 -J1A)

Tee Connector with Classic Series Barbs, 1/16" (1.6 mm) ID Tubing



T210 (-1 -2 -6005 -9 -J1A)

Tee Connector with 200 Series Barbs, 1/16" (1.6 mm) ID Tubing



T410 (-1 -2 -6005 -J1A)

Tee Connector with 400 Series Barbs, 1/16" (1.6 mm) ID Tubing



T20 (-1 -2 -6005 -J1A)

Tee Connector with Classic Series Barbs, 3/32" (2.4 mm) ID Tubing



T220 (-1 -2 -6005 -6012 -9 -J1A)

Tee Connector with 200 Series Barbs, 3/32" (2.4 mm) ID Tubing



TUBE FITTINGS:

Tee

T420 (-1 -2 -6005 -J1A)

Tee Connector with 400 Series Barbs, 3/32" (2.4 mm) ID Tubing



T30 (-1 -2 -6005 -9 -J1A)

Tee Connector with Classic Series Barbs, 1/8" (3.2 mm) ID Tubing



T230 (-1 -2 -5 -6005 -6012 -9 -J1A)

Tee Connector with 200 Series Barbs, 1/8" (3.2 mm) ID Tubing



T430 (-1 -2 -6005 -9002 -J1A)

Tee Connector with 400 Series Barbs, 1/8" (3.2 mm) ID Tubing



T013 (-6005)

Tee Connector with 500 Series Barbs, 1/8" (3.2 mm) ID Tubing



T40 (-1 -2 -6005 -J1A)

Tee Connector with Classic Series Barbs, 5/32" (4 mm) ID Tubing



T250 (-1 -2 -6005 -9 -J1A)

Tee Connector with 200 Series Barbs, 3/16" (4.8 mm) ID Tubing



T60 (-1 -2 -6005 -J1A)

Tee Connector with Classic Series Barbs, 1/4" (6.4 mm) ID Tubing



TUBE FITTINGS:
Tee

T055 (-6005)

Tee Connector with 500 Series Barbs, 1/4" (6.4 mm) ID Tubing



LARGE TUBE-TO-TUBE:
Tee



Reference part shown: **T070**

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

T070 (-0 -6005 -9)

Tee Connector with 500 Series Barbs, 3/8" (9.5 mm) ID Tubing

T670 (-1 -6005 -9 -J1A)

Tee Connector with 600 Series Barbs, 3/8" (9.5 mm) ID Tubing

T970 (-9)

Tee Connector with 900 Series Barbs, 3/8" (9.5 mm) ID Tubing

T080 (-0 -6005 -9)

Tee Connector with 500 Series Barbs, 1/2" (12.7 mm) ID Tubing

T680 (-6005 -9 -J1A)

Tee Connector with 600 Series Barbs, 1/2" (12.7 mm) ID Tubing

T975 (-9)

Tee Connector with 900 Series Barbs, 1/2" (12.7 mm) ID Tubing

T690 (-6005 -9)

Tee Connector with 600 Series Barbs, 5/8" (16.0 mm) ID Tubing

T100 (-0 -6005 -9)

Tee Connector with 500 Series Barbs, 3/4" (19.0 mm) ID Tubing

T6100 (-6005)

Tee Connector with 600 Series Barbs, 3/4" (19.0 mm) ID Tubing

T110 (-0 -6005)

Tee Connector with 500 Series Barbs, 1" (25.4 mm) ID Tubing

T6110 (-6005)

Tee Connector with 600 Series Barbs, 1" (25.4 mm) ID Tubing

TUBE FITTINGS:
Tee Reducer

T410/420 (-1 -2 -6005 -J1A)


Tee Reduction Connector with 400 Series Barbs, 1/16" (1.6 mm) and 3/32" (2.4 mm) ID Tubing



← 3/32"

T210/230 (-1 -2 -6005 -9 -J1A)


Tee Reduction Connector with 200 Series Barbs, 1/16" (1.6 mm) and 1/8" (3.2 mm) ID Tubing

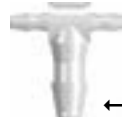

← 1/8"

T410/430 (-1 -2 -6005 -J1A)

Tee Reduction Connector with 400 Series Barbs, 1/16" (1.6 mm) and 1/8" (3.2 mm) ID Tubing

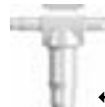

← 1/8"

T220/210 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



← 1/16"

T220R210 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



↑ 1/16"

T420/410 (-1 -2 -6005 -J1A)

Tee Reduction Connector with 400 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



← 1/16"

T420R410 (-1 -2 -6005 -J1A)

Tee Reduction Connector with 400 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



↑ 1/16"

T210/220 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



← 3/32"

T220/230 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 Series Barbs, 3/32" (2.4 mm) and 1/8" (3.2 mm) ID Tubing



← 1/8"

T420/430 (-1 -2 -6005 -J1A)

Tee Reduction Connector with 400 Series Barbs, 3/32" (2.4 mm) and 1/8" (3.2 mm) ID Tubing



T230/210 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing




T230R210 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing



T430/410 (-1 -2 -6005 -J1A)

Tee Reduction Connector with 400 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing

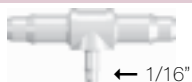


T430R410 (-1 -2 -6005 -J1A)

Tee Reduction Connector with 400 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing



T230/220 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing

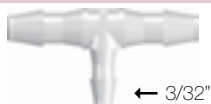


T230R220 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing




T430/420 (-1 -2 -3 -6005 -J1A)

Tee Reduction Connector with 400 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing

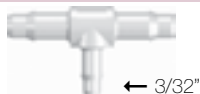


T430R420 (-1 -2 -6005 -J1A)

Tee Reduction Connector with 400 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing

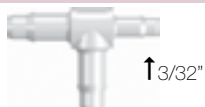


T360/220 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 and 300 Series Barbs, 1/4" (6.4 mm) and 3/32" (2.4 mm) ID Tubing

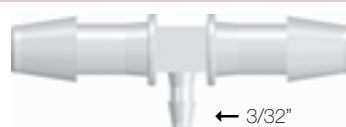


T360/230 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 and 300 Series Barbs, 1/4" (6.4 mm) and 1/8" (3.2 mm) ID Tubing

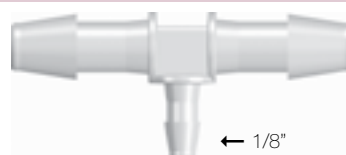

T360/250 (-1 -2 -6005 -9 -J1A)

Tee Reduction Connector with 200 and 300 Series Barbs, 1/4" (6.4 mm) and 3/16" (4.8 mm) ID Tubing

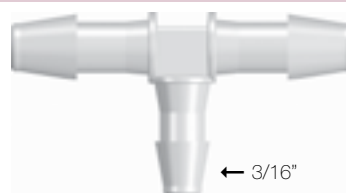

Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

LARGE TUBE FITTINGS:
Tee Reducer



Reference part shown: **T070/055**

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

T070/055 (-0 -6005 -9)

Tee Reduction Connector with 500 Series Barbs, 3/8" (9.5 mm) and 1/4" (6.4 mm) ID Tubing

T080/055 (-0 -6005 -9)

Tee Reduction Connector with 500 Series Barbs, 1/2" (12.7 mm) and 1/4" (6.4 mm) ID Tubing

T080/070 (-0 -6005 -9)

Tee Reduction Connector with 500 Series Barbs, 1/2" (12.7 mm) and 3/8" (9.5 mm) ID Tubing

T100/080 (-0 -6005 -9)

Tee Reduction Connector with 500 Series Barbs, 3/4" (19.0 mm) and 1/2" (12.7 mm) ID Tubing

T110/100 (-0 -6005)

Tee Reduction Connector with 500 Series Barbs, 1" (25.4 mm) and 3/4" (19.0 mm) ID Tubing

T670/655 (-6005)

Tee Reduction Connector with 600 Series Barbs, 3/8" x 2 (9.5 mm) and 1/4" (6.4 mm) ID Tubing

T680/655 (-6005)

Tee Reduction Connector with 600 Series Barbs, 1/2" x 2 (12.7 mm) and 1/4" (6.4 mm) ID Tubing

T680/670 (-6005)

Tee Reduction Connector with 600 Series Barbs, 1/2" x 2 (12.7 mm) and 3/8" (9.5 mm) ID Tubing

T6100/670 (-6005)

Tee Reduction Connector with 600 Series Barbs, 3/4" x 2 (19.0 mm) and 3/8" (9.5 mm) ID Tubing

T6100/680 (-6005)

Tee Reduction Connector with 600 Series Barbs, 3/4" x 2 (19.0 mm) and 1/2" (12.7 mm) ID Tubing

T6110/6100 (-6005)

Tee Reduction Connector with 600 Series Barbs, 1" x 2 (25.4 mm) and 3/4" (19.0 mm) ID Tubing

TUBE FITTINGS:
Y

Y210 (-1 -2 -6005 -9 -J1A)

Y Connector with 200 Series Barbs, 1/16" (1.6 mm) ID Tubing



Y410 (-1 -2 -6005 -J1A)

Y Connector with 400 Series Barbs, 1/16" (1.6 mm) ID Tubing



Y220 (-1 -2 -6005 -6012 -9 -J1A)

Y Connector with 200 Series Barbs, 3/32" (2.4 mm) ID Tubing



Y420 (-1 -2 -6005 -J1A)

Y Connector with 400 Series Barbs, 3/32" (2.4 mm) ID Tubing



Y230 (-1 -2 -3 -6005 -6012 -9 -9002 -9010 -J1A)

Y Connector with 200 Series Barbs, 1/8" (3.2 mm) ID Tubing



Y430 (-1 -2 -6005 -J1A)

Y Connector with 400 Series Barbs, 1/8" (3.2 mm) ID Tubing



Y013 (-6005)

Y Connector with 500 Series Barbs, 1/8" (3.2 mm) ID Tubing



Y240 (-1 -2 -6005 -6012 -9 -J1A)

Y Connector with 200 Series Barbs, 5/32" (4.0 mm) ID Tubing



Y250 (-1 -2 -6005 -9 -9002 -J1A)

Y Connector with 200 Series Barbs, 3/16" (4.8 mm) ID Tubing

**Y450 (-1 -2 -6005 -6012 -J1A)**

Y Connector with 400 Series Barbs, 3/16" (4.8 mm) ID Tubing

**Y260 (-1 -2 -6005 -9 -9002 -J1A)**

Y Connector with 200 Series Barbs, 1/4" (6.4 mm) ID Tubing

**Y460 (-1 -2 -6005 -9 -J1A)**

Y Connector with 400 Series Barbs, 1/4" (6.4 mm) ID Tubing

**Y055 (-6005)**

Y Connector with 500 Series Barbs, 1/4" (6.4 mm) ID Tubing

**LARGE TUBE FITTINGS:**

Y

**Reference part shown: Y070**

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

Y070 (-0 -6005 -9)

Y Connector with 500 Series Barbs, 3/8" (9.5 mm) ID Tubing

Y670 (-6005 -9 -J1A)

Y Connector with 600 Series Barbs, 3/8" (9.5 mm) ID Tubing

Y080 (-0 -6005 -9)

Y Connector with 500 Series Barbs, 1/2" (12.7 mm) ID Tubing

Y680 (-6005 -9 -J1A)

Y Connector with 600 Series Barbs, 1/2" (12.7 mm) ID Tubing

Y090 (-0 -6005 -J1A)

Y Connector with 500 Series Barbs, 5/8" (16.0 mm) ID Tubing

Y690 (-6005)

Y Connector with 600 Series Barbs, 5/8" (16.0 mm) ID Tubing

Y100 (-0 -6005 -J1A)

Y Connector with 500 Series Barbs, 3/4" (19.0 mm) ID Tubing

Y6100 (-6005)

Y Connector with 600 Series Barbs, 3/4" (19.0 mm) ID Tubing

Y110 (-0 -6005)

Y Connector with 500 Series Barbs, 1" (25.4 mm) ID Tubing

Y6110 (-6005)

Y Connector with 600 Series Barbs, 1" (25.4 mm) ID Tubing

LARGE TUBE FITTINGS:
Y Reducer



Reference part shown: **Y655/670**

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

Y630/655 (-6005)

Y Reduction Connector with 600 Series Barbs, 1/8" x 2 (3.2 mm) and 1/4" (6.4 mm) ID Tubing

Y655/670 (-6005)

Y Reduction Connector with 600 Series Barbs, 1/4" x 2 (6.4 mm) and 3/8" (9.5 mm) ID Tubing

Y670/680 (-6005)

Y Reduction Connector with 600 Series Barbs, 3/8" x 2 (9.5 mm) and 1/2" (12.7 mm) ID Tubing

Y6100/6110 (-6005)

Y Reduction Connector with 600 Series Barbs, 3/4" x 2 (19.0 mm) and 1" (25.4 mm) ID Tubing

TUBE FITTINGS:
Multiport

3PF210 (-1 -6005)

Three-Port F Style Manifold with 200 Series Barbs for 1/16" (1.6 mm) ID Tubing



3PF220 (-1 -6005)

Three-Port F Style Manifold with 200 Series Barbs for 3/32" (2.4 mm) ID Tubing



3PF230 (-1 -6005)


Three-Port F Style Manifold with 200 Series Barbs for 1/8" (3.2 mm) ID Tubing



TUBE FITTINGS:
Multiport


4PX210 (-1 -6005)

Four-Port Cross Style Manifold with 200 Series Barbs for 1/16" (1.6 mm) ID Tubing

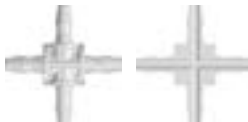


4PX220 (-1 -6005)

Four-Port Cross Style Manifold with 200 Series Barbs for 3/32" (2.4 mm) ID Tubing

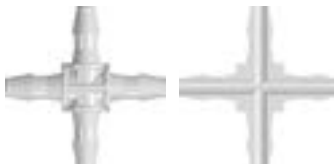


4PX230 (-1 -6005)

Four-Port Cross Style Manifold with 200 Series Barbs for 1/8" (3.2 mm) ID Tubing

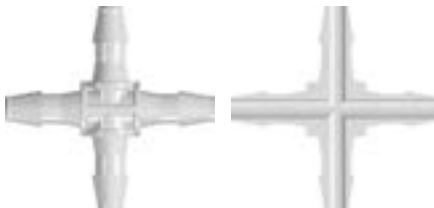

LARGE TUBE FITTINGS:
Multiport

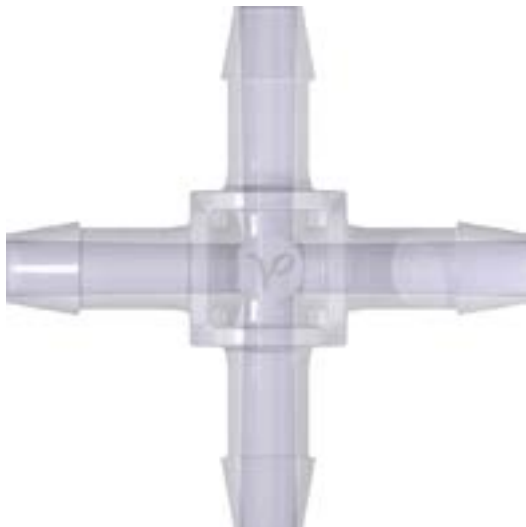

Reference part shown: **4PX670**

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

4PX655 (-6005)

Four-Port Cross Style Manifold with 600 Series Barbs for 1/4" (6.4 mm) ID Tubing

4PX670 (-6005)

Four-Port Cross Style Manifold with 600 Series Barbs for 3/8" (9.5 mm) ID Tubing

4PX680 (-6005)

Four-Port Cross Style Manifold with 600 Series Barbs for 1/2" (12.7 mm) ID Tubing



Reference part shown: 4PX670/655

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

4PX670/655 (-6005)

Four-Port Cross Style Manifold with 600 Series Barbs for 3/8" (9.5 mm) and 1/4" x 3 (6.4 mm) ID Tubing

4PX680/670 (-6005)

Four-Port Cross Style Manifold with 600 Series Barbs for 1/2" (12.7 mm) and 3/8" x 3 (9.5 mm) ID Tubing

TUBE FITTINGS:
Multiport

4PDL210 (-1 -6005)

Four-Port Double Elbow Style Manifold with 200 Series Barbs, 1/16" (1.6 mm) ID Tubing



4PDL220 (-1 -6005)

Four-Port Double Elbow Style Manifold with 200 Series Barbs, 3/32" (2.4 mm) ID Tubing



4PDL230 (-1 -6005)

Four-Port Double Elbow Style Manifold with 200 Series Barbs, 1/8" (3.2 mm) ID Tubing



4PLL210 (-1 -6005)

Four-Port Thru Flow Double Elbow Style Manifold with 200 Series Barbs, 1/16" (1.6 mm) ID Tubing



4PLL220 (-1 -6005 -9)

Four-Port Thru Flow Double Elbow Style Manifold with 200 Series Barbs, 3/32" (2.4 mm) ID Tubing




4PLL230 (-1 -2 -6005)

Four-Port Thru Flow Double Elbow Style Manifold with 200 Series Barbs, 1/8" (3.2 mm) ID Tubing




6PDT210 (-1 -6005)

Six-Port Double T Style Manifold with 200 Series Barbs for 1/16" (1.6 mm) ID Tubing

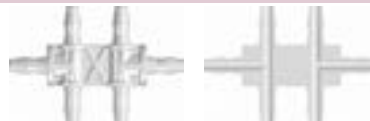


6PDT220 (-1 -6005)

Six-Port Double T Style Manifold with 200 Series Barbs for 3/32" (2.4 mm) ID Tubing

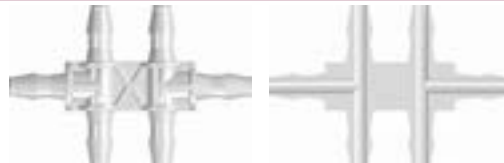

6PDT230 (-1 -6005)

Six-Port Double T Style Manifold with 200 Series Barbs for 1/8" (3.2 mm) ID Tubing

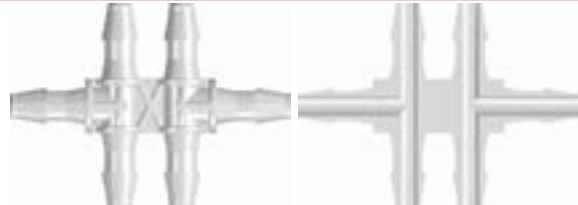


6PTT210 (-1 -6005)

Six-Port Thru Flow Double T Style Manifold with 200 Series Barbs for 1/16" (1.6 mm) ID Tubing




6PTT220 (-1 -6005)

Six-Port Thru Flow Double T Style Manifold with 200 Series Barbs for 3/32" (2.4 mm) ID Tubing

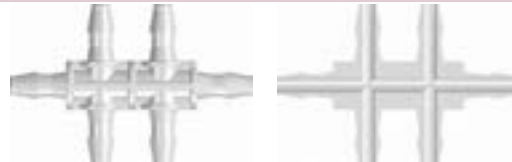

6PTT230 (-1 -6005)

Six-Port Thru Flow Double T Style Manifold with 200 Series Barbs for 1/8" (3.2 mm) ID Tubing

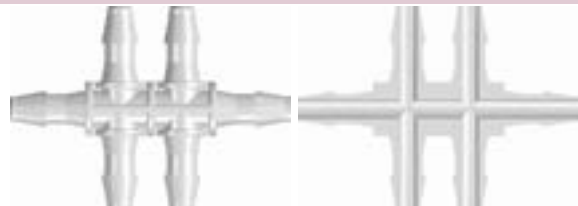

TUBE FITTINGS:
Press-In Plugs

PIP10 (-1 -2 -6005 -J1A)

Press-In Plug with Classic Series Barb, 1/16" (1.6 mm) ID Tubing



PIP210 (-1 -2 -6005 -J1A)

Press-In Plug with 200 Series Barb, 1/16" (1.6 mm) ID Tubing



PIP20 (-1 -2 -6005 -J1A)

Press-In Plug with Classic Series Barb, 3/32" (2.4 mm) ID Tubing



PIP220 (-1 -2 -6005 -6012 -J1A)

Press-In Plug with 200 Series Barb, 3/32" (2.4 mm) ID Tubing



PIP30 (-1 -2 -6005 -J1A)

Press-In Plug with Classic Series Barb, 1/8" (3.2 mm) ID Tubing



PIP230 (-1 -2 -6005 -6012 -J1A)

Press-In Plug with 200 Series Barb, 1/8" (3.2 mm) ID Tubing



PIP40 (-1 -2 -6005 -J1A)

Press-In Plug with Classic Series Barb, 5/32" (4.0 mm) ID Tubing



PIP240 (-1 -2 -6005 -J1A)

Press-In Plug with 200 Series Barb, 5/32" (4.0 mm) ID Tubing



PIP50 (-1 -2 -6005 -J1A)

Press-In Plug with Classic Series Barb, 3/16" (4.8 mm) ID Tubing




PIP60 (-1 -2 -6005 -J1A)

Press-In Plug with Classic Series Barb, 1/4" (6.4 mm) ID Tubing



LARGE TUBE FITTINGS:

Press-In Plugs

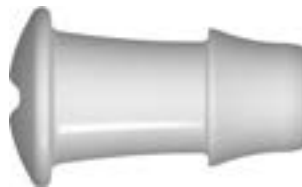

Reference part shown: PIP090

Due to size we are not able to display all large tube fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

PIP070 (-0 -6005)

Press-In Plug with 500 Series Barb, 3/8" (9.5 mm) ID Tubing

PIP670 (-6005)

Press-In Plug with 600 Series Barb, 3/8" (9.5 mm) ID Tubing

PIP080 (-6005)

Press-In Plug with 500 Series Barb, 1/2" (12.7 mm) ID tubing

PIP680 (-6005)

Press-In Plug with 600 Series Barb, 1/2" (12.7 mm) ID Tubing

PIP090 (-6005)

Press-In Plug with 500 Series Barb, 5/8" (16.0 mm) ID Tubing

PIP6100 (-6005)

Press-In Plug with 600 Series Barb, 3/4" (19.0 mm) ID Tubing

PIP6110 (-6005)

Press-In Plug with 600 Series Barb, 1" (25.4 mm) ID Tubing

TUBE FITTINGS: Ultrasonically Weldable Barb

U10 (-1 -2 -7 -9)

Ultrasonically Weldable Fitting to Classic Series Barb, 1/16" (1.6 mm) ID Tubing



U20 (-1 -2 -7)

Ultrasonically Weldable Fitting to Classic Series Barb, 3/32" (2.4 mm) ID Tubing



U30 (-1 -2 -7)

Ultrasonically Weldable Fitting to Classic Series Barb, 1/8" (3.2 mm) ID Tubing



TUBE FITTINGS PANEL MOUNTS: Straight Through and Reducer

PMK210 (-1 -6005)

Panel Mount Connector 10-32 UNF to 200 Series Barbs, 1/16" (1.6 mm) ID Tubing



PMK220 (-1 -6005)

Panel Mount Connector 10-32 UNF to 200 Series Barbs, 3/32" (2.4 mm) ID Tubing



PMS230 (-1 -6005)

Panel Mount Connector 1/4-28 UNF to 200 Series Barbs, 1/8" (3.2 mm) ID Tubing



PMK220-210 (-1 -6005)

Panel Mount Reduction Connector 10-32 UNF to 200 Series Barbs, 3/32" (2.4 mm) and 1/16" (1.6 mm) ID Tubing



PMS230-210 (-1 -6005 -J1A)

Panel Mount Reduction Connector 1/4-28 UNF to 200 Series Barbs, 1/8" (3.2 mm) and 1/16" (1.6 mm) ID Tubing



PMS230-220 (-1 -6005)

Panel Mount Reduction Connector 1/4-28 UNF to 200 Series Barbs, 1/8" (3.2 mm) and 3/32" (2.4 mm) ID Tubing



TUBE FITTINGS PANEL MOUNT ACCESSORIES: Lock Nuts and Color Coding Rings

CCR (-1 -2 -3 -4 -5 -06 -07)

Color Coded Ring (For use with PMS, PMK, FTLB or FTLB panel mount fittings)



CCLR (-1 -2 -3 -4 -5 -06 -07)

Color Coded Lock Ring (For use with PMS, PMK, FTLB or FTLB panel mount fittings)



PMKN (-XO)

Panel Mount Nut 10-32 UNF with 3/8" Hex (For use with PMK panel mount fittings)



PMSN (-XO)

Panel Mount Nut 1/4-28 UNF with 7/16" Hex (For use with PMS, FTLB, and FTLB panel mount fittings)





LNS (-1 -2 -3 -4 -5 -06 -6005 -07 -81 -J1A)

1/4-28 UNF Panel Mount Lock Nut (For use with PMS, FTLB or FTLB panel mount fittings)



Select the Correct Lock Nut for Your Panel Mount:

PART	Dimension A	Use Lock Nut
PMK Series	.193"	PMKN
PMS Series	.255"	PMSN or LNS

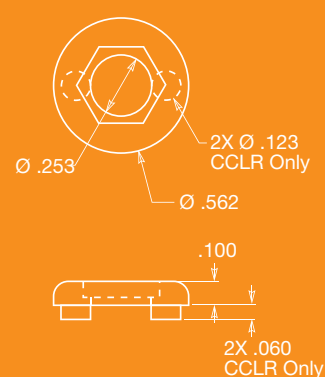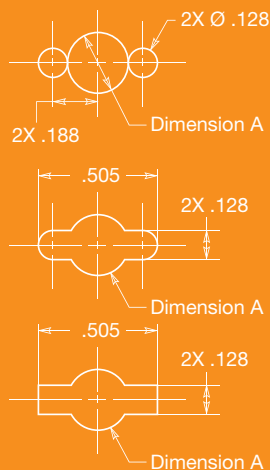

ISO 80369

Compliant Connectors

ISO 80369-3 (Enteral Feeding) Compliant Connectors PG-39



ISO 80369-5 (Blood Pressure) Compliant Connectors PG-39



What is ISO 80369

Hospitals and other healthcare facilities depend on a variety of catheters, tubing and syringes to deliver medications and other substances to patients through vascular, enteral, respiratory, epidural and intrathecal delivery systems. These delivery systems frequently employ fittings called luer connectors to link various system components. Unfortunately, because luer connectors are ubiquitous, easy-to-use, and compatible between different delivery systems, end users can inadvertently connect the incorrect systems together. These misconnections can cause medication or other fluids to be delivered through the wrong route, resulting in serious patient injuries and even death. In an effort to reduce the occurrence of these misconnections, the FDA is in collaboration with the International Organization for Standardization (ISO) and the Association for the Advancement of Medical Instrumentation (AAMI) to develop standards for the non-interconnectability of small bore medical connectors. A joint working group established by the International Organization for Standardization (ISO) and the International Electrotechnical Commission (IEC) leads this effort to develop a series of standards for incompatible connectors.

Nordson MEDICAL is well positioned to help you with your tubing connection needs. Our parts are designed, tested, qualified, and documented to be compliant with these new standards. Because Nordson MEDICAL's Value Plastics product line team has participated on the committees developing these standards, you can be assured of a full understanding of the testing and compliance requirements.

ISO80369-3

Enteral Feeding Connectors

ENTME050 (-8003)

Enteral Male Connector with 5.0 FR Tubing Bond-In Pocket in Natural ABS



ENTME064 (-8003)

Enteral Male Connector with 6.4 FR Tubing Bond-In Pocket in Natural ABS



ENTME078 (-8003)

Enteral Male Connector with 8.0 FR Tubing Bond-In Pocket in Natural ABS



ENTME098 (-8003)

Enteral Male Connector with 10 FR Tubing Bond-In Pocket in Natural ABS



ENTM230 (-6005)

Enteral Male Connector with Barb for 1/8" (3.2 mm) ID Tubing in Animal-Free Polypropylene



ENTFE098 (-8003)

Enteral Female Connector with 10 FR Tubing Bond-In Pocket in Natural ABS



ISO 80369-3 Accessories:

Tethered caps and plugs are not covered under the ISO 80369-3 standard, but are compatible with our enteral male and female connectors.

ENTFTC (-VP1)

Tethered Cap Fits Enteral Male Connector in Animal-Free Polyethylene



ENTMTP (-VP1)

Tethered Plug Fits ENT Series Female Connector in Animal-Free Polyethylene



ISO
80369

ISO80369-5

Blood Pressure Fittings

ISO 80369-5 COMPLIANT PARTS:

Bayonet Male Cuff Connectors

BPM013 (-1)

Blood Pressure Monitor Bayonet Connector to 500 Series Barb, 80369-5 Compliant, 1/8" (3.2 mm) ID Tubing



BPM240 (-1)

Blood Pressure Monitor Bayonet Connector to 200 Series Barb, 80369-5 Compliant, 5/32" (4.0 mm) ID Tubing



BPM025 (-1)

Blood Pressure Monitor Bayonet Connector to 500 Series Barb, 80369-5 Compliant, 5/32" (4.0 mm) ID Tubing



BPM035 (-1)

Blood Pressure Monitor Bayonet Connector to 500 Series Barb, 80369-5 Compliant, 3/16" (4.8 mm) ID Tubing



Extension Sets

ISO 80369-3 (ENFIT) compliant extension sets are also available and can be customized to meet your specific needs. We offer a variety of tubing diameter sizes, lengths, and materials, with or without tethered caps. Please contact Nordson MEDICAL to discuss your requirements.



NEURAXIAL

Connectors

Male Connectors

PG-41



Plugs

PG-41



Accessories

PG-41



NE Series

Neuraxial Connectors

NE Series connectors were designed and tested to be a reliable solution in neuraxial settings. Small and lightweight, these connectors are easy to clean and offer comfort and ease of use. The medical-grade material is clear with no colorants added to simplify regulatory management and also features excellent solvent bonding properties. The male connectors have an easy-to-grip hub, and the available caps and plugs seal off the connectors for irradiation. There are also vented caps and plugs for use with ethylene oxide (EtO) sterilization.

Features & Benefits:

- Small, lightweight, easy to clean
- Male connectors are ISO 80369-6 compliant
- Caps available vented and nonvented
- Easy-to-use bond-in tube pockets
- Clear materials; no colorants added
- Male connectors have easy-to-grip hub

Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!

www.nordsonmedical.com

NEME21-50

Neuraxial Male Connector with 2.1 mm Tubing Bond-In Pocket, Clear Acrylic



NEME21-57

Neuraxial Male Connector with 2.1 mm Tubing Bond-In Pocket, Yellow Acrylic



NEME31-50

Neuraxial Male Connector with 3.1 mm Tubing Bond-In Pocket, Clear Acrylic



NEMP-50

Neuraxial Male Plug, Fits Neuraxial Female Connector, Clear Acrylic



NEMP-57

Neuraxial Male Plug, Fits Neuraxial Female Connector, Yellow Acrylic



NEMP-CM064

Neuraxial Male Plug, Fits Neuraxial Female Connector, Yellow Tinted Polycarbonate



NEFCV-50

Neuraxial Female Vented Cap, Fits Neuraxial Male Connector, Clear Acrylic



NEFCV-CM064

Neuraxial Female Vented Cap, Fits Neuraxial Male Connector, Yellow Tinted Polycarbonate



NEMPV-50

Neuraxial Male Vented Plug, Fits Neuraxial Female Connector, Clear Acrylic



Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

QUICK CONNECT

Couplings



10 Series - 1/16" Flow Size
PG-42



20 Series - 1/8" Flow Size
PG-46



30 Series - 1/8" Flow Size
PG-54



33 Series - 1/8" Flow Size
PG-55



40 Series - 1/4" Flow Size
PG-56



42 Series - 3/32" Flow Size
PG-63



50 Series - 1/4" Flow Size
PG-63



55 Series - 1/8" - 1/4" Flow Size
PG-68



60 Series - 3/8" Flow Size
PG-68



65 Series - 1/2" Flow Size
PG-75



70 Series - 1/4" Flow Size
PG-75



100 Series - 1/16" Flow Size
PG-76



MQC Series
PG-78



RQ Series
PG-79



SBL Series
PG-80



XQ Series
PG-81



Our medical quick connect line is incredible -

Quick connect couplings allow you to quickly and easily connect or disconnect fluid transfer lines, without using tools and without disrupting tubing pressure. Sizes range from 1/16" (1.6 mm) to 1/2" (12.7 mm) ID tubing.

10Series

Quick Connect Couplings

10 SERIES - 1/16" FLOW SIZE: Female Open Flow Straight



Reference part shown: 10AC-S2-02

10AC-S2-01

LinkTech 10AC Series Straight Female, Non-Valved, 1/16" Hose Barb, Natural Acetal

10AC-S2-01BLK

LinkTech 10AC Series Straight Female, Non-Valved, 1/16" Hose Barb, Black, Acetal

10AC-S2-01NS

LinkTech 10AC Series Straight Female, Non-Valved, 1/16" Hose Barb, No Sleeve, Natural Acetal

10AC-S2-01NSBLK

LinkTech 10AC Series Straight Female, Non-Valved, 1/16" Hose Barb, No Sleeve, Black, Acetal

10AC-S2-02

LinkTech 10AC Series Straight Female, Non-Valved, 1/8" Hose Barb, Natural Acetal

10AC-S2-02BLK

LinkTech 10AC Series Straight Female, Non-Valved, 1/8" Hose Barb, Black, Acetal

10AC-S2-02NS

LinkTech 10AC Series Straight Female, Non-Valved, 1/8" Hose Barb, No Sleeve, Natural Acetal

10AC-S2-02NSBLK

LinkTech 10AC Series Straight Female, Non-Valved, 1/8" Hose Barb, No Sleeve, Black, Acetal

10AC-S2-M3

LinkTech 10AC Series Straight Female, Non-Valved, 3 mm Hose Barb, Natural Acetal

10 SERIES - 1/16" FLOW SIZE: Female Open Flow Straight

10AC-S2-M3BLK

LinkTech 10AC Series Straight Female, Non-Valved, 3 mm Hose Barb, Black, Acetal

10AC-S2-M3NS

LinkTech 10AC Series Straight Female, Non-Valved, 3 mm Hose Barb, No Sleeve, Natural Acetal

10AC-S2-M3NSBLK

LinkTech 10AC Series Straight Female, Non-Valved, 3 mm Hose Barb, No Sleeve, Black, Acetal

10PP-S2-01

LinkTech 10PP Series Straight Female, Non-Valved, 1/16" Hose Barb, Polypropylene

10PP-S2-01NS

LinkTech 10PP Series Straight Female, Non-Valved, 1/16" Hose Barb, No Sleeve, Polypropylene

10PP-S2-02

LinkTech 10PP Series Straight Female, Non-Valved, 1/8" Hose Barb, Polypropylene

10PP-S2-02NS

LinkTech 10PP Series Straight Female, Non-Valved, 1/8" Hose Barb, No Sleeve, Polypropylene

10PP-S2-M3

LinkTech 10PP Series Straight Female, Non-Valved, 3 mm Hose Barb, Polypropylene

10PP-S2-M3NS

LinkTech 10PP Series Straight Female, Non-Valved, 3 mm Hose Barb, No Sleeve, Polypropylene

10 SERIES - 1/16" FLOW SIZE: Male Open Flow Straight

10AC-PB2-02BLK

LinkTech 10AC Series Straight Male, Non-Valved, 1/8" Hose Barb, Black, Acetal, Buna-N O-ring

10AC-PB2-02NS

LinkTech 10AC Series Straight Male, Non-Valved, 1/8" Hose Barb, No Sleeve, Natural Acetal, Buna-N O-ring

10AC-PB2-02NSBLK

LinkTech 10AC Series Straight Male, Non-Valved, 1/8" Hose Barb, No Sleeve, Black, Acetal, Buna-N O-ring

10AC-PB2-M3

LinkTech 10AC Series Straight Male, Non-Valved, 3 mm Hose Barb, Natural Acetal, Buna-N O-ring

10AC-PB2-M3BLK

LinkTech 10AC Series Straight Male, Non-Valved, 3 mm Hose Barb, Black, Acetal, Buna-N O-ring

10AC-PB2-M3NS

LinkTech 10AC Series Straight Male, Non-Valved, 3 mm Hose Barb, No Sleeve, Natural Acetal, Buna-N O-ring

10AC-PB2-M3NSBLK

LinkTech 10AC Series Straight Male, Non-Valved, 3 mm Hose Barb, No Sleeve, Black, Acetal, Buna-N O-ring

10PP-PE2-01

LinkTech 10PP Series Straight Male, Non-Valved, 1/16" Hose Barb, Polypropylene, EPDM O-ring

10PP-PE2-01NS

LinkTech 10PP Series Straight Male, Non-Valved, 1/16" Hose Barb, No Sleeve, Polypropylene, EPDM O-ring

10PP-PE2-02

LinkTech 10PP Series Straight Male, Non-Valved, 1/8" Hose Barb, Polypropylene, EPDM O-ring

10PP-PE2-02NS

LinkTech 10PP Series Straight Male, Non-Valved, 1/8" Hose Barb, No Sleeve, Polypropylene, EPDM O-ring

10PP-PE2-M3

LinkTech 10PP Series Straight Male, Non-Valved, 3 mm Hose Barb, Polypropylene, EPDM O-ring

10PP-PE2-M3NS

LinkTech 10PP Series Straight Male, Non-Valved, 3 mm Hose Barb, No Sleeve, Polypropylene, EPDM O-ring

QUICK
CONNECTS

10 SERIES - 1/16" FLOW SIZE: Male Open Flow Straight



Reference part shown: 10AC-PB2-02

10AC-PB2-01

LinkTech 10AC Series Straight Male, Non-Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

10AC-PB2-01BLK

LinkTech 10AC Series Straight Male, Non-Valved, 1/16" Hose Barb, Black, Acetal, Buna-N O-ring

10AC-PB2-01NS

LinkTech 10AC Series Straight Male, Non-Valved, 1/16" Hose Barb, No Sleeve, Natural Acetal, Buna-N O-ring

10AC-PB2-01NSBLK

LinkTech 10AC Series Straight Male, Non-Valved, 1/16" Hose Barb, No Sleeve, Black, Acetal, Buna-N O-ring

10AC-PB2-02

LinkTech 10AC Series Straight Male, Non-Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

Did You Know?

Nordson MEDICAL has an interactive PDF available online at www.nordsonmedical.com. Access part drawings, reference charts, data sheets, and more from just one location!

Need More Details?

Check out the quick connect Datasheets at www.nordsonmedical.com or contact Sales at FMF-sales@nordsonmedical.com.

10 SERIES - 1/16" FLOW SIZE:
Female Open Flow Panel Mount



Reference part shown: 10AC-S3-02

10AC-S3-01
LinkTech 10AC Series Panel Mount Female, Non-Valved, 1/16" Hose Barb, Natural Acetal

10AC-S3-01BLK
LinkTech 10AC Series Panel Mount Female, Non-Valved, 1/16" Hose Barb, Black, Acetal

10AC-S3-02
LinkTech 10AC Series Panel Mount Female, Non-Valved, 1/8" Hose Barb, Natural Acetal

10AC-S3-02BLK
LinkTech 10AC Series Panel Mount Female, Non-Valved, 1/8" Hose Barb, Black, Acetal

10AC-S3-M3
LinkTech 10AC Series Panel Mount Female, Non-Valved, 3 mm Hose Barb, Natural Acetal

10AC-S3-M3BLK
LinkTech 10AC Series Panel Mount Female, Non-Valved, 3 mm Hose Barb, Black, Acetal

10PP-S3-01
LinkTech 10PP Series Panel Mount Female, Non-Valved, 1/16" Hose Barb, Polypropylene

10PP-S3-02
LinkTech 10PP Series Panel Mount Female, Non-Valved, 1/8" Hose Barb, Polypropylene

10PP-S3-M3
LinkTech 10PP Series Panel Mount Female, Non-Valved, 3 mm Hose Barb, Polypropylene

10 SERIES - 1/16" FLOW SIZE:
Male Open Flow Panel Mount



Reference part shown: 10AC-PB3-02

10AC-PB3-01
LinkTech 10AC Series Panel Mount Male, Non-Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

10AC-PB3-01BLK
LinkTech 10AC Series Panel Mount Male, Non-Valved, 1/16" Hose Barb, Black, Acetal, Buna-N O-ring

10AC-PB3-02
LinkTech 10AC Series Panel Mount Male, Non-Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

10AC-PB3-02BLK
LinkTech 10AC Series Panel Mount Male, Non-Valved, 1/8" Hose Barb, Black, Acetal, Buna-N O-ring

10 SERIES - 1/16" FLOW SIZE:
Male Open Flow Panel Mount

10AC-PB3-M3
LinkTech 10AC Series Panel Mount Male, Non-Valved, 3 mm Hose Barb, Natural Acetal, Buna-N O-ring

10AC-PB3-M3BLK
LinkTech 10AC Series Panel Mount Male, Non-Valved, 3 mm Hose Barb, Black, Acetal, Buna-N O-ring

10PP-PE3-01
LinkTech 10PP Series Panel Mount Male, Non-Valved, 1/16" Hose Barb, Polypropylene, EPDM O-ring

10PP-PE3-02
LinkTech 10PP Series Panel Mount Male, Non-Valved, 1/8" Hose Barb, Polypropylene, EPDM O-ring

10 SERIES - 1/16" FLOW SIZE:
Female Valved Straight



Reference part shown: 10ACV-SB2-02

10ACV-SB2-01
LinkTech 10AC Series Straight Female, Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

10ACV-SB2-01NS
LinkTech 10AC Series Straight Female, Valved, 1/16" Hose Barb, No Sleeve, Natural Acetal, Buna-N O-ring

10ACV-SB2-01NSBLK
LinkTech 10AC Series Straight Female, Valved, 1/16" Hose Barb, No Sleeve, Black, Acetal, Buna-N O-ring

10ACV-SB2-02
LinkTech 10AC Series Straight Female, Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

10ACV-SB2-02BLK
LinkTech 10AC Series Straight Female, Valved, 1/8" Hose Barb, Black, Acetal, Buna-N O-ring

10ACV-SB2-02NS
LinkTech 10AC Series Straight Female, Valved, 1/8" Hose Barb, No Sleeve, Natural Acetal, Buna-N O-ring

10ACV-SB2-02NSBLK
LinkTech 10AC Series Straight Female, Valved, 1/8" Hose Barb, No Sleeve, Black, Acetal, Buna-N O-ring

10ACV-SB2-M3
LinkTech 10AC Series Straight Female, Valved, 3 mm Hose Barb, Natural Acetal, Buna-N O-ring

10ACV-SB2-M3BLK
LinkTech 10AC Series Straight Female, Valved, 3 mm Hose Barb, Black, Acetal, Buna-N O-ring

10ACV-SB2-M3NS
LinkTech 10AC Series Straight Female, Valved, 3 mm Hose Barb, No Sleeve, Natural Acetal, Buna-N O-ring

10ACV-SB2-M3NSBLK
LinkTech 10AC Series Straight Female, Valved, 3 mm Hose Barb, No Sleeve, Black, Acetal, Buna-N O-ring



Reference part shown: 10ACV-SB3-02

10ACV-SB3-01

LinkTech 10AC Series Panel Mount Female, Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

10ACV-SB3-01BLK

LinkTech 10AC Series Panel Mount Female, Valved, 1/16" Hose Barb, Black, Acetal, Buna-N O-ring

10ACV-SB3-02

LinkTech 10AC Series Panel Mount Female, Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

10ACV-SB3-02BLK

LinkTech 10AC Series Panel Mount Female, Valved, 1/8" Hose Barb, Black, Acetal, Buna-N O-ring

10ACV-SB3-M3

LinkTech 10AC Series Panel Mount Female, Valved, 3 mm Hose Barb, Natural Acetal, Buna-N O-ring

10ACV-SB3-M3BLK

LinkTech 10AC Series Panel Mount Female, Valved, 3 mm Hose Barb, Black, Acetal, Buna-N O-ring



Reference part shown: 10AC2-02

10AC2-01

LinkTech 10AC Series Straight Coupling, Non-Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

10AC2-01BLK

LinkTech 10AC Series Straight Coupling, Non-Valved, 1/16" Hose Barb, Black, Acetal, Buna-N O-ring

10AC2-02

LinkTech 10AC Series Straight Coupling, Non-Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

10AC2-02BLK

LinkTech 10AC Series Straight Coupling, Non-Valved, 1/8" Hose Barb, Black, Acetal, Buna-N O-ring

10AC2-M3

LinkTech 10AC Series Straight Coupling, Non-Valved, 3 mm Hose Barb, Natural Acetal, Buna-N O-ring

10AC2-M3BLK

LinkTech 10AC Series Straight Coupling, Non-Valved, 3 mm Hose Barb, Black, Acetal, Buna-N O-ring

10ACV2-01

LinkTech 10AC Series Straight Coupling, Valved, 1/16" Hose Barb, Natural Acetal, O-ring

10ACV2-02

LinkTech 10AC Series Straight Coupling, Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

10ACV2-02BLK

LinkTech 10AC Series Straight Coupling, Valved, 1/8" Hose Barb, Black, Acetal, Buna-N O-ring

10ACV2-M3

LinkTech 10AC Series Straight Coupling, Valved, 3 mm Hose Barb, Natural Acetal, Buna-N O-ring

10ACV2-M3BLK

LinkTech 10AC Series Straight Coupling, Valved, 3 mm Hose Barb, Black, Acetal, Buna-N O-ring

10PP2-01

LinkTech 10PP Series Straight Coupling, Non-Valved, 1/16" Hose Barb, Polypropylene, EPDM O-ring

10PP2-02

LinkTech 10PP Series Straight Coupling, Non-Valved, 1/8" Hose Barb, Polypropylene, EPDM O-ring

10PP2-M3

LinkTech 10PP Series Straight Coupling, Non-Valved, 3 mm Hose Barb, Polypropylene, EPDM O-ring

Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

20Series

Quick Connect Couplings

20 SERIES - 1/8" FLOW SIZE:
Female Open Flow Straight



Reference part shown: **20AC-S2-04**

20AC-S1-02
LinkTech 20AC Series Male Thread Female, Non-Valved, 1/8" MNPT, Natural Acetal

20AC-S1-02BSPT
LinkTech 20AC Series Male Thread Female, Non-Valved, 1/8" BSPT, Natural Acetal

20AC-S1-04
LinkTech 20AC Series Male Thread Female, Non-Valved, 1/4" MNPT, Natural Acetal

20AC-S1-04BSPT
LinkTech 20AC Series Male Thread Female, Non-Valved, 1/4" BSPT, Natural Acetal

20AC-S2-01
LinkTech 20AC Series Straight Female, Non-Valved, 1/16" Hose Barb, Natural Acetal

20AC-S2-02
LinkTech 20AC Series Straight Female, Non-Valved, 1/8" Hose Barb, Natural Acetal

20AC-S2-03
LinkTech 20AC Series Straight Female, Non-Valved, 3/16" Hose Barb, Natural Acetal

20AC-S2-04
LinkTech 20AC Series Straight Female, Non-Valved, 1/4" Hose Barb, Natural Acetal

20AC-S9-04
LinkTech 20AC Series Straight Female, Non-Valved, 1/4" PTF, Natural Acetal

20AC-S9-025
LinkTech 20AC Series Straight Female, Non-Valved, 5/32" PTF, Natural Acetal

20CB-S1-02
LinkTech 20CB Series Male Thread Female, Non-Valved, 1/8" MNPT, Chrome Brass

20CB-S1-02BSPT
LinkTech 20CB Series Male Thread Female, Non-Valved, 1/8" BSPT, Chrome Brass

20CB-S1-04
LinkTech 20CB Series Male Thread Female, Non-Valved, 1/4" MNPT, Chrome Brass

20 SERIES - 1/8" FLOW SIZE:
Female Open Flow Straight

20CB-S1-04BSPT
LinkTech 20CB Series Male Thread Female, Non-Valved, 1/4" BSPT, Chrome Brass

20CB-S2-02
LinkTech 20CB Series Straight Female, Non-Valved, 1/8" Hose Barb, Chrome Brass

20CB-S2-03
LinkTech 20CB Series Straight Female, Non-Valved, 3/16" Hose Barb, Chrome Brass

20CB-S2-04
LinkTech 20CB Series Straight Female, Non-Valved, 1/4" Hose Barb, Chrome Brass

20CB-S9-04
LinkTech 20CB Series Straight Female, Non-Valved, 1/4" PTF, Chrome Brass

20CB-S9-025
LinkTech 20CB Series Straight Female, Non-Valved, 5/32" PTF, Chrome Brass

20PP-S1-02
LinkTech 20PP Series Male Thread Female, Non-Valved, 1/8" NPT, Molded Grey Color, Polypropylene

20PP-S1-02MALD
LinkTech 20PP Series Male Thread Female, Non-Valved, 1/8" NPT, Molded Almond Color, Polypropylene

20PP-S1-04
LinkTech 20PP Series Male Thread Female, Non-Valved, 1/4" NPT, Molded Grey Color, Polypropylene

20PP-S1-04MALD
LinkTech 20PP Series Male Thread Female, Non-Valved, 1/4" NPT, Molded Almond Color, Polypropylene

20PP-S2-01
LinkTech 20PP Series Straight Female, Non-Valved, 1/16" Hose Barb, Molded Grey Color, Polypropylene

20PP-S2-01MALD
LinkTech 20PP Series Straight Female, Non-Valved, 1/16" Hose Barb, Molded Almond Color, Polypropylene

20PP-S2-02
LinkTech 20PP Series Straight Female, Non-Valved, 1/8" Hose Barb, Molded Grey Color, Polypropylene

20PP-S2-02MALD
LinkTech 20PP Series Straight Female, Non-Valved, 1/8" Hose Barb, Molded Almond Color, Polypropylene

20PP-S2-04
LinkTech 20PP Series Straight Female, Non-Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene

20PP-S2-04MALD
LinkTech 20PP Series Straight Female, Non-Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene

20PP-S9-04
LinkTech 20PP Series Straight Female, Non-Valved, 1/4" PTF, Molded Grey Color, Polypropylene

20PP-S9-04MALD
LinkTech 20PP Series Straight Female, Non-Valved, 1/4" PTF, Molded Almond Color, Polypropylene

20PP-S9-025
LinkTech 20PP Series Straight Female, Non-Valved, 5/32" PTF, Molded Grey Color, Polypropylene



Reference part shown: 20AC-PB2-04

20AC-PB1-02

LinkTech 20AC Series Male Thread Male, Non-Valved, 1/8" MNPT, Natural Acetal, Buna-N O-ring

20AC-PB1-02BSPT

LinkTech 20AC Series Male Thread Male, Non-Valved, 1/8" BSPT, Natural Acetal, Buna-N O-ring

20AC-PB1-04

LinkTech 20AC Series Male Thread Male, Non-Valved, 1/4" MNPT, Natural Acetal, Buna-N O-ring

20AC-PB1-04BSPT

LinkTech 20AC Series Male Thread Male, Non-Valved, 1/4" BSPT, Natural Acetal, Buna-N O-ring

20AC-PB2-01

LinkTech 20AC Series Straight Male, Non-Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

20AC-PB2-02

LinkTech 20AC Series Straight Male, Non-Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

20AC-PB2-03

LinkTech 20AC Series Straight Male, Non-Valved, 3/16" Hose Barb, Natural Acetal, Buna-N O-ring

20AC-PB2-04

LinkTech 20AC Series Straight Male, Non-Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

20AC-PB9-04

LinkTech 20AC Series Straight Male, Non-Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

20AC-PB9-025

LinkTech 20AC Series Straight Male, Non-Valved, 5/32" PTF, Natural Acetal, Buna-N O-ring

20CB-PB1-02

LinkTech 20CB Series Male Thread Male, Non-Valved, 1/8" MNPT, Chrome Brass, Buna-N O-ring

20CB-PB1-02BSPT

LinkTech 20CB Series Male Thread Male, Non-Valved, 1/8" BSPT, Chrome Brass, Buna-N O-ring

20CB-PB1-04

LinkTech 20CB Series Male Thread Male, Non-Valved, 1/4" MNPT, Chrome Brass, Buna-N O-ring

20CB-PB1-04BSPT

LinkTech 20CB Series Male Thread Male, Non-Valved, 1/4" BSPT, Chrome Brass, Buna-N O-ring

20CB-PB2-02

LinkTech 20CB Series Straight Male, Non-Valved, 1/8" Hose Barb, Chrome Brass, Buna-N O-ring

20CB-PB2-03

LinkTech 20CB Series Straight Male, Non-Valved, 3/16" Hose Barb, Chrome Brass, Buna-N O-ring

20CB-PB2-04

LinkTech 20CB Series Straight Male, Non-Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

20CB-PB9-04

LinkTech 20CB Series Straight Male, Non-Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

20CB-PB9-025

LinkTech 20CB Series Straight Male, Non-Valved, 5/32" PTF, Chrome Brass, Buna-N O-ring

20CB-PB13-02

LinkTech 20CB Series Female Thread Male, Non-Valved, 1/8" FNPT, Chrome Brass, Buna-N O-ring

20CB-PB13-02BSPT

LinkTech 20CB Series Female Thread Male, Non-Valved, 1/8" BSPT, Chrome Brass, Buna-N O-ring

20CB-PB13-04

LinkTech 20CB Series Female Thread Male, Non-Valved, 1/4" FNPT, Chrome Brass, Buna-N O-ring

20CB-PB13-1032

LinkTech 20CB Series Female Thread Male, Non-Valved, 10/32", Chrome Brass, Buna-N O-ring

20PP-PE2-01

LinkTech 20PP Series Straight Male, Non-Valved, 1/16" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PP-PE2-01MALD

LinkTech 20PP Series Straight Male, Non-Valved, 1/16" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PP-PE2-02

LinkTech 20PP Series Straight Male, Non-Valved, 1/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PP-PE2-02MALD

LinkTech 20PP Series Straight Male, Non-Valved, 1/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PP-PE2-04

LinkTech 20PP Series Straight Male, Non-Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PP-PE2-04MALD

LinkTech 20PP Series Straight Male, Non-Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PPX-PE1-02

LinkTech 20PP Series Male Thread Male, Non-Valved, 1/8" NPT, Molded Grey Color, Polypropylene, EPDM O-ring

20PPX-PE1-02MALD

LinkTech 20PP Series Male Thread Male, Non-Valved, 1/8" NPT, Molded Almond Color, Polypropylene, EPDM O-ring

20PPX-PE9-04

LinkTech 20PP Series Straight Male, Non-Valved, 1/4" PTF, Molded Grey Color, Polypropylene, EPDM O-ring

20PPX-PE9-04MALD

LinkTech 20PP Series Straight Male, Non-Valved, 1/4" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

20 SERIES - 1/8" FLOW SIZE:
Female Open Flow Panel Mount



Reference part shown: 20AC-S3-04

20AC-S3-01

LinkTech 20AC Series Panel Mount Female, Non-Valved, 1/16" Hose Barb, Natural Acetal

20AC-S3-02

LinkTech 20AC Series Panel Mount Female, Non-Valved, 1/8" Hose Barb, Natural Acetal

20AC-S3-02GS

LinkTech 20AC Series Panel Mount Gasket Seal Female, Non-Valved, 1/8" Hose Barb, Natural Acetal

20AC-S3-03

LinkTech 20AC Series Panel Mount Female, Non-Valved, 3/16" Hose Barb, Natural Acetal

20AC-S3-03GS

LinkTech 20AC Series Panel Mount Gasket Seal Female, Non-Valved, 3/16" Hose Barb, Natural Acetal

20AC-S3-04

LinkTech 20AC Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Natural Acetal

20AC-S3-04GS

LinkTech 20AC Series Panel Mount Gasket Seal Female, Non-Valved, 1/4" Hose Barb, Natural Acetal

20AC-S8-04

LinkTech 20AC Series Panel Mount Female, Non-Valved, 1/4" PTF, Natural Acetal

20AC-S8-025

LinkTech 20AC Series Panel Mount Female, Non-Valved, 5/32" PTF, Natural Acetal

20AC-S8-025GS

LinkTech 20AC Series Panel Mount Gasket Seal Female, Non-Valved, 5/32" PTF, Natural Acetal

20AC-S19-1032

LinkTech 20AC Series Panel Mount Female, Female Thread, Non-Valved, 10/32", Natural Acetal

20CB-S3-02

LinkTech 20CB Series Panel Mount Female, Non-Valved, 1/8" Hose Barb, Chrome Brass

20CB-S3-03

LinkTech 20CB Series Panel Mount Female, Non-Valved, 3/16" Hose Barb, Chrome Brass

20CB-S3-04

LinkTech 20CB Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Chrome Brass

20CB-S8-04

LinkTech 20CB Series Panel Mount Female, Non-Valved, 1/4" PTF, Chrome Brass

20CB-S8-025

LinkTech 20CB Series Panel Mount Female, Non-Valved, 5/32" PTF, Chrome Brass

20CB-S12-02

LinkTech 20CB Series Panel Mount Female, Male Thread, Non-Valved, 1/8" MNPT, Chrome Brass

20CB-S19-1032

LinkTech 20CB Series Panel Mount Female, Female Thread, Non-Valved, 10/32", Chrome Brass

20PP-S3-01

LinkTech 20PP Series Panel Mount Female, Non-Valved, 1/16" Hose Barb, Molded Grey Color, Polypropylene

20PP-S3-02

LinkTech 20PP Series Panel Mount Female, Non-Valved, 1/8" Hose Barb, Molded Grey Color, Polypropylene

20PP-S3-02MALD

LinkTech 20PP Series Panel Mount Female, Non-Valved, 1/8" Hose Barb, Molded Almond Color, Polypropylene

20PP-S3-04

LinkTech 20PP Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene

20PP-S3-04MALD

LinkTech 20PP Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene

20PP-S8-04

LinkTech 20PP Series Panel Mount Female, Non-Valved, 1/4" PTF, Molded Grey Color, Polypropylene

20PP-S8-04MALD

LinkTech 20PP Series Panel Mount Female, Non-Valved, 1/4" PTF, Molded Almond Color, Polypropylene

20PP-S8-025

LinkTech 20PP Series Panel Mount Female, Non-Valved, 5/32" PTF, Molded Grey Color, Polypropylene

20PP-S8-025MALD

LinkTech 20PP Series Panel Mount Female, Non-Valved, 5/32" PTF, Molded Almond Color, Polypropylene

QUICK
CONNECTS

20 SERIES - 1/8" FLOW SIZE:
Male Open Flow Panel Mount



Reference part shown: 20AC-PB3-04

20AC-PB3-01

LinkTech 20AC Series Panel Mount Male, Non-Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

20AC-PB3-02

LinkTech 20AC Series Panel Mount Male, Non-Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

20 SERIES - 1/8" FLOW SIZE: Male Open Flow Panel Mount

20AC-PB3-03

LinkTech 20AC Series Panel Mount Male, Non-Valved, 3/16" Hose Barb, Natural Acetal, Buna-N O-ring

20AC-PB3-04

LinkTech 20AC Series Panel Mount Male, Non-Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

20AC-PB8-04

LinkTech 20AC Series Panel Mount Male, Non-Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

20CB-PB3-02

LinkTech 20CB Series Panel Mount Male, Non-Valved, 1/8" Hose Barb, Chrome Brass, Buna-N O-ring

20CB-PB3-03

LinkTech 20CB Series Panel Mount Male, Non-Valved, 3/16" Hose Barb, Chrome Brass, Buna-N O-ring

20CB-PB3-04

LinkTech 20CB Series Panel Mount Male, Non-Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

20CB-PB8-04

LinkTech 20CB Series Panel Mount Male, Non-Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

20CB-PB8-025

LinkTech 20CB Series Panel Mount Male, Non-Valved, 5/32" PTF, Chrome Brass, Buna-N O-ring

20 SERIES - 1/8" FLOW SIZE: Male Open Flow Elbow



Reference part shown: 20AC-PB4-04

20AC-PB4-02

LinkTech 20AC Series Elbow Male, Non-Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

20AC-PB4-04

LinkTech 20AC Series Elbow Male, Non-Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

20AC-PB10-04

LinkTech 20AC Series Elbow Male, Non-Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

20AC-PB10-025

LinkTech 20AC Series Elbow Male, Non-Valved, 5/32" PTF, Natural Acetal, Buna-N O-ring

20CB-PB4-02

LinkTech 20CB Series Elbow Male, Non-Valved, 1/8" Hose Barb, Chrome Brass, Buna-N O-ring

20CB-PB4-04

LinkTech 20CB Series Elbow Male, Non-Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

20CB-PB10-04

LinkTech 20CB Series Elbow Male, Non-Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

20PPX-PE4-02

LinkTech 20PP Series Elbow Male, Non-Valved, 1/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PPX-PE4-02MALD

LinkTech 20PP Series Elbow Male, Non-Valved, 1/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PPX-PE4-04

LinkTech 20PP Series Elbow Male, Non-Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PPX-PE4-04MALD

LinkTech 20PP Series Elbow Male, Non-Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PPX-PE10-04

LinkTech 20PP Series Elbow Male, Non-Valved, 1/4" PTF, Molded Grey Color, Polypropylene, EPDM O-ring

20PPX-PE10-04MALD

LinkTech 20PP Series Elbow Male, Non-Valved, 1/4" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

20 SERIES - 1/8" FLOW SIZE: Female Open Flow Elbow



Reference part shown: 20AC-S4-04

20AC-S4-02

LinkTech 20AC Series Elbow Female, Non-Valved, 1/8" Hose Barb, Natural Acetal

20AC-S4-04

LinkTech 20AC Series Elbow Female, Non-Valved, 1/4" Hose Barb, Natural Acetal

20AC-S10-04

LinkTech 20AC Series Elbow Female, Non-Valved, 1/4" PTF, Natural Acetal



Reference part shown: 20ACV-SB2-04

QUICK
CONNECTS

20ACV-SB1-02

LinkTech 20AC Series Male Thread Female, Valved, 1/8" MNPT,
Natural Acetal, Buna-N O-ring

20ACV-SB1-02BSPT

LinkTech 20AC Series Male Thread Female, Valved, 1/8" BSPT,
Natural Acetal, Buna-N O-ring

20ACV-SB1-02UT

LinkTech 20AC Series Male Thread Female, Valved, 1/8" Universal Thread,
Natural Acetal, Buna-N O-ring

20ACV-SB1-04

LinkTech 20AC Series Male Thread Female, Valved, 1/4" MNPT,
Natural Acetal, Buna-N O-ring

20ACV-SB1-04BSPT

LinkTech 20AC Series Male Thread Female, Valved, 1/4" BSPT,
Natural Acetal, Buna-N O-ring

20ACV-SB1-06

LinkTech 20AC Series Male Thread Female, Valved, 3/8" MNPT,
Natural Acetal, Buna-N O-ring

20ACV-SB2-01

LinkTech 20AC Series Straight Female, Valved, 1/16" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-SB2-02

LinkTech 20AC Series Straight Female, Valved, 1/8" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-SB2-03

LinkTech 20AC Series Straight Female, Valved, 3/16" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-SB2-04

LinkTech 20AC Series Straight Female, Valved, 1/4" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-SB9-04

LinkTech 20AC Series Straight Female, Valved, 1/4" PTF,
Natural Acetal, Buna-N O-ring

20ACV-SB9-025

LinkTech 20AC Series Straight Female, Valved, 5/32" PTF,
Natural Acetal, Buna-N O-ring

20CBV-SB1-02

LinkTech 20CB Series Male Thread Female, Valved, 1/8" MNPT,
Chrome Brass, Buna-N O-ring

20CBV-SB1-02BSPT

LinkTech 20CB Series Male Thread Female, Valved, 1/8" BSPT,
Chrome Brass, Buna-N O-ring

20CBV-SB1-04

LinkTech 20CB Series Male Thread Female, Valved, 1/4" MNPT,
Chrome Brass, Buna-N O-ring

20CBV-SB1-04BSPT

LinkTech 20CB Series Male Thread Female, Valved, 1/4" BSPT,
Chrome Brass, Buna-N O-ring

20CBV-SB2-02

LinkTech 20CB Series Straight Female, Valved, 1/8" Hose Barb,
Chrome Brass, Buna-N O-ring

20CBV-SB2-03

LinkTech 20CB Series Straight Female, Valved, 3/16" Hose Barb,
Chrome Brass, Buna-N O-ring

20CBV-SB2-04

LinkTech 20CB Series Straight Female, Valved, 1/4" Hose Barb,
Chrome Brass, Buna-N O-ring

20CBV-SB9-04

LinkTech 20CB Series Straight Female, Valved, 1/4" PTF,
Chrome Brass, Buna-N O-ring

20CBV-SB9-025

LinkTech 20CB Series Straight Female, Valved, 5/32" PTF,
Chrome Brass, Buna-N O-ring

20PPV-SE1-02

LinkTech 20PP Series Male Thread Female, Valved, 1/8" NPT,
Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE1-02MALD

LinkTech 20PP Series Male Thread Female, Valved, 1/8" NPT,
Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE1-04

LinkTech 20PP Series Male Thread Female, Valved, 1/4" NPT,
Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE1-04MALD

LinkTech 20PP Series Male Thread Female, Valved, 1/4" NPT,
Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE2-01

LinkTech 20PP Series Straight Female, Valved, 1/16" Hose Barb,
Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE2-01MALD

LinkTech 20PP Series Straight Female, Valved, 1/16" Hose Barb,
Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE2-02

LinkTech 20PP Series Straight Female, Valved, 1/8" Hose Barb,
Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE2-02MALD

LinkTech 20PP Series Straight Female, Valved, 1/8" Hose Barb,
Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE2-04

LinkTech 20PP Series Straight Female, Valved, 1/4" Hose Barb,
Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE2-04MALD

LinkTech 20PP Series Straight Female, Valved, 1/4" Hose Barb,
Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE9-04

LinkTech 20PP Series Straight Female, Valved, 1/4" PTF,
Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE9-04MALD

LinkTech 20PP Series Straight Female, Valved, 1/4" PTF,
Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE9-025MALD

LinkTech 20PP Series Straight Female, Valved, 5/32" PTF,
Molded Almond Color, Polypropylene, EPDM O-ring



Reference part shown: 20ACV-PB2-04

20ACV-PB1-02

LinkTech 20AC Series Male Thread Male, Valved, 1/8" MNPT, Natural Acetal, Buna-N O-ring

20ACV-PB1-02BSPT

LinkTech 20AC Series Male Thread Male, Valved, 1/8" BSPT, Natural Acetal, Buna-N O-ring

20ACV-PB1-04

LinkTech 20AC Series Male Thread Male, Valved, 1/4" MNPT, Natural Acetal, Buna-N O-ring

20ACV-PB1-04BSPT

LinkTech 20AC Series Male Thread Male, Valved, 1/4" BSPT, Natural Acetal, Buna-N O-ring

20ACV-PB2-01

LinkTech 20AC Series Straight Male, Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-PB2-02

LinkTech 20AC Series Straight Male, Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-PB2-03

LinkTech 20AC Series Straight Male, Valved, 3/16" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-PB2-04

LinkTech 20AC Series Straight Male, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-PB9-04

LinkTech 20AC Series Straight Male, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

20ACV-PB9-025

LinkTech 20AC Series Straight Male, Valved, 5/32" PTF, Natural Acetal, Buna-N O-ring

20CBV-PB1-02

LinkTech 20CB Series Male Thread Male, Valved, 1/8" MNPT, Chrome Brass, Buna-N O-ring

20CBV-PB1-02BSPT

LinkTech 20CB Series Male Thread Male, Valved, 1/8" BSPT, Chrome Brass, Buna-N O-ring

20CBV-PB1-04

LinkTech 20CB Series Male Thread Male, Valved, 1/4" MNPT, Chrome Brass, Buna-N O-ring

20CBV-PB1-04BSPT

LinkTech 20CB Series Male Thread Male, Valved, 1/4" BSPT, Chrome Brass, Buna-N O-ring

20CBV-PB2-02

LinkTech 20CB Series Straight Male, Valved, 1/8" Hose Barb, Chrome Brass, Buna-N O-ring

20CBV-PB2-03

LinkTech 20CB Series Straight Male, Valved, 3/16" Hose Barb, Chrome Brass, Buna-N O-ring

20CBV-PB2-04

LinkTech 20CB Series Straight Male, Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

20CBV-PB9-04

LinkTech 20CB Series Straight Male, Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

20CBV-PB9-025

LinkTech 20CB Series Straight Male, Valved, 5/32" PTF, Chrome Brass, Buna-N O-ring

20CBV-PB13-1032

LinkTech 20CB Series Female Thread Male, Valved, 10/32", Chrome Brass, Buna-N O-ring

20PPV-PE1-02

LinkTech 20PP Series Male Thread Male, Valved, 1/8" NPT, Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-PE1-02MALD

LinkTech 20PP Series Male Thread Male, Valved, 1/8" NPT, Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-PE2-01

LinkTech 20PP Series Straight Male, Valved, 1/16" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-PE2-01MALD

LinkTech 20PP Series Straight Male, Valved, 1/16" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-PE2-02

LinkTech 20PP Series Straight Male, Valved, 1/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-PE2-02MALD

LinkTech 20PP Series Straight Male, Valved, 1/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-PE2-04

LinkTech 20PP Series Straight Male, Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-PE2-04MALD

LinkTech 20PP Series Straight Male, Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-PE9-04

LinkTech 20PP Series Straight Male, Valved, 1/4" PTF, Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-PE9-04MALD

LinkTech 20PP Series Straight Male, Valved, 1/4" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

LinkTech Line Part Number Guide





Reference part shown: 20ACV-SB3-04

QUICK
CONNECTS

20ACV-SB3-01

LinkTech 20AC Series Panel Mount Female, Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-SB3-01GS

LinkTech 20AC Series Panel Mount Gasket Seal Female, Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-SB3-02

LinkTech 20AC Series Panel Mount Female, Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-SB3-02GS

LinkTech 20AC Series Panel Mount Gasket Seal Female, Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-SB3-03

LinkTech 20AC Series Panel Mount Female, Valved, 3/16" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-SB3-03GS

LinkTech 20AC Series Panel Mount Gasket Seal Female, Valved, 3/16" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-SB3-04

LinkTech 20AC Series Panel Mount Female, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-SB3-04GS

LinkTech 20AC Series Panel Mount Gasket Seal Female, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

20ACV-SB8-04

LinkTech 20AC Series Panel Mount Female, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

20ACV-SB8-04GS

LinkTech 20AC Series Panel Mount Gasket Seal Female, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

20ACV-SB8-025

LinkTech 20AC Series Panel Mount Female, Valved, 5/32" PTF, Natural Acetal, Buna-N O-ring

20ACV-SB19-1032

LinkTech 20AC Series Panel Mount Female, Female Thread, Valved, 10/32", Natural Acetal, Buna-N O-ring

20CBV-SB3-02

LinkTech 20CB Series Panel Mount Female, Valved, 1/8" Hose Barb, Chrome Brass, Buna-N O-ring

20CBV-SB3-03

LinkTech 20CB Series Panel Mount Female, Valved, 3/16" Hose Barb, Chrome Brass, Buna-N O-ring

20CBV-SB3-04

LinkTech 20CB Series Panel Mount Female, Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

20CBV-SB8-04

LinkTech 20CB Series Panel Mount Female, Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

20CBV-SB8-025

LinkTech 20CB Series Panel Mount Female, Valved, 5/32" PTF, Chrome Brass, Buna-N O-ring

20CBV-SB12-02

LinkTech 20CB Series Panel Mount Female, Male Thread, Valved, 1/8" MNPT, Chrome Brass, Buna-N O-ring

20CBV-SB19-1032

LinkTech 20CB Series Panel Mount Female, Female Thread, Valved, 10/32", Chrome Brass, Buna-N O-ring

20PPV-SE3-01

LinkTech 20PP Series Panel Mount Female, Valved, 1/16" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE3-01MALD

LinkTech 20PP Series Panel Mount Female, Valved, 1/16" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE3-02

LinkTech 20PP Series Panel Mount Female, Valved, 1/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE3-02MALD

LinkTech 20PP Series Panel Mount Female, Valved, 1/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE3-04

LinkTech 20PP Series Panel Mount Female, Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE3-04MALD

LinkTech 20PP Series Panel Mount Female, Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE8-04

LinkTech 20PP Series Panel Mount Female, Valved, 1/4" PTF, Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-SE8-04MALD

LinkTech 20PP Series Panel Mount Female, Valved, 1/4" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-SE8-025MALD

LinkTech 20PP Series Panel Mount Female, Valved, 5/32" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!

www.nordsonmedical.com



Reference part shown: 20ACV-PB3-04

20ACV-PB3-01

LinkTech 20AC Series Panel Mount Male, Valved, 1/16" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-PB3-02

LinkTech 20AC Series Panel Mount Male, Valved, 1/8" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-PB3-03

LinkTech 20AC Series Panel Mount Male, Valved, 3/16" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-PB3-04

LinkTech 20AC Series Panel Mount Male, Valved, 1/4" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-PB8-04

LinkTech 20AC Series Panel Mount Male, Valved, 1/4" PTF,
Natural Acetal, Buna-N O-ring

20CBV-PB3-02

LinkTech 20CB Series Panel Mount Male, Valved, 1/8" Hose Barb,
Chrome Brass, Buna-N O-ring

20CBV-PB3-03

LinkTech 20CB Series Panel Mount Male, Valved, 3/16" Hose Barb,
Chrome Brass, Buna-N O-ring

20CBV-PB3-04

LinkTech 20CB Series Panel Mount Male, Valved, 1/4" Hose Barb,
Chrome Brass, Buna-N O-ring

20CBV-PB8-04

LinkTech 20CB Series Panel Mount Male, Valved, 1/4" PTF,
Chrome Brass, Buna-N O-ring

20CBV-PB8-025

LinkTech 20CB Series Panel Mount Male, Valved, 5/32" PTF,
Chrome Brass, Buna-N O-ring

20 SERIES - 1/8" FLOW SIZE:
Female Valved Elbow



Reference part shown: 20ACV-SB4-04

20ACV-SB4-02

LinkTech 20AC Series Elbow Female, Valved, 1/8" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-SB4-04

LinkTech 20AC Series Elbow Female, Valved, 1/4" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-SB10-04

LinkTech 20AC Series Elbow Female, Valved, 1/4" PTF,
Natural Acetal, Buna-N O-ring



Reference part shown: 20ACV-PB4-04

20ACV-PB4-02

LinkTech 20AC Series Elbow Male, Valved, 1/8" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-PB4-04

LinkTech 20AC Series Elbow Male, Valved, 1/4" Hose Barb,
Natural Acetal, Buna-N O-ring

20ACV-PB10-04

LinkTech 20AC Series Elbow Male, Valved, 1/4" PTF,
Natural Acetal, Buna-N O-ring

20ACV-PB10-025

LinkTech 20AC Series Elbow Male, Valved, 5/32" PTF,
Natural Acetal, Buna-N O-ring

20CBV-PB4-02

LinkTech 20CB Series Elbow Male, Valved, 1/8" Hose Barb,
Chrome Brass, Buna-N O-ring

20CBV-PB4-04

LinkTech 20CB Series Elbow Male, Valved, 1/4" Hose Barb,
Chrome Brass, Buna-N O-ring

20CBV-PB10-04

LinkTech 20CB Series Elbow Male, Valved, 1/4" PTF,
Chrome Brass, Buna-N O-ring

20PPV-PE4-02

LinkTech 20PP Series Elbow Male, Valved, 1/8" Hose Barb,
Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-PE4-02MALD

LinkTech 20PP Series Elbow Male, Valved, 1/8" Hose Barb,
Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-PE4-04

LinkTech 20PP Series Elbow Male, Valved, 1/4" Hose Barb,
Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-PE4-04MALD

LinkTech 20PP Series Elbow Male, Valved, 1/4" Hose Barb,
Molded Almond Color, Polypropylene, EPDM O-ring

20PPV-PE10-04

LinkTech 20PP Series Elbow Male, Valved, 1/4" PTF,
Molded Grey Color, Polypropylene, EPDM O-ring

20PPV-PE10-04MALD

LinkTech 20PP Series Elbow Male, Valved, 1/4" PTF,
Molded Almond Color, Polypropylene, EPDM O-ring

30Series

Quick Connect Couplings

QUICK
CONNECTS

30 SERIES - 1/8" FLOW SIZE:
Female Open Flow Straight



Reference part shown: 30AC-S2-04

30AC-S2-01
LinkTech 30AC Series Straight Female, Non-Valved, 1/16" Hose Barb, Natural Acetal

30AC-S2-02
LinkTech 30AC Series Straight Female, Non-Valved, 1/8" Hose Barb, Natural Acetal

30AC-S2-03
LinkTech 30AC Series Straight Female, Non-Valved, 3/16" Hose Barb, Natural Acetal

30AC-S2-04
LinkTech 30AC Series Straight Female, Non-Valved, 1/4" Hose Barb, Natural Acetal

30AC-S9-025
LinkTech 30AC Series Straight Female, Non-Valved, 5/32" PTF, Natural Acetal

30 SERIES - 1/8" FLOW SIZE:
Male Open Flow Straight



Reference part shown: 30AC-PB2-02

30AC-PB2-02
LinkTech 30AC Series Straight Male, Non-Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

30ACX-PB1-02
LinkTech 30AC Series Male Thread Male, Non-Valved, 1/8" MNPT, Natural Acetal, Buna-N O-ring

30ACX-PB1-02BSPT
LinkTech 30AC Series Male Thread Male, Non-Valved, 1/8" BSPT, Natural Acetal, Buna-N O-ring

30ACX-PB1-04
LinkTech 30AC Series Male Thread Male, Non-Valved, 1/4" MNPT, Natural Acetal, Buna-N O-ring

30ACX-PB1-04BSPT
LinkTech 30AC Series Male Thread Male, Non-Valved, 1/4" BSPT, Natural Acetal, Buna-N O-ring

30ACX-PB9-04
LinkTech 30AC Series Straight Male, Non-Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring



Reference part shown: 30ACX-PB4-04

30ACX-PB4-02
LinkTech 30AC Series Elbow Male, Non-Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

30ACX-PB4-04
LinkTech 30AC Series Elbow Male, Non-Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

30ACX-PB10-04
LinkTech 30AC Series Elbow Male, Non-Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

30 SERIES - 1/8" FLOW SIZE:
Female Valved Straight



Reference part shown: 30ACV-SB2-04

30ACV-SB2-01
LinkTech 30AC Series Straight Female, Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-SB2-02
LinkTech 30AC Series Straight Female, Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-SB2-03
LinkTech 30AC Series Straight Female, Valved, 3/16" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-SB2-04
LinkTech 30AC Series Straight Female, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-SB9-04
LinkTech 30AC Series Straight Female, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

30ACV-SB9-025
LinkTech 30AC Series Straight Female, Valved, 5/32" PTF, Natural Acetal, Buna-N O-ring



Reference part shown: 30ACV-PB2-02

30ACV-PB1-02

LinkTech 30AC Series Male Thread Male, Valved, 1/8" MNPT, Natural Acetal, Buna-N O-ring

30ACV-PB1-02BSPT

LinkTech 30AC Series Male Thread Male, Valved, 1/8" BSPT, Natural Acetal, Buna-N O-ring

30ACV-PB1-04

LinkTech 30AC Series Male Thread Male, Valved, 1/4" MNPT, Natural Acetal, Buna-N O-ring

30ACV-PB1-04BSPT

LinkTech 30AC Series Male Thread Male, Valved, 1/4" BSPT, Natural Acetal, Buna-N O-ring

30ACV-PB2-01

LinkTech 30AC Series Straight Male, Valved, 1/16" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-PB2-02

LinkTech 30AC Series Straight Male, Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-PB2-03

LinkTech 30AC Series Straight Male, Valved, 3/16" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-PB2-04

LinkTech 30AC Series Straight Male, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-PB9-04

LinkTech 30AC Series Straight Male, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

30ACV-PB9-025

LinkTech 30AC Series Straight Male, Valved, 5/32" PTF, Natural Acetal, Buna-N O-ring



Reference part shown: 30ACV-PB4-04

30ACV-PB4-02

LinkTech 30AC Series Elbow Male, Valved, 1/8" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-PB4-04

LinkTech 30AC Series Elbow Male, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

30ACV-PB10-04

LinkTech 30AC Series Elbow Male, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

33Series

Quick Connect Couplings

33 SERIES - 1/8" FLOW SIZE:
Female Valved Straight



Reference part shown: 33PCV-SPS2-02

33PCV-SPS2-02

LinkTech 33PC Series Straight Female, Valved, Duckbill, 1/8" Hose Barb, Polycarbonate

LinkTech Line Part Number Guide

30	AC	V	-	P	B	1	-	04
Flow Size or Latch Type	Material	Valve	Config.	No O-Ring	Type	Connection Size		

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

40Series

Quick Connect Couplings

40 SERIES - 1/4" FLOW SIZE:
Female Open Flow Straight

QUICK
CONNECTS



Reference part shown: 40AC-S2-04

- 40AC-S1-04**
LinkTech 40AC Series Male Thread Female, Non-Valved, 1/4" MNPT, Natural Acetal
- 40AC-S1-04BSPT**
LinkTech 40AC Series Male Thread Female, Non-Valved, 1/4" BSPT, Natural Acetal
- 40AC-S1-06**
LinkTech 40AC Series Male Thread Female, Non-Valved, 3/8" MNPT, Natural Acetal
- 40AC-S1-06BSPT**
LinkTech 40AC Series Male Thread Female, Non-Valved, 3/8" BSPT, Natural Acetal
- 40AC-S2-04**
LinkTech 40AC Series Straight Female, Non-Valved, 1/4" Hose Barb, Natural Acetal
- 40AC-S2-05**
LinkTech 40AC Series Straight Female, Non-Valved, 5/16" Hose Barb, Natural Acetal
- 40AC-S2-06**
LinkTech 40AC Series Straight Female, Non-Valved, 3/8" Hose Barb, Natural Acetal
- 40AC-S9-04**
LinkTech 40AC Series Straight Female, Non-Valved, 1/4" PTF, Natural Acetal
- 40AC-S9-06**
LinkTech 40AC Series Straight Female, Non-Valved, 3/8" PTF, Natural Acetal
- 40AC-S17-04**
LinkTech 40AC Series Straight Female, Non-Valved, 1/4" OD Push-In, Natural Acetal
- 40CB-S1-04**
LinkTech 40CB Series Male Thread Female, Non-Valved, 1/4" MNPT, Chrome Brass
- 40CB-S1-04BSPT**
LinkTech 40CB Series Male Thread Female, Non-Valved, 1/4" BSPT, Chrome Brass
- 40CB-S1-06**
LinkTech 40CB Series Male Thread Female, Non-Valved, 3/8" MNPT, Chrome Brass

- 40CB-S1-06BSPT**
LinkTech 40CB Series Male Thread Female, Non-Valved, 3/8" BSPT, Chrome Brass
- 40CB-S2-04**
LinkTech 40CB Series Straight Female, Non-Valved, 1/4" Hose Barb, Chrome Brass
- 40CB-S2-06**
LinkTech 40CB Series Straight Female, Non-Valved, 3/8" Hose Barb, Chrome Brass
- 40CB-S9-04**
LinkTech 40CB Series Straight Female, Non-Valved, 1/4" PTF, Chrome Brass
- 40CB-S9-06**
LinkTech 40CB Series Straight Female, Non-Valved, 3/8" PTF, Chrome Brass
- 40CB-S13-04**
LinkTech 40CB Series Straight Female, Non-Valved, Female Thread 1/4" FNPT, Chrome Brass
- 40CB-S3557**
LinkTech 40CB Series Female Connector, Non-Valved, 15/16-27 Male Sink Thread with NBR Gasket, Chrome Brass
- 40CB-S3559**
LinkTech 40CB Series Female Connector, Non-Valved, 55/64-27 Female Sink Thread with NBR Gasket, Chrome Brass
- 40CB-S3562**
LinkTech 40CB Series Female Connector, Non-Valved, M24 X 1 Male Sink Thread with NBR Gasket, Chrome Brass
- 40CB-S3564**
LinkTech 40CB Series Female Connector, Non-Valved, M22 X 1 Female Sink Thread with NBR Gasket, Chrome Brass
- 40PP-S1-04**
LinkTech 40PP Series Male Thread Female, Non-Valved, 1/4" NPT, Molded Grey Color, Polypropylene
- 40PP-S1-04MALD**
LinkTech 40PP Series Male Thread Female, Non-Valved, 1/4" NPT, Molded Almond Color, Polypropylene
- 40PP-S1-06**
LinkTech 40PP Series Male Thread Female, Non-Valved, 3/8" NPT, Molded Grey Color, Polypropylene
- 40PP-S1-06MALD**
LinkTech 40PP Series Male Thread Female, Non-Valved, 3/8" NPT, Molded Almond Color, Polypropylene
- 40PP-S2-04**
LinkTech 40PP Series Straight Female, Non-Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene
- 40PP-S2-04MALD**
LinkTech 40PP Series Straight Female, Non-Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene
- 40PP-S2-06**
LinkTech 40PP Series Straight Female, Non-Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene
- 40PP-S2-06MALD**
LinkTech 40PP Series Straight Female, Non-Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene
- 40PP-S9-06**
LinkTech 40PP Series Straight Female, Non-Valved, 3/8" PTF, Molded Grey Color, Polypropylene
- 40PP-S9-06MALD**
LinkTech 40PP Series Straight Female, Non-Valved, 3/8" PTF, Molded Almond Color, Polypropylene



Reference part shown: 40AC-PB2-04

40AC-PB1-04

LinkTech 40AC Series Male Thread Male, Non-Valved, 1/4" MNPT, Natural Acetal, Buna-N O-ring

40AC-PB1-04BSPT

LinkTech 40AC Series Male Thread Male, Non-Valved, 1/4" BSPT, Natural Acetal, Buna-N O-ring

40AC-PB1-06

LinkTech 40AC Series Male Thread Male, Non-Valved, 3/8" MNPT, Natural Acetal, Buna-N O-ring

40AC-PB1-06BSPT

LinkTech 40AC Series Male Thread Male, Non-Valved, 3/8" BSPT, Natural Acetal, Buna-N O-ring

40AC-PB2-04

LinkTech 40AC Series Straight Male, Non-Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

40AC-PB2-05

LinkTech 40AC Series Straight Male, Non-Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

40AC-PB2-06

LinkTech 40AC Series Straight Male, Non-Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

40AC-PB6

LinkTech 40AC Series End Stop Male, Natural Acetal, Buna-N O-ring

40AC-PB9-04

LinkTech 40AC Series Straight Male, Non-Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

40AC-PB9-06

LinkTech 40AC Series Straight Male, Non-Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

40AC-PB17-04

LinkTech 40AC Series Straight Male, Non-Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring

40CB-PB1-04

LinkTech 40CB Series Male Thread Male, Non-Valved, 1/4" MNPT, Chrome Brass, Buna-N O-ring

40CB-PB1-04BSPT

LinkTech 40CB Series Male Thread Male, Non-Valved, 1/4" BSPT, Chrome Brass, Buna-N O-ring

40CB-PB1-06

LinkTech 40CB Series Male Thread Male, Non-Valved, 3/8" MNPT, Chrome Brass, Buna-N O-ring

40CB-PB1-06BSPT

LinkTech 40CB Series Male Thread Male, Non-Valved, 3/8" BSPT, Chrome Brass, Buna-N O-ring

40CB-PB2-04

LinkTech 40CB Series Straight Male, Non-Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

40CB-PB2-05

LinkTech 40CB Series Straight Male, Non-Valved, 5/16" Hose Barb, Chrome Brass, Buna-N O-ring

40CB-PB2-06

LinkTech 40CB Series Straight Male, Non-Valved, 3/8" Hose Barb, Chrome Brass, Buna-N O-ring

40CB-PB9-04

LinkTech 40CB Series Straight Male, Non-Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

40CB-PB9-06

LinkTech 40CB Series Straight Male, Non-Valved, 3/8" PTF, Chrome Brass, Buna-N O-ring

40CB-PB13-04

LinkTech 40CB Series Female Thread Male, Non-Valved, 1/4" FNPT, Chrome Brass, Buna-N O-ring

40PP-PE2-04

LinkTech 40PP Series Straight Male, Non-Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PP-PE2-04MALD

LinkTech 40PP Series Straight Male, Non-Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40PP-PE2-06

LinkTech 40PP Series Straight Male, Non-Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PP-PE2-06MALD

LinkTech 40PP Series Straight Male, Non-Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40PPX-PE1-04

LinkTech 40PP Series Male Thread Male, Non-Valved, 1/4" NPT, Molded Grey Color, Polypropylene, EPDM O-ring

40PPX-PE1-04MALD

LinkTech 40PP Series Male Thread Male, Non-Valved, 1/4" NPT, Molded Almond Color, Polypropylene, EPDM O-ring

40PPX-PE9-06

LinkTech 40PP Series Straight Male, Non-Valved, 3/8" PTF, Molded Grey Color, Polypropylene, EPDM O-ring

40PPX-PE9-06MALD

LinkTech 40PP Series Straight Male, Non-Valved, 3/8" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

Need More Details?

Check out the quick connect Datasheets at
www.nordsonmedical.com or contact Sales at
FMF-sales@nordsonmedical.com.

40 SERIES - 1/4" FLOW SIZE:
Female Open Flow Panel Mount



Reference part shown: **40AC-S3-06**

- 40AC-S3-04**
LinkTech 40AC Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Natural Acetal
- 40AC-S3-05**
LinkTech 40AC Series Panel Mount Female, Non-Valved, 5/16" Hose Barb, Natural Acetal
- 40AC-S3-06**
LinkTech 40AC Series Panel Mount Female, Non-Valved, 3/8" Hose Barb, Natural Acetal
- 40AC-S8-04**
LinkTech 40AC Series Panel Mount Female, Non-Valved, 1/4" PTF, Natural Acetal
- 40AC-S8-06**
LinkTech 40AC Series Panel Mount Female, Non-Valved, 3/8" PTF, Natural Acetal
- 40AC-S18-04**
LinkTech 40AC Series Panel Mount Female, Non-Valved, 1/4" OD Push-In, Natural Acetal
- 40CB-S3-04**
LinkTech 40CB Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Chrome Brass
- 40CB-S3-05**
LinkTech 40CB Series Panel Mount Female, Non-Valved, 5/16" Hose Barb, Chrome Brass
- 40CB-S3-06**
LinkTech 40CB Series Panel Mount Female, Non-Valved, 3/8" Hose Barb, Chrome Brass
- 40CB-S8-04**
LinkTech 40CB Series Panel Mount Female, Non-Valved, 1/4" PTF, Chrome Brass
- 40CB-S8-06**
LinkTech 40CB Series Panel Mount Female, Non-Valved, 3/8" PTF, Chrome Brass
- 40CB-S12-04**
LinkTech 40CB Series Panel Mount Female, Male Thread, Non-Valved, 1/4" MNPT, Chrome Brass
- 40PP-S3-04**
LinkTech 40PP Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene
- 40PP-S3-04MALD**
LinkTech 40PP Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene
- 40PP-S3-06**
LinkTech 40PP Series Panel Mount Female, Non-Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene
- 40PP-S3-06MALD**
LinkTech 40PP Series Panel Mount Female, Non-Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene
- 40PP-S8-06**
LinkTech 40PP Series Panel Mount Female, Non-Valved, 3/8" PTF, Molded Grey Color, Polypropylene

40 SERIES - 1/4" FLOW SIZE:
Male Open Flow Panel Mount



Reference part shown: **40AC-PB3-04**

- 40AC-PB3-04**
LinkTech 40AC Series Panel Mount Male, Non-Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring
- 40AC-PB3-05**
LinkTech 40AC Series Panel Mount Male, Non-Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring
- 40AC-PB3-06**
LinkTech 40AC Series Panel Mount Male, Non-Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring
- 40AC-PB8-04**
LinkTech 40AC Series Panel Mount Male, Non-Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring
- 40AC-PB8-06**
LinkTech 40AC Series Panel Mount Male, Non-Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring
- 40AC-PB18-04**
LinkTech 40AC Series Panel Mount Male, Non-Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring
- 40CB-PB3-04**
LinkTech 40CB Series Panel Mount Male, Non-Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring
- 40CB-PB3-05**
LinkTech 40CB Series Panel Mount Male, Non-Valved, 5/16" Hose Barb, Chrome Brass, Buna-N O-ring
- 40CB-PB3-06**
LinkTech 40CB Series Panel Mount Male, Non-Valved, 3/8" Hose Barb, Chrome Brass, Buna-N O-ring
- 40CB-PB8-04**
LinkTech 40CB Series Panel Mount Male, Non-Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring
- 40CB-PB8-06**
LinkTech 40CB Series Panel Mount Male, Non-Valved, 3/8" PTF, Chrome Brass, Buna-N O-ring
- 40PP-PE3-04**
LinkTech 40PP Series Panel Mount Male, Non-Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring
- 40PP-PE3-04MALD**
LinkTech 40PP Series Panel Mount Male, Non-Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring
- 40PP-PE3-06**
LinkTech 40PP Series Panel Mount Male, Non-Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring
- 40PP-PE3-06MALD**
LinkTech 40PP Series Panel Mount Male, Non-Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring



Reference part shown: 40AC-PB4-06

40AC-PB4-04

LinkTech 40AC Series Elbow Male, Non-Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

40AC-PB4-06

LinkTech 40AC Series Elbow Male, Non-Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

40AC-PB10-06

LinkTech 40AC Series Elbow Male, Non-Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

40ACX-PB10-04

LinkTech 40AC Series Elbow Male, Non-Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

40CB-PB4-04

LinkTech 40CB Series Elbow Male, Non-Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

40CB-PB4-06

LinkTech 40CB Series Elbow Male, Non-Valved, 3/8" Hose Barb, Chrome Brass, Buna-N O-ring

40CB-PB10-04

LinkTech 40CB Series Elbow Male, Non-Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

40PPX-PE4-04

LinkTech 40PP Series Elbow Male, Non-Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PPX-PE4-04MALD

LinkTech 40PP Series Elbow Male, Non-Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40PPX-PE4-06

LinkTech 40PP Series Elbow Male, Non-Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PPX-PE4-06MALD

LinkTech 40PP Series Elbow Male, Non-Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring



Reference part shown: 40ACV-SB2-04

40ACV-SB1-04

LinkTech 40AC Series Male Thread Female, Valved, 1/4" MNPT, Natural Acetal, Buna-N O-ring

40ACV-SB1-04BSPT

LinkTech 40AC Series Male Thread Female, Valved, 1/4" BSPT, Natural Acetal, Buna-N O-ring

40ACV-SB1-06

LinkTech 40AC Series Male Thread Female, Valved, 3/8" MNPT, Natural Acetal, Buna-N O-ring

40ACV-SB1-06BSPT

LinkTech 40AC Series Male Thread Female, Valved, 3/8" BSPT, Natural Acetal, Buna-N O-ring

40ACV-SB2-04

LinkTech 40AC Series Straight Female, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

40ACV-SB2-05

LinkTech 40AC Series Straight Female, Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

40ACV-SB2-06

LinkTech 40AC Series Straight Female, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

40ACV-SB9-04

LinkTech 40AC Series Straight Female, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

40ACV-SB9-06

LinkTech 40AC Series Straight Female, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

40ACV-SB17-04

LinkTech 40AC Series Straight Female, Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring

40CB-S35271

LinkTech 40CB Series Compact Female Connector, Valved, M24 X 1 Male Sink Thread with NBR Gasket, Chrome Brass

40CBV-SB1-04

LinkTech 40CB Series Male Thread Female, Valved, 1/4" MNPT, Chrome Brass, Buna-N O-ring

40CBV-SB1-04BSPT

LinkTech 40CB Series Male Thread Female, Valved, 1/4" BSPT, Chrome Brass, Buna-N O-ring

40CBV-SB1-06

LinkTech 40CB Series Male Thread Female, Valved, 3/8" MNPT, Chrome Brass, Buna-N O-ring

40CBV-SB1-06BSPT

LinkTech 40CB Series Male Thread Female, Valved, 3/8" BSPT, Chrome Brass, Buna-N O-ring

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!
www.nordsonmedical.com

40 SERIES - 1/4" FLOW SIZE:
Female Valved Straight

40CBV-SB2-04

LinkTech 40CB Series Straight Female, Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

40CBV-SB2-05

LinkTech 40CB Series Straight Female, Valved, 5/16" Hose Barb, Chrome Brass, Buna-N O-ring

40CBV-SB2-06

LinkTech 40CB Series Straight Female, Valved, 3/8" Hose Barb, Chrome Brass, Buna-N O-ring

40CBV-SB9-04

LinkTech 40CB Series Straight Female, Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

40CBV-SB9-06

LinkTech 40CB Series Straight Female, Valved, 3/8" PTF, Chrome Brass, Buna-N O-ring

40CBV-SB13-04

LinkTech 40CB Series Female Thread Female, Valved, 1/4" FNPT, Chrome Brass, Buna-N O-ring

40PPV-SE1-04

LinkTech 40PP Series Male Thread Female, Valved, 1/4" NPT, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-SE1-04MALD

LinkTech 40PP Series Male Thread Female, Valved, 1/4" NPT, Molded Almond Color, Polypropylene, EPDM O-ring

40PPV-SE1-06

LinkTech 40PP Series Male Thread Female, Valved, 3/8" NPT, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-SE1-06MALD

LinkTech 40PP Series Male Thread Female, Valved, 3/8" NPT, Molded Almond Color, Polypropylene, EPDM O-ring

40PPV-SE2-04

LinkTech 40PP Series Straight Female, Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-SE2-04MALD

LinkTech 40PP Series Straight Female, Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40PPV-SE2-06

LinkTech 40PP Series Straight Female, Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-SE2-06MALD

LinkTech 40PP Series Straight Female, Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40PPV-SE9-06

LinkTech 40PP Series Straight Female, Valved, 3/8" PTF, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-SE9-06MALD

LinkTech 40PP Series Straight Female, Valved, 3/8" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

40 SERIES - 1/4" FLOW SIZE:
Male Valved Straight



Reference part shown: 40ACV-PB2-04

40ACV-PB1-04BSPT

LinkTech 40AC Series Male Thread Male, Valved, 1/4" BSPT, Natural Acetal, Buna-N O-ring

40ACV-PB1-06

LinkTech 40AC Series Male Thread Male, Valved, 3/8" MNPT, Natural Acetal, Buna-N O-ring

40ACV-PB1-06BSPT

LinkTech 40AC Series Male Thread Male, Valved, 3/8" BSPT, Natural Acetal, Buna-N O-ring

40ACV-PB2-04

LinkTech 40AC Series Straight Male, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

40ACV-PB2-05

LinkTech 40AC Series Straight Male, Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

40ACV-PB2-06

LinkTech 40AC Series Straight Male, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

40ACV-PB9-04

LinkTech 40AC Series Straight Male, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

40ACV-PB9-06

LinkTech 40AC Series Straight Male, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

40ACV-PB17-04

LinkTech 40AC Series Straight Male, Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring

40CBV-PB1-04

LinkTech 40CB Series Male Thread Male, Valved, 1/4" MNPT, Chrome Brass, Buna-N O-ring

40CBV-PB1-04BSPT

LinkTech 40CB Series Male Thread Male, Valved, 1/4" BSPT, Chrome Brass, Buna-N O-ring

40CBV-PB1-06

LinkTech 40CB Series Male Thread Male, Valved, 3/8" MNPT, Chrome Brass, Buna-N O-ring

40CBV-PB1-06BSPT

LinkTech 40CB Series Male Thread Male, Valved, 3/8" BSPT, Chrome Brass, Buna-N O-ring

40CBV-PB2-04

LinkTech 40CB Series Straight Male, Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

40CBV-PB2-05

LinkTech 40CB Series Straight Male, Valved, 5/16" Hose Barb, Chrome Brass, Buna-N O-ring

40 SERIES - 1/4" FLOW SIZE:

Male Valved Straight

40CBV-PB2-06

LinkTech 40CB Series Straight Male, Valved, 3/8" Hose Barb, Chrome Brass, Buna-N O-ring

40CBV-PB9-04

LinkTech 40CB Series Straight Male, Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

40CBV-PB9-06

LinkTech 40CB Series Straight Male, Valved, 3/8" PTF, Chrome Brass, Buna-N O-ring

40CBV-PB13-04

LinkTech 40CB Series Female Thread Male, Valved, 1/4" FNPT, Chrome Brass, Buna-N O-ring

40PPV-PE1-04

LinkTech 40PP Series Male Thread Male, Valved, 1/4" NPT, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-PE1-04MALD

LinkTech 40PP Series Male Thread Male, Valved, 1/4" NPT, Molded Almond Color, Polypropylene, EPDM O-ring

40PPV-PE2-04

LinkTech 40PP Series Straight Male, Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-PE2-04MALD

LinkTech 40PP Series Straight Male, Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40PPV-PE2-06

LinkTech 40PP Series Straight Male, Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-PE2-06MALD

LinkTech 40PP Series Straight Male, Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40PPV-PE9-06

LinkTech 40PP Series Straight Male, Valved, 3/8" PTF, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-PE9-06MALD

LinkTech 40PP Series Straight Male, Valved, 3/8" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

40 SERIES - 1/4" FLOW SIZE:

Female Valved Panel Mount



Reference part shown: 40ACV-SB3-06

40ACV-SB3-04

LinkTech 40AC Series Panel Mount Female, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

40ACV-SB3-05

LinkTech 40AC Series Panel Mount Female, Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

40 SERIES - 1/4" FLOW SIZE:

Female Valved Panel Mount

40ACV-SB3-06

LinkTech 40AC Series Panel Mount Female, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

40ACV-SB8-04

LinkTech 40AC Series Panel Mount Female, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

40ACV-SB8-06

LinkTech 40AC Series Panel Mount Female, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

40ACV-SB18-04

LinkTech 40AC Series Panel Mount Female, Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring

40CBV-SB3-04

LinkTech 40CB Series Panel Mount Female, Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring

40CBV-SB3-05

LinkTech 40CB Series Panel Mount Female, Valved, 5/16" Hose Barb, Chrome Brass, Buna-N O-ring

40CBV-SB3-06

LinkTech 40CB Series Panel Mount Female, Valved, 3/8" Hose Barb, Chrome Brass, Buna-N O-ring

40CBV-SB8-04

LinkTech 40CB Series Panel Mount Female, Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring

40CBV-SB8-06

LinkTech 40CB Series Panel Mount Female, Valved, 3/8" PTF, Chrome Brass, Buna-N O-ring

40CBV-SB12-04

LinkTech 40CB Series Panel Mount Female, Male Thread, Valved, 1/4" MNPT, Chrome Brass, Buna-N O-ring

40PPV-SE3-04

LinkTech 40PP Series Panel Mount Female, Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-SE3-04MALD

LinkTech 40PP Series Panel Mount Female, Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40PPV-SE3-06

LinkTech 40PP Series Panel Mount Female, Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-SE3-06MALD

LinkTech 40PP Series Panel Mount Female, Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40PPV-SE8-06

LinkTech 40PP Series Panel Mount Female, Valved, 3/8" PTF, Molded Grey Color, Polypropylene, EPDM O-ring

40PPV-SE8-06MALD

LinkTech 40PP Series Panel Mount Female, Valved, 3/8" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

40 SERIES - 1/4" FLOW SIZE:
Male Valved Panel Mount



Reference part shown: 40ACV-PB3-06

- 40ACV-PB3-04**
LinkTech 40AC Series Panel Mount Male, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring
- 40ACV-PB3-05**
LinkTech 40AC Series Panel Mount Male, Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring
- 40ACV-PB3-06**
LinkTech 40AC Series Panel Mount Male, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring
- 40ACV-PB8-04**
LinkTech 40AC Series Panel Mount Male, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring
- 40ACV-PB8-06**
LinkTech 40AC Series Panel Mount Male, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring
- 40ACV-PB18-04**
LinkTech 40AC Series Panel Mount Male, Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring
- 40CBV-PB3-04**
LinkTech 40CB Series Panel Mount Male, Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring
- 40CBV-PB3-05**
LinkTech 40CB Series Panel Mount Male, Valved, 5/16" Hose Barb, Chrome Brass, Buna-N O-ring
- 40CBV-PB3-06**
LinkTech 40CB Series Panel Mount Male, Valved, 3/8" Hose Barb, Chrome Brass, Buna-N O-ring
- 40CBV-PB8-04**
LinkTech 40CB Series Panel Mount Male, Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring
- 40CBV-PB8-06**
LinkTech 40CB Series Panel Mount Male, Valved, 3/8" PTF, Chrome Brass, Buna-N O-ring
- 40PPV-PE3-04**
LinkTech 40PP Series Panel Mount Male, Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring
- 40PPV-PE3-04MALD**
LinkTech 40PP Series Panel Mount Male, Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring
- 40PPV-PE3-06**
LinkTech 40PP Series Panel Mount Male, Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring
- 40PPV-PE3-06MALD**
LinkTech 40PP Series Panel Mount Male, Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring

40 SERIES - 1/4" FLOW SIZE:
Male Valved Elbow



Reference part shown: 40ACV-PB4-06

- 40ACV-PB4-04**
LinkTech 40AC Series Elbow Male, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring
- 40ACV-PB4-06**
LinkTech 40AC Series Elbow Male, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring
- 40ACV-PB10-04**
LinkTech 40AC Series Elbow Male, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring
- 40ACV-PB10-06**
LinkTech 40AC Series Elbow Male, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring
- 40CBV-PB4-04**
LinkTech 40CB Series Elbow Male, Valved, 1/4" Hose Barb, Chrome Brass, Buna-N O-ring
- 40CBV-PB4-06**
LinkTech 40CB Series Elbow Male, Valved, 3/8" Hose Barb, Chrome Brass, Buna-N O-ring
- 40CBV-PB10-04**
LinkTech 40CB Series Elbow Male, Valved, 1/4" PTF, Chrome Brass, Buna-N O-ring
- 40CBV-PB10-06**
LinkTech 40CB Series Elbow Male, Valved, 3/8" PTF, Chrome Brass, Buna-N O-ring
- 40PPV-PE4-04**
LinkTech 40PP Series Elbow Male, Valved, 1/4" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring
- 40PPV-PE4-04MALD**
LinkTech 40PP Series Elbow Male, Valved, 1/4" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring
- 40PPV-PE4-06**
LinkTech 40PP Series Elbow Male, Valved, 3/8" Hose Barb, Molded Grey Color, Polypropylene, EPDM O-ring
- 40PPV-PE4-06MALD**
LinkTech 40PP Series Elbow Male, Valved, 3/8" Hose Barb, Molded Almond Color, Polypropylene, EPDM O-ring
- 40PPV-PE10-06MALD**
LinkTech 40PP Series Elbow Male, Valved, 3/8" PTF, Molded Almond Color, Polypropylene, EPDM O-ring

42Series

Quick Connect Couplings

42 SERIES - 3/32" FLOW SIZE:
Male Open Flow Straight



Reference part shown: 42AB-PB2-02

42AB-PB2-02

LinkTech 42AB Series Straight Male, Non-Valved, 1/8" Hose Barb, ABS, Buna-N O-ring

42 SERIES - 3/32" FLOW SIZE:
Female Open Flow Panel Mount



Reference part shown: 42AB-S3-02

42AB-S3-02

LinkTech 42AB Series Panel Mount Female, Non-Valved, 1/8" Hose Barb, ABS

Need More Details?

Check out the quick connect Datasheets at
www.nordsonmedical.com or contact Sales at
FMF-sales@nordsonmedical.com.

50Series

Quick Connect Couplings

50 SERIES - 1/4" FLOW SIZE:
Female Open Flow Straight



Reference part shown: 50AC-S2-06LB

50AC-S1-04

LinkTech 50AC Series Male Thread Female, Non-Valved, 1/4" MNPT, Natural Acetal

50AC-S1-06

LinkTech 50AC Series Male Thread Female, Non-Valved, 3/8" MNPT, Natural Acetal

50AC-S2-04

LinkTech 50AC Series Straight Female, Non-Valved, 1/4" Hose Barb, Natural Acetal

50AC-S2-04LB

LinkTech 50AC Series Straight Female, Non-Valved, 1/4" Hose Barb, Long Body, Natural Acetal

50AC-S2-05LB

LinkTech 50AC Series Straight Female, Non-Valved, 5/16" Hose Barb, Long Body, Natural Acetal

50AC-S2-06

LinkTech 50AC Series Straight Female, Non-Valved, 3/8" Hose Barb, Natural Acetal

50AC-S2-06LB

LinkTech 50AC Series Straight Female, Non-Valved, 3/8" Hose Barb, Long Body, Natural Acetal

50AC-S9-04LB

LinkTech 50AC Series Straight Female, Non-Valved, 1/4" PTF, Long Body, Natural Acetal

50AC-S9-06LB

LinkTech 50AC Series Straight Female, Non-Valved, 3/8" PTF, Long Body, Natural Acetal

50AC-S17-04

LinkTech 50AC Series Straight Female, Non-Valved, 1/4" OD Push-In, Natural Acetal

50AC-S35107

LinkTech 50AC Series Male Thread Female, Non-Valved, 1/4" MNPT with Latch Guard, Natural Acetal

QUICK
CONNECTS

50 SERIES - 1/4" FLOW SIZE:

Female Open Flow Straight

50BAC-S1-04

LinkTech 50BAC Series Male Thread Female, Non-Valved, 1/4" MNPT, Natural Acetal

50BAC-S2-04

LinkTech 50BAC Series Straight Female, Non-valved, 1/4" Hose Barb, Natural Acetal

50BAC-S2-05LB

LinkTech 50BAC Series Straight Female, Non-Valved, 5/16" Hose Barb, Long Body, Natural Acetal

50BAC-S2-06

LinkTech 50BAC Series Straight Female, Non-Valved, 3/8" Hose Barb, Natural Acetal

50BAC-S2-06LB

LinkTech 50BAC Series Straight Female, Non-Valved, 3/8" Hose Barb, Long Body, Natural Acetal

50BAC-S9-06LB

LinkTech 50BAC Series Straight Female, Non-Valved, 3/8" PTF, Long Body, Natural Acetal

50GP-S1-04

LinkTech 50GP Series Male Thread Female, Non-Valved, 1/4" MNPT, Glass-Filled Polypropylene

50GP-S5-04

LinkTech 50GP Series Straight Female, Non-Valved, 1/4" Hose Barb, Glass-Filled Polypropylene

50GP-S5-06

LinkTech 50GP Series Straight Female, Non-Valved, 3/8" Hose Barb, Glass-Filled Polypropylene

50GP-S5-06LB

LinkTech 50GP Series Straight Female, Non-Valved, 3/8" Hose Barb, Long Body, Glass-Filled Polypropylene

50 SERIES - 1/4" FLOW SIZE:

Male Open Flow Straight



Reference part shown: 50AC-PB2-06

50AC-PB1-06

LinkTech 50AC Series Male Thread Male, Non-Valved, 3/8" MNPT, Natural Acetal, Buna-N O-ring

50AC-PB2-04

LinkTech 50AC Series Straight Male, Non-Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

50AC-PB2-05

LinkTech 50AC Series Straight Male, Non-Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

50AC-PB2-06

LinkTech 50AC Series Straight Male, Non-Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

50AC-PB9-04

LinkTech 50AC Series Straight Male, Non-Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

50AC-PB9-06

LinkTech 50AC Series Straight Male, Non-Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

50AC-PB17-04

LinkTech 50AC Series Straight Male, Non-Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring

50GP-PF1-04

LinkTech 50GP Series Male Thread Male, Non-Valved, 1/4" MNPT, Glass-Filled Polypropylene, Food Grade EPDM O-ring

50GP-PF5-04

LinkTech 50GP Series Straight Male, Non-Valved, 1/4" Hose Barb, Glass-Filled Polypropylene, Food Grade EPDM O-ring

50GP-PF5-06

LinkTech 50GP Series Straight Male, Non-Valved, 3/8" Hose Barb, Glass-Filled Polypropylene, Food Grade EPDM O-ring

50GP-PF6

LinkTech 50GP Series End Stop Male, Glass-Filled Polypropylene, Food Grade EPDM O-ring

50PPX-PE1-04

LinkTech 50PP Series Male Thread Male, Non-Valved, 1/4" MNPT, Polypropylene, EPDM O-ring

50PPX-PE1-06

LinkTech 50PP Series Male Thread Male, Non-Valved, 3/8" MNPT, Polypropylene, EPDM O-ring

50PPX-PE2-04

LinkTech 50PP Series Straight Male, Non-Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring


50PPX-PE2-06

LinkTech 50PP Series Straight Male, Non-Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!
www.nordsonmedical.com

LinkTech Line Part Number Guide

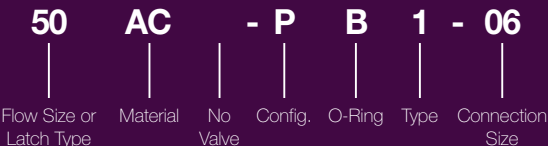



Reference part shown: 50AC-S3-04

50AC-S3-04

LinkTech 50AC Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Natural Acetal

50AC-S3-05

LinkTech 50AC Series Panel Mount Female, Non-Valved, 5/16" Hose Barb, Natural Acetal

50AC-S3-06

LinkTech 50AC Series Panel Mount Female, Non-Valved, 3/8" Hose Barb, Natural Acetal

50AC-S8-04

LinkTech 50AC Series Panel Mount Female, Non-Valved, 1/4" PTF, Natural Acetal

50AC-S8-06

LinkTech 50AC Series Panel Mount Female, Non-Valved, 3/8" PTF, Natural Acetal

50AC-S18-04

LinkTech 50AC Series Panel Mount Female, Non-Valved, 1/4" OD Push-In, Natural Acetal

50BAC-S3-04

LinkTech 50BAC Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Natural Acetal

50BAC-S3-05

LinkTech 50BAC Series Panel Mount Female, Non-Valved, 5/16" Hose Barb, Natural Acetal

50BAC-S8-04

LinkTech 50BAC Series Panel Mount Female, Non-Valved, 1/4" PTF, Natural Acetal

50BAC-S8-06

LinkTech 50BAC Series Panel Mount Female, Non-Valved, 3/8" PTF, Natural Acetal



Reference part shown: 50GP-S7-06

50GP-S7-06

LinkTech 50GP Series Elbow Female, Non-Valved, 3/8" Hose Barb, Glass-Filled Polypropylene

50 SERIES - 1/4" FLOW SIZE:
Male Open Flow Elbow



Reference part shown: 50AC-PB4-06

50AC-PB4-06

LinkTech 50AC Series Elbow Male, Non-Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

50GP-PF7-06

LinkTech 50GP Series Elbow Male, Non-Valved, 3/8" Hose Barb, Glass-Filled Polypropylene, Food Grade EPDM O-ring

50PPX-PE4-04

LinkTech 50PP Series Elbow Male, Non-Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

50PPX-PE4-06

LinkTech 50PP Series Elbow Male, Non-Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

50 SERIES - 1/4" FLOW SIZE:
Female Valved Straight



Reference part shown: 50ACV-SB2-06

50ACV-S35107

LinkTech 50AC Series Male Thread Female, Valved, 1/4" MNPT with Latch Guard, Natural Acetal, Buna-N O-ring

50ACV-SB1-06

LinkTech 50AC Series Male Thread Female, Valved, 3/8" MNPT, Natural Acetal, Buna-N O-ring

50ACV-SB2-04

LinkTech 50AC Series Straight Female, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-SB2-05

LinkTech 50AC Series Straight Female, Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-SB2-06

LinkTech 50AC Series Straight Female, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-SB9-04

LinkTech 50AC Series Straight Female, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

50ACV-SB9-06

LinkTech 50AC Series Straight Female, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

50ACV-SB17-04

LinkTech 50AC Series Straight Female, Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring

50BACV-SB1-04

LinkTech 50BAC Series Male Thread Female, Valved, 1/4" MNPT, Natural Acetal, Buna-N O-ring

50BACV-SB1-06

LinkTech 50BAC Series Male Thread Female, Valved, 3/8" MNPT, Natural Acetal, Buna-N O-ring

50BACV-SB2-04

LinkTech 50BAC Series Straight Female, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

50BACV-SB2-05

LinkTech 50BAC Series Straight Female, Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

50BACV-SB2-06

LinkTech 50BAC Series Straight Female, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

50BACV-SB9-04

LinkTech 50BAC Series Straight Female, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

50GPV-SF5-06

LinkTech 50GP Series Straight Female, Valved, 3/8" Hose Barb, Glass-Filled Polypropylene, Food Grade EPDM O-ring

50PPV-SE2-04

LinkTech 50PP Series Straight Female, Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

50PPV-SE2-06

LinkTech 50PP Series Straight Female, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

50 SERIES - 1/4" FLOW SIZE:
Male Valved Straight



Reference part shown: 50ACV-PB2-06

50ACV-PB1-04

LinkTech 50AC Series Male Thread Male, Valved, 1/4" MNPT, Natural Acetal, Buna-N O-ring

50ACV-PB1-06

LinkTech 50AC Series Male Thread Male, Valved, 3/8" MNPT, Natural Acetal, Buna-N O-ring

50ACV-PB2-04

LinkTech 50AC Series Straight Male, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-PB2-05

LinkTech 50AC Series Straight Male, Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-PB2-06

LinkTech 50AC Series Straight Male, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-PB9-04

LinkTech 50AC Series Straight Male, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

50ACV-PB9-06

LinkTech 50AC Series Straight Male, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

50ACV-PB17-04

LinkTech 50AC Series Straight Male, Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring

50PPV-PE1-04

LinkTech 50PP Series Male Thread Male, Valved, 1/4" MNPT, Polypropylene, EPDM O-ring

50PPV-PE1-06

LinkTech 50PP Series Male Thread Male, Valved, 3/8" MNPT, Polypropylene, EPDM O-ring

50PPV-PE2-04

LinkTech 50PP Series Straight Male, Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

50PPV-PE2-06

LinkTech 50PP Series Straight Male, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.



Reference part shown: 50ACV-SB3-06

50ACV-SB3-04

LinkTech 50AC Series Panel Mount Female, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-SB3-05

LinkTech 50AC Series Panel Mount Female, Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-SB3-06

LinkTech 50AC Series Panel Mount Female, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-SB8-04

LinkTech 50AC Series Panel Mount Female, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

50ACV-SB8-06

LinkTech 50AC Series Panel Mount Female, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

50ACV-SB18-04

LinkTech 50AC Series Panel Mount Female, Valved, 1/4" OD Push-In, Natural Acetal, Buna-N O-ring

50BACV-SB3-04

LinkTech 50BAC Series Panel Mount Female, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

50BACV-SB3-05

LinkTech 50BAC Series Panel Mount Female, Valved, 5/16" Hose Barb, Natural Acetal, Buna-N O-ring

50BACV-SB3-06

LinkTech 50BAC Series Panel Mount Female, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

50BACV-SB8-04

LinkTech 50BAC Series Panel Mount Female, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

50BACV-SB8-06

LinkTech 50BAC Series Panel Mount Female, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

50PPV-SE3-04

LinkTech 50PP Series Panel Mount Female, Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

50PPV-SE3-06

LinkTech 50PP Series Panel Mount Female, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring



Reference part shown: 50ACV-PB4-06

50ACV-PB4-04

LinkTech 50AC Series Elbow Male, Valved, 1/4" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-PB4-06

LinkTech 50AC Series Elbow Male, Valved, 3/8" Hose Barb, Natural Acetal, Buna-N O-ring

50ACV-PB10-04

LinkTech 50AC Series Elbow Male, Valved, 1/4" PTF, Natural Acetal, Buna-N O-ring

50ACV-PB10-06

LinkTech 50AC Series Elbow Male, Valved, 3/8" PTF, Natural Acetal, Buna-N O-ring

50PPV-PE4-04

LinkTech 50PP Series Elbow Male, Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

50PPV-PE4-06

LinkTech 50PP Series Elbow Male, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

55Series

Quick Connect Couplings

55 SERIES - 1/8" - 1/4" FLOW SIZE:
Female Valved Straight



Reference part shown: 55GPV-SE2-04

55GPV-SE1-04
LinkTech 55GP Series Male Thread Female, Valved, 1/4 NPT,
Glass-Filled Polypropylene, EPDM O-ring

55GPV-SE2-04
LinkTech 55GP Series Straight Female, Valved, 1/4" Hose Barb,
Glass-Filled Polypropylene, EPDM O-ring

55GPV-SE2-06
LinkTech 55GP Series Straight Female, Valved, 3/8" Hose Barb,
Glass-Filled Polypropylene, EPDM O-ring

55 SERIES - 1/8" - 1/4" FLOW SIZE:
Male Valved Straight



Reference part shown: 55GPV-PE2-04

55GPV-PE1-04
LinkTech 55GP Series Male Thread Male, Valved, 1/4 NPT,
Glass-Filled Polypropylene, EPDM O-ring

55GPV-PE2-04
LinkTech 55GP Series Straight Male, Valved, 1/4" Hose Barb,
Glass-Filled Polypropylene, EPDM O-ring

55GPV-PE2-06
LinkTech 55GP Series Straight Male, Valved, 3/8" Hose Barb,
Glass-Filled Polypropylene, EPDM O-ring

60Series

Quick Connect Couplings

60 SERIES - 3/8" FLOW SIZE:
Female Open Flow Straight



Reference part shown: 60PP-S2-08

60PP-S1-06
LinkTech 60PP Series Male Thread Female, Non-Valved, 3/8" MNPT,
Polypropylene

60PP-S1-08
LinkTech 60PP Series Male Thread Female, Non-Valved, 1/2" MNPT,
Polypropylene

60PP-S1-08BLK
LinkTech 60PP Series Male Thread Female, Non-Valved, 1/2" MNPT,
Molded Black Color, Polypropylene

60PP-S1-12
LinkTech 60PP Series Male Thread Female, Non-Valved, 3/4" MNPT,
Polypropylene

60PP-S1-12BLK
LinkTech 60PP Series Male Thread Female, Non-Valved, 3/4" MNPT,
Molded Black Color, Polypropylene

60PP-S2-04
LinkTech 60PP Series Straight Female, Non-Valved, 1/4" Hose Barb,
Polypropylene

60PP-S2-04BLK
LinkTech 60PP Series Straight Female, Non-Valved, 1/4" Hose Barb,
Molded Black Color, Polypropylene

60PP-S2-06
LinkTech 60PP Series Straight Female, Non-Valved, 3/8" Hose Barb,
Polypropylene

60PP-S2-08
LinkTech 60PP Series Straight Female, Non-Valved, 1/2" Hose Barb,
Polypropylene

60PP-S2-08BLK
LinkTech 60PP Series Straight Female, Non-Valved, 1/2" Hose Barb,
Molded Black Color, Polypropylene

60PP-S2-12
LinkTech 60PP Series Straight Female, Non-Valved, 3/4" Hose Barb,
Polypropylene

QUICK
CONNECTS

60 SERIES - 3/8" FLOW SIZE: Female Open Flow Straight

60PS-S1-06

LinkTech 60PS Series Male Thread Female, Non-Valved, 3/8" MNPT, Polysulfone

60PS-S1-06BLK

LinkTech 60PS Series Male Thread Female, Non-Valved, 3/8" MNPT, Molded Black Color, Polysulfone

60PS-S1-08

LinkTech 60PS Series Male Thread Female, Non-Valved, 1/2" MNPT, Polysulfone

60PS-S1-08BLK

LinkTech 60PS Series Male Thread Female, Non-Valved, 1/2" MNPT, Molded Black Color, Polysulfone

60PS-S2-06

LinkTech 60PS Series Straight Female, Non-Valved, 3/8" Hose Barb, Polysulfone

60PS-S2-06BLK

LinkTech 60PS Series Straight Female, Non-Valved, 3/8" Hose Barb, Molded Black Color, Polysulfone

60PS-S2-08

LinkTech 60PS Series Straight Female, Non-Valved, 1/2" Hose Barb, Polysulfone

60PS-S2-08BLK

LinkTech 60PS Series Straight Female, Non-Valved, 1/2" Hose Barb, Molded Black Color, Polysulfone

60PS-S2-12

LinkTech 60PS Series Straight Female, Non-Valved, 3/4" Hose Barb, Polysulfone

60PS-S2-12BLK

LinkTech 60PS Series Straight Female, Non-Valved, 3/4" Hose Barb, Molded Black Color, Polysulfone

60 SERIES - 3/8" FLOW SIZE: Male Open Flow Straight



Reference part shown: 60PP-PE2-04

60PP-PE1-06

LinkTech 60PP Series Male Thread Male, Non-Valved, 3/8" MNPT, Polypropylene, EPDM O-ring

60PP-PE1-08

LinkTech 60PP Series Male Thread Male, Non-Valved, 1/2" MNPT, Polypropylene, EPDM O-ring

60PP-PE1-08BLK

LinkTech 60PP Series Male Thread Male, Non-Valved, 1/2" MNPT, Molded Black Color, Polypropylene, EPDM O-ring

60PP-PE1-12

LinkTech 60PP Series Male Thread Male, Non-Valved, 3/4" MNPT, Polypropylene, EPDM O-ring

60 SERIES - 3/8" FLOW SIZE: Male Open Flow Straight

60PP-PE1-12BLK

LinkTech 60PP Series Male Thread Male, Non-Valved, 3/4" MNPT, Molded Black Color, Polypropylene, EPDM O-ring

60PP-PE2-04

LinkTech 60PP Series Straight Male, Non-Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

60PP-PE2-04BLK

LinkTech 60PP Series Straight Male, Non-Valved, 1/4" Hose Barb, Molded Black Color, Polypropylene, EPDM O-ring

60PP-PE2-06

LinkTech 60PP Series Straight Male, Non-Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

60PP-PE2-08

LinkTech 60PP Series Straight Male, Non-Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PP-PE2-08BLK

LinkTech 60PP Series Straight Male, Non-Valved, 1/2" Hose Barb, Molded Black Color, Polypropylene, EPDM O-ring

60PP-PE2-12

LinkTech 60PP Series Straight Male, Non-Valved, 3/4" Hose Barb, Polypropylene, EPDM O-ring

60PP-PE2-12BLK

LinkTech 60PP Series Straight Male, Non-Valved, 3/4" Hose Barb, Molded Black Color, Polypropylene, EPDM O-ring

60PS-PE1-06

LinkTech 60PS Series Male Thread Male, Non-Valved, 3/8" MNPT, Polysulfone, EPDM O-ring

60PS-PE1-06BLK

LinkTech 60PS Series Male Thread Male, Non-Valved, 3/8" MNPT, Molded Black Color, Polysulfone, EPDM O-ring

60PS-PE1-08

LinkTech 60PS Series Male Thread Male, Non-Valved, 1/2" MNPT, Polysulfone, EPDM O-ring

60PS-PE1-08BLK

LinkTech 60PS Series Male Thread Male, Non-Valved, 1/2" MNPT, Molded Black Color, Polysulfone, EPDM O-ring

60PS-PE2-06

LinkTech 60PS Series Straight Male, Non-Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

60PS-PE2-06BLK

LinkTech 60PS Series Straight Male, Non-Valved, 3/8" Hose Barb, Molded Black Color, Polysulfone, EPDM O-ring

60PS-PE2-08

LinkTech 60PS Series Straight Male, Non-Valved, 1/2" Hose Barb, Polysulfone, EPDM O-ring

60PS-PE2-08BLK

LinkTech 60PS Series Straight Male, Non-Valved, 1/2" Hose Barb, Molded Black Color, Polysulfone, EPDM O-ring

60PS-PE2-12

LinkTech 60PS Series Straight Male, Non-Valved, 3/4" Hose Barb, Polysulfone, EPDM O-ring

60PS-PE2-12BLK

LinkTech 60PS Series Straight Male, Non-Valved, 3/4" Hose Barb, Molded Black Color, Polysulfone, EPDM O-ring

60 SERIES - 3/8" FLOW SIZE: Female Open Flow Panel Mount



Reference part shown: 60PP-S3-08

60PP-S3-04

LinkTech 60PP Series Panel Mount Female, Non-Valved, 1/4" Hose Barb, Polypropylene

60PP-S3-06

LinkTech 60PP Series Panel Mount Female, Non-Valved, 3/8" Hose Barb, Polypropylene

60PP-S3-08

LinkTech 60PP Series Panel Mount Female, Non-Valved, 1/2" Hose Barb, Polypropylene

60PP-S3-12

LinkTech 60PP Series Panel Mount Female, Non-Valved, 3/4" Hose Barb, Polypropylene

60PP-S22-06

LinkTech 60PP Series Panel Mount Female, Non-Valved, 3/8" Hose Barb Elbow, Polypropylene

60PP-S22-08

LinkTech 60PP Series Panel Mount Female, Non-Valved, 1/2" Hose Barb Elbow, Polypropylene

60PS-S3-06

LinkTech 60PS Series Panel Mount Female, Non-Valved, 3/8" Hose Barb, Polysulfone

60PS-S3-06BLK

LinkTech 60PS Series Panel Mount Female, Non-Valved, 3/8" Hose Barb, Molded Black Color, Polysulfone

60PS-S3-08

LinkTech 60PS Series Panel Mount Female, Non-Valved, 1/2" Hose Barb, Polysulfone

60PS-S3-08BLK

LinkTech 60PS Series Panel Mount Female, Non-Valved, 1/2" Hose Barb, Molded Black Color, Polysulfone

60PS-S3-12

LinkTech 60PS Series Panel Mount Female, Non-Valved, 3/4" Hose Barb, Polysulfone

60 SERIES - 3/8" FLOW SIZE: Male Open Flow Panel Mount



Reference part shown: 60PP-PE3-06

60PP-PE3-04

LinkTech 60PP Series Panel Mount Male, Non-Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

60PP-PE3-06

LinkTech 60PP Series Panel Mount Male, Non-Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

60PP-PE3-08

LinkTech 60PP Series Panel Mount Male, Non-Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PP-PE3-08BLK

LinkTech 60PP Series Panel Mount Male, Non-Valved, 1/2" Hose Barb, Molded Black Color, Polypropylene, EPDM O-ring

60PP-PE3-12

LinkTech 60PP Series Panel Mount Male, Non-Valved, 3/4" Hose Barb, Polypropylene, EPDM O-ring

60PS-PE3-06

LinkTech 60PS Series Panel Mount Male, Non-Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

60PS-PE3-08

LinkTech 60PS Series Panel Mount Male, Non-Valved, 1/2" Hose Barb, Polysulfone, EPDM O-ring

60PS-PE3-12

LinkTech 60PS Series Panel Mount Male, Non-Valved, 3/4" Hose Barb, Polysulfone, EPDM O-ring

Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

LinkTech Line Part Number Guide

60	PP	-	S	-	3	-	08
Flow Size or Latch Type	Material	No Valve	Config.	No O-Ring	Type	Connection Size	



Reference part shown: 60PP-PE22-08

Reference part shown: 60PPV-SE2-08

60PP-PE4-06

LinkTech 60PP Series Elbow Male, Non-Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

60PPV-SE1-06

LinkTech 60PP Series Male Thread Female, Valved, 3/8" MNPT, Polypropylene, EPDM O-ring

60PP-PE4-08

LinkTech 60PP Series Elbow Male, Non-Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PPV-SE1-08

LinkTech 60PP Series Male Thread Female, Valved, 1/2" MNPT, Polypropylene, EPDM O-ring

60PP-PE22-08

LinkTech 60PP Series Panel Mount Elbow Male, Non-Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PPV-SE1-08BLK

LinkTech 60PP Series Male Thread Female, Valved, 1/2" MNPT, Molded Black Color, Polypropylene, EPDM O-ring

60PS-PE4-06

LinkTech 60PS Series Elbow Male, Non-Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

60PPV-SE1-12

LinkTech 60PP Series Male Thread Female, Valved, 3/4" MNPT, Polypropylene, EPDM O-ring

60PS-PE4-08

LinkTech 60PS Series Elbow Male, Non-Valved, 1/2" Hose Barb, Polysulfone, EPDM O-ring

60PPV-SE2-04

LinkTech 60PP Series Straight Female, Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

60PS-PE4-08BLK

LinkTech 60PS Series Elbow Male, Non-Valved, 1/2" Hose Barb, Molded Black Color, Polysulfone, EPDM O-ring

60PPV-SE2-06

LinkTech 60PP Series Straight Female, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

60PS-PE22-06

LinkTech 60PS Series Panel Mount Elbow Male, Non-Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

60PPV-SE2-06BLK

LinkTech 60PP Series Straight Female, Valved, 3/8" Hose Barb, Molded Black Color, Polypropylene, EPDM O-ring

60PS-PE22-08

LinkTech 60PS Series Panel Mount Elbow Male, Non-Valved, 1/8" Hose Barb, Polysulfone, EPDM O-ring

60PPV-SE2-08

LinkTech 60PP Series Straight Female, Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PPV-SE2-08BLK

LinkTech 60PP Series Straight Female, Valved, 1/2" Hose Barb, Molded Black Color, Polypropylene, EPDM O-ring

60PPV-SE2-12

LinkTech 60PP Series Straight Female, Valved, 3/4" Hose Barb, Polypropylene, EPDM O-ring

60PSV-SE1-06

LinkTech 60PS Series Male Thread Female, Valved, 3/8" MNPT, Polysulfone, EPDM O-ring

60PSV-SE1-08

LinkTech 60PS Series Male Thread Female, Valved, 1/2" MNPT, Polysulfone, EPDM O-ring

60PSV-SE1-08BLK

LinkTech 60PS Series Male Thread Female, Valved, 1/2" MNPT, Molded Black Color, Polysulfone, EPDM O-ring

60PSV-SE2-06

LinkTech 60PS Series Straight Female, Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

60 SERIES - 3/8" FLOW SIZE:

Female Valved Straight

60PSV-SE2-06BLK

LinkTech 60PS Series Straight Female, Valved, 3/8" Hose Barb,
Molded Black Color, Polysulfone, EPDM O-ring

60PSV-SE2-08

LinkTech 60PS Series Straight Female, Valved, 1/2" Hose Barb,
Polysulfone, EPDM O-ring

60PSV-SE2-08BLK

LinkTech 60PS Series Straight Female, Valved, 1/2" Hose Barb,
Molded Black Color, Polysulfone, EPDM O-ring

60PSV-SE2-12

LinkTech 60PS Series Straight Female, Valved, 3/4" Hose Barb,
Polysulfone, EPDM O-ring

60PSV-SE2-12BLK

LinkTech 60PS Series Straight Female, Valved, 3/4" Hose Barb,
Molded Black Color, Polysulfone, EPDM O-ring

60 SERIES - 3/8" FLOW SIZE:

Male Valved Straight

Reference part shown: **60PPV-PE2-08****60PPV-PE1-06**

LinkTech 60PP Series Male Thread Male, Valved, 3/8" MNPT,
Polypropylene, EPDM O-ring

60PPV-PE1-06BLK

LinkTech 60PP Series Male Thread Male, Valved, 3/8" MNPT,
Molded Black Color, Polypropylene, EPDM O-ring

60PPV-PE1-08

LinkTech 60PP Series Male Thread Male, Valved, 1/2" MNPT,
Polypropylene, EPDM O-ring

60PPV-PE1-08BLK

LinkTech 60PP Series Male Thread Male, Valved, 1/2" MNPT,
Molded Black Color, Polypropylene, EPDM O-ring

60PPV-PE1-12

LinkTech 60PP Series Male Thread Male, Valved, 3/4" MNPT,
Polypropylene, EPDM O-ring

60PPV-PE2-04

LinkTech 60PP Series Straight Male, Valved, 1/4" Hose Barb,
Polypropylene, EPDM O-ring

60PPV-PE2-04BLK

LinkTech 60PP Series Straight Male, Valved, 1/4" Hose Barb,
Molded Black Color, Polypropylene, EPDM O-ring

60PPV-PE2-06

LinkTech 60PP Series Straight Male, Valved, 3/8" Hose Barb,
Polypropylene, EPDM O-ring

60 SERIES - 3/8" FLOW SIZE:

Male Valved Straight

60PPV-PE2-06BLK

LinkTech 60PP Series Straight Male, Valved, 3/8" Hose Barb,
Molded Black Color, Polypropylene, EPDM O-ring

60PPV-PE2-08

LinkTech 60PP Series Straight Male, Valved, 1/2" Hose Barb,
Polypropylene, EPDM O-ring

60PPV-PE2-08BLK

LinkTech 60PP Series Straight Male, Valved, 1/2" Hose Barb,
Molded Black Color, Polypropylene, EPDM O-ring

60PPV-PE2-12

LinkTech 60PP Series Straight Male, Valved, 3/4" Hose Barb,
Polypropylene, EPDM O-ring

60PSV-PE1-06

LinkTech 60PS Series Male Thread Male, Valved, 3/8" MNPT,
Polysulfone, EPDM O-ring

60PSV-PE1-06BLK

LinkTech 60PS Series Male Thread Male, Valved, 3/8" MNPT,
Molded Black Color, Polysulfone, EPDM O-ring

60PSV-PE1-08

LinkTech 60PS Series Male Thread Male, Valved, 1/2" MNPT,
Polysulfone, EPDM O-ring

60PSV-PE1-08BLK

LinkTech 60PS Series Male Thread Male, Valved, 1/2" MNPT,
Molded Black Color, Polysulfone, EPDM O-ring

60PSV-PE2-06

LinkTech 60PS Series Straight Male, Valved, 3/8" Hose Barb,
Polysulfone, EPDM O-ring

60PSV-PE2-06BLK

LinkTech 60PS Series Straight Male, Valved, 3/8" Hose Barb,
Molded Black Color, Polysulfone, EPDM O-ring

60PSV-PE2-08

LinkTech 60PS Series Straight Male, Valved, 1/2" Hose Barb,
Polysulfone, EPDM O-ring

60PSV-PE2-08BLK

LinkTech 60PS Series Straight Male, Valved, 1/2" Hose Barb,
Molded Black Color, Polysulfone, EPDM O-ring

60PSV-PE2-12

LinkTech 60PS Series Straight Male, Valved, 3/4" Hose Barb,
Polysulfone, EPDM O-ring

60PSV-PE2-12BLK

LinkTech 60PS Series Straight Male, Valved, 3/4" Hose Barb,
Molded Black Color, Polysulfone, EPDM O-ring

Need it Fast?

You can order any of our 9,600+ parts online now with no
order \$ minimums, 24 hours a day, 7 days a week!

www.nordsonmedical.com



Reference part shown: 60PPV-SE3-08

60PPV-SE3-04

LinkTech 60PP Series Panel Mount Female, Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

60PPV-SE3-06

LinkTech 60PP Series Panel Mount Female, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

60PPV-SE3-08

LinkTech 60PP Series Panel Mount Female, Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PPV-SE3-12

LinkTech 60PP Series Panel Mount Female, Valved, 3/4" Hose Barb, Polypropylene, EPDM O-ring

60PSV-SE3-06

LinkTech 60PS Series Panel Mount Female, Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

60PSV-SE3-06BLK

LinkTech 60PS Series Panel Mount Female, Valved, 3/8" Hose Barb, Molded Black Color, Polysulfone, EPDM O-ring

60PSV-SE3-08

LinkTech 60PS Series Panel Mount Female, Valved, 1/2" Hose Barb, Polysulfone, EPDM O-ring

60PSV-SE3-12

LinkTech 60PS Series Panel Mount Female, Valved, 3/4" Hose Barb, Polysulfone, EPDM O-ring

60PSV-SE3-12BLK

LinkTech 60PS Series Panel Mount Female, Valved, 3/4" Hose Barb, Molded Black Color, Polysulfone, EPDM O-ring



Reference part shown: 60PPV-PE3-08

60PPV-PE3-04

LinkTech 60PP Series Panel Mount Male, Valved, 1/4" Hose Barb, Polypropylene, EPDM O-ring

60PPV-PE3-06

LinkTech 60PP Series Panel Mount Male, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

60PPV-PE3-08

LinkTech 60PP Series Panel Mount Male, Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PPV-PE3-12

LinkTech 60PP Series Panel Mount Male, Valved, 3/4" Hose Barb, Polypropylene, EPDM O-ring

60PSV-PE3-06

LinkTech 60PS Series Panel Mount Male, Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

60PSV-PE3-06BLK

LinkTech 60PS Series Panel Mount Male, Valved, 3/8" Hose Barb, Molded Black Color, Polysulfone, EPDM O-ring

60PSV-PE3-08

LinkTech 60PS Series Panel Mount Male, Valved, 1/2" Hose Barb, Polysulfone, EPDM O-ring

60PSV-PE3-12

LinkTech 60PS Series Panel Mount Male, Valved, 3/4" Hose Barb, Polysulfone, EPDM O-ring

Did You Know?

Nordson MEDICAL has an interactive PDF available online at www.nordsonmedical.com. Access part drawings, reference charts, data sheets, and more from just one location!

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.


Reference part shown: 60PPV-SE22-08
60PPV-SE22-06

LinkTech 60PP Series Panel Mount Elbow Female, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

60PPV-SE22-08

LinkTech 60PP Series Panel Mount Elbow Female, Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PPV-SE22-08BLK

LinkTech 60PP Series Panel Mount Elbow Female, Valved, 1/2" Hose Barb, Molded Black Color, Polypropylene, EPDM O-ring

60PSV-SE22-06

LinkTech 60PS Series Panel Mount Elbow Female, Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

60PSV-SE22-08

LinkTech 60PS Series Panel Mount Elbow Female, Valved, 1/2" Hose Barb, Polysulfone, EPDM O-ring


Reference part shown: 60PPV-PE4-08
60PPV-PE4-06

LinkTech 60PP Series Elbow Male, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

60PPV-PE4-08

LinkTech 60PP Series Elbow Male, Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PPV-PE22-06

LinkTech 60PP Series Panel Mount Elbow Male, Valved, 3/8" Hose Barb, Polypropylene, EPDM O-ring

60PPV-PE22-08

LinkTech 60PP Series Panel Mount Elbow Male, Valved, 1/2" Hose Barb, Polypropylene, EPDM O-ring

60PPV-PE22-08BLK

LinkTech 60PP Series Panel Mount Elbow Male, Valved, 1/2" Hose Barb, Molded Black Color, Polypropylene, EPDM O-ring

60PSV-PE4-06

LinkTech 60PS Series Elbow Male, Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

60PSV-PE4-08

LinkTech 60PS Series Elbow Male, Valved, 1/2" Hose Barb, Polysulfone, EPDM O-ring

60PSV-PE22-06

LinkTech 60PS Series Panel Mount Elbow Male, Valved, 3/8" Hose Barb, Polysulfone, EPDM O-ring

60PSV-PE22-06BLK

LinkTech 60PS Series Panel Mount Elbow Male, Valved, 3/8" Hose Barb, Molded Black Color, Polysulfone, EPDM O-ring

60PSV-PE22-08

LinkTech 60PS Series Panel Mount Elbow Male, Valved, 1/2" Hose Barb, Polysulfone, EPDM O-ring

Technical Assistance

 Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

Part Drawings Online!

 Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.
www.nordsonmedical.com
Need Samples for Evaluation?

 You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.
+1 (970) 267-5200 Toll Free: (888) 404-5837

65Series

Quick Connect Couplings

65 SERIES - 1/2" FLOW SIZE:

Female Open Flow Straight



Reference part shown: 65GP-S2-12

65GP-S2-12

LinkTech 65GP Series Straight Female, Non-Valved, 3/4" Hose Barb, Glass-Filled Polypropylene

65 SERIES - 1/2" FLOW SIZE:

Male Open Flow Straight



Reference part shown: 65GP-PF1-08

65GP-PF1-08

LinkTech 65GP Series Male Thread Male, Non-Valved, 1/2" MNPT, Glass-Filled Polypropylene, Food Grade EPDM O-ring

65GP-PF2-08

LinkTech 65GP Series Straight Male, Non-Valved, 1/2" Hose Barb, Glass-Filled Polypropylene, Food Grade EPDM O-ring

70Series

Quick Connect Couplings

70 SERIES - 1/4" FLOW SIZE:

Female Valved Straight



Reference part shown: 70GPV-SE2-04

70GPV-SE2-04

LinkTech 70 Series Straight Female, Valved, 1/4" Hose Barb, Glass-Filled Polypropylene, Non-Exposed EPDM O-ring Seals

70GPV-SE2-06

LinkTech 70 Series Straight Female, Valved, 3/8" Hose Barb, Glass-Filled Polypropylene, Non-Exposed EPDM O-ring Seals

70 SERIES - 1/4" FLOW SIZE:

Male Valved Straight



Reference part shown: 70GPV-PE2-04

70GPV-PE2-04

LinkTech 70 Series Straight Male, Valved, 1/4" Hose Barb, Glass-Filled Polypropylene, Non-Exposed EPDM O-ring Seals

70GPV-PE2-06

LinkTech 70 Series Straight Male, Valved, 3/8" Hose Barb, Glass-Filled Polypropylene, Non-Exposed EPDM O-ring Seals

QUICK
CONNECTS

100Series

Quick Connect Couplings

100 SERIES - 1/16" FLOW SIZE:
Female Open Flow Straight



Reference part shown: 100AC-S2-02

100AC-S2-01

LinkTech 100AC Series Straight Female, Non-Valved, 1/16" Hose Barb, Acetal

100AC-S2-02

LinkTech 100AC Series Straight Female, Non-Valved, 1/8" Hose Barb, Acetal

100AC-S2-M3

LinkTech 100AC Series Straight Female, Non-Valved, 3 mm Hose Barb, Acetal

100PP-S2-01

LinkTech 100PP Series Straight Female, Non-Valved, 1/16" Hose Barb, Polypropylene

100PP-S2-02

LinkTech 100PP Series Straight Female, Non-Valved, 1/8" Hose Barb, Polypropylene

100PP-S2-M3

LinkTech 100PP Series Straight Female, Non-Valved, 3 mm Hose Barb, Polypropylene

100-S2

LinkTech 100 Series Female Housing without Inserts, Acetal



Reference part shown: 100AC-PB2-02

100AC-PB2-01

LinkTech 100AC Series Straight Male, Non-Valved, 1/16" Hose Barb, Acetal, Buna-N O-ring

100AC-PB2-02

LinkTech 100AC Series Straight Male, Non-Valved, 1/8" Hose Barb, Acetal, Buna-N O-ring

100-P2

LinkTech 100 Series Male Housing without Inserts, Acetal

100PP-PE2-01

LinkTech 100PP Series Straight Male, Non-Valved, 1/16" Hose Barb, Polypropylene, EPDM O-ring

100PP-PE2-02

LinkTech 100PP Series Straight Male, Non-Valved, 1/8" Hose Barb, Polypropylene, EPDM O-ring

100 SERIES - 1/16" FLOW SIZE:
Male Open Flow Panel Mount



Reference part shown: 100-PM3

100-PM3

LinkTech 100 Series Panel Mount Adapter and Nut, 7/16", Acetal

Need More Details?

Check out the quick connect Datasheets at
www.nordsonmedical.com or contact Sales at
FMF-sales@nordsonmedical.com.



Reference part shown: 100ACV-SB2-02

100ACV-SB2-01

LinkTech 100AC Series Straight Female, Valved, 1/16" Hose Barb, Acetal, Buna-N O-ring

100ACV-SB2-02

LinkTech 100AC Series Straight Female, Valved, 1/8" Hose Barb, Acetal, Buna-N O-ring



Reference part shown: 100AC2-02

100AC2-01

LinkTech 100AC Series Straight Coupling, Non-Valved, 1/16" Hose Barb, Acetal, Buna-N O-ring

100AC2-02

LinkTech 100AC Series Straight Coupling, Non-Valved, 1/8" Hose Barb, Acetal, Buna-N O-ring

100AC2-M3

LinkTech 100AC Series Straight Coupling, Non-Valved, 3 mm Hose Barb, Acetal, Buna-N O-ring

100ACV2-01

LinkTech 100AC Series Straight Coupling, Valved, 1/16" Hose Barb, Acetal, Buna-N O-ring

100ACV2-02

LinkTech 100AC Series Straight Coupling, Valved, 1/8" Hose Barb, Acetal, Buna-N O-ring

100ACV2-M3

LinkTech 100AC Series Straight Coupling, Valved, 3 mm Hose Barb, Acetal, Buna-N O-ring

100PP2-01

LinkTech 100PP Series Straight Coupling, Non-Valved, 1/16" Hose Barb, Polypropylene, EPDM O-ring

100PP2-02

LinkTech 100PP Series Straight Coupling, Non-Valved, 1/8" Hose Barb, Polypropylene, EPDM O-ring

100PP2-M3

LinkTech 100PP Series Straight Coupling, Non-Valved, 3 mm Hose Barb, Polypropylene, EPDM O-ring

Did You Know?

Nordson MEDICAL has an interactive PDF available online at www.nordsonmedical.com. Access part drawings, reference charts, data sheets, and more from just one location!

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

LinkTech Line Part Number Guide

100	AC	V - S	B	2 - 01
Flow Size or Latch Type	Material	Valve Config.	O-Ring	Type Connection Size

MQC Series

Quick Connect Couplings

MQC SERIES: Female Open Flow Straight



Reference part shown: MQCF230

MQCF230 (-1 -2 -3 -4 -5 -8)

Dual Latch Quick Connect Female 200 Series Barb 1/8" (3.2 mm) Tubing ID

MQCF240 (-1 -2 -3 -4 -5 -8)

Dual Latch Quick Connect Female 200 Series Barb 5/32" (4.0 mm) Tubing ID

MQCF250 (-1 -2 -3 -4 -5 -8)

Dual Latch Quick Connect Female 200 Series Barb 3/16" (4.8 mm) Tubing ID

MQC SERIES: Male Open Flow Straight



Reference part shown: MQCRM230

MQCRM230 (-001 -002 -003 -004 -005 -08B)

Nylon, Buna-N O-Ring

Dual Latch Quick Connect Male 200 Series Barb 1/8" (3.2 mm) Tubing ID

MQCRM240 (-001 -002 -003 -004 -005 -08B)

Nylon, Buna-N O-Ring

Dual Latch Quick Connect Male 200 Series Barb 5/32" (4.0 mm) Tubing ID

MQCRM250 (-001 -002 -003 -004 -005 -08B)

Nylon, Buna-N O-Ring

Dual Latch Quick Connect Male 200 Series Barb 3/16" (4.8 mm) Tubing ID



Reference part shown: MQCMF230

MQCMF230-901

Bag of 50 Dual Latch Quick Connect Male and Female Sets, 200 Series Barb 1/8" (3.2 mm) Tubing ID, Blue Female and White Male

MQCMF240-901

Bag of 50 Dual Latch Quick Connect Male and Female Sets, 200 Series Barb 5/32" (4.0 mm) Tubing ID, Blue Female and White Male

MQCMF250-901

Bag of 50 Dual Latch Quick Connect Male and Female Sets, 200 Series Barb 3/16" (4.8 mm) Tubing ID, Blue Female and White Male

Looking for Something Different?

Most of our best product ideas have come from you, our customers. While we pride ourselves on our extensive product offering, we recognize that your needs are ever-changing and that our products, too, must adapt to your needs.

If there's a new product you think Nordson MEDICAL should make, please send us an Email at designcenter@nordsonmedical.com describing your product concept in as much detail as you desire. We will contact you with any questions, and who knows - the next catalog edition could include products YOU told us to build.

See page 10 for more information on the Design Center!



RQSeries

Quick Connect Couplings

RQ SERIES:

Female Open Flow Straight



Reference part shown: RQXFL670-9024

RQCF655 (-9024 -001)**RQCF655 (-6005 -001)**

Female Open Flow, 600 Series Barb, 1/4" (6.4 mm) Tubing ID

RQCF670 (-9024 -001)**RQCF670 (-6005 -001)**

Female Open Flow, 600 Series Barb, 3/8" (9.5 mm) Tubing ID

RQCFL655 (-9024 -001)**RQCFL655 (-6005 -001)**

Female Open Flow, 600 Series Barb, 1/4" (6.4 mm) Tubing ID, Locking Feature

RQCFL670 (-9024 -001)**RQCFL670 (-6005 -001)**

Female Open Flow, 600 Series Barb, 3/8" (9.5 mm) Tubing ID, Locking Feature

RQXF670 (-9024 -001)**RQXF670 (-6005 -001)**

Female Open Flow, 600 Series Barb, 3/8" (9.5 mm) Tubing ID

RQXF680 (-9024 -001)**RQXF680 (-6005 -001)**

Female Open Flow, 600 Series Barb, 1/2" (12.7 mm) Tubing ID

RQXFL670 (-9024 -001)**RQXFL670 (-6005 -001)**

Female Open Flow, 600 Series Barb, 3/8" (9.5 mm) Tubing ID, Locking Feature

RQXFL680 (-9024 -001)**RQXFL680 (-6005 -001)**

Female Open Flow, 600 Series Barb, 1/2" (12.7 mm) Tubing ID, Locking Feature

RQ SERIES:

Male Open Flow Straight



Reference part shown: RQXM670-9024

RQCM655 (-9024 -001) Platinum-Cured Silicone O-Ring**RQCM655 (-6005 -001) Platinum-Cured Silicone O-Ring**

Male Open Flow, 600 Series Barb, 1/4" (6.4 mm) Tubing ID

RQCM670 (-9024 -001) Platinum-Cured Silicone O-Ring**RQCM670 (-6005 -001) Platinum-Cured Silicone O-Ring**

Male Open Flow, 600 Series Barb, 3/8" (9.5 mm) Tubing ID

RQXM670 (-9024 -001) Platinum-Cured Silicone O-Ring**RQXM670 (-6005 -001) Platinum-Cured Silicone O-Ring**

Male Open Flow, 600 Series Barb, 3/8" (9.5 mm) Tubing ID

RQXM680 (-9024 -001) Platinum-Cured Silicone O-Ring**RQXM680 (-6005 -001) Platinum-Cured Silicone O-Ring**

Male Open Flow, 600 Series Barb, 1/2" (12.7 mm) Tubing ID

RQ SERIES:

Caps



Reference part shown: RQXFLC-9024

RQCFC (-9024 -001)**RQCFC (-6005 -001)**

Female Cap

RQCFLC (-9024 -001)**RQCFLC (-6005 -001)**

Female Cap, Locking Feature

RQXFC (-9024 -001)**RQXFC (-6005 -001)**

Female Cap

RQXFLC (-9024 -001)**RQXFLC (-6005 -001)**

Female Cap, Locking Feature

RQ SERIES:

Plugs



Reference part shown: RQXMP-9024

RQCMP (-9024 -001) Platinum-Cured Silicone O-Ring**RQCMP (-6005 -001) Platinum-Cured Silicone O-Ring**

Male Plug

RQXMP (-9024 -001) Platinum-Cured Silicone O-Ring**RQXMP (-6005 -001) Platinum-Cured Silicone O-Ring**

Male Plug

SBL Series

Quick Connect Couplings

QUICK
CONNECTS

SBL SERIES:
Female Open Flow Straight



Reference part shown: **SBLF035**

SBLF013 (-1 -2 -8B -CM030)
Side Latch Quick Connect Female 500 Series Barb, 1/8" (3.2 mm) Tubing ID

SBLF035 (-1 -2 -8B -CM030)
Side Latch Quick Connect Female 500 Series Barb, 3/16" (4.8 mm) Tubing ID

SBL SERIES:
Male Open Flow Straight



Reference part shown: **SBLM013**

SBLM013 (-001 -002 -004 -005) Nylon, Buna-N O-Ring
Side Latch Quick Connect Male 500 Series Barb, 1/8" (3.2 mm) Tubing ID

SBLM035 (-001 -002 -004 -005) Nylon, Buna-N O-Ring
Side Latch Quick Connect Male 500 Series Barb, 3/16" (4.8 mm) Tubing ID



Reference part shown: **SBLMFBG035-003**

SBLMFBG013 (-002)
Bag of 50 Side Latch Quick Connect Male and Female Sets, 500 Series Barb 1/8" (3.2 mm) Tubing ID, Blue Male and Grey Female Nylon

SBLMFBG035 (-004)
Bag of 50 Side Latch Quick Connect Male and Female Sets, 500 Series Barb 3/16" (4.8 mm) Tubing ID, Blue Male and Grey Female Nylon

SBLMFBG013 (-001)
Bag of 50 Side Latch Quick Connect Male and Female Sets, 500 Series Barb 1/8" (3.2 mm) Tubing ID, Grey Male and Blue Female Nylon

SBLMFBG035 (-003)
Bag of 50 Side Latch Quick Connect Male and Female Sets, 500 Series Barb 3/16" (4.8 mm) Tubing ID, Grey Male and Blue Female Nylon

Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.
+1 (970) 267-5200 Toll Free: (888) 404-5837

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

XQSeries

Quick Connect Couplings



Reference part shown: XQSBF755

XQSBF755 (-1006)

XQ Series Female Open Flow Panel Mount Connector with 11/16-24 UNEF Thread, 700 Series Barb, 1/4" (6.4 mm) Tubing ID, Nickel-Plated Brass Hex Nut (included)

XQSBF770 (-1006)

XQ Series Female Open Flow Panel Mount Connector with 11/16-24 UNEF Thread, 700 Series Barb, 3/8" (9.5 mm) Tubing ID, Nickel-Plated Brass Hex Nut (included)

QUICK
CONNECTS

XQ SERIES - 1/4" FLOW SIZE:
Female Open Flow Straight



Reference part shown: XQSF755

XQSF755 (-1006)

XQ Series Female Open Flow Connector, 700 Series Barb, 1/4" (6.4 mm) Tubing ID

XQSF770 (-1006 -1007)

XQ Series Female Open Flow Connector, 700 Series Barb, 3/8" (9.5 mm) Tubing ID

XQ SERIES - 1/4" FLOW SIZE:
Male Open Flow Panel Mount



Reference part shown: XQCBM755

XQCBM755 (-1006-B)

XQ Series Male Open Flow Panel Mount Connector with 11/16-24 UNEF Thread, 700 Series Barb, 1/4" (6.4 mm) Tubing ID, Nickel-Plated Brass Hex Nut (included)

XQCBM770 (-1006-B)

XQ Series Male Open Flow Panel Mount Connector with 11/16-24 UNEF Thread, 700 Series Barb, 3/8" (9.5 mm) Tubing ID, Nickel-Plated Brass Hex Nut (included)

XQ SERIES - 1/4" FLOW SIZE:
Male Open Flow Straight



Reference part shown: XQCM755

XQCM755 (-1006-B)

XQ Series Male Open Flow Connector, 700 Series Barb, 1/4" (6.4 mm) Tubing ID

XQCM770 (-1006-B)

XQ Series Male Open Flow Connector, 700 Series Barb, 3/8" (9.5 mm) Tubing ID

XQ SERIES - 1/4" FLOW SIZE:
Male Open Flow Elbow



Reference part shown: XQCLM755

XQCLM755 (-1006-B)

XQ Series Male Open Flow Elbow Connector, 700 Series Barb, 1/4" (6.4 mm) Tubing ID

XQCLM770 (-1006-B)

XQ Series Male Open Flow Elbow Connector, 700 Series Barb, 3/8" (9.5 mm) Tubing ID

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!

www.nordsonmedical.com

XQ SERIES - 1/4" FLOW SIZE:
Female Valved Straight



Reference part shown: **XQSVF755**

XQSVF755 (-1006-B -1007-B -6006-EP)
XQ Series Female Valved Connector, 700 Series Barb,
1/4" (6.4 mm) Tubing ID

XQSVF765 (-1006-B -1007-B)
XQ Series Female Valved Connector, 700 Series Barb,
5/16" (8.0 mm) Tubing ID

XQSVF770 (-1006-B -1007-B -6006-EP)
XQ Series Female Valved Connector, 700 Series Barb,
3/8" (9.5 mm) Tubing ID

XQ SERIES - 1/4" FLOW SIZE:
Male Valved Panel Mount



Reference part shown: **XQCVBM755**

XQCVBM755 (-1006-B)
XQ Series Male Valved Panel Mount Connector with 11/16-24 UNEF
Thread, 700 Series Barb, 1/4" (6.4 mm) Tubing ID, Nickel-Plated Brass
Hex Nut (included)

XQCVBM770 (-1006-B)
XQ Series Male Valved Panel Mount Connector with 11/16-24 UNEF
Thread, 700 Series Barb, 3/8" (9.5 mm) Tubing ID, Nickel-Plated Brass
Hex Nut (included)

XQ SERIES - 1/4" FLOW SIZE:
Male Valved Straight



Reference part shown: **XQCVM755**

XQCVM755 (-1006-B -6006-EP)
XQ Series Male Valved Connector, 700 Series Barb, 1/4" (6.4 mm) Tubing ID

XQCVM765 (-1006-B -1007-B)
XQ Series Male Valved Connector, 700 Series Barb, 5/16" (7.9 mm) Tubing ID

XQCVM770 (-1006-B -1007-B -6006-EP)
XQ Series Male Valved Connector, 700 Series Barb, 3/8" (9.5 mm) Tubing ID

XQ SERIES - 1/4" FLOW SIZE:
Male Valved Elbow



Reference part shown: **XQCVLM755**

XQCVLM755 (-1006-B)
XQ Series Elbow Male Valved Connector, 700 Series Barbs,
1/4" (6.4 mm) Tubing ID

XQ SERIES - 1/4" FLOW SIZE:
Female Valved Panel Mount



Reference part shown: **XQSVBF755**

XQSVBF755 (-1006-B)
XQ Series Female Valved Panel Mount Connector with 11/16-24
UNEF Thread, 700 Series Barb, 1/4" (6.4 mm) Tubing ID, Nickel-Plated
Brass Hex Nut (included)

XQSVBF770 (-1006-B)
XQ Series Female Valved Panel Mount Connector with 11/16-24
UNEF Thread, 700 Series Barb, 3/8" (9.5 mm) Tubing ID, Nickel-Plated
Brass Hex Nut (included)

XQ SERIES - 1/4" FLOW SIZE:
Panel Mount Nut

PM4N (Nickel-Plated Brass)
XQ Series 11/16-24 UNEF Panel Mount Nut



CUSTOM

Solutions

Nordson MEDICAL has developed hundreds of custom connection solutions for challenging applications. Our engineering staff is ready to work with your designers to solve application problems beyond the scope of our standard product offering. Whether it is a simple o-ring, color change, or a completely new and unique design, Nordson MEDICAL offers its customers the best possible solution. To submit a new product request, email us at designcenter@nordsonmedical.com.

CUSTOM SOLUTIONS: Multi-port Couplings



CUSTOM SOLUTIONS: Unique Valved Couplings



CUSTOM SOLUTIONS: Value Added Services



CUSTOM SOLUTIONS: Molded Color Couplings



CUSTOM SOLUTIONS: Non-Spill Capabilities



CUSTOM SOLUTIONS: Non-Standard, High Performance Materials



CUSTOM SOLUTIONS: Custom Mounting Options



CHECK

Valves

Anti-Siphon Check Valves

PG-84




High Flow Check Valves

PG-85



One-Way Check Valves

PG-85

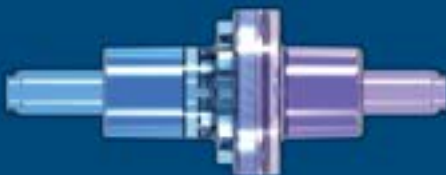

Barbed Check Valves

PG-86



Double Check Valves

PG-87



SCV Series

Check Valves

SCV Series products are available in barbed, anti-siphon or one-way check valve configurations in a variety of materials and end fittings. Barbed check valves reduce the risk of contamination and feature more precise, lower crack pressures and higher flow rates. They are ideal for medical device designers focused on infusion, drainage and irrigation applications. All the valves feature a low priming volume, ease of removing air bubbles during priming and operate in any spatial orientation. All materials of construction satisfy USP Class VI and ISO 10993 criteria.

Features & Benefits:

- Materials meet USP Class VI and ISO 10993 criteria
- DEHP- and latex-free silicone diaphragms
- All products are suitable for EtO or gamma sterilization
- Several connection configurations and materials
- Air bubbles completely removable during priming
- Operate in any spatial orientation
- Normally closed valve positioning for optimum flow control
- Available with the lowest cracking pressure on the market
- Air entrapment readily alleviated during priming

SCV SERIES:

Anti-Siphon Check Valves

SCV06018 (Clear and White MABS w/Silicone Diaphragm)

Check Valve, Socket Fitment for 0.095" (2.4 mm) OD Tubing



Cracking Pressure: $\leq .100$ psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 110 ml/min
Flow Direction: Left to right

SCV06246 (Clear Transparent MABS and Blue Transparent SAN w/Silicone Diaphragm)

Check Valve, Pocket for .122" (3.1 mm) OD Tubing to Male Locking Luer



Cracking Pressure: 2.9 ± 0.725 psig
Max. Back Pressure: 304.5 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV06248 (Blue SAN and Clear MABS w/Silicone Diaphragm)

Check Valve, Pocket for .161" (4.0 mm) OD Tubing to Male Locking Luer



Cracking Pressure: 2.176-3.626 psig
Max. Back Pressure: 304.579 psig
max 30 sec.
Flow Rate: 170 ml/min
Flow Direction: Left to right

SCV06249 (Blue Transparent and Transparent MABS w/Silicone Diaphragm)

Check Valve, Pocket for .165" (4.2 mm) OD Tubing to Male Locking Luer



Cracking Pressure: 2.9 +/- 0.725 psig
Max. Back Pressure: 304.5 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV06252 (Clear Transparent SAN and Blue Transparent SAN w/Silicone Diaphragm)

Check Valve, Female Locking Luer to Male Locking Luer



Cracking Pressure: 2.9 +/- 0.725 psig
Max. Back Pressure: 304.5 psig
Flow Rate with Glucose: ≥ 150 ml/min
Flow Direction: Left to right

SCV06264 (Clear and Blue MABS w/Silicone Diaphragm)

Check Valve, Pocket for .161" (4.0 mm) OD Tubing to .161" (4.0 mm) OD Tubing



Cracking Pressure: 1.450-4.351 psig
Max. Back Pressure: 304.579 psig
Flow Rate: max 200 ml/min
Flow Direction: Left to right

SCV06272 (Transparent SAN and Blue Transparent MABS w/Silicone Diaphragm)

Check Valve, Pocket for .161" (4.1 mm) OD Tubing to Female Locking Luer



Cracking Pressure: 2.9 +/- 0.725 psig
Max. Back Pressure: 304.5 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV06308 (Clear Transparent Polycarbonate and Green Polycarbonate w/Silicone Diaphragm)

Check Valve, Female Locking Luer to Pocket for .098" (2.5 mm) OD Tubing



Cracking Pressure: 2.9 +/- 0.725 psig
Max. Back Pressure: 304.5 psig
Flow Rate with Glucose: ≥ 150 ml/min
Flow Direction: Left to right

SCV06329 (Clear and Blue SAN w/Silicone Diaphragm)

Check Valve, Female Locking Luer to Linden



Cracking Pressure: 1.450-4.351 psig
Max. Back Pressure: 304.579 psig
max 30 sec.
Flow Rate: 200 ml/min
Flow Direction: Left to right

SCV06397 (Clear MABS w/Silicone Diaphragm)

High Flow Check Valve, Pocket for 8.8 mm OD Tubing to Pocket for 8.8 mm OD Tubing



Cracking Pressure: < .29 psig
Max. Back Pressure: 14,500 psig
Flow Rate with Water at 100 mbar:
> 1000 ml/min
Flow Direction: Left to right

SCV06398 (Clear MABS w/Silicone Diaphragm)

High Flow Check Valve, Pocket for 6.8 mm OD Tubing to Pocket for 6.8 mm OD Tubing



Cracking Pressure: < .29 psig
Max. Back Pressure: 14,500 psig
Flow Rate with Water at 100 mbar:
> 1000 ml/min
Flow Direction: Left to right

SCV06265 (Blue and Clear Polycarbonate w/Silicone Diaphragm)

Check Valve, Autoclavable, Socket for .122" (3.1 mm) ID x .165" (4.2 mm) OD Tubing x 2



Cracking Pressure: ≤ .087 psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV08053 (White SAN w/Silicone Diaphragm w/Natural LDPE)

Check Valve, Female Locking Luer to Male Locking Luer w/Cap



Cracking Pressure: ≤ .174 psig
Max. Back Pressure: 87.023 psig
Flow Rate: ≥ 90 ml/min
Flow Direction: Left to right

SCV10090 (Clear Polycarbonate w/Silicone Diaphragm)

Check Valve, Female Locking Luer to Pocket for .118" (3.0 mm) OD Tubing



Cracking Pressure: ≤ .174 psig
Max. Back Pressure: 87.023 psig
Flow Rate: ≥ 90 ml/min
Flow Direction: Left to right

SCV17090 (Clear Polycarbonate w/Silicone Diaphragm)

Check Valve, Female Locking Luer to Socket for .118" (3.0 mm) ID x .161" (4.1 mm) OD Tubing



Cracking Pressure: ≤ .174 psig
Max. Back Pressure: 87.023 psig
Flow Rate: ≥ 90 ml/min
Flow Direction: Left to right

SCV21053 (Clear and White MABS w/Silicone Diaphragm)

Check Valve, Socket for .118" (3.0 mm) ID x .161" (4.1 mm) OD Tubing x 2



Cracking Pressure: ≤ .174 psig
Max. Back Pressure: 87.023 psig
Flow Rate: ≥ 90 ml/min
Flow Direction: Left to right

SCV22050 (Clear SAN and Clear MABS w/Silicone Diaphragm)

Check Valve, Socket for .118" (3.0 mm) ID x .161" (4.1 mm) OD to Male Locking Luer



Cracking Pressure: ≤ .174 psig
Max. Back Pressure: 87.023 psig
Flow Rate: ≥ 90 ml/min
Flow Direction: Left to right

SCV SERIES:
One-Way Check Valves

SCV23050 (Clear SAN w/Silicone Diaphragm)

Check Valve, Female Locking Luer to Male Locking Luer



Cracking Pressure: ≤ .174 psig
Max. Back Pressure: 87.023 psig
Flow Rate: ≥ 90 ml/min
Flow Direction: Left to right

VPS5401001N (ABS Housing with Silicone Disk)

Check Valve, Male Slip Luer, 420 mL/min Flow Rate at 1M Head, 30 psi Back Pressure, 1/4 psi Inlet Cracking Pressure



VPS5401020N (Polycarbonate with Silicone Disk)

Check Valve, Male/Female Luer, 129 mL/min Flow Rate at 100 mmHg, 45 psi Back Pressure, 1-5 psi Inlet Cracking Pressure



VPS5401036SN (Polycarbonate with Silicone Disk)

Reflux Check Valve, Male/Female Luer, 350 mL/min Flow Rate at 1 M Head, 45 psi Back Pressure, Opens at 284 mmHg



VPS5401068N (Polycarbonate with Silicone Disk)

Check Valve, Male/Female Luer, 420 mL/min Flow Rate at 1 M Head 30 psi Back Pressure, 1/4 psi Inlet Cracking Pressure



SCV SERIES:
Barbed Check Valves

SCV67FL/220 (SAN and MABS w/Silicone Diaphragm)

Check Valve, Female Locking Luer to Barb for 3/32" (2.4 mm) ID Tubing



Cracking Pressure: ≤ .087 psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV67FL/230 (SAN and MABS w/Silicone Diaphragm)

Check Valve, Female Locking Luer to Barb for 1/8" (3.2 mm) ID Tubing



Cracking Pressure: ≤ .087 psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV67220 (MABS w/Silicone Diaphragm)

Check Valve, Socket for 3/32" (2.4 mm) x 3/32" (2.4 mm) ID Tubing



Cracking Pressure: ≤ .087 psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV67220/230 (MABS w/Silicone Diaphragm)

Check Valve, Socket for 3/32" (2.4 mm) x 1/8" (3.2 mm) ID Tubing



Cracking Pressure: ≤ .087 psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV67220/ML (SAN and MABS w/Silicone Diaphragm)

Check Valve, Socket for 3/32" (2.4 mm) ID Tubing to Male Locking Luer



Cracking Pressure: ≤ .087 psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV67230 (MABS w/Silicone Diaphragm)

Check Valve, Socket for 1/8" (3.2 mm) x 1/8" (3.2 mm) ID Tubing



Cracking Pressure: ≤ .087 psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV67230/220 (MABS w/Silicone Diaphragm)

Check Valve, Socket for 1/8" (3.2 mm) x 3/32" (2.4 mm) ID Tubing



Cracking Pressure: ≤ .087 psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

SCV67230/ML (SAN and MABS w/Silicone Diaphragm)

Check Valve, Socket for 1/8" (3.2 mm) ID Tubing to Male Locking Luer



Cracking Pressure: ≤ .087 psig
Max. Back Pressure: 73 psig
Flow Rate: ≥ 150 ml/min
Flow Direction: Left to right

CHECK
VALVES

DCV114 (-001) Radiation-Stable Polycarbonate

Double Check Valve, Tubing Pocket Chimney Port

Cracking Pressure: 2-5 psi



DCV115 (-001) Radiation-Stable Polycarbonate

Double Check Valve, Tubing Pocket Chimney Port, Rotating Luer Lock Exit Port

Cracking Pressure: 2-5 psi



DCV SERIES:

Double Check Valves, Female Luer Chimney Port

DCV118 (-001) Radiation-Stable Polycarbonate

Double Check Valve, Female Luer Chimney Port

Cracking Pressure: 2-5 psi



DCV125 (-001) Radiation-Stable Polycarbonate

Double Check Valve, Female Luer Chimney Port, Rotating Luer Lock Exit Port

Cracking Pressure: 2-5 psi



DCV Series


Check Valves

Nordson MEDICAL's unique DCV Series bi-directional double check valve is designed to transfer clean fluid from a supply vessel to a use site. The fluid supply vessel is connected to the valve supply port (chimney) using a luer connection or tubing. The fluid is withdrawn from the supply vessel by a syringe or other device connected to the aspiration port. When the syringe is compressed the fluid is then transferred, through the exit port, to the use site without adulterating the fluid. DCV Series bi-directional valves are produced in polycarbonate resin using a silicone diaphragm, with two available inlet ports connecting to the female luer aspiration and male luer discharge ports.


Features & Benefits:

- Bi-directional double check valve
- No need for disconnection to fill syringe
- Satisfies USP Class VI criteria
- Low cracking pressure
- Clean fluid transfer from supply to site
- High flow rate

Step 1: Clean fluid drawn into syringe

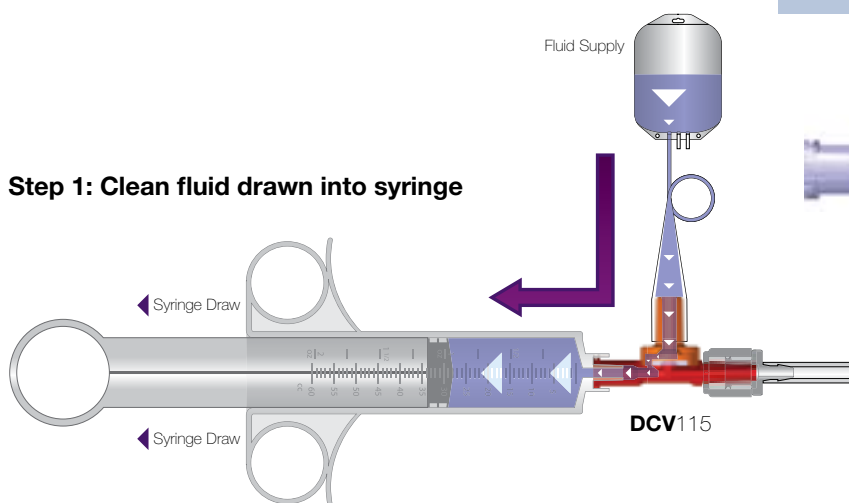

Step 2: Clean fluid expelled to use site

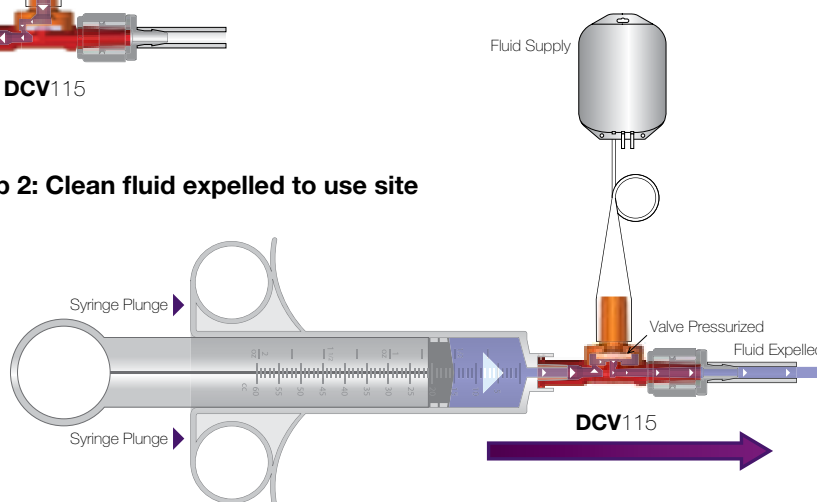

NEEDLEFREE

Swabable Valves



Vial Caps
PG-88

Tube End Valves
PG-89



Swabable Bag Valves
PG-89

Valved Male Luers
PG-89



Swabable Break-Off Valves
PG-89

4-Way Swabable Stopcock
PG-89



Swabable 'Y' Valves
PG-89

Swabable Straight Valves
PG-89



Bondable Valves
PG-89

Swabable T-Port Valves
PG-89



NV Series

Needlefree Swabable Valves

Nordson MEDICAL now offers cost efficient needlefree swabable valves featuring higher flow rates, low priming volumes, and a straight through design for unobstructed flow. The valves don't require caps, and the septum surface can be easily disinfected. Designed for the injection and aspiration of fluids, particularly for IV applications, cardiology, and catheters, the hermetic septum seal reduces the risk of biofilm formation and contamination. Available in widely used polycarbonate and copolyester materials, these valves are compatible with all industry standard luer syringes and connects, validating them as an efficient component solution for the OEM market.

Features & Benefits:

- Valve stems and bodies securely mate with all standard luer syringes and connectors
- Several configurations that are easy to customize
- Normally closed system to reduce the risk of infection
- Low priming volume and high flow rates
- Straight through design for unobstructed flow
- Gamma radiation stable and EtO compatible
- Not made with natural rubber latex or DEHP

The NV Series products are shipped bulk, non-sterile and are single patient use medical device components requiring further processing (e.g. assembly, packaging, sterilization) before clinical use. The buyer is responsible for determining effects of processing/multiple usage on these components, appropriateness of the component in the final application, and pre/post shelf life.

NV SERIES:

Vial Caps



Reference part shown:

NVF13VPC

NVF13VPC (Polycarbonate)

13 mm Swabable Vial Cap, Female Luer to 13 mm Vial Cap

NVF20VPC (Polycarbonate)

20 mm Swabable Vial Cap, Female Luer to 20 mm Vial Cap

NV SERIES: Tube End Valves



Reference part shown:
NVF20PC

NVF20PC (Polycarbonate)

Tube End Valve, Female Luer to 2.0 mm Tube Socket

NVF25PC (Polycarbonate)

Tube End Valve, Female Luer to 2.5 mm Tube Socket

NVF28PC (Polycarbonate)

Tube End Valve, Female Luer to 2.8 mm Tube Socket

NVF37PC (Polycarbonate)

Tube End Valve, Female Luer to 3.7 mm Tube Socket

NVF41PC (Polycarbonate)

Tube End Valve, Female Luer to 4.1 mm Tube Socket

NV SERIES: Swabable Bag Valves

NVFBGPC (Polycarbonate)

Swabable Bag Valve, Female Luer to Bag Access Port



NV SERIES: Valved Male Luers



Reference
part shown:
NVM31S-001

NVM31S-001 (Polycarbonate)

Valved Male Luer to .125" Barb for 3.1 mm ID Tubing, Swabable Stem

NVM41S-001 (Polycarbonate)

Valved Male Luer to 4.0 / 4.1 mm Bondable Socket

NV SERIES: Swabable Break-Off Valves



Reference part shown:
NVFBV60PC

NVFBV60PC (Polycarbonate)

Swabable Break-Off Valve, Female Luer to Break Off 6.0 mm Port

NVFBV66PC (Polycarbonate)

Swabable Break-Off Valve, Female Luer to Break Off 6.6 mm Port

NV SERIES: 4-Way Swabable Stopcock

NVFSCMLLPC (Polycarbonate)

4-Way Lever Swabable Stopcock, Female Luer Stopcock to Male Luer Lock



NV SERIES: Swabable 'Y' Valves



Reference part shown:
NV FY36PC

NV FY36PC (Polycarbonate)

Swabable 'Y' Valve, Female Luer Y Port to 3.6 mm Tube Sockets

NV FY42S-001 (Polycarbonate)

Swabable 'Y' Valve, Female Luer Y Port to 4.2 mm Bondable Sockets

NV FY401CP (Copolyester)

Swabable 'Y' Valve, Female Luer Y Port to 4.01 mm Tube Sockets

NV FY401PC (Polycarbonate)

Swabable 'Y' Valve, Female Luer Y Port to 4.01 mm Tube Sockets

NV SERIES: Swabable Straight Valves



Reference part shown:
NV FM LLS-001

NV FM L L-001 (Polycarbonate)

Luer Activated Straight Valve, Female Luer to Male Luer Lock

NV FM L LCP (Copolyester)

Swabable Straight Valve, Female Luer to Male Luer Lock

NV FM L LPC (Polycarbonate)

Swabable Straight Valve, Female Luer to Male Luer Lock

NV FM L L S-001 (Polycarbonate)

Swabable Straight Valve, Female Luer to Male Luer Lock

NV SERIES: Bondable Valves



Reference part shown:
NV FM L PC

NV FM L PC (Polycarbonate)

Bondable Valve, Female Luer to Male Luer Slip

NV FM L PCB (Polycarbonate)

Bondable Valve, Female Luer to Male Luer Slip, Blue Cap

NV SERIES: Swabable T-Port Valves

NV FT PC (Polycarbonate)

Swabable T-Port Valve, Female Luer T Port to 6.6 mm Tube Sockets



STOPCOCK

FITTINGS

High Pressure Stopcocks
PG-90



Medium Pressure Stopcocks
PG-91



Low Pressure Stopcocks
PG-91



High Flow Stopcocks
PG-93



DIVERT™

Stopcocks

Diverter™ stopcocks are designed to easily control the flow of a liquid or gas through the use of an on/off diverter, and are available in three basic configurations – 1-way, 3-way, and 4-way. These chemically resistant stopcocks bring efficiency, convenience, and practicality to your tubing applications.

Features & Benefits:

- 1-way, 3-way, and 4-way configurations
- Chemically resistant
- Not made with DEHP
- Diverter arrows provide a clear indication of flow
- Custom solutions available, including manifolds, ganged manifolds, and tubing assemblies
- Standard flow and high flow options
- Polycarbonate, PEI, nylon, and HDPE materials

Polycarbonate is known to experience environmental stress cracking when exposed to solvents. Please contact Nordson MEDICAL technical support if polycarbonate is planning on being used in a solvent bonding application. Looking for a stress-free connection? Ask our technical experts at 1-888-404-5837 how annealing and other secondary processes can alleviate the risk of environmental stress cracking in polycarbonate materials.

DIVERT STOPCOCKS: High Pressure Stopcocks

FCHPFFMR-003

3-Way High Pressure Stopcock, Female Luer to Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body and Lock Ring, 4-Arrow Gold Color Acetal Diverter



DIVERT STOPCOCKS: Medium Pressure Stopcocks

FCHPFFMR-001

3-Way High Pressure Stopcock, Female Luer to Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body and Lock Ring, 4-Arrow Gold Color Acetal Diverter



FCHPFFMR-002

3-Way High Pressure Stopcock, Female Luer to Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body and Lock Ring, 4-Arrow Grey Color Acetal Diverter



FCHPFFMR-002

1-Way High Pressure Stopcock, Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body and Lock Ring, 2-Arrow Gold Color Acetal Diverter



FCHPFFMR-003

1-Way High Pressure Stopcock, Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body and Lock Ring, 2-Arrow Gold Acetal Diverter



FCHPFFMR-005

1-Way High Pressure Stopcock, Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body and Lock Ring, 2-Arrow Grey Color Acetal Diverter



DIVERT STOPCOCKS: Low Pressure Stopcocks

FC4FFMR-001

4-Way Stopcock, Female Luer to Female Luer to Male Luer with Luer Lock Ring, Natural Nylon Body, Tinted Polycarbonate Lock Ring, 1-Arrow White HDPE Diverter



FCFFM-002

3-Way Stopcock, Female Luer to Female Luer to Bondable Socket, Clear Polycarbonate Body, 1-Arrow Red Color HDPE Diverter



FCFFM-003

3-Way Stopcock, Female Luer to Female Luer to Bondable Socket, Clear Polycarbonate Body, 1-Arrow Natural HDPE Diverter



FCFFM-008

3-Way Stopcock, Female Luer to Female Luer to Male Luer, Clear Polycarbonate Body, 1-Arrow Natural HDPE Diverter



FCFFM-009

3-Way Stopcock, Female Luer to Female Luer to Male Luer, Tinted Polycarbonate Body, 1-Arrow White HDPE Diverter



DIVERT STOPCOCKS:
Low Pressure Stopcocks

FCFFM-010

3-Way Stopcock, Female Luer to Female Luer to Male Luer, Clear Polycarbonate Body, 1-Arrow Natural HDPE Diverter



FCFFM-013

3-Way Stopcock, Female Luer to Female Luer to Male Luer, Tinted Polycarbonate Body, 1-Arrow Natural HDPE Diverter



FCFFMR-002

3-Way Stopcock, Female Luer to Female Luer to Male Luer with Luer Lock Ring, Tinted Polycarbonate Body and Lock Ring, 1-Arrow White HDPE Diverter



FCFFMR-003

3-Way Stopcock, Female Luer to Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body and Lock Ring, 1-Arrow Natural HDPE Diverter



FCFFMR-004

3-Way Stopcock, Female Luer to Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body and Lock Ring, 1-Arrow Natural HDPE Diverter



FCFM-001

1-Way Stopcock, Female Luer to Male Luer, Clear Polycarbonate Body, 2-Arrow Natural HDPE Diverter



FCFM-002

1-Way Stopcock, Female Luer to Male Luer, Clear Polycarbonate Body, 2-Arrow Blue HDPE Diverter



FCFM-006

1-Way Stopcock, Female Luer to Male Luer, Tinted Polycarbonate Body, 2-Arrow Natural HDPE Diverter



FCFM-010

1-Way Stopcock, Female Luer to Male Luer, Tinted Polycarbonate Body, 2-Arrow Natural POM Diverter



FCFML-2A-001

1-Way Stopcock, Female Luer to Male Luer with Positive Stop, Clear Polycarbonate Body, 2-Arrow Blue HDPE Diverter



STOPCOCK
FITTINGS

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

DIVERT STOPCOCKS:

Low Pressure Stopcocks

FCFMR-001

1-Way Non-Sterile Stopcock, Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body, 2-Arrow Natural HDPE Diverter

**FCFMR-2A-001**

1-Way Stopcock, Female Luer to Male Luer with Luer Lock Ring, Tinted Polycarbonate Body, 2-Arrow Natural HDPE Diverter

**FCFMR-2A-004**

1-Way Stopcock, Female Luer to Male Luer with Luer Lock Ring, Clear Polycarbonate Body, 2-Arrow Natural HDPE Diverter

**DIVERT STOPCOCKS:**

High Flow Stopcocks

FC4LB3F-001

4-Way Stopcock, Female Luer on All Ports, Large Bore, Tinted Polycarbonate Body, 1-Arrow White HDPE Diverter

**FCLBFFM-001**

3-Way High Flow Stopcock, Female Luer to Female Luer to Male Luer, Large Bore, Natural Nylon Body, 1-Arrow White HDPE Diverter

**FCLBFFM-1A-003**

3-Way High Flow Stopcock, Female Luer to Female Luer to Male Luer, Large Bore, Clear Polycarbonate Body, 1-Arrow Natural HDPE Diverter

**FCLBFFM-003**

3-Way High Flow Stopcock, Female Luer to Female Luer to Bondable Socket, Large Bore, Clear Polycarbonate Body, 1-Arrow Natural HDPE Diverter

**FCLBFFMR-1A-003**

3-Way High Flow Stopcock, Female Luer to Female Luer to Male Luer with Luer Lock Ring, Large Bore, Clear Polycarbonate Body and Lock Ring, 1-Arrow Natural HDPE Diverter

**FCLBFMR-001**

1-Way High Pressure Stopcock, Female Luer to Male Luer with Luer Lock Ring, Large Bore, Clear Polycarbonate Body and Lock Ring, 2-Arrow Natural HDPE Diverter

STOPCOCK
FITTINGS**Need Samples for Evaluation?**

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

BIOPROCESS

Fittings

Bag Ports
PG-94



Opposing Barb Ports
PG-98



BioClamp®
PG-99



BioValve®
PG-100



PharmaLok™ Clamps
PG-101



PharmaLok™ Tube Clamps
PG-104



Spaulding™ Aseptic Disconnects
PG-105



CYLINDRAFlow™ Manifolds
PG-107



Sanitary Fittings
PG-108



Instrument Tees
PG-111



Gaskets
PG-112



BAG Ports

Nordson MEDICAL's new high tech series of bag ports offers a variety of benefits new to the biopharmaceutical industry through a unique design. Value Pharma™ resin holds credentials that many manufacturers deem fundamental, such as animal derivative free qualification, USP Class VI certification, and the guarantee of heat sealing and conformance with polyethylene bags. The resin then goes even further to provide a strong, 24-hour resistance to alcohol stress cracking. RENOLIT SOLMED Granuflex 4301 resin demonstrates lower leachable levels and high purity.

Their unique design sets these ports apart from existing solutions. The most notable of the bag ports' specifications are the parabolic lead and the alignment ribs. These two features, respectively, reduce the amount of time it takes to drain a bag by up to 24 percent and improve flow dynamics through the port.

Features & Benefits:

Compliant to ASME BPE Ferrule Specification

- Works with leading industry single-use bags

Barb Sizes Consist of Range From 1/8" (3.2 mm) to 1" (25.4 mm)

- Works with the frequently utilized tubing sizes found in the pharmaceutical industry

Sanitary TC Ports Available in 1-1/2" (38.1 mm) and 3" (76.2 mm)

- Allows the option to connect directly to sanitary TC fittings

Opposing Barb Ports Available in 3/16" (4.8 mm) to 3/4" (19.0 mm) Hose Barb

- Ideal for dip tubes

Parabolic Lead

- Reduces the time it takes to drain a bag by 17 to 24%
- Similar results for both fluid and powder applications

Alignment Ribs

- Improves fluid dynamics of flow through the port

Value Pharma™ Resin

- Ensures heat sealing with polyethylene bags
- Includes a USP Class VI certification
- Animal derivative free
- Resistant to alcohol stress cracking

RENOLIT SOLMED Granuflex 4301 Resin

- Demonstrates lower leachable levels and high purity

BAG PORTS: 600 Series Barb Style

SFBP630 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene

Bag Port with 600 Series Barb, 1/8" (3.2 mm) Tubing ID



SFBP655 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene

Bag Port with 600 Series Barb, 1/4" (6.4 mm) Tubing ID



SFBP670 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene

Bag Port with 600 Series Barb, 3/8" (9.5 mm) Tubing ID



BAG PORTS: 600 Series Barb Style

SFBP680 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene

Bag Port with 600 Series Barb, 1/2" (12.7 mm) Tubing ID



SFBP6100 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene

Bag Port with 600 Series Barb, 3/4" (19.0 mm) Tubing ID



SFBP6110 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene

Bag Port with 600 Series Barb, 1" (25.4 mm) Tubing ID



Need More Details?

Check out all of our Bag Ports at www.nordsonmedical.com
or contact Sales at FMF-sales@nordsonmedical.com.

Need it Fast?

You can order any of our 9,600+ parts online now with no
order \$ minimums, 24 hours a day, 7 days a week!
www.nordsonmedical.com

BAG PORTS:

Ferrule Fitting to 1-1/2" TC Connection

SFBPF1 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene

Bag Port with Ferrule Fitting to 1-1/2" TC Connection



Innovative Parabolic Lead &
Alignment Ribs

BAG PORTS:

4" TC Connection

SFBPF4 (-VP1) Value Pharma™ Resin, Bag Port with Ferrule Fitting to 4" TC Connection



BAG PORTS:

3" TC Connection

SFBPF3 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene

Bag Port with Ferrule Fitting to 3" TC Connection



BIOPROCESSING
FITTINGS

Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.
www.nordsonmedical.com

SFM3C (-6005)

Flat TC Style End Cap, 3", Works with SFBPF3



SFM4C (-6005)

Flat TC Style End Cap, 4" (101.6 mm) Flange



Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.
+1 (970) 267-5200 Toll Free: (888) 404-5837

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

OPPOSING

Barb Ports

Nordson MEDICAL now offers the biomedical market opposing barb ports. The larger ports can be used in filling applications, and they are also ideally suited to work with sensors.

Features & Benefits:

- Easy to heat seal into film material
- Can be used in filling applications
- Port plates featuring Value Pharma™ or Granuflex resin
- Works great with silicone and TPE style tubing
- Better tube retention
- 600 Series barb sizes
- Animal-free polyethylene
- Materials meet USP Class VI criteria and are ADF compliant
- Fits tubing sizes from 3/16" to 3/4"

OPPOSING BARB PORTS:
600 Series Barb Style

SFBP635/635 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene
Dip Tube (Opposable Barbs) with Two 600 Series Barbs, 3/16" (4.8 mm) Tubing ID




SFBP655/655 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene
Dip Tube (Opposable Barbs) with Two 600 Series Barbs, 1/4" (6.4 mm) Tubing ID




SFBP670/670 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene
Dip Tube (Opposable Barbs) with Two 600 Series Barbs, 3/8" (9.5 mm) Tubing ID



SFBP680/680 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene
Dip Tube (Opposable Barbs) with Two 600 Series Barbs, 1/2" (12.7 mm) Tubing ID

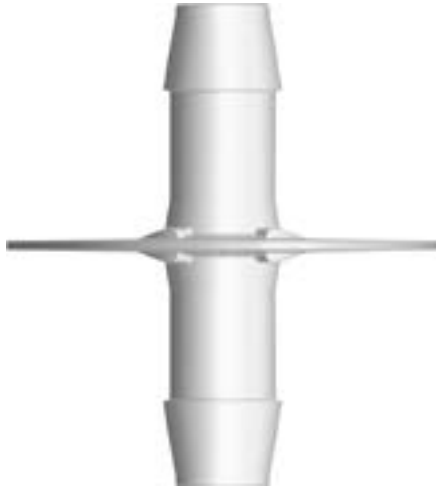

SFBP6100/6100 (-VP1 -CM056) Value Pharma™ Resin, Granuflex Polyethylene
Dip Tube (Opposable Barbs) with Two 600 Series Barbs, 3/4" (19.0 mm) Tubing ID

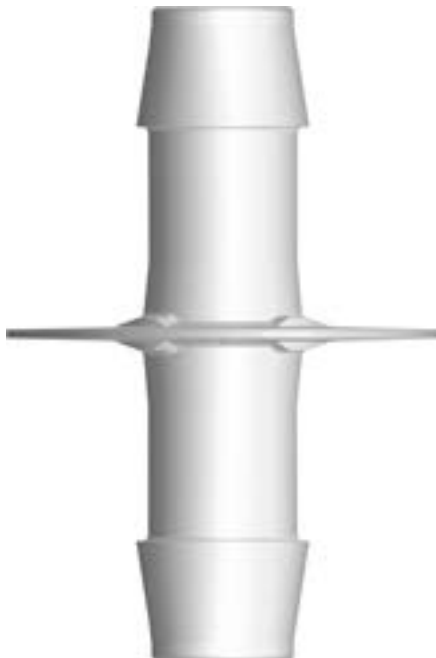

BIOPROCESSING
FITTINGS

BIOClamp®

BioClamp® plastic quick-release clamps are specifically designed to meet the needs of busy pharmaceutical and bioprocessing laboratories, offering full flexibility in single-use applications. The clamp is lightweight and easy for technicians to install and provides for quick and convenient assembly release.

Features & Benefits:

- Can be SIP repeatedly
- Can be autoclaved repeatedly at 134°C for 5 minutes
- Can be subjected to 40 kiloGrays without weakening
- Reduced distortion on polymeric components
- More compatible with current manufacturing than steel
 - Greater flexibility in operation
 - Single-use disposal
- Applications
 - Bioprocessing
 - Single-use technology (SUS)
 - Drug discovery and development
 - Fluid transport
 - Pharmaceutical scale up
 - Sterile bag systems

Specifications:

Materials: Glass reinforced Nylon 66

Max. Operating Pressure: 150 psig for 1/2" to 4" clamps & 90 psig for 6" and 8" clamps

Max. Cumulative Gamma Dose: 40 kiloGrays

USP Class VI Certified Materials



Reference part shown: VNG150WHT

Due to size we are not able to display all large BioClamp® fittings. Please see www.nordsonmedical.com for part reference. Minimum package quantities vary by tube IDs, please verify minimum package sizes when placing orders.

VNG075WHT (Glass Reinforced White Nylon)

Clamp, Mini Flange, 1/2" (12.7 mm) to 3/4" (19.0 mm)

VNG150WHT (Glass Reinforced White Nylon)

Clamp, Maxi Flange, 1" (25.4 mm) to 1.5" (38.1 mm)

VNG200WHT (Glass Reinforced White Nylon)

Clamp, Maxi Flange, 2" (50.8 mm)

VNG250WHT (Glass Reinforced White Nylon)

Clamp, Maxi Flange, 2.5" (63.5 mm)

VNG300WHT (Glass Reinforced White Nylon)

Clamp, Maxi Flange, 3" (76.2 mm)

VNG400WHT (Glass Reinforced White Nylon)

Clamp, Maxi Flange, 4" (101.6 mm)

VNG600WH (Glass Reinforced White Nylon)

Clamp, Maxi Flange, 6" (152.4 mm)

VNG800WH (Glass Reinforced White Nylon)

Clamp, Maxi Flange, 8" (203.2 mm)

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!
www.nordsonmedical.com

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

BIOValve®

The BioValve® pinch valve is a precision restriction flow controller and shut-off valve for silicone tubing used in bioprocessing and pharmaceutical manufacturing applications.

Features & Benefits:

- Can be SIP repeatedly
- Can be autoclaved repeatedly at 134°C for 5 minutes
- Can be subjected to 40 kiloGrays without weakening
- Can be calibrated to minimize flow path turbulence
- Can be used with one hand
- Thread pitch is calibrated to 2 mm per turn

Specifications:


Materials: Glass reinforced Nylon 66

Tubing Sizes: 1/4" to 1" OD

Max. Operating Pressure: 150 psig

Max. Cumulative Gamma Dose: 40 kiloGrays

USP Class VI Certified Materials

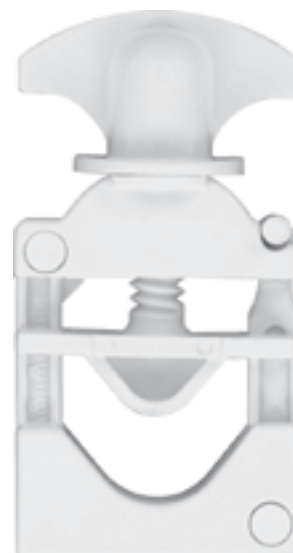

VBV1000NW (Glass Reinforced White Nylon)

Pinch Valve, Maxi, 1/2" (12.7 mm) to 1" (25.4 mm) Tubing OD



Looking for Something Different?

Most of our best product ideas have come from you, our customers. While we pride ourselves on our extensive product offering, we recognize that your needs are ever-changing and that our products, too, must adapt to your needs.

If there's a new product you think Nordson MEDICAL should make, please send us an Email at designcenter@nordsonmedical.com describing your product concept in as much detail as you desire. We will contact you with any questions, and who knows - the next catalog edition could include products YOU told us to build.

See page 10 for more information on the Design Center!



PHARMATMLok

Clamps

The biopharmaceutical industry has made a strong move towards single-use technology, so Nordson MEDICAL has developed ground-breaking clamp designs that simplify the assembly process. Effortlessly and securely attach sanitary fittings to bag ports and filters with PharmaLok™ Easy-Open Clamps, and close off flow with PharmaLok™ Tube Clamps.

Nordson MEDICAL's PharmaLok™ clamps are ASME BPE compliant. Single-use and easy-open clamp sizes include .5" to .75" and 1" to 1.5".

PharmaLok™ clamps provide a significant advancement for the industry, as a number of existing solutions either require tools for closure or consist of multiple components, which turns what should be a simple operation into a difficult and awkward process.

Easily operated with one hand, PharmaLok's patented design allows effortless attachment of fittings to bag ports and filters. The user simply rests the fitting in the PharmaLok's regulation-sized slot, lightly grips the surrounding handles, and listens for the sound of clicks to confirm a secure connection.

Single-Use Clamps & Easy-Open Clamps

Features & Benefits:

- ASME BPE 018 DT-7-1 compliant & DIN compliant
- Easy to use, one-handed closure - no tools needed
- Repeatable use with constantly equal pressure on the gasket
- Tamper evident and zip tie compatible
- 360° compression
- Connects both plastic and stainless steel fittings
- SIP/CIP compatible
- Applications
 - Bioprocessing
 - Single-use technology (SUS)
 - Drug discovery and development
 - Fluid transport
 - Pharmaceutical scale up
 - Sterile bag systems

Specifications:

Materials: Zytel® FGFE5171 glass-filled nylon certified to USP Class VI & ISO 10993 and available in white, blue, or orange

Temperature & Sterilization Specifications:

Normal operating temperature -94°F (-70°C) to 140°F (60°C)


Cryogenic freezing to -112°F (-80°C)

Autoclavable up to 273°F (134°C) for a maximum of 30 min.

Gamma radiation up to 50 kGy


Maximum Pressure: 80 psi (conforms to ASME BPE-2014 PM-4.2.2)

Sizes: Fitting .5" (12.7 mm) to .75" (19.0 mm) and 1" (25.4 mm) to 1.5" (31.8 mm)

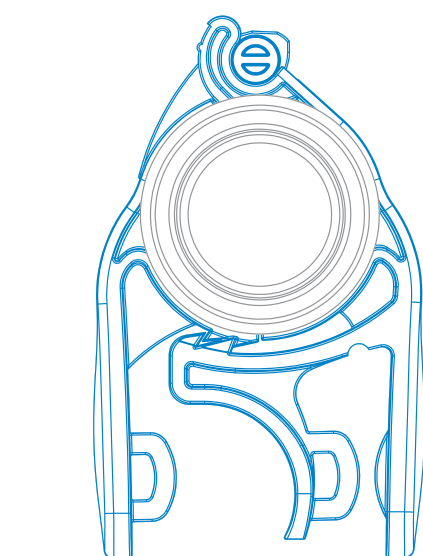

**Easy,
2 Click
Closure**

360° Compression

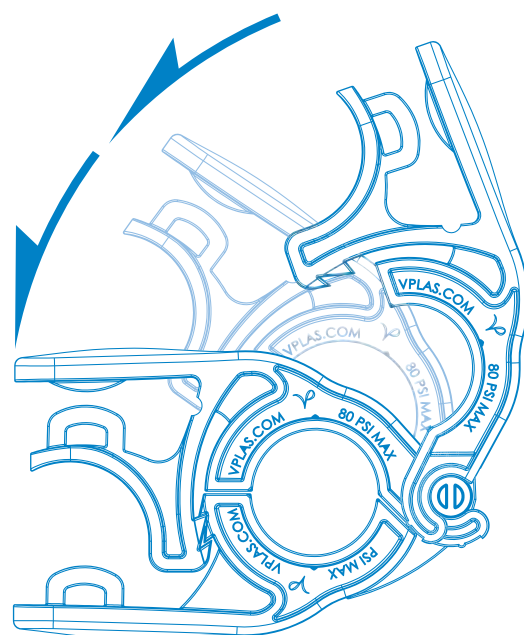

PHARMALOK SERIES:
Easy-Open Clamps

CPAX150 (-U71)

PharmaLok™ Easy-Open Clamp Series, Hygienic Union TC Version
1" (25.4 mm) to 1.5" (38.1 mm), Glass Reinforced White Nylon



CPAN75 (-U71)

PharmaLok™ Easy-Open Clamp Series, Hygienic Union TC Version
.5" (12.7 mm) to .75" (19.0 mm), Glass Reinforced White Nylon



CPAX150 (-U75)

PharmaLok™ Easy-Open Clamp Series, Hygienic Union TC Version
1" (25.4 mm) to 1.5" (38.1 mm), Glass Reinforced Blue Nylon



CPAN75 (-U75)

PharmaLok™ Easy-Open Clamp Series, Hygienic Union TC Version
.5" (12.7 mm) to .75" (19.0 mm), Glass Reinforced Blue Nylon



CPAX150 (-U706)

PharmaLok™ Easy-Open Clamp Series, Hygienic Union TC Version
1" (25.4 mm) to 1.5" (38.1 mm), Glass Reinforced Orange Nylon



CPAN75 (-U706)

PharmaLok™ Easy-Open Clamp Series, Hygienic Union TC Version
.5" (12.7 mm) to .75" (19.0 mm), Glass Reinforced Orange Nylon



BIOPROCESSING
FITTINGS

CPX150 (-U71)

Single-Use Series Clamp, TC Version 1" (25.4 mm) to 1.5" (38.1 mm)

CPN75 (-U71)

Single-Use Series Clamp, TC Version .5" (12.7 mm) to .75" (19.0 mm)



These clamps are specifically designed to work with Nordson MEDICAL's SF Series sanitary fitting Mini and Maxi flanges, as well as related caps and gaskets.

Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!

www.nordsonmedical.com

Looking for Something Different?

Most of our best product ideas have come from you, our customers. While we pride ourselves on our extensive product offering, we recognize that your needs are ever-changing and that our products, too, must adapt to your needs.

If there's a new product you think Nordson MEDICAL should make, please send us an Email at designcenter@nordsonmedical.com describing your product concept in as much detail as you desire. We will contact you with any questions, and who knows - the next catalog edition could include products YOU told us to build.

See page 10 for more information on the Design Center!



Tube Clamps

Features & Benefits:

- Designed for use with silicone and TPE tubing
- Easy to use, one-handed closure - no tools needed
- Tamper evident and zip tie compatible

Specifications:

Materials: Zytel® FGFE5171 glass-filled nylon certified to USP Class VI & ISO 10993 and available in white

Temperature & Sterilization Specifications:

Normal operating temperature -94°F (-70°C) to 140°F (60°C)

Cryogenic freezing to -112°F (-80°C)

Autoclavable up to 273°F (134°C) for a maximum of 30 min.

Gamma radiation up to 50 kGy

Maximum Pressure: Up to 150 psi

Thin-Walled Tube Clamp Size: Wall thickness 1/32" (0.79 mm) to 3/32" (2.38 mm), Max OD 1/2" (12.7 mm)

Tube Clamp Size: Wall thickness 1/8" (3.2 mm) to 1/4" (6.4 mm), Max OD 1" (25.4 mm)

Large Tube Clamp Size: Wall thickness 1/8" (3.2 mm) to 1/4" (6.4 mm), Max OD 2.5" (63.5 mm)

CPT90 (-U71 -U706)

PharmaLok™ Small and Thin-Walled Tube Clamp, Wall Thickness 1/32" (0.79 mm) to 3/32" (2.4 mm), Max OD 1/2" (12.7 mm), Glass-Filled White Nylon



CPT100 (-U71 -U706)

PharmaLok™ Tube Clamp, Wall Thickness 1/8" (3.2 mm) to 1/4" (6.4 mm), Max OD 1" (25.4 mm), Glass-Filled White Nylon



CPT110 (-U71 -U706)

PharmaLok™ Large Tube Clamp, Wall Thickness 1/8" (3.2 mm) to 1/4" (6.4 mm), Max OD 2.5" (63.5 mm), Glass-Filled White Nylon



BIOPROCESSING
FITTINGS

SPAULDING™

Aseptic Disconnects

Providing a leak-free connection, the Spaulding™ Series Aseptic Disconnects eliminate the need for pinch clamps and tube welders. Smooth disconnects are achievable with the easy-to-use design, and valuable fluids are protected by valves that automatically shut off flow when the fittings are disconnected. Available in USP Class VI, ADCF Polysulfone, validation and test results confirm their suitability in pre-sterilized systems and are also gamma radiation stable. These connectors are interchangeable with similar plastic connectors.

Features & Benefits:

- Valves that automatically shut off flow at disconnection
- Interchangeable with other plastic & metal connectors
- Plugs and caps available to maintain fluid path integrity
- Made with BSE/TSE-free materials
- Validation and test results confirm their suitability in pre-sterilized systems
- Materials satisfy USP Class VI criteria
- Fittings meet industry requirements
- Aseptic disconnect operation

Specifications:

- **Materials:** Body - USP Class VI, ADCF tinted polysulfone, o-rings - platinum-cured silicone, springs - 316 stainless steel
- **Sterilization:** Gamma radiation stable up to 50kGy, autoclavable to 270°F (132°C) for a maximum of 60 minutes
- **Fits Tubing IDs:** 1/4" (6.4 mm), 3/8" (9.5 mm), 1/2" (12.7 mm)
- **Operating Temperature:** -40°F (-40°C) to 280°F (138°C)
- **Leak Resistance:** 125 psi (8.62 bar)

SDAVF655NPS (-001)

Spaulding™ Aseptic Disconnect, Female Valved Connector, 600 Series Barb, 1/4" (6.4 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



SDAM655NPS (-001)

Spaulding™ Aseptic Disconnect, Male Open Flow Connector, 600 Series Barb, 1/4" (6.4 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



SDAVM655NPS (-001)

Spaulding™ Aseptic Disconnect, Male Valved Connector, 600 Series Barb, 1/4" (6.4 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



Need More Details?

Check out all of our bioprocessing components at www.nordsonmedical.com or contact Sales at FMF-sales@nordsonmedical.com.

SPAULDING SERIES:
Aseptic Disconnects

SDAVF670NPS (-001)

Spaulding™ Aseptic Disconnect, Female Valved Connector, 600 Series Barb, 3/8" (9.5 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



SDAVF680NPS (-001)

Spaulding™ Aseptic Disconnect, Female Valved Connector, 600 Series Barb, 1/2" (12.7 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



SDAM670NPS (-001)

Spaulding™ Aseptic Disconnect, Male Open Flow Connector, 600 Series Barb, 3/8" (9.5 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



SDAM680NPS (-001)

Spaulding™ Aseptic Disconnect, Male Open Flow Connector, 600 Series Barb, 1/2" (12.7 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



SDAVM670NPS (-001)

Spaulding™ Aseptic Disconnect, Male Valved Connector, 600 Series Barb, 3/8" (9.5 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



SDAVM680NPS (-001)

Spaulding™ Aseptic Disconnect, Male Valved Connector, 600 Series Barb, 1/2" (12.7 mm) ID Tubing, Animal-Free Polysulfone, Silicone O-Ring



SDAFCNPS (-001)

Spaulding™ Aseptic Disconnect, Female Cap Connector, Animal-Free Polysulfone, Silicone O-Ring



SDAMPNPS (-001)

Spaulding™ Aseptic Disconnect, Male Plug Connector, Animal-Free Polysulfone, Silicone O-Ring



BIOPROCESSING
FITTINGS

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!
www.nordsonmedical.com

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

CYLINDRAFlow™

Manifolds

Nordson MEDICAL is proud to introduce the first all-in-one manifold for pharmaceutical and single-use biopharmaceutical fluid transfer applications. The revolutionary design of this patent-pending, compact fitting not only saves space, but also simplifies use and assembly. With five ports available in one component, gone are the many leak points found in traditional stopcock and Y connector manifolds. Ideal for downstream processing, multiple bags can easily be filled from one large bag. Flow is consistent as it travels around the inner body of the part, maintaining a steady pressure. The large flow selector is effortless to turn, and the lock-in-place feature gives you assurance that you are utilizing the desired port.

CYLINDRAFlow: Manifolds

FMD5P670-001

CYLINDRAFlow™ Manifold Director, 600 Series Barbs, 3/8" (9.5 mm) Flow, 1 Input, 5 Outputs, 15 psi Liquid Rated, -9024 Polycarbonate Material

Features & Benefits:

- Can accommodate up to 5 fluid lines
- Flow restricted to one port at a time
- No pressure drop as fluid moves around & through the ports
- Intuitive push-to-turn flow selector
- Easier to use and fewer leak points than stopcock and Y connector manifolds
- Materials meet USP Class VI & ISO 10993 criteria
- Gamma radiation stable
- Can be daisy chained for additional ports
- Arrows to provide a clear indication of how to turn the flow selector and which port is in use
- Etched data code for easy access to product resources and material information


Specifications:

Materials: Body & flow selector - gamma stable polycarbonate, internal spring - 316 stainless steel, gasket - platinum-cured silicone

Maximum Pressure: 15psi (1.0 bar) liquid rated

Fits Tubing IDs: 3/8" (9.5 mm)

Barb Size: 600 series

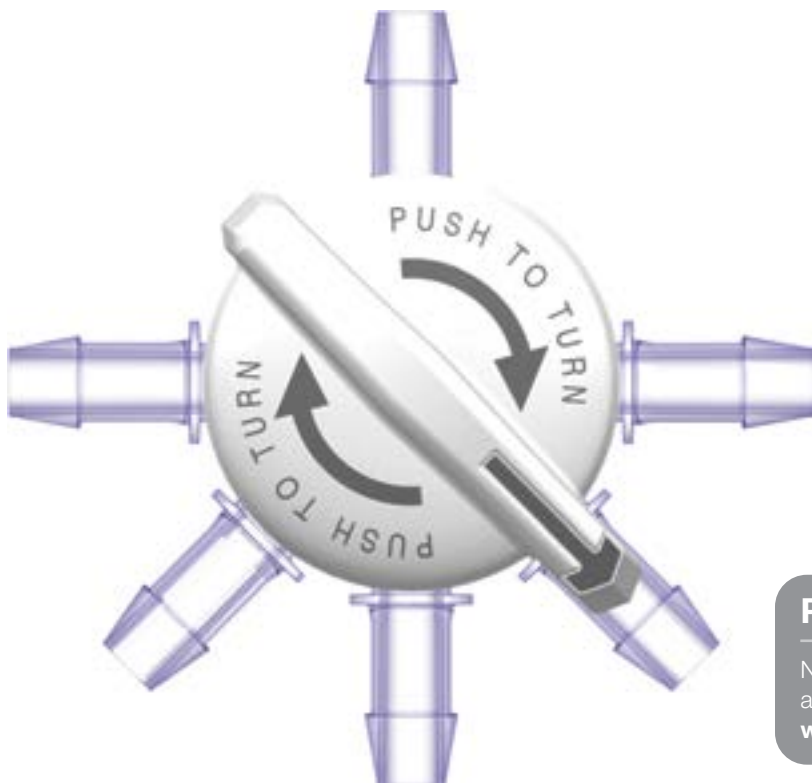

Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

SF Series

Sanitary Fittings

The SF Series sanitary fittings are based on design concepts calling for a heavy-duty part with superior side load strength, while providing superior flow-through characteristics, making them the connectors of choice for single-use bioprocesses. Pharmaceutical professionals find these plastic fittings offer all the benefits of disposable parts while offering performance capabilities not found in earlier versions of flanged fittings.

Features & Benefits:

- Animal-free materials: No BSE documentation required
- Radiused flow surfaces: Smooth uninterrupted flow
- "Mini" & "Maxi" sizes: Easy to retrofit current operations, works with commercially available biopharm clamps
- No Particulate build-up: Interior surface controlled to 15 Ra microinch
- Reducer sizes available

Specifications:

Materials: Animal-Free Polypropylene, Natural PVDF

Tubing ID Sizes:

- 1/8, 1/4", 3/8", 1/2", 3/4", 1"
- Reducers 1" x 3/4", 1-1/2" x 3/4", 1-1/2" x 1", 2" x 3/4", 2" x 1-1/2"

Leak Resistance: Tested at 80 psi air under water for 90 seconds with no leaks

Tubing Removal Force:

- 1/4" PP with C-Flex® shore A tubing: 21 lbs.
- 1/2" PP with C-Flex® shore A tubing: 89 lbs.
- 1" PP with C-Flex® shore A tubing: 121 lbs.

Interior Surface Finish: 15 Ra Microinch

SFMN630 (-6005 -J1A)

Sanitary Fitting, Mini Flange 3/4" to 600 Series Barb, 1/8" (3.2 mm) Tubing ID



SFMN655 (-6005 -J1A)

Sanitary Fitting, Mini Flange 3/4" to 600 Series Barb, 1/4" (6.4 mm) Tubing ID



SFMN670 (-6005 -J1A)

Sanitary Fitting, Mini Flange 3/4" to 600 Series Barb, 3/8" (9.5 mm) Tubing ID



SFMN680 (-6005 -J1A)

Sanitary Fitting, Mini Flange 3/4" to 600 Series Barb, 1/2" (12.7 mm) Tubing ID



SFMN6100 (-6005 -J1A)

Sanitary Fitting, Mini Flange 3/4" to 600 Series Barb, 3/4" (19.0 mm) Tubing ID



SFMNC (-6005 -J1A)

Sanitary Fitting Mini Cap 3/4"



Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

SFMX655 (-6005 -J1A)

Sanitary Fitting, Maxi Flange 1-1/2" to 600 Series Barb, 1/4" (6.4 mm) Tubing ID



SFMX6100 (-6005 -J1A)

Sanitary Fitting, Maxi Flange 1-1/2" to 600 Series Barb, 3/4" (19.0 mm) Tubing ID



SFMX670 (-6005 -J1A)

Sanitary Fitting, Maxi Flange 1-1/2" to 600 Series Barb, 3/8" (9.5 mm) Tubing ID



SFMX6110 (-6005 -J1A)

Sanitary Fitting, Maxi Flange 1-1/2" to 600 Series Barb, 1.0" (25.4 mm) Tubing ID



SFMX680 (-6005 -J1A)

Sanitary Fitting, Maxi Flange 1-1/2" to 600 Series Barb, 1/2" (12.7 mm) Tubing ID



SFMXC (-6005 -J1A)

Sanitary Fitting Maxi Cap 1-1/2"



SFM100/SFM075 (-6005 -J1A)

Sanitary Fitting, Concentric Reducer 1" (25.4 mm) Landish TC x .75" (19.0 mm) TC



SFM200/SFM075 (-6005 -J1A)

Sanitary Fitting, Concentric Reducer 2" (50.8 mm) Landish TC x .75" (19.0 mm) TC



SFM150/SFM075 (-6005 -J1A)

Sanitary Fitting, Concentric Reducer 1.5" (38.1 mm) Landish TC x .75" (19.0 mm) TC



SFM200/SFM150 (-6005 -J1A)

Sanitary Fitting, Concentric Reducer 2" (50.8 mm) Landish TC x 1.5" (38.1 mm) TC



SFM150/SFM100 (-6005 -J1A)

Sanitary Fitting, Concentric Reducer 1.5" (38.1 mm) Landish TC x 1" (25.4 mm) TC



BIOPROCESSING
FITTINGS

SF Series

Instrument Tees

SF Series Instrument Tees work well with existing pharmaceutical instrumentation and are available in 1/2" and 1" tubing IDs. Critical to a successful single-use implementation, the fittings work for pressure testing, fluid transfer, and filtration applications. The configuration is a standard 1.5" TC style interface on the instrument side and tee style for the fluid path. The fittings also utilize Nordson MEDICAL's popular 600 series barb, which is ideally suited for silicone and TPE tubing.

Features & Benefits:

- Works with existing pharmaceutical instrumentation
- Single-use technology
- ASME BPE compliant interface on the instrument side
- Tee style for the fluid path
- 600 series barb
- Common industry materials

Specifications:

Materials: Animal-Free Polyethylene

Tubing Sizes: 1/2" and 1" tubing IDs

Configuration: 1.5" TC style

INSTRUMENT TEES:

Maxi Flange

SFMXT680 (-VP1) Value Pharma™ Resin

Sanitary Fitting Instrument Tee, Maxi Flange to 600 Series Barb, 1/2" (12.7 mm) Tubing ID



INSTRUMENT TEES:

Maxi Flange

SFMXT6110 (-VP1) Value Pharma™ Resin

Sanitary Fitting Instrument Tee, Maxi Flange to 600 Series Barb, 1" (25.4 mm) Tubing ID



Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

GASKETS

Material options include Class VI compliant PTFE/EPDM and platinum-cured silicone in gasket sizes that fit both the 1-1/2" (Maxi) and the 3/4" (Mini) Flange Series. Keeping customer requirements in mind, these components were manufactured for applications requiring the highest sanitary standards and chemical purity.

GMN450 (-PX) Platinum-Cured Silicone
Molded Gasket for Mini Flange SFMN Series Fittings



G42MPGR-TF075 (-E) PTFE
Molded Tuf-Flex® Gasket for Mini Flange SFMN Series Fittings



GMX455 (-PX) Platinum-Cured Silicone
Molded Gasket for Maxi Flange SFMX Series Fittings



G40MPGR-TF150 (-E) PTFE
Molded Tuf-Flex® Gasket for Maxi Flange SFMX Series Fittings



G424TRX-0200 Platinum-Cured Silicone
Molded Gasket for 2" Flange Fittings



G424TRX-0300 Platinum-Cured Silicone
Molded Gasket for 3" Flange Fittings



BIOPROCESSING
FITTINGS

Need More Details?

Check out all of our bioprocessing components at www.nordsonmedical.com or contact Sales at FMF-sales@nordsonmedical.com.

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

GM4749TRX-0600 Platinum-Cured Silicone
Molded Gasket for 6" Flange Fittings




GASKETS:
Platinum-Cured Silicone Gaskets

G40RX-PX-0250 Platinum-Cured Silicone
Molded Gasket, 2.5"



GASKETS:
Platinum-Cured Silicone Gaskets

G40RX-PX-0400 Platinum-Cured Silicone
Molded Gasket, 4"

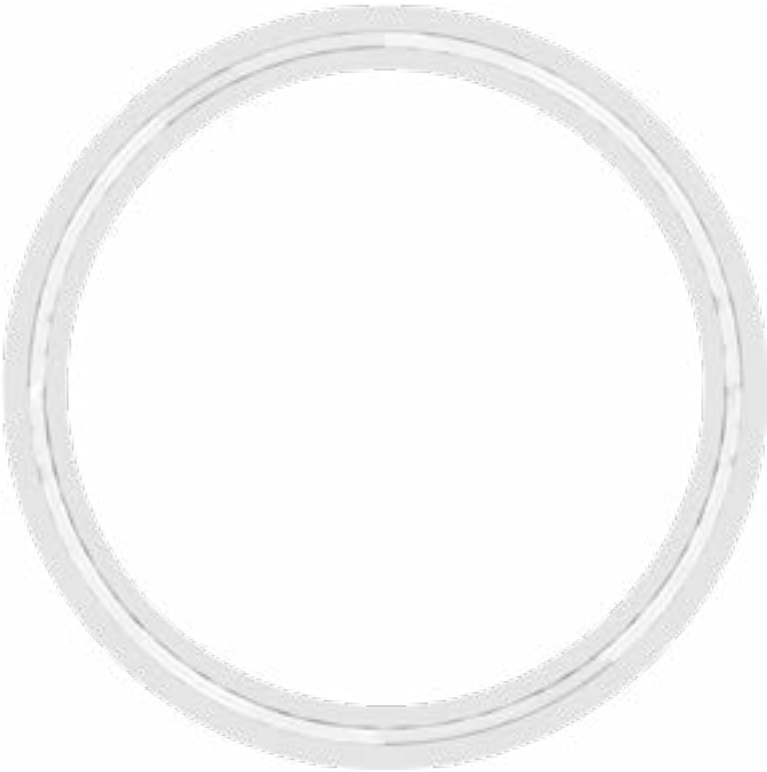

BIOPROCESSING
FITTINGS

BLOOD PRESSURE

Fittings



BPF Series Female Bayonet Style Connectors PG-114



Bayonet Male Cuff Connectors PG-115



Bayonet Panel Mount Connectors PG-115



Dinamap® Monitor Adapters PG-116



Dinamap® Cuff Connectors PG-116



BPL Series 1/4 Turn Cuff Connectors PG-117



BPF Series

Blood Pressure Fittings

The BPF Series features an ergonomic design for simple one-handed operation. Easy to connect and disconnect to existing tubing, many locations are replacing their metal bayonet connectors with these light weight connectors. The BPF Series also works well with the new ISO 80369-5 limb cuff inflation standard compliant parts. The standard defines the bayonet as the only approved design.

BPF SERIES:

Female Bayonet Style Quick Connect

BPF220 (-001) Black ABS Housing

Bayonet-Style Female Quick Connect Fitting, 200 Series Barb, 3/32" (2.4 mm) ID Tubing



works with...

BPM013	BPMB230
BPM240	BPMB240
BPM025	BPMB250
BPM035	BPAM
BPMMTLL	

BPF230 (-001) Black ABS Housing

Bayonet-Style Female Quick Connect Fitting, 200 Series Barb, 1/8" (3.2 mm) ID Tubing



works with...

BPM013	BPMB230
BPM240	BPMB240
BPM025	BPMB250
BPM035	BPAM
BPMMTLL	

BPF240 (-001) Black ABS Housing

Bayonet-Style Female Quick Connect Fitting, 200 Series Barb, 5/32" (4.0 mm) ID Tubing



works with...

BPM013	BPMB230
BPM240	BPMB240
BPM025	BPMB250
BPM035	BPAM
BPMMTLL	

BPF250 (-001) Black ABS Housing

Bayonet-Style Female Quick Connect Fitting, 200 Series Barb, 3/16" (4.8 mm) ID Tubing



works with...

BPM013	BPMB230
BPM240	BPMB240
BPM025	BPMB250
BPM035	BPAM
BPMMTLL	

BLOOD PRESSURE: Bayonet Male Cuff Connectors

BPM013 (-1)

Blood Pressure Monitor Bayonet Connector to 500 Series Barb, 80369-5 Compliant, 1/8" (3.2 mm) ID Tubing



BPM240 (-1)

Blood Pressure Monitor Bayonet Connector to 200 Series Barb, 80369-5 Compliant, 5/32" (4.0 mm) ID Tubing



BPM025 (-1)

Blood Pressure Monitor Bayonet Connector to 500 Series Barb, 80369-5 Compliant, 5/32" (4.0 mm) ID Tubing



BPM035 (-1)

Blood Pressure Monitor Bayonet Connector to 500 Series Barb, 80369-5 Compliant, 3/16" (4.8 mm) ID Tubing



BPMMTLL (-1 -5)

Blood Pressure Monitor Bayonet Connector to Male Luer Integral Lock Ring



works with...



BLOOD PRESSURE: Bayonet Panel Mount Connectors

BPMB230 (-1)

Blood Pressure Monitor Bayonet Panel Mount Connector to 200 Series Barb, 1/4-28 UNF Thread, 1/8" (3.2 mm) ID Tubing



works with...



BPMB240 (-1)

Blood Pressure Monitor Bayonet Panel Mount Connector to 200 Series Barb, 5/16-24 UNF Thread, 5/32" (4.0 mm) ID Tubing



works with...



BPMB250 (-1)

Blood Pressure Monitor Bayonet Panel Mount Connector to 200 Series Barb, 3/8-24 UNF Thread, 3/16" (4.8 mm) ID Tubing



works with...



BLOOD PRESSURE: Panel Mount Nuts

PMSN (-XO)

Panel Mount Nut 1/4-28 UNF with 7/16" Hex (For use with BPMB230)



PM3N (-XO)

Panel Mount Nut 5/16-24 UNF Zinc-Plated Class 2 Steel Nut (For use with BPMB240)



PM2N (-XO)

Panel Mount Nut 3/8-24 UNF Zinc-Plated Class 2 Steel Nut (For use with BPMB250)



Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

BLOOD PRESSURE:
Dinamap® Monitor Adapters

BPAM (-001) Black Nylon over White Nylon

Blood Pressure Monitor Bayonet Connector to Internal 5/16-24 UNF Thread and O-Ring



BPA013 (-001) Black Nylon over White Nylon

Blood Pressure Adapter with Internal 5/16-24 UNF Thread and O-Ring to 500 Series Barb, 1/8" (3.2 mm) ID Tubing



BPA025 (-001) Black Nylon over White Nylon

Blood Pressure Adapter with Internal 5/16-24 UNF Thread and O-Ring to 500 Series Barb, 5/32" (4.0 mm) ID Tubing



BPA035 (-001) Black Nylon over White Nylon

Blood Pressure Adapter with Internal 5/16-24 UNF Thread and O-Ring to 500 Series Barb, 3/16" (4.8 mm) ID Tubing



BPA055 (-001) Black Nylon over White Nylon

Blood Pressure Adapter with Internal 5/16-24 UNF Thread and O-Ring to 500 Series Barb, 1/4" (6.4 mm) ID Tubing



BLOOD PRESSURE:
Dinamap® Cuff Connectors

BPCL220 (-10 -81)

Blood Pressure Adapter with 5/16-24 UNF Thread and 6% Taper to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



BPC220 (-10)

Blood Pressure Adapter with 5/16-24 UNF Thread to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



BPCL230 (-10 -81)

Blood Pressure Adapter with 5/16-24 UNF Thread and 6% Internal Taper to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



BPC230 (-10)

Blood Pressure Adapter with 5/16-24 UNF Thread to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



BPCL240 (-10 -81)

Blood Pressure Adapter with 5/16-24 UNF Thread and 6% Internal Taper to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



BPC240 (-10)

Blood Pressure Adapter with 5/16-24 UNF Thread to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



BLOOD
PRESSURE

Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.
www.nordsonmedical.com

BPCL025 (-10)

Blood Pressure Adapter with 5/16-24 UNF Thread 6% Internal Taper to 500 Series Barb, 5/32" (4.0 mm) ID Tubing



BPC025 (-10)

Blood Pressure Adapter with 5/16-24 UNF Thread to 500 Series Barb, 5/32" (4.0 mm) ID Tubing



BPCL250 (-10 -81)

Blood Pressure Adapter with 5/16-24 UNF Thread and 6% Internal Taper to 200 Series Barb, 3/16" (4.8 mm) ID Tubing



BPC250 (-10)

Blood Pressure Adapter with 5/16-24 UNF Thread to 200 Series Barb, 3/16" (4.8 mm) ID Tubing



CF229 (-1)

Neonatal Blood Pressure Slip Fit (non luer) with 200 Series Barb, 3/32" (2.4 mm) ID Tubing



CF230 (-1)

Neonatal Blood Pressure Slip Fit (non luer) with 200 Series Barb, 5/32" (4.0 mm) ID Tubing



BPL Series

Blood Pressure Fittings

The BPL Series quick connect fittings offers superior ergonomics - a requirement for many designers in their effort to accommodate end users. The BPL's ease-of-use stems largely from a winged male fitting and Nordson MEDICAL's patented 10-point grip, making the quarter-turn connection simple and easy to accomplish. A solid snap, that can be both heard and felt, will then reassure the end user that a secure fit has been made.

Another unique characteristic of the BPL Series is the true tube stop on the rotating half of the two connectors. This allows tubing to be fully assembled onto the barb without its edge pushing against the body of the connector.

BPLF230 (-8012)

Quarter-turn Female Blood Pressure Cuff Connector to 200 Series Barb, 1/8" (3.2 mm) Tubing ID



BPLF240 (-8012)

Quarter-turn Female Blood Pressure Cuff Connector to 200 Series Barb, 5/32" (4.0 mm) Tubing ID



BPLM230 (-001) White ABS, White Acetal, Buna-N O-Ring

Quarter-turn Male Blood Pressure Cuff Connector to 200 Series Barb, 1/8" (3.2 mm) Tubing ID



BPLM240 (-001) White ABS, White Acetal, Buna-N O-Ring

Quarter-turn Male Blood Pressure Cuff Connector to 200 Series Barb, 5/32" (4.0 mm) Tubing ID



Did You Know?

Nordson MEDICAL has an interactive PDF available online at www.nordsonmedical.com. Access part drawings, reference charts, data sheets, and more from just one location!

THREADED

Fittings

10-32 UNF Special Tapered Thread **PG-119** 1/4-28 UNF **PG-122**




1/4-28 UNF Bottom Sealing **PG-123** Metric - M6x1 **PG-124**



Metric - M5x.8 **PG-124** 1/16-27 NPT **PG-125**



1/8-27 NPT **PG-126** 1/4-18 NPT Reducers **PG-127**

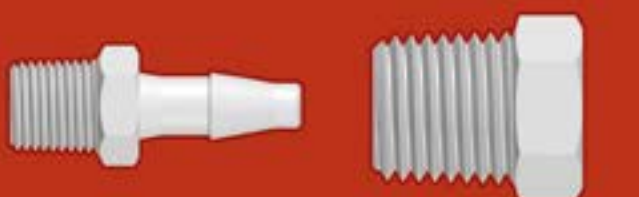

The X Thread

Nordson MEDICAL X Series components have a tapered 10-32 thread that functions like a pipe thread. In comparison to standard 10-32 UNF threaded fittings, the tapered design works better with oversized ports and ports which are not perpendicular to the face of the component.

The 120 degree chamfer under the hex provides a seal when the fitting is tightened into a port with a countersink as shown below. It also seals on the thread due to interference on the pitch diameter. When using an elbow or tee, the tapered thread tightens radially and facilitates directional alignment of the barb.


Other Thread Series

K, B, C and D Series fittings, designed for use with 10-32 UNF ports, are straight-threaded components. These fittings seal against the 120 degree chamfer under the hex on the fitting and do not tighten radially as they are assembled. S Series fittings are straight threaded components for 1/4-28 UNF ports, M5 Series are for M5 x .8 ports and M6 Series are for M6 x 1 ports.

Part numbers for fittings with pipe threads are noted below. Those for use with 1/16-27 NPT ports have part numbers prefixed with 16; those for 1/8-27 NPT ports have part numbers prefixed with 18; and those for 1/4-18 NPT ports have part numbers prefixed with 14.

The features shown below must be incorporated into the port when using our B, C, D, X, K, S, M5 and M6 Series fittings. Failure to do so will jeopardize the performance of the seal between the fitting and port. For applications in which this seal is highly critical, we recommend that you contact a sealant specialist.

When working with threaded parts, Nordson MEDICAL strongly recommends that you test sample parts in your application.

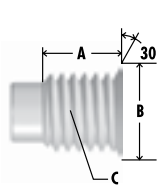

Nordson MEDICAL's Thread Series	Dimension A Thread Length	Clearance Needed to Achieve Chamfer Seal	Dimension B Countersink Diameter	Feature C Internal Screw Thread Series
X	.200"	.250"	.210"-.230"	10-32 UNF 2B
K, B, C, D	.156"	.200"	.210"-.230"	10-32 UNF 2B
S	.204"	.250"	.270"-.290"	1/4-28 UNF 2B
M5	.192"	.250"	.210"-.230"	M5 X .8
M6	.235"	.375"	.270"-.290"	M6 X 1

K10 (-1 -2 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to Classic Series Barb, 1/16" (1.6 mm) ID Tubing



X10 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread with 1/4" Hex to Classic Series Barb, 1/16" (1.6 mm) ID Tubing



K210 (-1 -2 -3 -5 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



X210 (-1 -2 -3 -5 -6005 -6012 -J1A)

10-32 Special Tapered Thread with 1/4" Hex to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



K410 (-1 -2 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to 400 Series Barb, 1/16" (1.6 mm) ID Tubing



X410 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread with 1/4" Hex to 400 Series Barb, 1/16" (1.6 mm) ID Tubing



D (-1 -2 -3 -4 -5 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to Classic Series Barb, 3/32" (2.4 mm) ID Tubing



K220 (-1 -2 -6005 -6012 -J1A)

10-32 UNF Thread with 1/4" Hex to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



X220 (-1 -2 -3 -6005 -J1A)

10-32 Special Tapered Thread with 1/4" Hex to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



K420 (-1 -2 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to 400 Series Barb, 3/32" (2.4 mm) ID Tubing



X420 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread with 1/4" Hex to 400 Series Barb, 3/32" (2.4 mm) ID Tubing



B (-1 -2 -3 -4 -5 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to Classic Series Barb, 1/8" (3.2 mm) ID Tubing



C (-1 -2 -3 -4 -5 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to Classic Series Barb, Short Neck, 1/8" (3.2 mm) ID Tubing



K230 (-1 -2 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



X230 (-1 -2 -3 -5 -6005 -6012 -J1A)

10-32 Special Tapered Thread with 1/4" Hex to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



K430 (-1 -2 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to 400 Series Barb, 1/8" (3.2 mm) ID Tubing



X430 (-1 -2 -3 -6005 -J1A)

10-32 Special Tapered Thread with 1/4" Hex to 400 Series Barb, 1/8" (3.2 mm) ID Tubing



K240 (-1 -2 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



X240 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread with 1/4" Hex to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



K445 (-1 -2 -6005 -J1A)

10-32 UNF Thread with 1/4" Hex to 400 Series Barb, .170" (4.3 mm) ID Tubing



X445 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread with 1/4" Hex to 400 Series Barb, .170" (4.3 mm) ID Tubing



THREAD:
10-32 UNF, Special Tapered Thread Elbow

KL210 (-1 -2 -6005 -J1A)

10-32 UNF Thread Elbow with 1/4" Hex to 200 Series Barb,
1/16" (1.6 mm) ID Tubing



XL210 (-1 -2 -6005 -6012 -J1A)

10-32 Special Tapered Thread Elbow with 1/4" Hex to 200 Series Barb,
1/16" (1.6 mm) ID Tubing



KL410 (-1 -2 -6005 -J1A)

10-32 UNF Thread Elbow with 1/4" Hex to 400 Series Barb,
1/16" (1.6 mm) ID Tubing



XL410 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread Elbow with 1/4" Hex to 400 Series Barb,
1/16" (1.6 mm) ID Tubing



KL220 (-1 -2 -6005 -J1A)

10-32 UNF Thread Elbow with 1/4" Hex to 200 Series Barb,
3/32" (2.4 mm) ID Tubing



XL220 (-1 -2 -6005 -6012 -J1A)

10-32 Special Tapered Thread Elbow with 1/4" Hex to 200 Series Barb,
3/32" (2.4 mm) ID Tubing



KL420 (-1 -2 -6005 -J1A)

10-32 UNF Thread Elbow with 1/4" Hex to 400 Series Barb,
3/32" (2.4 mm) ID Tubing



XL420 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread Elbow with 1/4" Hex to 400 Series Barb,
3/32" (2.4 mm) ID Tubing



KL230 (-1 -2 -6005 -J1A)

10-32 UNF Thread Elbow with 1/4" Hex to 200 Series Barb,
1/8" (3.2 mm) ID Tubing



XL230 (-1 -2 -6005 -6012 -J1A)

10-32 Special Tapered Thread Elbow with 1/4" Hex to 200 Series Barb,
1/8" (3.2 mm) ID Tubing



KL430 (-1 -2 -6005 -J1A)

10-32 UNF Thread Elbow with 1/4" Hex to 400 Series Barb,
1/8" (3.2 mm) ID Tubing



XL430 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread Elbow with 1/4" Hex to 400 Series Barb,
1/8" (3.2 mm) ID Tubing



Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

KT210 (-1 -2 -6005 -J1A)

10-32 UNF Thread Tee with 1/4" Hex to 200 Series Barb,
1/16" (1.6 mm) ID Tubing



KT420 (-1 -2 -6005 -J1A)

10-32 UNF Thread Tee with 1/4" Hex to 400 Series Barb,
3/32" (2.4 mm) ID Tubing



XT210 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread Tee with 1/4" Hex to 200 Series Barb,
1/16" (1.6 mm) ID Tubing



XT420 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread Tee with 1/4" Hex to 400 Series Barb,
3/32" (2.4 mm) ID Tubing



KT410 (-1 -2 -6005 -J1A)

10-32 UNF Thread Tee with 1/4" Hex to 400 Series Barb,
1/16" (1.6 mm) ID Tubing



KT230 (-1 -2 -6005 -J1A)

10-32 UNF Thread Tee with 1/4" Hex to 200 Series Barb,
1/8" (3.2 mm) ID Tubing



XT410 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread Tee with 1/4" Hex to 400 Series Barb,
1/16" (1.6 mm) ID Tubing



XT230 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread Tee with 1/4" Hex to 200 Series Barb,
1/8" (3.2 mm) ID Tubing



KT220 (-1 -2 -6005 -J1A)

10-32 UNF Thread Tee with 1/4" Hex to 200 Series Barb,
3/32" (2.4 mm) ID Tubing



KT430 (-1 -2 -6005 -J1A)

10-32 UNF Thread Tee with 1/4" Hex to 400 Series Barb,
1/8" (3.2 mm) ID Tubing



XT220 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread Tee with 1/4" Hex to 200 Series Barb,
3/32" (2.4 mm) ID Tubing



XT430 (-1 -2 -6005 -J1A)

10-32 Special Tapered Thread Tee with 1/4" Hex to 400 Series Barb,
1/8" (3.2 mm) ID Tubing



Part Drawings Online!

Nordson MEDICAL provides dimensional drawings
available for download for all of its fittings.

www.nordsonmedical.com

THREAD:
10-32 UNF, Special Tapered Thread Couplers and Plugs

KN (-1 -2 -6005 -J1A)
10-32 UNF Thread Coupler with 1/4" Hex



XN (-1 -2 -6005 -J1A)
10-32 Special Tapered Thread Coupler with 1/4" Hex



KP (-1 -2 -6005 -6012 -J1A)
10-32 UNF Thread Plug with 1/4" Hex



XP (-1 -2 -6005 -J1A)
10-32 Special Tapered Thread Plug with 1/4" Hex



THREAD:
1/4-28 UNF

S10 (-1 -2 -6005 -J1A)
1/4-28 UNF Thread with 5/16" Hex to Classic Series Barb,
1/16" (1.6 mm) ID Tubing



S210 (-1 -2 -6005 -J1A)
1/4-28 UNF Thread with 5/16" Hex to 200 Series Barb,
1/16" (1.6 mm) ID Tubing



S410 (-1 -2 -6005 -J1A)
1/4-28 UNF Thread with 5/16" Hex to 400 Series Barb,
1/16" (1.6 mm) ID Tubing



S20 (-1 -2 -6005 -J1A)
1/4-28 UNF Thread with 5/16" Hex to Classic Series Barb,
3/32" (2.4 mm) ID Tubing



S220 (-1 -2 -6005 -J1A)
1/4-28 UNF Thread with 5/16" Hex to 200 Series Barb,
3/32" (2.4 mm) ID Tubing



S420 (-1 -2 -6005 -J1A)
1/4-28 UNF Thread with 5/16" Hex to 400 Series Barb,
3/32" (2.4 mm) ID Tubing



S30 (-1 -2 -6005 -J1A)
1/4-28 UNF Thread with 5/16" Hex to Classic Series Barb,
1/8" (3.2 mm) ID Tubing



S230 (-1 -2 -6005 -6012 -J1A)
1/4-28 UNF Thread with 5/16" Hex to 200 Series Barb,
1/8" (3.2 mm) ID Tubing



Looking for Something Different?

Most of our best product ideas have come from you, our customers. While we pride ourselves on our extensive product offering, we recognize that your needs are ever-changing and that our products, too, must adapt to your needs.

If there's a new product you think Nordson MEDICAL should make, please send us an Email at designcenter@nordsonmedical.com describing your product concept in as much detail as you desire. We will contact you with any questions, and who knows - the next catalog edition could include products YOU told us to build.

See page 10 for more information on the Design Center!



THREAD:
1/4-28 UNF**S430 (-1 -2 -6005 -J1A)**1/4-28 UNF Thread with 5/16" Hex to 400 Series Barb,
1/8" (3.2 mm) ID Tubing**S40 (-1 -2 -6005 -J1A)**1/4-28 UNF Thread with 5/16" Hex to Classic Series Barb,
5/32" (4.0 mm) ID Tubing**S240 (-1 -2 -6005 -J1A)**1/4-28 UNF Thread with 5/16" Hex to 200 Series Barb,
5/32" (4.0 mm) ID Tubing**S445 (-1 -2 -6005 -J1A)**1/4-28 UNF Thread with 5/16" Hex to 400 Series Barb,
.170" (4.3 mm) ID Tubing**S50 (-1 -2 -6005 -J1A)**1/4-28 UNF Thread with 5/16" Hex to Classic Series Barb,
3/16" (4.8 mm) ID Tubing**S250 (-1 -2 -6005 -J1A)**1/4-28 UNF Thread with 5/16" Hex to 200 Series Barb,
3/16" (4.8 mm) ID Tubing**S450 (-1 -2 -6005 -J1A)**1/4-28 UNF Thread with 5/16" Hex to 400 Series Barb,
3/16" (4.8 mm) ID Tubing**SP (-1 -2 -6005 -J1A)**

1/4-28 UNF Thread Plug with 5/16" Hex

**THREAD:**
1/4-28 UNF Bottom Sealing**ABR004-6005 UNDER (-1 -2 -3 -4 -5)**1/4-28 UNF Bottom Sealing Rotating Thread with 5/16" Hex to
500 Series Barb, 1/16" (1.6 mm) ID Tubing*Recommended torque:*
Polypropylene - 1.5-2 in. lbs.
Kynar® - 2-4 in. lbs.**ABR004-J1A UNDER (-1 -2 -3 -4 -5 -J1A)**1/4-28 UNF Bottom Sealing Rotating Thread with 5/16" Hex to
500 Series Barb, 1/16" (1.6 mm) ID Tubing*Recommended torque:*
Polypropylene - 1.5-2 in. lbs.
Kynar® - 2-4 in. lbs.**ABR007-6005 UNDER (-1 -2 -3 -4 -5)**1/4-28 UNF Bottom Sealing Rotating Thread with 5/16" Hex to
500 Series Barb, 3/32" (2.4 mm) ID Tubing*Recommended torque:*
Polypropylene - 1.5-2 in. lbs.
Kynar® - 2-4 in. lbs.**ABR007-J1A UNDER (-1 -2 -3 -4 -5 -J1A)**1/4-28 UNF Bottom Sealing Rotating Thread with 5/16" Hex to
500 Series Barb, 3/32" (2.4 mm) ID Tubing*Recommended torque:*
Polypropylene - 1.5-2 in. lbs.
Kynar® - 2-4 in. lbs.**ABR013-6005 UNDER (-1 -2 -3 -4 -5 -6012)**1/4-28 UNF Bottom Sealing Rotating Thread with 5/16" Hex to
500 Series Barb, 1/8" (3.2 mm) ID Tubing*Recommended torque:*
Polypropylene - 1.5-2 in. lbs.
Kynar® - 2-4 in. lbs.**ABR013-J1A UNDER (-1 -2 -3 -4 -5 -J1A)**1/4-28 UNF Bottom Sealing Rotating Thread with 5/16" Hex to
500 Series Barb, 1/8" (3.2 mm) ID Tubing*Recommended torque:*
Polypropylene - 1.5-2 in. lbs.
Kynar® - 2-4 in. lbs.**TECS116 (-6005 -J1A)**1/4-28 UNF Bottom Sealing Thread Tube End Coupler with 5/16" Hex,
1/16" (1.6 mm) OD Tubing**TECS118 (-6005 -J1A)**1/4-28 UNF Bottom Sealing Thread Tube End Coupler with 5/16" Hex,
1/8" (3.2 mm) OD Tubing

THREAD:
M6x1

M6210 (-1 -2 -6005)
M6x1 Thread with 5/16" Hex to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



M6410 (-1 -2 -6005 -J1A)
M6x1 Thread with 5/16" Hex to 400 Series Barb, 1/16" (1.6 mm) ID Tubing



M6220 (-1 -2 -6005 -J1A)
M6x1 Thread with 5/16" Hex to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



M6420 (-1 -2 -6005 -J1A)
M6x1 Thread with 5/16" Hex to 400 Series Barb, 3/32" (2.4 mm) ID Tubing



M6230 (-1 -2 -6005 -J1A)
M6x1 Thread with 5/16" Hex to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



M6430 (-1 -2 -6005 -J1A)
M6x1 Thread with 5/16" Hex to 400 Series Barb, 1/8" (3.2 mm) ID Tubing



M6240 (-1 -2 -6005 -J1A)
M6x1 Thread with 5/16" Hex to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



M6445 (-1 -2 -6005 -J1A)
M6x1 Thread with 5/16" Hex to 400 Series Barb, .170" (4.3 mm) ID Tubing



THREAD:
M6x1

M6250 (-1 -2 -6005 -J1A)
M6x1 Thread with 5/16" Hex to 200 Series Barb, 3/16" (4.8 mm) ID Tubing



M6450 (-1 -2 -6005 -J1A)
M6x1 Thread with 5/16" Hex to 400 Series Barb, 3/16" (4.8 mm) ID Tubing



M6P (-1 -2 -6005 -J1A)
M6x1 Thread Plug with 5/16" Hex



THREAD:
M5x.8

M5210 (-1 -2 -6005)
M5x.8 Thread with 1/4" Hex to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



M5220 (-1 -2 -6005)
M5x.8 Thread with 1/4" Hex to 200 Series Barb, 3/32" (2.4 mm) ID Tubing



M5230 (-1 -2 -6005)
M5x.8 Thread with 1/4" Hex to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



M5240 (-1 -2 -6005)
M5x.8 Thread with 1/4" Hex to 200 Series Barb, 5/32" (4.0 mm) ID Tubing



Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.
www.nordsonmedical.com

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

1610 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to Classic Series Barb,
1/16" (1.6 mm) ID Tubing



16430 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 400 Series Barb,
1/8" (3.2 mm) ID Tubing



16210 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 200 Series Barb,
1/16" (1.6 mm) ID Tubing



1640 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to Classic Series Barb,
5/32" (4.0 mm) ID Tubing



16410 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 400 Series Barb,
1/16" (1.6 mm) ID Tubing



16240 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 200 Series Barb,
5/32" (4.0 mm) ID Tubing



1620 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to Classic Series Barb,
3/32" (2.4 mm) ID Tubing



16445 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 400 Series Barb,
.170" (4.3 mm) ID Tubing



16220 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 200 Series Barb,
3/32" (2.4 mm) ID Tubing



1650 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to Classic Series Barb,
3/16" (4.8 mm) ID Tubing



16420 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 400 Series Barb,
3/32" (2.4 mm) ID Tubing



16250 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 200 Series Barb,
3/16" (4.8 mm) ID Tubing



1630 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to Classic Series Barb,
1/8" (3.2 mm) ID Tubing



16450 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 400 Series Barb,
3/16" (4.8 mm) ID Tubing



16230 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread with 5/16" Hex to 200 Series Barb,
1/8" (3.2 mm) ID Tubing



1600 (-1 -2 -6005 -J1A)

1/16-27 NPT Thread Plug with 5/16" Hex



THREAD:
1/8-27 NPT

1810 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to Classic Series Barb,
1/16" (1.6 mm) ID Tubing



18230 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to 200 Series Barb,
1/8" (3.2 mm) ID Tubing



18210 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to 200 Series Barb,
1/16" (1.6 mm) ID Tubing



18430 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to 400 Series Barb,
1/8" (3.2 mm) ID Tubing



18410 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to 400 Series Barb,
1/16" (1.6 mm) ID Tubing



1840 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to Classic Series Barb,
5/32" (4.0 mm) ID Tubing



1820 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to Classic Series Barb,
3/32" (2.4 mm) ID Tubing



18240 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to 200 Series Barb,
5/32" (4.0 mm) ID Tubing



18220 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to 200 Series Barb,
3/32" (2.4 mm) ID Tubing



18445 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to 400 Series Barb,
1/70" (4.3 mm) ID Tubing



18420 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to 400 Series Barb,
3/32" (2.4 mm) ID Tubing



1850 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to Classic Series Barb,
3/16" (4.8 mm) ID Tubing



1830 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread with 7/16" Hex to Classic Series Barb,
1/8" (3.2 mm) ID Tubing



1860 (-1 -2 -6005 -6012 -J1A)

1/8-27 NPT Thread with 7/16" Hex to Classic Series Barb,
1/4" (6.4 mm) ID Tubing



THREADED
FITTINGS

THREAD: 1/8-27 NPT Elbows

18L230 (-1 -2 -6005 -6012 -J1A)

1/8-27 NPT Thread Elbow with 7/16" Hex to 200 Series Barb,
1/8" (3.2 mm) ID Tubing



18L250 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread Elbow with 7/16" Hex to 200 Series Barb,
3/16" (4.8 mm) ID Tubing



18L260 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread Elbow with 7/16" Hex to 200 Series Barb,
1/4" (6.4 mm) ID Tubing



THREAD: 1/8-27 NPT Tees

18T250 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread Tee with 7/16" Hex to 200 Series Barb,
3/16" (4.8 mm) ID Tubing



18T260 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread Tee with 7/16" Hex to 200 Series Barb,
1/4" (6.4 mm) ID Tubing



THREAD: 1/8-27 NPT Plugs

1800 (-1 -2 -6005 -J1A)

1/8-27 NPT Thread Plug with 7/16" Hex



THREAD: 1/8-27 NPT Reducers

1801 (-1 -2 -6005 -J1A)

Reducer with 1/8-27 NPT Male Thread and 7/16" Hex to
10-32 UNF Female Thread



1802 (-1 -2 -6005 -J1A)


Reducer with 1/8-27 NPT Male Thread and 7/16" Hex to
1/4-28 UNF Female Thread



THREAD: 1/4-18 NPT Reducers


1401 (-1 -2 -6005 -J1A)

Reducer with 1/4-18 NPT Male Thread and 9/16" Hex to
10-32 UNF Female Thread

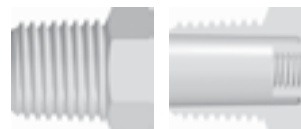

1402 (-1 -2 -6005 -J1A)

Reducer with 1/4-18 NPT Male Thread and 9/16" Hex to
1/4-28 UNF Female Thread

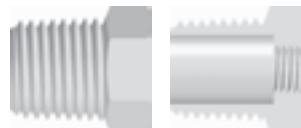

Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

BUILD-A-PART®

Fittings



Junction Blocks

PG-129



Accessories

PG-130



Barbs

PG-130



Threads

PG-131



Bonding Suggestions

Prior to actually bonding, press together a model of the manifold to use as a reference.

Do not mix plastic resins. Consult your adhesive specialist if you need to bond dissimilar materials.

These components are designed for bonding with any one of a number of solvents, but other adhesive systems may be used. When solvent bonding, use a very small gauge needle (typically 27 or 30) for best solvent control. Consult your adhesive specialist for guidance on other systems.

Apply solvent sparingly. Most solvents are low viscosity and readily wick into the bondline. Excessive solvent may compromise the aesthetics and the function of the assembly.


Press the mating components together leaving a .01-.03" gap between them for the solvent. Place the components in an orientation such that they may be rotated ¼ turn to their final position after applying the solvent.

After applying solvent, wait a few seconds, then rotate the components to their final position. When making a long manifold with bondable junction blocks it may be helpful to push the assembly against a flat surface to align the blocks before the solvent has set.

Affordable Unique Solutions:

Design Your New Component with Build-A-Part® Bondable Components!

Component designers no longer need to be constrained by the limitations of injection molding tool design. The exact fitting, in the optimal layout for the process or product concept, is now available with Nordson MEDICAL's line of Build-a-Part® bondable components, designed for use with gases or liquids.

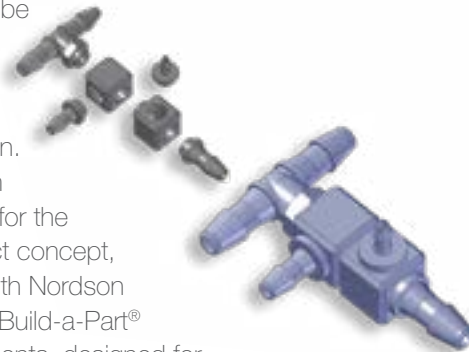

The Build-a-Part® system is based on combining interchangeable components, joining junction blocks and luers, hose barbs or threaded features to create the desired configuration. The features are bonded to the junction blocks using commonly available adhesives. Junction blocks and features may be bonded together to create multi-port, multi-axis products that receive and redirect fluid flow in whatever direction the designer desires. The only limitation is the designer's imagination.


114 component options are available from stock in two primary materials. Options include ISO 594-1 compliant male and female luers and hose barbs in 1/16" (1.6 mm), 3/32" (2.4 mm) and 1/8" (3.2 mm) in straight, elbow and branch tee configurations. Also available are male threads including 10-32 UNF, 1/16-27 NPT, 1/8-27 NPT, 1/4-28 UNF and M6x1 versions, and mounting devices.

The complete range is available in two primary resins: white ABS or translucent polysulfone. For ease of testing concepts, Nordson MEDICAL offers designer kits in both resins. Each kit contains hundreds of components for prototype purposes, enabling designs to be created and modified prior to production commitment.

BLOCK	NUMBER	MATERIAL
Build-A-Part® Junction Blocks. See Configuration Table for specifications.		
	BDJ01-40	Polysulfone
	BDJ01-81	Lace White ABS
	BDJ02-40	Polysulfone
	BDJ02-81	Lace White ABS
	BDJ03-40	Polysulfone
	BDJ03-81	Lace White ABS
	BDJ04-40	Polysulfone
	BDJ04-81	Lace White ABS
	BDJ05-40	Polysulfone
	BDJ05-81	Lace White ABS
	BDJ06-40	Polysulfone
	BDJ06-81	Lace White ABS
	BDJ07-40	Polysulfone
	BDJ07-81	Lace White ABS
	BDJ08-40	Polysulfone
	BDJ08-81	Lace White ABS
	BDJ09-40	Polysulfone
	BDJ09-81	Lace White ABS
	BDJ10-40	Polysulfone
	BDJ10-81	Lace White ABS
	BDJ11-40	Polysulfone
	BDJ11-81	Lace White ABS
	BDJ12-40	Polysulfone
	BDJ12-81	Lace White ABS
	BDJ13-40	Polysulfone
	BDJ13-81	Lace White ABS
	BDJ14-40*	Polysulfone
	BDJ14-81*	Lace White ABS
	BDJ15-40*	Polysulfone
	BDJ15-81*	Lace White ABS

*Check valve junction block. For use with RB01 and SB01

BLOCK	NUMBER	MATERIAL
Build-A-Part® Junction Blocks with Mounting Port. See Configuration Table for specifications.		
	BDJ16-40*	Polysulfone
	BDJ16-81*	Lace White ABS
	BDJ17-40*	Polysulfone
	BDJ17-81*	Lace White ABS
	BDJ18-40	Polysulfone
	BDJ18-81	Lace White ABS
	BDJ19-40	Polysulfone
	BDJ19-81	Lace White ABS
	BDJ20-40	Polysulfone
	BDJ20-81	Lace White ABS
	BDJ21-40	Polysulfone
	BDJ21-81	Lace White ABS
	BDJ22-40	Polysulfone
	BDJ22-81	Lace White ABS

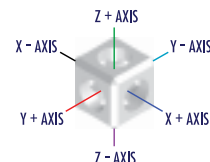

Build-A-Part® Junction Block Configuration Table

BLOCK	X+ AXIS	X- AXIS	Y+ AXIS	Y- AXIS	Z+ AXIS	Z- AXIS
BDJ01	Female	Female				
BDJ02			Female		Female	
BDJ03			Female		Female	Female
BDJ04	Female	Female	Female	Female		
BDJ05	Female		Female		Female	
BDJ06	Female		Female		Female	Female
BDJ07	Female	Female	Female	Female	Female	
BDJ08	Female	Female	Female	Female	Female	Female
BDJ09	Male	Female	Female		Female	Female
BDJ10	Male	Female	Female	Female	Female	Female
BDJ11	Male	Female	Female	Female		
BDJ12	Male	Female	Female			
BDJ13	Male	Female	Female		Female	
BDJ14	Male	Female				
BDJ15	Female	Female				
BDJ16	Female	Female		Female ^a		
BDJ17	Male	Female	Female ^a			
BDJ18	Male ^b	Female	Female		Female	Female
BDJ19	Male ^b	Female	Female	Female	Female	Female
BDJ20	Male ^b	Female	Female	Female		
BDJ21	Male ^b	Female	Female			
BDJ22	Male ^b	Female	Female		Female	

^aFor non-flow manifolding ^bFor use with BDMD

RB01 (-X2) SB01 (-X0)

Check Valve Ball 5/32" Diameter



BDMRP (-40 -81)

Bondable Plug



BDMD (-40 -81)

Bondable Mounting Device



BDMFTLL (-40 -81 -9)

Bondable Female Luer Thread Style




BDMMTL (-40 -81)

Bondable Male Luer (May be used with separate rotating lock ring; RMLLR)



RMLLR (-1 -2 -3 -4 -5 -06 -07 -6005 -9 -9002 -J1A)

Rotating Male Luer Lock Ring (For use with BDMMTL)

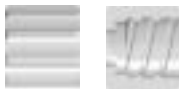

BDMR210 (-40 -81 -9 -9002)

Bondable 200 Series Barb, 1/16" (1.6 mm) ID Tubing



BDMR410 (-40 -81)

Bondable 400 Series Barb, 1/16" (1.6 mm) ID Tubing



BDMR220 (-40 -81)

Bondable 200 Series Barb, 3/32" (2.4 mm) ID Tubing



BDMR420 (-40 -81)

Bondable 400 Series Barb, 3/32" (2.4 mm) ID Tubing



BDMR230 (-40 -81)

Bondable 200 Series Barb, 1/8" (3.2 mm) ID Tubing



BDMR430 (-40 -81)

Bondable 400 Series Barb, 1/8" (3.2 mm) ID Tubing



BDMR240 (-40 -6005 -81)

Bondable 200 Series Barb, 5/32" (4.0 mm) ID Tubing



BDMR445 (-40 -81)

Bondable 400 Series Barb, .170" (4.3 mm) ID Tubing



BDML210 (-40 -81 -9)

Bondable Elbow 200 Series Barb, 1/16" (1.6 mm) ID Tubing



BDML410 (-40 -81)

Bondable Elbow 400 Series Barb, 1/16" (1.6 mm) ID Tubing



BDML220 (-40 -81 -9)

Bondable Elbow 200 Series Barb, 3/32" (2.4 mm) ID Tubing



BDML420 (-40 -81)

Bondable Elbow 400 Series Barb, 3/32" (2.4 mm) ID Tubing



BDML230 (-40 -81)

Bondable Elbow 200 Series Barb, 1/8" (3.2 mm) ID Tubing



BDML430 (-40 -81)

Bondable Elbow 400 Series Barb, 1/8" (3.2 mm) ID Tubing



BDMT210 (-40 -81)

Bondable Tee 200 Series Barb, 1/16" (1.6 mm) ID Tubing



BDMT410 (-40 -81)

Bondable Tee 400 Series Barb, 1/16" (1.6 mm) ID Tubing



BDMT220 (-40 -81)

Bondable Tee 200 Series Barb, 3/32" (2.4 mm) ID Tubing



BDMT420 (-40 -81)

Bondable Tee 400 Series Barb, 3/32" (2.4 mm) ID Tubing



BDMT230 (-40 -81)

Bondable Tee 200 Series Barb, 1/8" (3.2 mm) ID Tubing



BDMT430 (-40 -81)

Bondable Tee 400 Series Barb, 1/8" (3.2 mm) ID Tubing



BDMK (-40 -81)

Bondable 10-32 UNF Thread with 1/4" Hex



BDMX (-40 -81)

Bondable 10-32 Special Tapered Thread with 1/4" Hex



BDMS (-81)

Bondable 1/4-28 UNF Thread with 5/16" Hex



BDMM6 (-40 -81)

Bondable M6x1 Thread with 5/16" Hex



BDM16 (-40 -81)

Bondable 1/16-27 NPT Thread with 5/16" Hex



BDM18 (-40 -81)

Bondable 1/8-27 NPT Thread with 7/16" Hex



BDF16 (-40 -81)

Bondable 1/16-27 NPT Thread with 5/16" Hex



BDF18 (-40 -81)

Bondable 1/8-27 NPT Thread with 7/16" Hex



BDFM6 (-40 -81)

Bondable M6x1 Thread with 5/16" Hex



BDFS (-40 -81)

Bondable 1/4-28 UNF Thread with 5/16" Hex



Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!
www.nordsonmedical.com

TUBING

Medical Grade



Linear Low Density Polyethylene Tubing **PG-132** Thermoplastic Elastomer Tubing (TPE) **PG-132**



Platinum-Cured Silicone Tubing **PG-133** Pump Grade Platinum-Cured Silicone Tubing **PG-133**



Polyurethane Tubing **PG-133** Anti-Microbial Polyvinyl Chloride (PVC) Tubing **PG-133**



TUBING:

Linear Low Density Polyethylene

92 Durometer: This tubing has excellent mechanical and chemical resistance properties. Please use caution if your application involves possible bending or kinking the tubing. It is a semi-rigid material and will be permanently deformed if kinked or bent too tight.

Resin Grade & Technical Information: This tubing is extruded from Equistar® LLDPE, which is in compliance with all U.S. Food & Drug Administration regulations for food grade use under Title 21 CFR. It can be sterilized by gamma radiation or Ethylene Oxide (EtO) gas. Autoclaving is not recommended.

DIMENSION (ID X OD X LGTH)	NUMBER
1/16" X 1/10" X 100'	LLDPE02-1019C
1/16" X 1/10" X 1000'	LLDPE02-1019M
1/16" X 1/8" X 100'	LLDPE02-1031C
1/16" X 1/8" X 1000'	LLDPE02-1031M
3/32" X 5/32" X 100'	LLDPE02-2031C
3/32" X 5/32" X 1000'	LLDPE02-2031M
1/8" X 3/16" X 100'	LLDPE02-3031C
1/8" X 3/16" X 1000'	LLDPE02-3031M

TUBING:

Thermoplastic Elastomer (TPE)

65 Shore A: Our Biopharmaceutical-grade TPE tubing is sterile weldable, heat sealable, moldable, and offers excellent life in peristaltic pumps. It is ideal for use in single-use processes for sampling, filling, and storage tubing set assemblies.

Resin Grade & Technical Information: This tubing is made from FDA-approved ingredients, certified free of silicone oils and animal derived ingredients, and is sterilizable by autoclave or gamma irradiation. It meets various ISO and USP standards, including Class VI, as well as European Pharmacopoeia 3.2.2.1. It also offers excellent low absorption and permeability as compared to silicone tubing.

DIMENSION (ID X OD X LGTH)	NUMBER
1/8" X 1/4" X 50'	AF7500077
1/8" X 1/4" X 500'	AF7500077D
3/16" X 5/16" X 50'	AF7500155
1/4" X 3/8" X 100'	AF7500285C
1/4" X 7/16" X 50'	AF7500311
3/8" X 5/8" X 50'	AF7500506
1/2" X 3/4" X 50'	AF7500584
3/4" X 1" X 50'	AF7500753

50 Shore A: Extruded from a low-volatile grade, platinum-cured compound, this tubing is ideal for critical applications in pharmaceutical, biotech, medical, food, and beverage industries.

Resin Grade & Technical Information: This tubing is extruded and packaged in a Class 7 (Class 10,000) ISO-certified clean room. It has undergone extensive physical, chemical, and biological testing and meets USP Class VI, FDA CFR 177.2600, ISO 10993, European Pharmacopoeia 3.1.9, and 3-A standards. It is certified free of animal-derived ingredients and is certified by the National Sanitation Foundation for use in food equipment materials (NSF-51). It offers a temperature range use of -100° F (-73.3° C) to 400° F (204.4° C) and is sterilizable by autoclave or gamma irradiation.

DIMENSION (ID X OD X LGTH)	NUMBER
1/16" X 1/8" X 50'	SILC7001098
5/64" X 9/64" X 50'	SILC7001280
1/16" X 3/16" X 500'	SILC7001422
1/8" X 1/4" X 500'	SILC7001424
1/8" X 1/4" X 50'	SILC7001657
5/32" X 9/32" X 100'	SILC7001839
1/4" X 3/8" X 50'	SILC7002203
3/8" X 5/8" X 50'	SILC7002593
1/2" X 3/4" X 50'	SILC7002775
3/4" X 1" X 50'	SILC7003139

TUBING:
Pump Grade Platinum-Cured Silicone

DIMENSION (ID X OD X LGTH)	NUMBER
1/4" X 7/16" X 300'	SILPG7021041
3/16" X 3/8" X 250'	SILPG7021042
1/8" X 1/4" X 500'	SILPG7021043
1/16" X 3/16" X 500'	SILPG7021044

TUBING:
Polyurethane

80 Durometer: This tubing is extremely kink resistant and highly resistant to abrasion and scuffing. It is widely accepted as a premium material for medical device use because of its flexible nature and easy sterilizability.

Resin Grade & Technical Information: This tubing is extruded from Lubrizol's® compound Pellethane PUR80A-MED, which has passed USP Class VI testing. It can be sterilized by gamma radiation or Ethylene Oxide (EtO) gas. Please do not autoclave this tubing, as it may alter the chemical nature of the resin.

DIMENSION (ID X OD X LGTH)	NUMBER
1/16" X 1/8" X 100'	UR00-1031C
1/16" X 1/8" X 1000'	UR00-1031M
3/32" X 3/16" X 100'	UR00-2047C
3/32" X 3/16" X 1000'	UR00-2047M
1/8" X 1/4" X 100'	UR00-3062C
1/8" X 1/4" X 1000'	UR00-3062M

Ag-47 Clear Anti-Microbial: CLEARFLO Ag-47 is suitable for a wide variety of food contact applications and designed to protect your food or beverage product from harmful bacteria. It's manufactured from a clear flexible Shore A74 PVC compound with the addition of antimicrobial protection. The antimicrobial technology is especially effective where there is infrequent use of fluids in warm conditions.

Resin Grade & Technical Information: CLEARFLO Ag-47 offers antibacterial performance to ISO 22196:2011. Unlike some products CLEARFLO Ag-47 is 100% protected with the effective additive Silver-Magnesium-Aluminum-Phosphorus Complex. This is present throughout the tube and not only on the interior, offering you protection on all surfaces. The active ingredients in CLEARFLO Ag-47 will have no detrimental effect on yeasts used for brewing and will have no negative impact on the taste or flavor of products passed through it.

DIMENSION (ID X OD X LGTH)	NUMBER
3/16" X 5/16" X 100'	PVCA1190035
3/16" X 7/16" X 100'	PVCA1190049
1/4" X 3/8" X 100'	PVCA1190112
1/4" X 1/2" X 100'	PVCA1190126
5/16" X 7/16" X 100'	PVCA1190189
5/16" X 9/16" X 100'	PVCA1190203
3/8" X 1/2" X 100'	PVCA1190266
3/8" X 5/8" X 100'	PVCA1190280
1/2" X 5/8" X 100'	PVCA1190343
1/2" X 3/4" X 100'	PVCA1190357

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

OTHER

Components

Tapered Seal
PG-134



Clamps
PG-137



Injection Sites
PG-137



Spikes
PG-137



TSC Series

Tapered Seal

With TSC Series conical seal fittings, Nordson MEDICAL expands its range of products for small-bore connector applications in medical devices. TSC Series products offer the reliability, convenience and ease of use of luer fittings in a design that cannot be interconnected with conventional luers and avoid tubing misconnections. These luer alternatives are available in sizes to fit 1/8" (3.2 mm) and 3/16" (4.8 mm) ID flexible tubing.

Features & Benefits:

- Ergonomic grip
- Conical (tapered) seal
- 15% higher flow capacity than conventional luers
- Polypropylene resin
- Cannot interconnect with conventional luers taper fittings
- 500 Series single barbs

Need Samples for Evaluation?

You can call to request most of our standard product in quantities of 5, or less, for free application evaluation.

+1 (970) 267-5200 Toll Free: (888) 404-5837

Need it Fast?

You can order any of our 9,600+ parts online now with no order \$ minimums, 24 hours a day, 7 days a week!

www.nordsonmedical.com

SMALL BORE CONNECTORS:

Male

TSC2M210 (-6005)

Male Taper Seal Non-Luer Connector to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



TSC2M013 (-6005)

Male Taper Seal Non-Luer Connector to 500 Series Barb, 1/8" (3.2 mm) ID Tubing



TSC2M035 (-6005)

Male Taper Seal Non-Luer Connector to 500 Series Barb, 3/16" (4.8 mm) ID Tubing



SMALL BORE CONNECTORS:

Female

TSC2F210 (-6005)

Female Taper Seal Non-Luer Connector to 200 Series Barb, 1/16" (1.6 mm) ID Tubing



TSC2F007 (-6005)

Female Taper Seal Non-Luer Connector to 500 Series Barb, 3/32" (2.4 mm) ID Tubing



TSC2F013 (-6005)

Female Taper Seal Non-Luer Connector to 500 Series Barb, 1/8" (3.2 mm) ID Tubing



TSC2F035 (-6005)

Female Taper Seal Non-Luer Connector to 500 Series Barb, 3/16" (4.8 mm) ID Tubing



TSC2FB230 (-6005)

Female Taper Seal Non-Luer with 1/4-28 UNF Panel Mount to 200 Series Barb, 1/8" (3.2 mm) ID Tubing



SMALL BORE CONNECTORS:

Plugs and Caps

TSC2MP (-6005) 1,000 piece minimum

Male Taper Seal Non-Luer Plug with Tether Recess



Looking for Something Different?

Most of our best product ideas have come from you, our customers. While we pride ourselves on our extensive product offering, we recognize that your needs are ever-changing and that our products, too, must adapt to your needs.

If there's a new product you think Nordson MEDICAL should make, please send us an Email at designcenter@nordsonmedical.com describing your product concept in as much detail as you desire. We will contact you with any questions, and who knows - the next catalog edition could include products YOU told us to build.

See page 10 for more information on the Design Center!




LARGE BORE CONNECTORS:

Male


MLB035 (-1 -8012 -9)

Male Large Bore Non-Luer Connector to 500 Series Barb, 3/16" (4.8 mm) ID Tubing

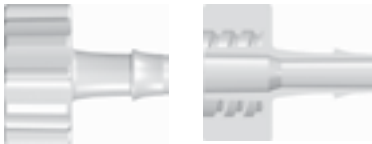


MLB055 (-1 -8012 -9)

Male Large Bore Non-Luer Connector to 500 Series Barb, 1/4" (6.4 mm) ID Tubing

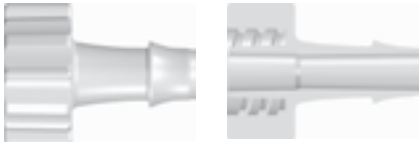

MLB065 (-1 -8012 -9)

Male Large Bore Non-Luer Connector to 500 Series Barb, 5/16" (8.0 mm) ID Tubing

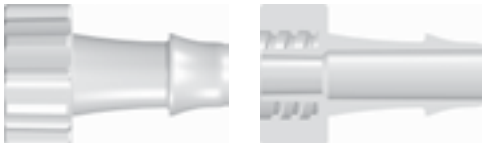

MLBC7 (-9002)

Male Large Bore Non-Luer Connector with Integral Lock Ring to Male Slip Luer (May be used with separate rotating lock ring; FSLLR)



MLBC8 (-9002)

Male Large Bore Non-Luer Connector with Integral Lock Ring to Female Thread-Style Locking Luer



LARGE BORE CONNECTORS:

Female


FLB035 (-1 -8012 -9)

Female Large Bore Non-Luer Connector to 500 Series Barb, 3/16" (4.8 mm) ID Tubing




FLB055 (-1 -2 -8012 -9)

Female Large Bore Non-Luer Connector to 500 Series Barb, 1/4" (6.4 mm) ID Tubing

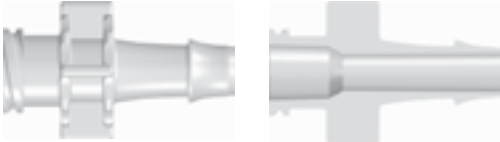

FLB065 (-1 -8012 -9)

Female Large Bore Non-Luer Connector to 500 Series Barb, 5/16" (8.0 mm) ID Tubing

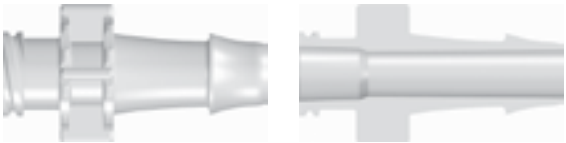

Technical Assistance

Trying to choose the correct part? Need more material information than is published in this catalog? Call Nordson MEDICAL Technical Support **(970) 267-5200**.

FLBP (-1 -8012 -9)

Female Large Bore Non-Luer Cap



MLBP (-1 -2 -8012 -9 -9002)

Male Large Bore Non-Luer Plug



VPM2930601N (White ABS)

Tube Clamp with Split, 1/8" (3.2 mm) OD Tubing



VPS3901074N (Polyisoprene with ABS Acrylic Body)

IN-4000 Intermittent Injection Site to Male Luer Lock



VPM1060200N (Clear Low Density Polyethylene)

Long Spike Guard Closed at Tip



VPM0480201N (White Polyethylene)

Luer Tip Syringe Cap



VPB0850400N (Clear ABS)

Female Luer Thread Style Coupler



VPM1080303N (Blue Polyethylene)

Male Luer Integral Lock Ring Plug



VPS5069200N (White ABS Spike with Clear LDPE Guard)

Bag Spike with Vented Guard (Replaces VPS5069700N)



Part Drawings Online!

Nordson MEDICAL provides dimensional drawings available for download for all of its fittings.

www.nordsonmedical.com

Custom Assembly

We produce standard components to fit your current assemblies and devices, but we also provide custom solutions by applying a variety of assembly options, such as joining fittings, connectors, filters, and other components to tubing. From systems using a single fitting to complex configurations utilizing various components, we can design complete assemblies to your requirements.

TECHNICAL

Information





Barb Dimensions PG-138

Barb Performance PG-140



Quick Connect Performance Ratings PG-141

BPF Series Performance PG-145

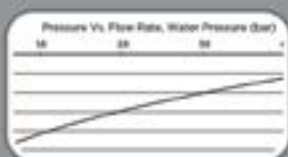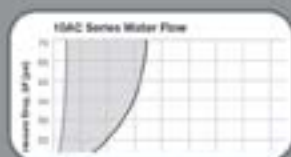

Panel Mount Instructions PG-145

Resin Information PG-146



CHEMICAL	%	TEMP °C	TEMP °F	APL
Acetic Acid	10	23	73	Satisfactory
Acetone	100	30	86	Satisfactory
Acrylonitrile	100	34	93	Satisfactory
Acrylonitrile	100	34	93	Satisfactory
Alk	100	82	180	Satisfactory
Ammonia, Liquid	100	34	93	Satisfactory
Ammonia, Gaseous	100	34	93	Satisfactory

Trademarks PG-149

Index PG-150



TECHNICAL: Barb Dimensions

The **CLASSIC SERIES** is characterized by a longer barb that extends the sealing surface of the fitting. Tubing is expanded from 42% to 69% above its nominal ID.



Barb Size	Barb OD	Barb ID	Tube Internal Diameter		
			Dec.	Frac.	Metric
10	0.105	0.046	0.063	1/16"	1.6 mm
20	0.145	0.060	0.094	3/32"	2.4 mm
30	0.187	0.078	0.125	1/8"	3.2 mm
40	0.234	0.100	0.156	5/32"	4.0 mm
50	0.280	0.130	0.188	3/16"	4.8 mm
60	0.355	0.187	0.250	1/4"	6.4 mm

The **200 SERIES** offers outstanding performance characteristics when used with flexible tubing. Tubing is expanded by 50% above its nominal ID, permitting high-pressure capabilities.



Barb Size	Barb OD	Barb ID	Tube Internal Diameter		
			Dec.	Frac.	Metric
210	0.094	0.047	0.063	1/16"	1.6 mm
220	0.141	0.070	0.094	3/32"	2.4 mm
230	0.188	0.094	0.125	1/8"	3.2 mm
240	0.234	0.117	0.156	5/32"	4.0 mm
250	0.282	0.141	0.188	3/16"	4.8 mm
260	0.375	0.188	0.250	1/4"	6.4 mm

The **300 SERIES** offers outstanding performance characteristics when used with flexible tubing. Tubing is expanded by 35% above its nominal ID, permitting high-pressure capabilities in the larger tube sizes.



Barb Size	Barb OD	Barb ID	Tube Internal Diameter		
			Dec.	Frac.	Metric
360	0.338	0.192	0.250	1/4"	6.4 mm

*Variations of barb IDs may exist. Please visit www.nordsonmedical.com to learn more about these alternatives. The data presented here is for reference only. It was compiled to provide our customers with a means of comparing the characteristics of components at the time of publication. The particular conditions of your use and application of our products are beyond our control. Thus, it is imperative that you test our products in your specific application to determine their suitability. All information is provided without implied or express warranty or guarantee by Nordson MEDICAL or other manufacturers. None of the information provided constitutes a recommendation or endorsement of any kind by Nordson MEDICAL.

The **400 SERIES** provides easy assembly while maintaining excellent pressure and tensile capabilities. Originally created for semi-rigid tubing, this series is preferred for applications where tensile/pressure properties are not highly critical, or when ease of assembly is desired. The 400 Series expands the tubing 25% above its nominal ID.



Barb Size	Barb OD	Barb ID	Tube Internal Diameter		
			Dec.	Frac.	Metric
410	0.078	0.046	0.063	1/16"	1.6 mm
420	0.117	0.068	0.094	3/32"	2.4 mm
430	0.156	0.092	0.125	1/8"	3.2 mm
440	0.206	0.125	0.156	5/32"	4.0 mm
445	0.213	0.125	0.170	-	4.3 mm
450	0.234	0.137	0.188	3/16"	4.8 mm
460	0.313	0.184	0.250	1/4"	6.4 mm

The **500 SERIES** profile is composed of blended radii, giving it the distinct look. Each barb is designed to accommodate a range of tube sizes and tube types. Consideration should be given to flow restriction, pressure capability, and the force required for assembly of the tube.



Barb Size	Barb OD	Barb ID*	Tube Internal Diameter		
			Dec.	Frac.	Metric
004	0.102	0.060	0.063	1/16"	1.6 mm
007	0.129	0.076	0.094	3/32"	2.4 mm
013	0.164	0.096	0.125	1/8"	3.2 mm
025	0.208	0.122	0.156	5/32"	4.0 mm
035	0.264	0.155	0.188	3/16"	4.8 mm
055	0.335	0.197	0.250	1/4"	6.4 mm
065	0.425	0.250	0.312	5/16"	8.0 mm
N070	0.483	0.262	0.375	3/8"	9.5 mm
L070	0.483	0.251	0.375	3/8"	9.5 mm
T070	0.483	0.290	0.375	3/8"	9.5 mm
Y070	0.483	0.290	0.375	3/8"	9.5 mm
N080	0.613	0.333	0.500	1/2"	12.7 mm
L080	0.613	0.319	0.500	1/2"	12.7 mm
T080	0.613	0.356	0.500	1/2"	12.7 mm
Y080	0.613	0.362	0.500	1/2"	12.7 mm
N090	0.718	0.390	0.625	5/8"	16.0 mm
L090	0.718	0.374	0.625	5/8"	16.0 mm
T090	0.718	0.374	0.625	5/8"	16.0 mm
Y090	0.718	0.390	0.625	5/8"	16.0 mm
N100	0.853	0.465	0.750	3/4"	19.0 mm
L100	0.853	0.443	0.750	3/4"	19.0 mm
T100	0.853	0.447	0.750	3/4"	19.0 mm
Y100	0.853	0.492	0.750	3/4"	19.0 mm
N110	1.125	0.623	1.000	1"	25.4 mm
L110	1.125	0.600	1.000	1"	25.4 mm
T110	1.125	0.600	1.000	1"	25.4 mm
Y110	1.125	0.623	1.000	1"	25.4 mm

The **600 SERIES** barbs are designed for use with silicone, C-Flex® and PVC tubing sizes commonly used in biopharm and other applications. Tubing is expanded from 18% to 34% above its nominal ID to maintain a low profile and ease of assembly, and the straight barb stem facilitates use of recommended cable tie retainers or hose clamps.



Barb Size	Barb OD	Barb ID*	Tube Internal Diameter		
			Dec.	Frac.	Metric
630	0.188	0.085	0.125	1/8"	3.2 mm
655	0.335	0.181	0.250	1/4"	6.4 mm
670	0.490	0.296	0.375	3/8"	9.5 mm
680	0.625	0.403	0.500	1/2"	12.7 mm
690	0.762	0.499	0.625	5/8"	16.0 mm
6100	0.900	0.615	0.750	3/4"	19.0 mm
6110	1.180	0.804	1.000	1"	25.4 mm

The **700 SERIES** barb is designed for use with Nordson MEDICAL's XQ Series quick connect fittings. This unique single barb is available in sizes to fit the 1/4" (6.4 mm), 5/16" (7.9 mm) and 3/8" (9.5 mm) ID tubing sizes and materials commonly used in disposable and reusable medical equipment. Reinforced PVC, unreinforced PVC and polyurethane tubing materials work well with this barb style.



Barb Size	Barb OD	Barb ID	Tube Internal Diameter		
			Dec.	Frac.	Metric
755	0.329	0.204	0.250	1/4"	6.4 mm
765	0.398	0.245	0.312	5/16"	7.9 mm
770	0.455	0.302	0.375	3/8"	9.5 mm

The **900 SERIES** offers outstanding performance characteristics when used with flexible tubing. Tubing is expanded from 8% to 17% above its nominal ID. This barb provides high flow efficiency and excellent tensile strength. A hose clamp is recommended for this series.



Barb Size	Barb OD	Barb ID	Tube Internal Diameter		
			Dec.	Frac.	Metric
970	0.437	0.269	0.375	3/8"	9.5 mm
975	0.573	0.358	0.500	1/2"	12.7 mm


Tensile pull and burst tests were performed by an independent laboratory under controlled conditions to provide comparative performance data for the various barb styles and sizes offered by Nordson MEDICAL. Graph numbers 1 through 8 present the data from each of the tests for each barb size and each barb style connected to PVC tube.

The values shown are the means of a statistically significant sample quantity of each configuration. Nordson MEDICAL barbed fittings specifically designed for use with the tubing size indicated were used in these tests. No oversized barbs, solvents, adhesives, or secondary fasteners of any kind were used.


Tensile Pull Test: Tensile properties were measured to determine the force required to remove the barbed fitting from the tubing. Tensile specimens were mounted axially in a holding fixture and tested at a rate of 20 inches per minute. Failure mode was tubing separating from the barb.

Hydrostatic Burst Test: Hydrostatic burst testing was performed in general accordance with ASTM D1599. Failure mode was either tubing separating from the barb, the tubing bursting, or a combination of the two.


Performance of Barb to Tube Connections - ■ Tensile Pull Test (lbs) - ■ Hydrostatic Burst Test (psi)

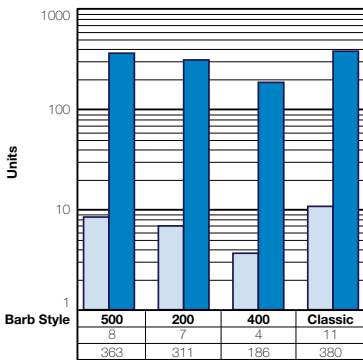


Graph 1: 1/16" Barb Performance

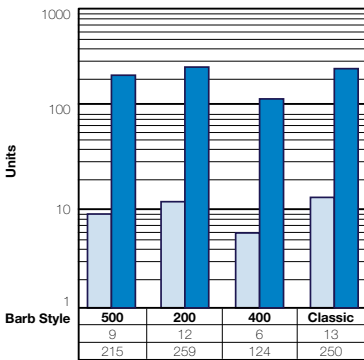


Graph 2: 3/32" Barb Performance

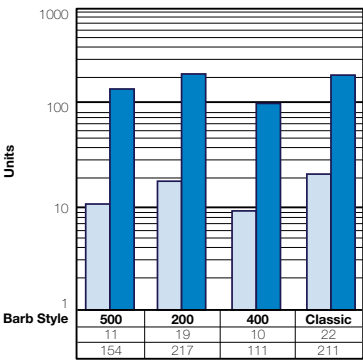


Graph 3: 1/8" Barb Performance

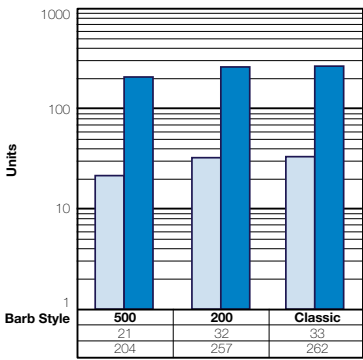



Graph 4: 5/32" Barb Performance

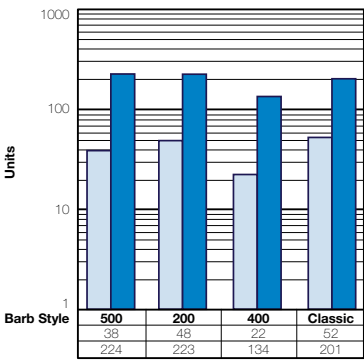

Graph 5: 3/16" Barb Performance




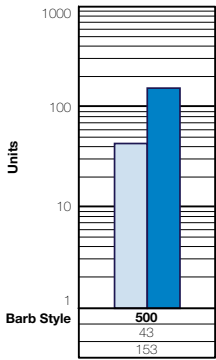


Graph 6: 1/4" Barb Performance






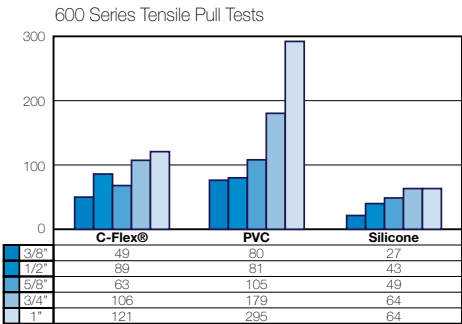
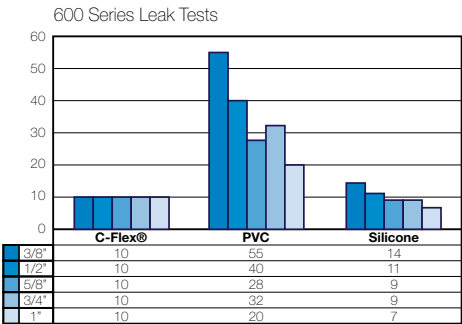


The data presented here is for reference only. It was compiled to provide our customers with a means of comparing the characteristics of components at the time of publication. The particular conditions of your use and application of our products are beyond our control. Thus, it is imperative that you test our products in your specific application to determine their suitability. All information is provided without implied or express warranty or guarantee by Nordson MEDICAL or other manufacturers. None of the information provided constitutes a recommendation or endorsement of any kind by Nordson MEDICAL.



10 SERIES - 1/16" Flow Size

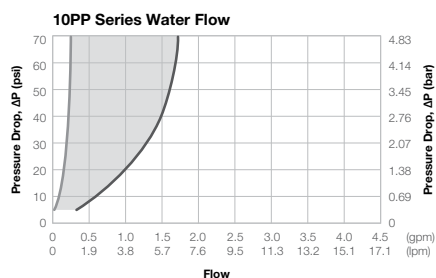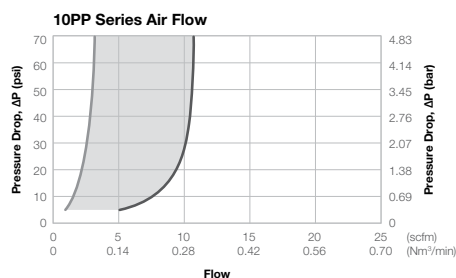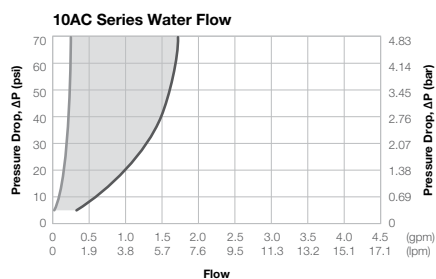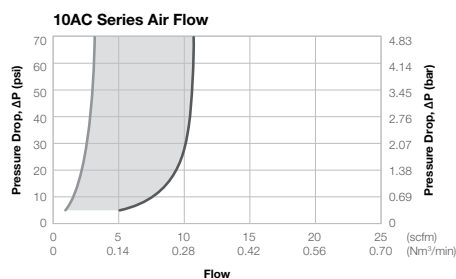








20 SERIES - 1/8" Flow Size

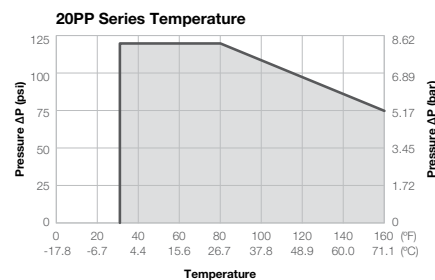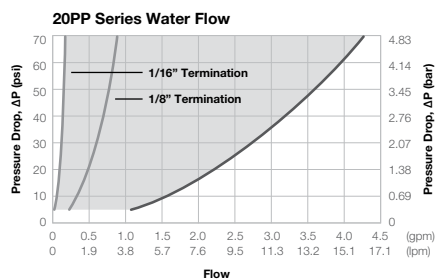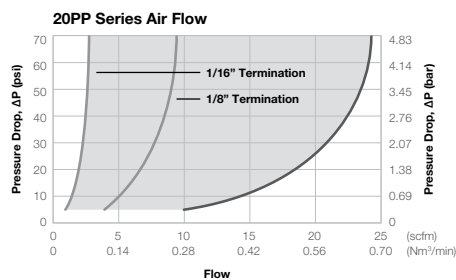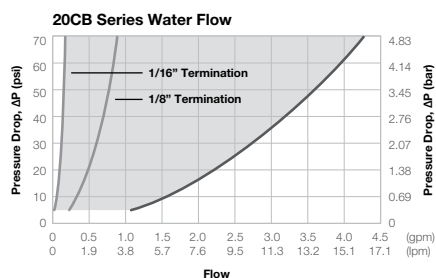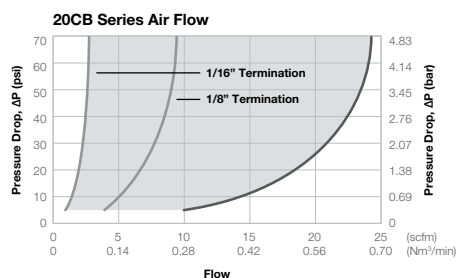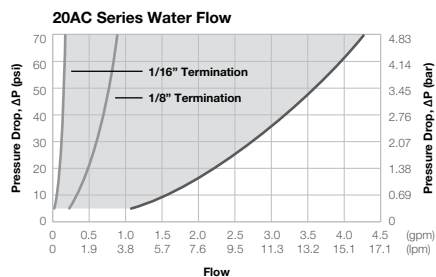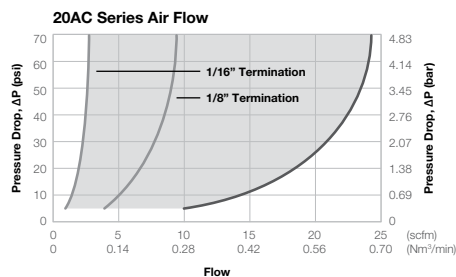



30 SERIES - 1/8" Flow Size

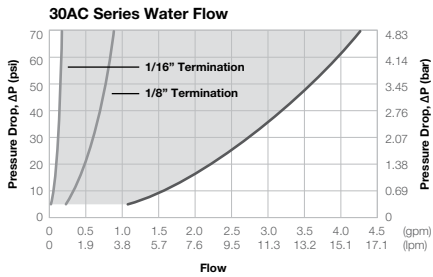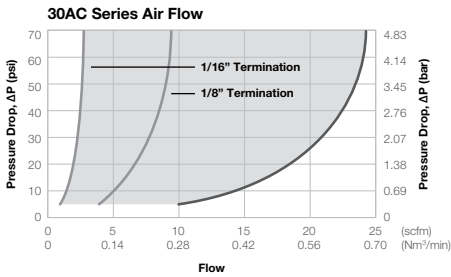









40 SERIES - 1/4" Flow Size

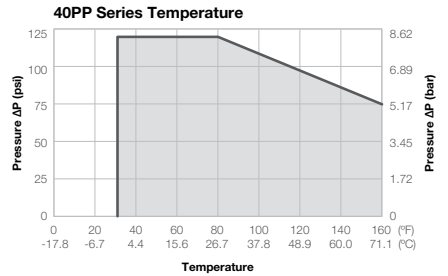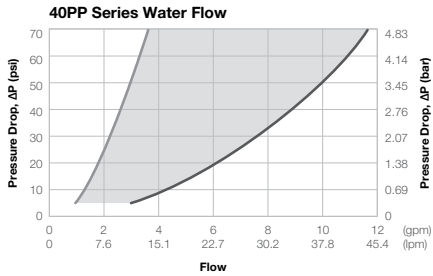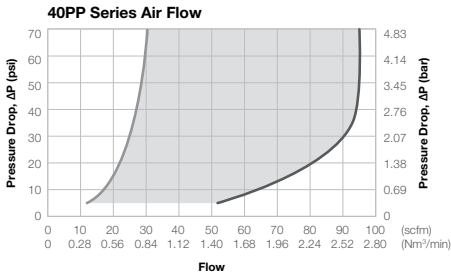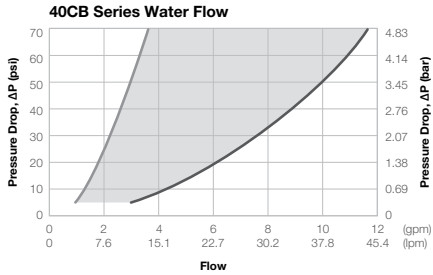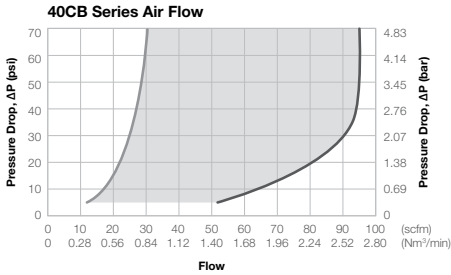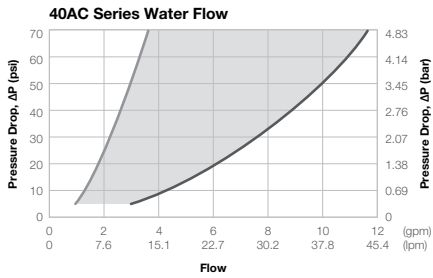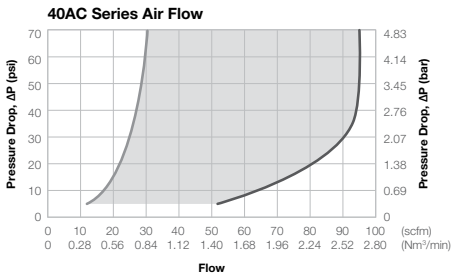



42 SERIES - 3/32" Flow Size

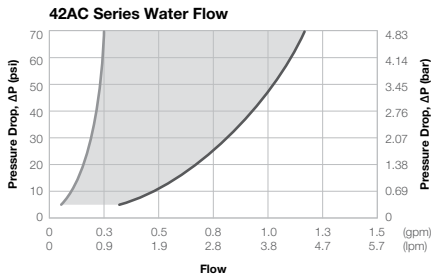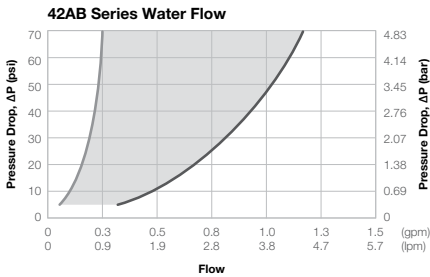

50 SERIES - 1/4" Flow Size

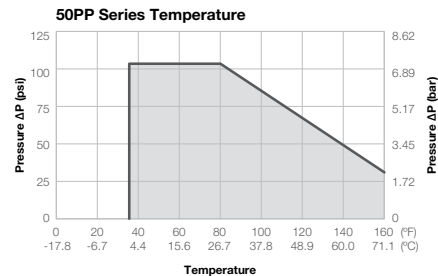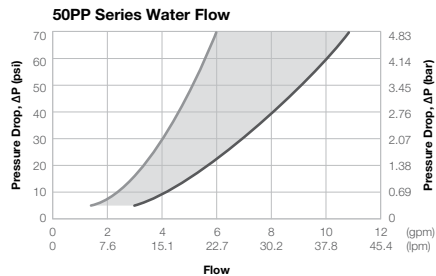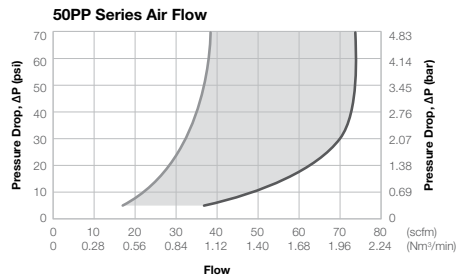
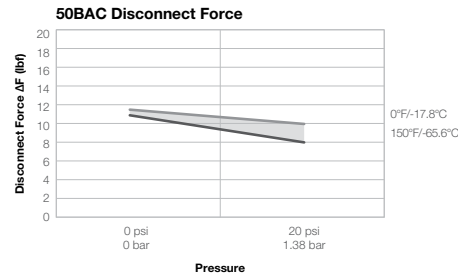
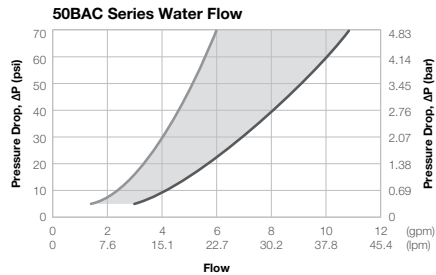
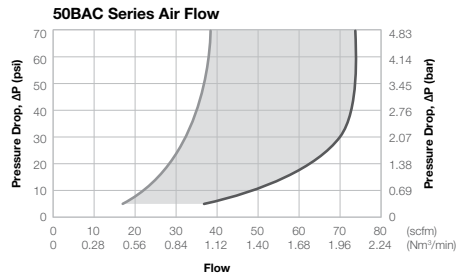
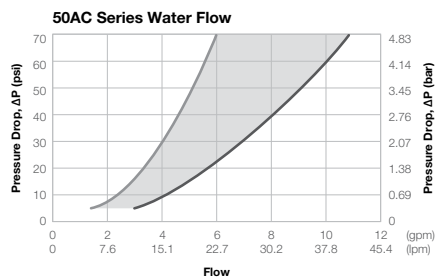
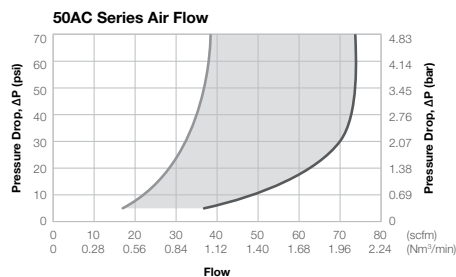


55 SERIES - 1/8" - 1/4" Flow Size

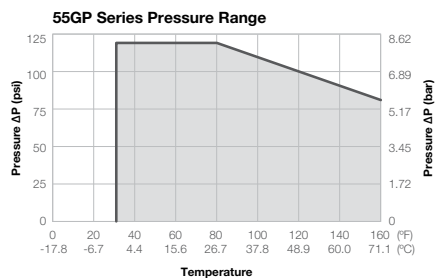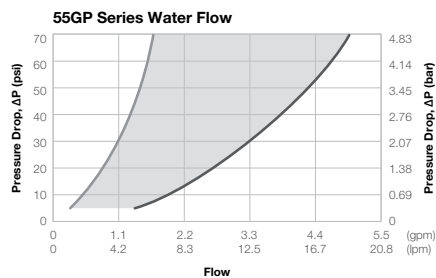


The data presented here is for reference only. It was compiled to provide our customers with a means of comparing the characteristics of components at the time of publication. The particular conditions of your use and application of our products are beyond our control. Thus, it is imperative that you test our products in your specific application to determine their suitability. All information is provided without implied or express warranty or guarantee by Nordson MEDICAL or other manufacturers. None of the information provided constitutes a recommendation or endorsement of any kind by Nordson MEDICAL.


60 SERIES - 3/8" Flow Size

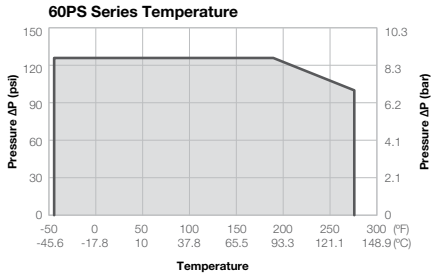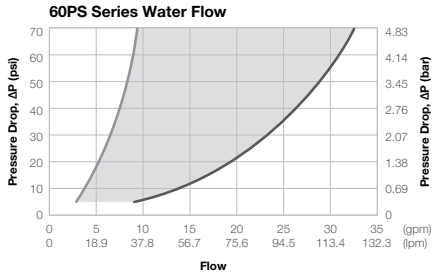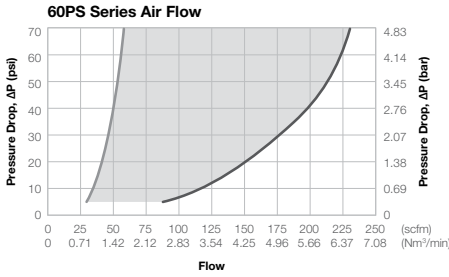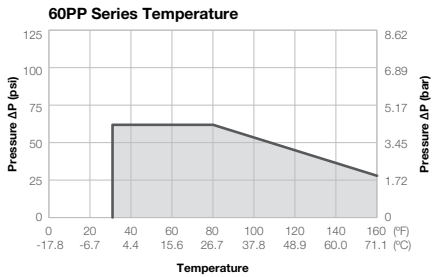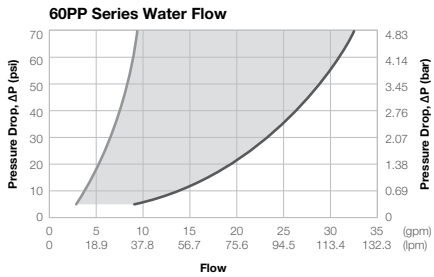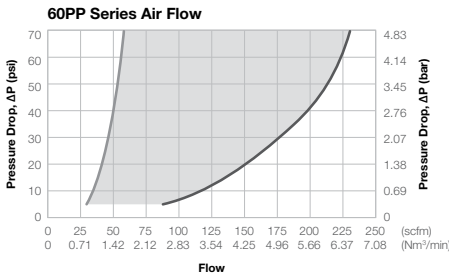





65 SERIES - 1/2" Flow Size

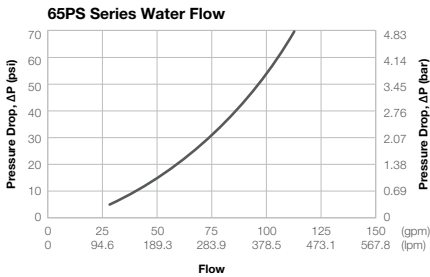

70 SERIES - 1/4" Flow Size

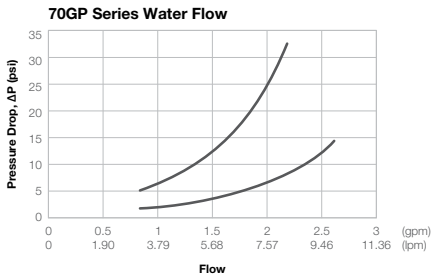

100 SERIES - 1/16" Flow Size

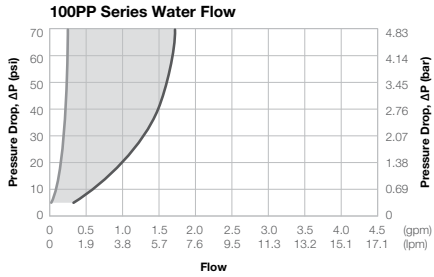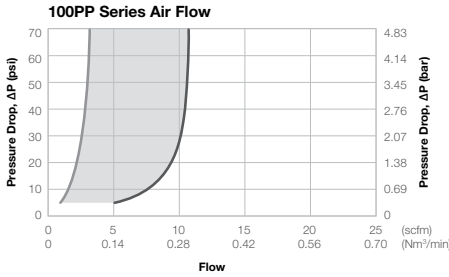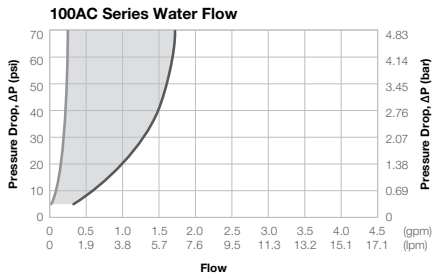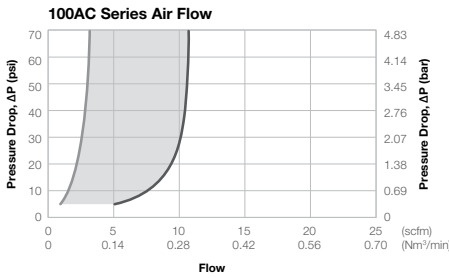

Maximum cycles _____ 50,000 connect-disconnect cycles


Maximum leak resistance _____ < .007 psi / 0.05 kPa / sec.

Maximum axial tensile resistance _____ 20 lbf.

Minimum Fitting ID _____ .113" / 2.9 mm

Maximum working pressure _____ 60 psi, 4.14 bar

Minimum Cv _____ .39

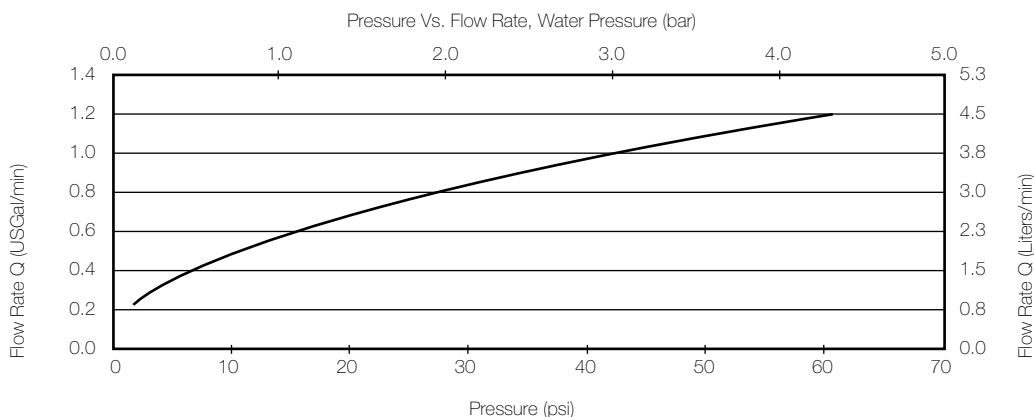

TECHNICAL: Panel Mount Instructions


Panel-mounted female luers, in both FTLLB (thread style) and FTLB (lug style) series, are shown beginning on page 20. When used with our CCLR lock rings or CCR color coding rings and LNS lock nuts (see pages 21, 37), color coding of connections is easy.

The information provided here will help you mount the panel mount luers. Please note that they are intended for panel thickness of 1/16" to 1/8" (1.5 to 3 mm).

FTLLB series shown with MLSL luer, FSLLR lock ring, CCLR color coding lock ring, and LNS lock nut.

Threaded Panel Openings: Tap the panel with threads from tube fittings, panel mounts or BPMBs and use a CCR color coded ring. In this case, the fitting tightened into the panel prevents fitting rotation. Use an LNS lock nut for additional security; tighten snug plus 1/4 turn only.

Non-Threaded Panel Openings: Use any one of the three panel patterns shown below with a CCLR lock ring and LNS lock nut. Together, the panel opening and CCLR prevent the panel mount fitting from rotating. Tighten LNS snug plus 1/4 turn only.

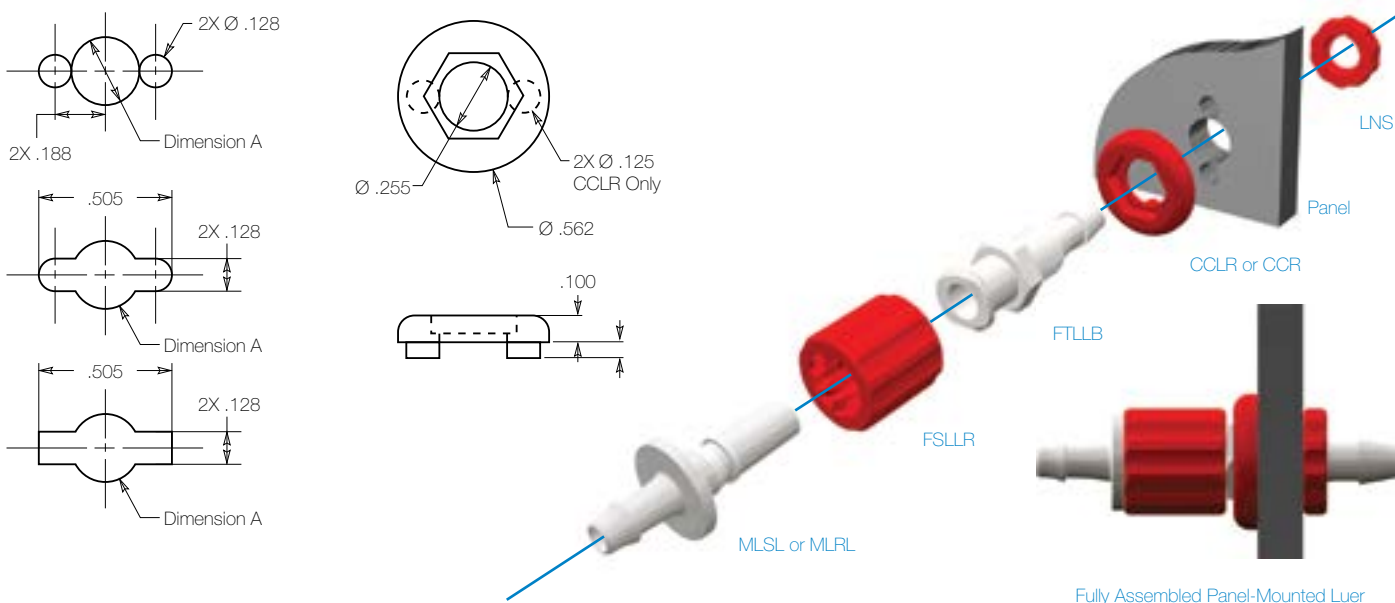

CHEMICAL	%	TEMP C	TEMP F	ABS	ACETAL	PVDF	NYLON	POLYCARBONATE	POLYPROPYLENE	POLYSULFONE
Acetic Acid	5%	23 C	73 F	Excellent	Unsatisfactory	Excellent	Satisfactory	Satisfactory	Excellent	Excellent
Acetone	100%	50 C	122 F	Unsatisfactory	Unsatisfactory	Unsatisfactory	Satisfactory	Unsatisfactory	Excellent	Unsatisfactory
Acetophenone	100%	24 C	75 F	Satisfactory	—	Unsatisfactory	Excellent	—	Satisfactory	—
Acetylene	100%	24 C	75 F	—	—	Excellent	Excellent	—	—	—
Air	100%	82 C	180 F	Excellent	Excellent	Excellent	Excellent	—	—	Excellent
Ammonia, Liquid	100%	24 C	75 F	Satisfactory	—	Unsatisfactory	Satisfactory	Unsatisfactory	Excellent	Excellent
Ammonium Hydroxide	10%	23 C	73 F	Satisfactory	—	Excellent	Excellent	Unsatisfactory	Excellent	Excellent
Ammonium Hydroxide	10%	70 C	158 F	Unsatisfactory	—	Excellent	Unsatisfactory	Unsatisfactory	Excellent	Excellent
Barium Sulfide	100%	24 C	75 F	Excellent	Excellent	Excellent	Satisfactory	—	Excellent	—
Benzene	100%	23 C	73 F	Satisfactory	Excellent	Excellent	Excellent	Unsatisfactory	Satisfactory	Unsatisfactory
Bleach	100%	23 C	73 F	—	Unsatisfactory	Excellent	Satisfactory	Unsatisfactory	Satisfactory	Excellent
Boric Acid	7%	35 C	95 F	Excellent	Excellent	Excellent	Unsatisfactory	Excellent	Excellent	—
Calcium Carbonate	Sat. sol.	24 C	75 F	Excellent	—	Excellent	—	—	Excellent	—
Carbon Dioxide	100%	50 C	122 F	Satisfactory	Excellent	Excellent	Excellent	—	Excellent	—
Carbon Monoxide	100%	50 C	122 F	Satisfactory	—	Excellent	Excellent	—	Excellent	—
Carbon Tetrachloride	100%	50 C	122 F	Unsatisfactory	—	Excellent	Excellent	Unsatisfactory	Unsatisfactory	Unsatisfactory
Chlorine Water	Dilute	23 C	73 F	Unsatisfactory	—	Excellent	Satisfactory	Unsatisfactory	Unsatisfactory	Satisfactory
Chlorine Water	Concen.	23 C	73 F	Unsatisfactory	—	Excellent	Unsatisfactory	Unsatisfactory	Unsatisfactory	Unsatisfactory
Chlorobenzene	100%	23 C	73 F	Satisfactory	—	Excellent	Excellent	Unsatisfactory	Unsatisfactory	Unsatisfactory
Chlorofluorocarbon 11	100%	24 C	75 F	—	—	Excellent	Excellent	Satisfactory	—	Excellent
Chloroform	100%	23 C	73 F	Unsatisfactory	Satisfactory	Excellent	Satisfactory	Unsatisfactory	Unsatisfactory	Unsatisfactory
Cyclohexanone	100%	24 C	75 F	Unsatisfactory	—	Excellent	Excellent	Unsatisfactory	Satisfactory	Unsatisfactory
Dichlorethylene	100%	23 C	73 F	—	—	Excellent	Satisfactory	—	Excellent	Unsatisfactory
Ethanol	95%	50 C	122 F	Satisfactory	—	Excellent	Excellent	Satisfactory	Excellent	Satisfactory
Ethyl Acetate	95%	50 C	122 F	Satisfactory	—	Unsatisfactory	Excellent	Unsatisfactory	Satisfactory	Unsatisfactory
Ethylene Glycol	100%	23 C	73 F	Excellent	Satisfactory	Excellent	Excellent	Satisfactory	Excellent	Excellent
Ethylene Oxide	100%	24 C	75 F	Unsatisfactory	—	Excellent	Satisfactory	Satisfactory	Satisfactory	Excellent
Ethylene Oxide	100%	79 C	175 F	Unsatisfactory	—	Excellent	Unsatisfactory	Satisfactory	Unsatisfactory	Excellent
Fatty Acids	—	—	—	—	Excellent	Excellent	—	Satisfactory	Excellent	—
Fluorine	100%	23 C	73 F	Unsatisfactory	—	Excellent	Unsatisfactory	—	—	—
Formaldehyde	37%	24 C	75 F	Unsatisfactory	Excellent	Excellent	—	Unsatisfactory	Excellent	Unsatisfactory
Gasoline	100%	85 C	185 F	Excellent	Satisfactory	Excellent	Excellent	Unsatisfactory	Satisfactory	Satisfactory
Glucose	Concen.	24 C	75 F	Excellent	—	Excellent	—	—	Excellent	—
Glycerin	100%	24 C	75 F	Excellent	—	Excellent	—	Excellent	Excellent	Excellent
Hydrochloric Acid	2%	23 C	73 F	Excellent	Satisfactory	Excellent	Excellent	Excellent	Excellent	Excellent
Hydrochloric Acid	10%	25 C	77 F	Excellent	Satisfactory	Excellent	Unsatisfactory	Excellent	Excellent	Excellent
Hydrofluoric Acid	10%	23 C	73 F	Satisfactory	—	Excellent	Unsatisfactory	—	Excellent	Excellent
Hydrogen Peroxide	1%	24 C	75 F	Excellent	Unsatisfactory	Excellent	Satisfactory	Excellent	Excellent	Excellent
Hydrogen Peroxide	5%	43 C	110 F	Satisfactory	Unsatisfactory	Excellent	Unsatisfactory	Excellent	Satisfactory	Excellent
Isopropanol	70%	23 C	73 F	—	Excellent	Excellent	Excellent	Call*	Excellent	Satisfactory
Kerosene	100%	85 C	185 F	Satisfactory	—	Excellent	Excellent	Satisfactory	Satisfactory	Satisfactory
Methyl Ethyl Ketone	100%	50 C	122 F	Unsatisfactory	—	Unsatisfactory	Excellent	Unsatisfactory	Satisfactory	Unsatisfactory
Methylene Chloride	100%	23 C	73 F	Unsatisfactory	—	Excellent	Satisfactory	Unsatisfactory	Excellent	Unsatisfactory
Methanol	100%	23 C	73 F	Unsatisfactory	—	Excellent	Excellent	Satisfactory	Excellent	Satisfactory
Nitric Acid	10%	23 C	73 F	Satisfactory	—	Excellent	Unsatisfactory	Unsatisfactory	Excellent	Satisfactory
Oxygen	100%	24 C	75 F	—	—	Excellent	Satisfactory	—	—	—
Ozone	100%	43 C	110 F	Satisfactory	—	Satisfactory	Unsatisfactory	Unsatisfactory	—	—
Phenol	90%	23 C	73 F	Unsatisfactory	Unsatisfactory	Excellent	Unsatisfactory	—	Excellent	—
Phosphoric Acid	5%	98 C	208 F	Satisfactory	Unsatisfactory	Excellent	Unsatisfactory	Unsatisfactory	Excellent	Excellent
Propane	100%	23 C	73 F	Satisfactory	—	Excellent	Excellent	—	—	—
Sodium Bicarbonate	Concen.	24 C	75 F	Excellent	Excellent	Excellent	Excellent	—	Excellent	Excellent
Sodium Chloride	10%	23 C	73 F	Excellent	—	Excellent	Excellent	—	Excellent	Excellent
Sodium Chloride	Sat. sol.	24 C	75 F	Excellent	—	Excellent	Excellent	—	Excellent	Excellent
Sodium Hydroxide	10%	70 C	158 F	Satisfactory	—	Excellent	Satisfactory	—	Excellent	Excellent
Steam	—	120 C	248 F	Unsatisfactory	—	Excellent	Unsatisfactory	Unsatisfactory	Satisfactory	Excellent
Sulfuric Acid	30%	23 C	73 F	Satisfactory	Unsatisfactory	Excellent	Unsatisfactory	Excellent	Excellent	Excellent
Tetrahydrofuran	100%	23 C	73 F	Satisfactory	—	Unsatisfactory	Excellent	—	Unsatisfactory	—
Toluene	100%	50 C	122 F	Satisfactory	Unsatisfactory	Excellent	Excellent	Unsatisfactory	Unsatisfactory	Unsatisfactory
Trichloroethylene	100%	23 C	73 F	Satisfactory	—	Excellent	Satisfactory	Unsatisfactory	Unsatisfactory	Unsatisfactory
Water	100%	79 C	175 F	Excellent	—	Excellent	Excellent	Unsatisfactory	Excellent	Excellent

* Polycarbonate is known to experience environmental stress cracking when exposed to solvents. Please contact Nordson MEDICAL technical support if polycarbonate is planning on being used in a solvent bonding application. Looking for a stress-free connection? Ask our technical experts at 1-888-404-5837 how annealing and other secondary processes can alleviate the risk of environmental stress cracking in polycarbonate materials.

The data presented here is for reference only. It was compiled to provide our customers with a means of comparing the characteristics of components at the time of publication. The particular conditions of your use and application of our products are beyond our control. Thus, it is imperative that you test our products in your specific application to determine their suitability. All information is provided without implied or express warranty or guarantee by Nordson MEDICAL or other manufacturers. None of the information provided constitutes a recommendation or endorsement of any kind by Nordson MEDICAL.

Material - Manufacturer & Grade		VP Code	Regulatory Compliance of Base Resin*		RoHS Compliant	BPA Free	DEHP Free	Animal Derivative Free*
ABS (Lace White) - INEOS Lustran 348-000000		-81		USP Class VI Certified; ISO 10993 Part I Compliant	Yes	Yes	Yes	No
ABS, MABS (Methyl Methacrylate) - BASF Terluc 2802 HD		SCV Series		USP Class VI Certified	Yes	Yes	Yes	No
ABS, MABS (Methyl Methacrylate) - BASF Terluc 2802 TR		-8003		General Purpose	Yes	Yes	Yes	No
ABS (Sno White) - INEOS Lustran 348-012002		-8012		USP Class VI Certified; ISO 10993 Part I Compliant	Yes	Yes	Yes	No
Acetal - DuPont Delrin 500P		-10		21 CFR 177.2480 FDA Approved	Yes	Yes	Yes	No
Acetal Copolymer - BASF Ultratorm N 23020 003		-1006, -1007		21 CFR 177.2470 FDA Approved	Yes	Yes	Yes	No
Acrylic - Cyro Cyrolite Med 2		-50		USP Class VI Certified; ISO 10993 Compliant; 21 CFR 177.1010 FDA Approved	Yes	Yes	Yes	Yes
Nylon (PA 66) - DuPont Zytel 101F		-0, -1, -2, -3, -4, -5, -06, -07, -08B, -8, -8B, -CM030		21 CFR 177.1500 FDA Approved	Yes	Yes	Yes	No
Nylon (PA 66) - DuPont Zytel SC310 NC010		-U1, -U5, -U01, -U05		USP Class VI Certified; ISO 10993 Compliant; 21 CFR 177.1500 FDA Approved	Yes	Yes	Yes	No
Nylon (Glass Reinforced) - DuPont Zytel FGF5171		-U71, -U75		USP Class VI Certified; ISO 10993 Compliant; 21 CFR 177.1500 FDA Approved	Yes	Yes	Yes	No
Nylon (Glass Reinforced) - DuPont Zytel 70G33L		-7, -72		21 CFR 177.1500 FDA Approved	Yes	Yes	Yes	No
Polycarbonate (Clear) - Bayer Makrolon 2558 550115		-9		USP Class VI Certified; ISO 10993 Compliant; 21 CFR 177.1580 FDA Approved	Yes	No	Yes	No
Polycarbonate (Clear) - Bayer Makrolon Rx 1805 451118		-9010, SCV Series		USP Class VI Certified; ISO 10993 Compliant	Yes	No	Yes	Yes
Polycarbonate (Tinted) - Sabic Lexan HPS1		-9024		USP Class VI Certified; ISO 10993 Compliant	Yes	No	Yes	No
Polycarbonate (Tinted) - Styron Calibre Megarad 2081-15		-9002, DCV Series		ISO 10993 Compliant	Yes	No	Yes	No
Polyethylene - LyondellBasell Purell PE 1840 H		SCV Series		Please call for details	Yes	Yes	Yes	Yes
Polyethylene - Trilliant HC5420-0002 LL Natural Resin		-VP1		USP Class VI Certified	Yes	Yes	Yes	Yes
Polypropylene - Flint Hills P5M6K-080		-6005, -6006		USP Class VI Certified; 21 CFR 177.1520 FDA Approved	Yes	Yes	Yes	Yes
Polysulfone - Solvay Udel P1700 CL2611		-40		USP Class VI Certified; ISO 10993 Compliant; 21 CFR 177.1655 FDA Approved	Yes	Yes	Yes	Yes
PVDF - Arkema Kynar 1000 HD		-J1A		USP Class VI Certified; 21 CFR 177.2510 FDA Approved	Yes	Yes	Yes	Yes
SAN - BASF 388S TR 77770		SCV Series		ISO 10993 Compliant	Yes	Yes	Yes	No
Silicone - Wacker Elastosil LR 3003		SCV Series		USP Class VI Certified	Yes	Yes	Yes	Yes

*As determined by resin manufacturer

Material - Manufacturer & Grade	VP Code	Gamma Radiation Sterilization*	Ethylene Oxide Sterilization	Autoclave Sterilization
ABS (Lace White) - INEOS Lustran 348-000000	-81	Compatible up to 5 Mrads	Good	Not suitable due to low heat deflection temperatures
ABS, MABS (Methyl Methacrylate) - BASF Terluc 2802 TR/HD	-8003, SCV Series	Withstands up to 5 Mrad doses without yellowing	Good	Poor; components may distort due to low heat deflection temperatures
ABS (Sno White) - INEOS Lustran 348-012002	-8012	Compatible up to 5 Mrads	Good	Not suitable due to low heat deflection temperatures
Acetal - DuPont Delfin 500P	-10	Compatible up to 1 Mrad	Excellent	Very good
Acetal Copolymer - BASF Ultralorm N 2320 003	-1006, -1007	Compatible up to 1.5 Mrads	Excellent	Very good
Acrylic - Cyro Cyrolite Med 2	-50	Very good up to commonly used doses (6 Mrads)	Excellent	Not recommended
Nylon (PA 66) - DuPont Zytel 101F	-0, -1, -2, -3, -4, -5, -06, -07	Very good; may discolor to brownish hue	Very good; some susceptibility to oxidizing agents	Very good; product may swell slightly due to water absorption
Nylon (PA 66) - DuPont Zytel SC310 NC010	-U1, U5, -U01, -U05	Very good; may discolor to brownish hue	Very good; some susceptibility to oxidizing agents	Very good; product may swell slightly due to water absorption
Nylon (Glass Reinforced) - DuPont Zytel FGF5171	-U71, -U75	Very good; may discolor to brownish hue	Very good; some susceptibility to oxidizing agents	Excellent
Nylon (Glass Reinforced) - DuPont Zytel 70G33L	-72	Very good; may discolor to brownish hue	Very good; some susceptibility to oxidizing agents	Excellent
Polycarbonate (Clear) - Bayer Makrolon 2558 5501 15	-9	Compatible up to 10 Mrads with minor loss of physical properties; will discolor to yellow-green hue	Highly compatible	Poor; may craze or stress crack due to molding stresses
Polycarbonate (Clear) - Bayer Makrolon FK 1805 4511 18	-9010, SCV Series	Compatible up to 10 Mrads with minor loss of physical properties; may discolor slightly dependent upon dose/number of cycles	Highly compatible	Poor; may craze or stress crack due to molding stresses
Polycarbonate (Tinted) - Sabic Lexan HPS1	-9024	Highly compatible with reduced color shift when compared to unstabilized polycarbonate resins	Highly compatible	Limited to 121°C for 20 minutes due to low heat deflection temperatures
Polycarbonate (Tinted) - Syron Calibre Megarad 2081-15	-9002, DCV Series	Excellent up to 10 Mrads with minor loss of physical properties; light violet hue turns clear upon sterilization	Highly compatible	Not recommended
Polyethylene - LyondellBasell Purell PE 1840 H	SCV Series	Compatible up to 5 Mrads	Highly compatible	Not recommended
Polyethylene - Trilliant HC5420-0002 LL Natural Resin	-VP1	Highly compatible	Highly compatible	Not recommended
Polypropylene - Flint Hills P5M6K-080	-6005, -6006, -6012	Very good to commonly used sterilization doses (6 Mrads)	Highly compatible	Limited to 121°C for 20 minutes due to low heat deflection temperatures
Polysulfone - Solvay Udel P1 700 QL2611	-40	Highly compatible; will discolor to brownish hue	Excellent	Excellent
PVDF - Arkema Kyrar 1000 HD	-J1A	Highly compatible; will discolor to brownish hue	Excellent	Highly compatible
SAN - BASF 388S TR 77770	SCV Series	Compatible up to 5 Mrads	Good	Not suitable due to low heat deflection temperatures
Silicone - Wacker Elastosil LR 3003	SCV Series	Compatible up to 5 Mrads	Highly compatible	Very good

*1 Megarad (Mrad) = 10 kiloGrays (kGy)

TECHNICAL:

Common Resin Sterilization Stability

Property	ASTM Test Method	PVDF	Polycarb	Radiation Stable PC	Polypro	Nylona	GF Nylona	ABS	MABS	Polysulfone	Acetal
Tensile Strength	D638	5,400	10,000	9,800	—	12,000	27,000	4,800	—	—	—
Elongation at Break, %	D638	100-400	120	150	375	300	3	50	12	50-100	40
Tensile Yield Strength, psi	D638	6,500	9,200	9,800	4,400	12,000	—	5,700	6,960	10,200	9,900
Flexural Modulus, kpsi	D790	360	330	350	170	410	1,300	330	—	390	450
Rockwell Hardness D	D785	—	R118	R118	R80-102	R121	—	R107	—	—	R122
Rockwell Hardness M	—	—	M75	M72	—	M79	M101	—	—	M69	M89
Shore Hardness D	D2240	D76-80	—	—	—	—	—	—	—	—	—
Coef. of Linear 10-5 in/in/°F	D696	7.6	3.9	3.8	—	4	1.3	5.5	—	3.1	—
Deflection Temp°F @264psi	D648	244	266	252	—	194	480	173	194	345	216
Deflection Temp°F @66psi	D648	—	273	—	190	455	—	—	201	358	332
Water Absorption, %, 24 h	D570	0.015	0.15	0.15	0	1.2	0.7	0.25-0.40	0.7	0.3	0.4
Izod Impact, (Notched) ft-lb/in	D256	3.1	16.0	14.0	0.5	1.0	2.2	6.0	1.31	1.3	—

The data presented here is for reference only. It was compiled to provide our customers with a means of comparing the characteristics of components at the time of publication. The particular conditions of your use and application of our products are beyond our control. Thus, it is imperative that you test our products in your specific application to determine their suitability. All information is provided without implied or express warranty or guarantee by Nordson MEDICAL or other manufacturers. None of the information provided constitutes a recommendation or endorsement of any kind by Nordson MEDICAL.

TECHNICAL: Trademarks

Trademark 1,428,380



Corporate Trademark. The stylized design associated with Value Plastics® products and services.

To protect our intellectual property and our customers, it is important to be aware of the patents and trademarks that apply to our products. Look for the Value Plastics logo (V) to ensure you do not mistakenly use a poor quality counterfeit of a Value Plastics® fitting, thereby jeopardizing the integrity and quality of your product.

BioClamp® is a registered trademark of BioPure Technology Ltd.
 BioValve® is a registered trademark of BioPure Technology Ltd.
 C-Flex® is a registered trademark of Consolidated Polymer Technologies, Inc.
 CYROLITE® Protect is a registered trademark of Evonik Corporation
 Dinamap® is a registered trademark of GE Healthcare
 Equistar® is a registered trademark of Equistar Chemicals
 Kynar® is a registered trademark of Atochem Atofina North America, Inc.
 Lubrizol® is a registered trademark of The Lubrizol Corporation
 Teknor® is a registered trademark of Teknor Apex Company
 Tuf-Flex® is a registered trademark of Rubber Fab Technologies Group

One or more components contained in this catalog may be covered by B. Braun OEM Division patents and trademarks. Patent numbers: 4,246,932; 4,535,820; 4,683,916.

©2020 Nordson MEDICAL. All rights reserved. No part of this catalog may be reproduced without the express consent of Nordson MEDICAL. All terms are current as of April, 2020 and are subject to change without notice.

Trademark 1,887,178

Value Plastics®

Corporate Trademark.

Trademark 3,838,623

BUILD-A-PART® Trademark for our non-metal bondable fittings

OTHER TRADEMARKS

1,827,467
 1,817,921
 1,817,922
 2,655,263
 2,616,170
 2,535,232
 2,535,234

3PF210.....	34	10ACV2-01.....	45	20AC-PB8-04.....	49
3PF220.....	34	10ACV2-02.....	45	20AC-PB9-04.....	47
3PF230.....	34	10ACV2-02BLK.....	45	20AC-PB9-025.....	47
4PDL210.....	35	10ACV2-M3.....	45	20AC-PB10-04.....	49
4PDL220.....	35	10ACV2-M3BLK.....	45	20AC-PB10-025.....	49
4PDL230.....	35	10ACV-SB2-01.....	44	20AC-S1-02.....	46
4PLL210.....	35	10ACV-SB2-01NS.....	44	20AC-S1-02BSPT.....	46
4PLL220.....	35	10ACV-SB2-01NSBLK.....	44	20AC-S1-04.....	46
4PLL230.....	35	10ACV-SB2-02.....	44	20AC-S1-04BSPT.....	46
4PX210.....	34	10ACV-SB2-02BLK.....	44	20AC-S2-01.....	46
4PX220.....	34	10ACV-SB2-02NS.....	44	20AC-S2-02.....	46
4PX230.....	34	10ACV-SB2-02NSBLK.....	44	20AC-S2-03.....	46
4PX655.....	34	10ACV-SB2-M3.....	44	20AC-S2-04.....	46
4PX670.....	34	10ACV-SB2-M3BLK.....	44	20AC-S3-01.....	48
4PX670/655.....	35	10ACV-SB2-M3NS.....	44	20AC-S3-02.....	48
4PX680.....	34	10ACV-SB2-M3NSBLK.....	44	20AC-S3-02GS.....	48
4PX680/670.....	35	10ACV-SB3-01.....	45	20AC-S3-03.....	48
6PDT210.....	35	10ACV-SB3-01BLK.....	45	20AC-S3-03GS.....	48
6PDT220.....	35	10ACV-SB3-02.....	45	20AC-S3-04.....	48
6PDT230.....	35	10ACV-SB3-02BLK.....	45	20AC-S3-04GS.....	48
6PTT210.....	35	10ACV-SB3-M3.....	45	20AC-S4-02.....	49
6PTT220.....	35	10ACV-SB3-M3BLK.....	45	20AC-S4-04.....	49
6PTT230.....	35	10PP2-01.....	45	20AC-S8-04.....	48
10AC2-01.....	45	10PP2-02.....	45	20AC-S8-025.....	48
10AC2-01BLK.....	45	10PP2-M3.....	45	20AC-S8-025GS.....	48
10AC2-02.....	45	10PP-PE2-01.....	43	20AC-S9-04.....	46
10AC2-02BLK.....	45	10PP-PE2-01NS.....	43	20AC-S9-025.....	46
10AC2-M3.....	45	10PP-PE2-02.....	43	20AC-S10-04.....	49
10AC2-M3BLK.....	45	10PP-PE2-02NS.....	43	20AC-S19-1032.....	48
10AC-PB2-01.....	43	10PP-PE2-M3.....	43	20ACV-PB1-02.....	51
10AC-PB2-01BLK.....	43	10PP-PE2-M3NS.....	43	20ACV-PB1-02BSPT.....	51
10AC-PB2-01NS.....	43	10PP-PE3-01.....	44	20ACV-PB1-04.....	51
10AC-PB2-01NSBLK.....	43	10PP-PE3-02.....	44	20ACV-PB1-04BSPT.....	51
10AC-PB2-02.....	43	10PP-S2-01.....	43	20ACV-PB2-01.....	51
10AC-PB2-02BLK.....	43	10PP-S2-01NS.....	43	20ACV-PB2-02.....	51
10AC-PB2-02NS.....	43	10PP-S2-02.....	43	20ACV-PB2-03.....	51
10AC-PB2-02NSBLK.....	43	10PP-S2-02NS.....	43	20ACV-PB2-04.....	51
10AC-PB2-M3.....	43	10PP-S2-M3.....	43	20ACV-PB3-01.....	53
10AC-PB2-M3BLK.....	43	10PP-S2-M3NS.....	43	20ACV-PB3-02.....	53
10AC-PB2-M3NS.....	43	10PP-S3-01.....	44	20ACV-PB3-03.....	53
10AC-PB2-M3NSBLK.....	43	10PP-S3-02.....	44	20ACV-PB3-04.....	53
10AC-PB3-01.....	44	10PP-S3-M3.....	44	20ACV-PB4-02.....	53
10AC-PB3-01BLK.....	44	14MTLL.....	19	20ACV-PB4-04.....	53
10AC-PB3-02.....	44	18FTLL.....	20	20ACV-PB8-04.....	53
10AC-PB3-02BLK.....	44	18L230.....	127	20ACV-PB9-04.....	51
10AC-PB3-M3.....	44	18L250.....	127	20ACV-PB9-025.....	51
10AC-PB3-M3BLK.....	44	18L260.....	127	20ACV-PB10-04.....	53
10AC-S2-01.....	42	18LFTLL.....	20	20ACV-PB10-025.....	53
10AC-S2-01BLK.....	42	18MTLL.....	19	20ACV-SB1-02.....	50
10AC-S2-01NS.....	42	18T250.....	127	20ACV-SB1-02BSPT.....	50
10AC-S2-01NSBLK.....	42	18T260.....	127	20ACV-SB1-02UT.....	50
10AC-S2-02.....	42	20AC-PB1-02.....	47	20ACV-SB1-04.....	50
10AC-S2-02BLK.....	42	20AC-PB1-02BSPT.....	47	20ACV-SB1-04BSPT.....	50
10AC-S2-02NS.....	42	20AC-PB1-04.....	47	20ACV-SB1-06.....	50
10AC-S2-02NSBLK.....	42	20AC-PB1-04BSPT.....	47	20ACV-SB2-01.....	50
10AC-S2-M3.....	42	20AC-PB2-01.....	47	20ACV-SB2-02.....	50
10AC-S2-M3BLK.....	43	20AC-PB2-02.....	47	20ACV-SB2-03.....	50
10AC-S2-M3NS.....	43	20AC-PB2-03.....	47	20ACV-SB2-04.....	50
10AC-S2-M3NSBLK.....	43	20AC-PB2-04.....	47	20ACV-SB3-01.....	52
10AC-S3-01.....	44	20AC-PB3-01.....	48	20ACV-SB3-01GS.....	52
10AC-S3-01BLK.....	44	20AC-PB3-02.....	48	20ACV-SB3-02.....	52
10AC-S3-02.....	44	20AC-PB3-03.....	49	20ACV-SB3-02GS.....	52
10AC-S3-02BLK.....	44	20AC-PB3-04.....	49	20ACV-SB3-03.....	52
10AC-S3-M3.....	44	20AC-PB4-02.....	49	20ACV-SB3-03GS.....	52
10AC-S3-M3BLK.....	44	20AC-PB4-04.....	49	20ACV-SB3-04.....	52

20ACV-SB3-04GS	52	20CBV-SB1-02	50	20PPV-SE2-01MALD	50
20ACV-SB4-02	53	20CBV-SB1-02BSPT	50	20PPV-SE2-02	50
20ACV-SB4-04	53	20CBV-SB1-04	50	20PPV-SE2-02MALD	50
20ACV-SB8-04	52	20CBV-SB1-04BSPT	50	20PPV-SE2-04	50
20ACV-SB8-04GS	52	20CBV-SB2-02	50	20PPV-SE2-04MALD	50
20ACV-SB8-025	52	20CBV-SB2-03	50	20PPV-SE3-01	52
20ACV-SB9-04	50	20CBV-SB2-04	50	20PPV-SE3-01MALD	52
20ACV-SB9-025	50	20CBV-SB3-02	52	20PPV-SE3-02	52
20ACV-SB10-04	53	20CBV-SB3-03	52	20PPV-SE3-02MALD	52
20ACV-SB19-1032	52	20CBV-SB3-04	52	20PPV-SE3-04	52
20CB-PB1-02	47	20CBV-SB8-04	52	20PPV-SE3-04MALD	52
20CB-PB1-02BSPT	47	20CBV-SB8-025	52	20PPV-SE8-04	52
20CB-PB1-04	47	20CBV-SB9-04	50	20PPV-SE8-04MALD	52
20CB-PB1-04BSPT	47	20CBV-SB9-025	50	20PPV-SE8-025MALD	52
20CB-PB2-02	47	20CBV-SB12-02	52	20PPV-SE9-04	50
20CB-PB2-03	47	20CBV-SB19-1032	52	20PPV-SE9-04MALD	50
20CB-PB2-04	47	20PP-PE2-01	47	20PPV-SE9-025MALD	50
20CB-PB3-02	49	20PP-PE2-01MALD	47	20PPX-PE1-02	47
20CB-PB3-03	49	20PP-PE2-02	47	20PPX-PE1-02MALD	47
20CB-PB3-04	49	20PP-PE2-02MALD	47	20PPX-PE4-02	49
20CB-PB4-02	49	20PP-PE2-04	47	20PPX-PE4-02MALD	49
20CB-PB4-04	49	20PP-PE2-04MALD	47	20PPX-PE4-04	49
20CB-PB8-04	49	20PP-S1-02	46	20PPX-PE4-04MALD	49
20CB-PB8-025	49	20PP-S1-02MALD	46	20PPX-PE9-04	47
20CB-PB9-04	47	20PP-S1-04	46	20PPX-PE9-04MALD	47
20CB-PB9-025	47	20PP-S1-04MALD	46	20PPX-PE10-04	49
20CB-PB10-04	49	20PP-S2-01	46	20PPX-PE10-04MALD	49
20CB-PB13-02	47	20PP-S2-01MALD	46	30AC-PB2-02	54
20CB-PB13-02BSPT	47	20PP-S2-02	46	30AC-S2-01	54
20CB-PB13-04	47	20PP-S2-02MALD	46	30AC-S2-02	54
20CB-PB13-1032	47	20PP-S2-04	46	30AC-S2-03	54
20CB-S1-02	46	20PP-S2-04MALD	46	30AC-S2-04	54
20CB-S1-02BSPT	46	20PP-S3-01	48	30AC-S9-025	54
20CB-S1-04	46	20PP-S3-02	48	30ACV-PB1-02	55
20CB-S1-04BSPT	46	20PP-S3-02MALD	48	30ACV-PB1-02BSPT	55
20CB-S2-02	46	20PP-S3-04	48	30ACV-PB1-04	55
20CB-S2-03	46	20PP-S3-04MALD	48	30ACV-PB1-04BSPT	55
20CB-S2-04	46	20PP-S8-04	48	30ACV-PB2-01	55
20CB-S3-02	48	20PP-S8-04MALD	48	30ACV-PB2-02	55
20CB-S3-03	48	20PP-S8-025	48	30ACV-PB2-03	55
20CB-S3-04	48	20PP-S8-025MALD	48	30ACV-PB2-04	55
20CB-S8-04	48	20PP-S9-04	46	30ACV-PB4-02	55
20CB-S8-025	48	20PP-S9-04MALD	46	30ACV-PB4-04	55
20CB-S9-04	46	20PP-S9-025	46	30ACV-PB9-04	55
20CB-S9-025	46	20PPV-PE1-02	51	30ACV-PB9-025	55
20CB-S12-02	48	20PPV-PE1-02MALD	51	30ACV-PB10-04	55
20CB-S19-1032	48	20PPV-PE2-01	51	30ACV-SB2-01	54
20CBV-PB1-02	51	20PPV-PE2-01MALD	51	30ACV-SB2-02	54
20CBV-PB1-02BSPT	51	20PPV-PE2-02	51	30ACV-SB2-03	54
20CBV-PB1-04	51	20PPV-PE2-02MALD	51	30ACV-SB2-04	54
20CBV-PB1-04BSPT	51	20PPV-PE2-04	51	30ACV-SB9-04	54
20CBV-PB2-02	51	20PPV-PE2-04MALD	51	30ACV-SB9-025	54
20CBV-PB2-03	51	20PPV-PE4-02	53	30ACX-PB1-02	54
20CBV-PB2-04	51	20PPV-PE4-02MALD	53	30ACX-PB1-02BSPT	54
20CBV-PB3-02	53	20PPV-PE4-04	53	30ACX-PB1-04	54
20CBV-PB3-03	53	20PPV-PE4-04MALD	53	30ACX-PB1-04BSPT	54
20CBV-PB3-04	53	20PPV-PE9-04	51	30ACX-PB4-02	54
20CBV-PB4-02	53	20PPV-PE9-04MALD	51	30ACX-PB4-04	54
20CBV-PB4-04	53	20PPV-PE10-04	53	30ACX-PB9-04	54
20CBV-PB8-04	53	20PPV-PE10-04MALD	53	30ACX-PB10-04	54
20CBV-PB8-025	53	20PPV-SE1-02	50	33PCV-SPS2-02	55
20CBV-PB9-04	51	20PPV-SE1-02MALD	50	40AC-PB1-04	57
20CBV-PB9-025	51	20PPV-SE1-04	50	40AC-PB1-04BSPT	57
20CBV-PB10-04	53	20PPV-SE1-04MALD	50	40AC-PB1-06	57
20CBV-PB13-1032	51	20PPV-SE2-01	50	40AC-PB1-06BSPT	57

40AC-PB2-04	57	40ACV-SB17-04	59	40CBV-SB2-05	60
40AC-PB2-05	57	40ACV-SB18-04	61	40CBV-SB2-06	60
40AC-PB2-06	57	40ACX-PB10-04	59	40CBV-SB3-04	61
40AC-PB3-04	58	40CB-PB1-04	57	40CBV-SB3-05	61
40AC-PB3-05	58	40CB-PB1-04BSPT	57	40CBV-SB3-06	61
40AC-PB3-06	58	40CB-PB1-06	57	40CBV-SB8-04	61
40AC-PB4-04	59	40CB-PB1-06BSPT	57	40CBV-SB8-06	61
40AC-PB4-06	59	40CB-PB2-04	57	40CBV-SB9-04	60
40AC-PB6	57	40CB-PB2-05	57	40CBV-SB9-06	60
40AC-PB8-04	58	40CB-PB2-06	57	40CBV-SB12-04	61
40AC-PB8-06	58	40CB-PB3-04	58	40CBV-SB13-04	60
40AC-PB9-04	57	40CB-PB3-05	58	40PP-PE2-04	57
40AC-PB9-06	57	40CB-PB3-06	58	40PP-PE2-04MALD	57
40AC-PB10-06	59	40CB-PB4-04	59	40PP-PE2-06	57
40AC-PB17-04	57	40CB-PB4-06	59	40PP-PE2-06MALD	57
40AC-PB18-04	58	40CB-PB8-04	58	40PP-PE3-04	58
40AC-S1-04	56	40CB-PB8-06	58	40PP-PE3-04MALD	58
40AC-S1-04BSPT	56	40CB-PB9-04	57	40PP-PE3-06	58
40AC-S1-06	56	40CB-PB9-06	57	40PP-PE3-06MALD	58
40AC-S1-06BSPT	56	40CB-PB10-04	59	40PP-S1-04	56
40AC-S2-04	56	40CB-PB13-04	57	40PP-S1-04MALD	56
40AC-S2-05	56	40CB-S1-04	56	40PP-S1-06	56
40AC-S2-06	56	40CB-S1-04BSPT	56	40PP-S1-06MALD	56
40AC-S3-04	58	40CB-S1-06	56	40PP-S2-04	56
40AC-S3-05	58	40CB-S1-06BSPT	56	40PP-S2-04MALD	56
40AC-S3-06	58	40CB-S2-04	56	40PP-S2-06	56
40AC-S8-04	58	40CB-S2-06	56	40PP-S2-06MALD	56
40AC-S8-06	58	40CB-S3-04	58	40PP-S3-04	58
40AC-S9-04	56	40CB-S3-05	58	40PP-S3-04MALD	58
40AC-S9-06	56	40CB-S3-06	58	40PP-S3-06	58
40AC-S17-04	56	40CB-S8-04	58	40PP-S3-06MALD	58
40AC-S18-04	58	40CB-S8-06	58	40PP-S8-06	58
40ACV-PB1-04BSPT	60	40CB-S9-04	56	40PP-S9-06	56
40ACV-PB1-06	60	40CB-S9-06	56	40PP-S9-06MALD	56
40ACV-PB1-06BSPT	60	40CB-S12-04	58	40PPV-PE1-04	61
40ACV-PB2-04	60	40CB-S13-04	56	40PPV-PE1-04MALD	61
40ACV-PB2-05	60	40CB-S3557	56	40PPV-PE2-04	61
40ACV-PB2-06	60	40CB-S3559	56	40PPV-PE2-04MALD	61
40ACV-PB3-04	62	40CB-S3562	56	40PPV-PE2-06	61
40ACV-PB3-05	62	40CB-S3564	56	40PPV-PE2-06MALD	61
40ACV-PB3-06	62	40CB-S35271	59	40PPV-PE3-04	62
40ACV-PB4-04	62	40CBV-PB1-04	60	40PPV-PE3-04MALD	62
40ACV-PB4-06	62	40CBV-PB1-04BSPT	60	40PPV-PE3-06	62
40ACV-PB8-04	62	40CBV-PB1-06	60	40PPV-PE3-06MALD	62
40ACV-PB8-06	62	40CBV-PB1-06BSPT	60	40PPV-PE4-04	62
40ACV-PB9-04	60	40CBV-PB2-04	60	40PPV-PE4-04MALD	62
40ACV-PB9-06	60	40CBV-PB2-05	60	40PPV-PE4-06	62
40ACV-PB10-04	62	40CBV-PB2-06	61	40PPV-PE4-06MALD	62
40ACV-PB10-06	62	40CBV-PB3-04	62	40PPV-PE9-06	61
40ACV-PB17-04	60	40CBV-PB3-05	62	40PPV-PE9-06MALD	61
40ACV-PB18-04	62	40CBV-PB3-06	62	40PPV-PE10-06MALD	62
40ACV-SB1-04	59	40CBV-PB4-04	62	40PPV-SE1-04	60
40ACV-SB1-04BSPT	59	40CBV-PB4-06	62	40PPV-SE1-04MALD	60
40ACV-SB1-06	59	40CBV-PB8-04	62	40PPV-SE1-06	60
40ACV-SB1-06BSPT	59	40CBV-PB8-06	62	40PPV-SE1-06MALD	60
40ACV-SB2-04	59	40CBV-PB9-04	61	40PPV-SE2-04	60
40ACV-SB2-05	59	40CBV-PB9-06	61	40PPV-SE2-04MALD	60
40ACV-SB2-06	59	40CBV-PB10-04	62	40PPV-SE2-06	60
40ACV-SB3-04	61	40CBV-PB10-06	62	40PPV-SE2-06MALD	60
40ACV-SB3-05	61	40CBV-PB13-04	61	40PPV-SE3-04	61
40ACV-SB3-06	61	40CBV-SB1-04	59	40PPV-SE3-04MALD	61
40ACV-SB8-04	61	40CBV-SB1-04BSPT	59	40PPV-SE3-06	61
40ACV-SB8-06	61	40CBV-SB1-06	59	40PPV-SE3-06MALD	61
40ACV-SB9-04	59	40CBV-SB1-06BSPT	59	40PPV-SE8-06	61
40ACV-SB9-06	59	40CBV-SB2-04	60	40PPV-SE8-06MALD	61

40PPV-SE9-06.....	60	50BAC-S2-05LB.....	64	60PP-PE3-06.....	70
40PPV-SE9-06MALD	60	50BAC-S2-06.....	64	60PP-PE3-08.....	70
40PPX-PE1-04.....	57	50BAC-S2-06LB.....	64	60PP-PE3-08BLK	70
40PPX-PE1-04MALD	57	50BAC-S3-04.....	65	60PP-PE3-12.....	70
40PPX-PE4-04.....	59	50BAC-S3-05.....	65	60PP-PE4-06.....	71
40PPX-PE4-04MALD	59	50BAC-S8-04.....	65	60PP-PE4-08.....	71
40PPX-PE4-06.....	59	50BAC-S8-06.....	65	60PP-PE22-08.....	71
40PPX-PE4-06MALD	59	50BAC-S9-06LB.....	64	60PP-S1-06	68
40PPX-PE9-06.....	57	50BACV-SB1-04.....	66	60PP-S1-08.....	68
40PPX-PE9-06MALD	57	50BACV-SB1-06.....	66	60PP-S1-08BLK	68
42AB-PB2-02	63	50BACV-SB2-04.....	66	60PP-S1-12	68
42AB-S3-02	63	50BACV-SB2-05.....	66	60PP-S1-12BLK	68
50AC-PB1-06.....	64	50BACV-SB2-06.....	66	60PP-S2-04	68
50AC-PB2-04.....	64	50BACV-SB3-04.....	67	60PP-S2-04BLK	68
50AC-PB2-05	64	50BACV-SB3-05.....	67	60PP-S2-06	68
50AC-PB2-06.....	64	50BACV-SB3-06.....	67	60PP-S2-08	68
50AC-PB4-06.....	65	50BACV-SB8-04.....	67	60PP-S2-08BLK	68
50AC-PB9-04.....	64	50BACV-SB8-06.....	67	60PP-S2-12	68
50AC-PB9-06.....	64	50BACV-SB9-04.....	66	60PP-S3-04.....	70
50AC-PB17-04.....	64	50GP-PF1-04.....	64	60PP-S3-06	70
50AC-S1-04.....	63	50GP-PF5-04.....	64	60PP-S3-08	70
50AC-S1-06.....	63	50GP-PF5-06.....	64	60PP-S3-12	70
50AC-S2-04.....	63	50GP-PF6.....	64	60PP-S22-06	70
50AC-S2-04LB.....	63	50GP-PF7-06.....	65	60PP-S22-08	70
50AC-S2-05LB.....	63	50GP-S1-04.....	64	60PPV-PE1-06.....	72
50AC-S2-06.....	63	50GP-S5-04.....	64	60PPV-PE1-06BLK	72
50AC-S2-06LB.....	63	50GP-S5-06.....	64	60PPV-PE1-08.....	72
50AC-S3-04.....	65	50GP-S5-06LB.....	64	60PPV-PE1-08BLK	72
50AC-S3-05.....	65	50GP-S7-06.....	65	60PPV-PE1-12.....	72
50AC-S3-06.....	65	50GPV-SF5-06.....	66	60PPV-PE2-04.....	72
50AC-S8-04.....	65	50PPV-PE1-04.....	66	60PPV-PE2-04BLK	72
50AC-S8-06.....	65	50PPV-PE1-06.....	66	60PPV-PE2-06.....	72
50AC-S9-04LB.....	63	50PPV-PE2-04.....	66	60PPV-PE2-06BLK	72
50AC-S9-06LB.....	63	50PPV-PE2-06.....	66	60PPV-PE2-08.....	72
50AC-S17-04.....	63	50PPV-PE4-04.....	67	60PPV-PE2-08BLK	72
50AC-S18-04.....	65	50PPV-PE4-06.....	67	60PPV-PE2-12.....	72
50AC-S35107	63	50PPV-SE2-04.....	66	60PPV-PE3-04.....	73
50ACV-PB1-04.....	66	50PPV-SE2-06.....	66	60PPV-PE3-06.....	73
50ACV-PB1-06.....	66	50PPV-SE3-04.....	67	60PPV-PE3-08.....	73
50ACV-PB2-04.....	66	50PPV-SE3-06.....	67	60PPV-PE3-12.....	73
50ACV-PB2-05.....	66	50PPX-PE1-04.....	64	60PPV-PE4-06.....	74
50ACV-PB2-06.....	66	50PPX-PE1-06.....	64	60PPV-PE4-08.....	74
50ACV-PB4-04.....	67	50PPX-PE2-04.....	64	60PPV-PE22-06.....	74
50ACV-PB4-06.....	67	50PPX-PE2-06.....	64	60PPV-PE22-08.....	74
50ACV-PB9-04.....	66	50PPX-PE4-04.....	65	60PPV-PE22-08BLK	74
50ACV-PB9-06.....	66	50PPX-PE4-06.....	65	60PPV-SE1-06.....	71
50ACV-PB10-04.....	67	55GPV-PE1-04.....	68	60PPV-SE1-08.....	71
50ACV-PB10-06.....	67	55GPV-PE2-04.....	68	60PPV-SE1-08BLK	71
50ACV-PB17-04.....	66	55GPV-PE2-06.....	68	60PPV-SE1-12.....	71
50ACV-S35107	66	55GPV-SE1-04.....	68	60PPV-SE2-04.....	71
50ACV-SB1-06.....	66	55GPV-SE2-04.....	68	60PPV-SE2-06.....	71
50ACV-SB2-04.....	66	55GPV-SE2-06.....	68	60PPV-SE2-06BLK	71
50ACV-SB2-05.....	66	60PP-PE1-06.....	69	60PPV-SE2-08.....	71
50ACV-SB2-06.....	66	60PP-PE1-08.....	69	60PPV-SE2-08BLK	71
50ACV-SB3-04.....	67	60PP-PE1-08BLK	69	60PPV-SE2-12.....	71
50ACV-SB3-05.....	67	60PP-PE1-12.....	69	60PPV-SE3-04.....	73
50ACV-SB3-06.....	67	60PP-PE1-12BLK	69	60PPV-SE3-06.....	73
50ACV-SB8-04.....	67	60PP-PE2-04.....	69	60PPV-SE3-08.....	73
50ACV-SB8-06.....	67	60PP-PE2-04BLK	69	60PPV-SE3-12.....	73
50ACV-SB9-04.....	66	60PP-PE2-06.....	69	60PPV-SE22-06.....	74
50ACV-SB9-06.....	66	60PP-PE2-08.....	69	60PPV-SE22-08.....	74
50ACV-SB17-04.....	66	60PP-PE2-08BLK	69	60PPV-SE22-08BLK	74
50ACV-SB18-04.....	67	60PP-PE2-12.....	69	60PS-PE1-06	69
50BAC-S1-04.....	64	60PP-PE2-12BLK	69	60PS-PE1-06BLK	69
50BAC-S2-04.....	64	60PP-PE3-04.....	70	60PS-PE1-08.....	69

60PS-PE1-08BLK	69	65GP-PF1-08	75	16250.....	125
60PS-PE2-06	69	65GP-PF2-08	75	16410.....	125
60PS-PE2-06BLK	69	65GP-S2-12	75	16420.....	125
60PS-PE2-08	69	70GPV-PE2-04	75	16430.....	125
60PS-PE2-08BLK	69	70GPV-PE2-06	75	16445.....	125
60PS-PE2-12	69	70GPV-SE2-04	75	16450.....	125
60PS-PE2-12BLK	69	70GPV-SE2-06	75	18210.....	126
60PS-PE3-06	70	100AC2-01.....	77	18220.....	126
60PS-PE3-08	70	100AC2-02.....	77	18230.....	126
60PS-PE3-12	70	100AC2-M3.....	77	18240.....	126
60PS-PE4-06	71	100AC-PB2-01	76	18410.....	126
60PS-PE4-08	71	100AC-PB2-02	76	18420.....	126
60PS-PE4-08BLK	71	100AC-S2-01	76	18430.....	126
60PS-PE22-06	71	100AC-S2-02	76	18445.....	126
60PS-PE22-08	71	100AC-S2-M3	76	AA	23
60PS-S1-06	69	100ACV2-01.....	77	ABR004	123
60PS-S1-06BLK	69	100ACV2-02.....	77	ABR007	123
60PS-S1-08	69	100ACV2-M3.....	77	ABR013.....	123
60PS-S1-08BLK	69	100ACV-SB2-01	77	ABSM L	19
60PS-S2-06	69	100ACV-SB2-02	77	AC	24
60PS-S2-06BLK	69	100-P2	76	AD	24
60PS-S2-08	69	100-PM3	76	AF7500077.....	132
60PS-S2-08BLK	69	100PP2-01	77	AF7500077D	132
60PS-S2-12	69	100PP2-02	77	AF7500155.....	132
60PS-S2-12BLK	69	100PP2-M3	77	AF7500285C.....	132
60PS-S3-06	70	100PP-PE2-01.....	76	AF7500311.....	132
60PS-S3-06BLK	70	100PP-PE2-02.....	76	AF7500506.....	132
60PS-S3-08	70	100PP-S2-01	76	AF7500584.....	132
60PS-S3-08BLK	70	100PP-S2-02	76	AF7500753.....	132
60PS-S3-12	70	100PP-S2-M3.....	76	B	119
60PSV-PE1-06.....	72	100-S2	76	Bag Ports.....	94
60PSV-PE1-06BLK	72	1040.....	25	Barb Dimensions	138
60PSV-PE1-08.....	72	1050.....	25	Barb Performance	140
60PSV-PE1-08BLK	72	1401.....	127	BDF16.....	131
60PSV-PE2-06.....	72	1402.....	127	BDF18.....	131
60PSV-PE2-06BLK	72	1600.....	125	BDFM6	131
60PSV-PE2-08.....	72	1610.....	125	BDFS.....	131
60PSV-PE2-08BLK	72	1620.....	125	BDJ01	129
60PSV-PE2-12.....	72	1630.....	125	BDJ02	129
60PSV-PE2-12BLK	72	1640.....	125	BDJ03	129
60PSV-PE3-06.....	73	1650.....	125	BDJ04	129
60PSV-PE3-06BLK	73	1800.....	127	BDJ05	129
60PSV-PE3-08.....	73	1801.....	127	BDJ06	129
60PSV-PE3-12.....	73	1802.....	127	BDJ07	129
60PSV-PE4-06.....	74	1810.....	126	BDJ08	129
60PSV-PE4-08.....	74	1820.....	126	BDJ09	129
60PSV-PE22-06.....	74	1830.....	126	BDJ10	129
60PSV-PE22-06BLK	74	1840.....	126	BDJ11	129
60PSV-PE22-08.....	74	1850.....	126	BDJ12	129
60PSV-SE1-06	71	1860.....	126	BDJ13	129
60PSV-SE1-08	71	2040.....	25	BDJ14	129
60PSV-SE1-08BLK	71	2050.....	25	BDJ15	129
60PSV-SE2-06	71	3040.....	25	BDJ16	129
60PSV-SE2-06BLK	72	3050.....	25	BDJ17	129
60PSV-SE2-08	72	3060.....	25	BDJ18	129
60PSV-SE2-08BLK	72	4040.....	23	BDJ19	129
60PSV-SE2-12	72	4050.....	25	BDJ20	129
60PSV-SE2-12BLK	72	4060.....	25	BDJ21	129
60PSV-SE3-06	73	5050.....	23	BDJ22	129
60PSV-SE3-06BLK	73	5060.....	26	BDM16.....	131
60PSV-SE3-08	73	6060.....	23	BDM18.....	131
60PSV-SE3-12	73	16210.....	125	BDMD.....	130
60PSV-SE3-12BLK	73	16220.....	125	BDMFTLL.....	130
60PSV-SE22-06	74	16230.....	125	BDMK.....	131
60PSV-SE22-08	74	16240.....	125	BDML210	130

BDML220.....	130	CC.....	23	FCLBFMR-001.....	93
BDML230.....	130	CCLR.....	21, 37	FLB035.....	136
BDML410.....	130	CCR.....	21, 37	FLB055.....	136
BDML420.....	130	CF229.....	117	FLB065.....	136
BDML430.....	130	CF230.....	117	FLBP.....	137
BDMM6.....	131	CF494.....	28	FMD5P670-001.....	107
BDMMTL.....	130	Check Valves.....	84	FSLLR.....	14, 17
BDMR210.....	130	Chemical Resistance of Resin.....	146	FTL10.....	16
BDMR220.....	130	CPAN75.....	102	FTL20.....	16
BDMR230.....	130	CPAX150.....	102	FTL30.....	16
BDMR240.....	130	CPN75.....	103	FTL10.....	16
BDMR410.....	130	CPT90.....	104	FTL220.....	16
BDMR420.....	130	CPT100.....	104	FTL230.....	17
BDMR430.....	130	CPT110.....	104	FTL240.....	17
BDMR445.....	130	CPX150.....	103	FTLB210.....	21
BDMRP.....	130	Custom Quick Connect Solutions.....	83	FTLB220.....	21
BDMS.....	131	CYLINDRAFlow™ Manifolds.....	107	FTLB230.....	21
BDMT210.....	131	D.....	119	FTLC.....	17
BDMT220.....	131	DC.....	24	FTLE.....	17
BDMT230.....	131	DCV114.....	87	FTLL004.....	15
BDMT410.....	131	DCV115.....	87	FTLL007.....	15
BDMT420.....	131	DCV118.....	87	FTLL013.....	15
BDMT430.....	131	DCV125.....	87	FTLL025.....	16
BDMX.....	131	DD.....	23	FTLL035.....	16
BioClamp®.....	99	Design Center.....	10	FTLL055.....	16
BioProcess Fittings.....	94	Diver™ Stopcocks.....	90	FTLL210.....	15
BioValve®.....	100	Double Check Valves.....	87	FTLL220.....	15
Blood Pressure Fittings.....	114	Elbow Connectors.....	27	FTLL230.....	15
BPA013.....	116	ENTFE098.....	39	FTLL240.....	16
BPA025.....	116	ENTFTC.....	39	FTLL250.....	16
BPA035.....	116	ENTM230.....	39	FTLL360.....	16
BPA055.....	116	ENTME050.....	39	FTLL410.....	15
BPAM.....	116	ENTME064.....	39	FTLL420.....	15
BPC025.....	117	ENTME078.....	39	FTLL430.....	15
BPC220.....	116	ENTME098.....	39	FTLL445.....	16
BPC230.....	116	ENTMTP.....	39	FTLL450.....	16
BPC240.....	116	FC4FFMR-001.....	91	FTLLB210.....	20
BPC250.....	117	FC4LB3F-001.....	93	FTLLB220.....	20
BPCL025.....	117	FCFFM-002.....	91	FTLLB230.....	20
BPCL220.....	116	FCFFM-003.....	91	FTLLBMLRL.....	20
BPCL230.....	116	FCFFM-008.....	91	FTLLC.....	17
BPCL240.....	116	FCFFM-009.....	91	FTLLE.....	17
BPCL250.....	117	FCFFM-010.....	92	FTLLP.....	18
BPF220.....	114	FCFFM-013.....	92	FTLLSB004.....	16
BPF230.....	114	FCFFMR-002.....	92	FTLLSB007.....	16
BPF240.....	114	FCFFMR-003.....	92	FTLLSB013.....	16
BPF250.....	114	FCFFMR-004.....	92	FTLP.....	18
BPF Series Properties & Performance.....	145	FCFM-001.....	92	FTLT.....	17
BPLF230.....	117	FCFM-002.....	92	G40MPGR-TF150.....	112
BPLF240.....	117	FCFM-006.....	92	G40RX-PX-0250.....	113
BPLM230.....	117	FCFM-010.....	92	G40RX-PX-0400.....	113
BPLM240.....	117	FCFML-2A-001.....	92	G42MPGR-TF075.....	112
BPM013.....	39, 115	FCFMR-001.....	93	G424TRX-0200.....	112
BPM025.....	39, 115	FCFMR-2A-001.....	93	G424TRX-0300.....	112
BPM035.....	39, 115	FCFMR-2A-004.....	93	Gaskets.....	112
BPM240.....	39, 115	FCHPFFMR-001.....	91	GM4749TRX-0600.....	113
BPMB230.....	115	FCHPFFMR-002.....	91	GMN450.....	112
BPMB240.....	115	FCHPFFMR-003.....	90	GMX455.....	112
BPMB250.....	115	FCHPFFMR-002.....	91	Instrument Tees.....	111
BPMMTL.....	115	FCHPFFMR-003.....	91	ISO 80369.....	38
BSFTLL.....	20	FCHPFFMR-005.....	91	ISO 80369-3 Compliant Connectors.....	39
BSMLRL.....	19	FCLBFFM-001.....	93	ISO 80369-5 Compliant Connectors.....	39
Build-a-Part® Fittings.....	128	FCLBFFM-1A-003.....	93	K10.....	119
C.....	119	FCLBFFM-003.....	93	K210.....	119
CA091.....	18	FCLBFFMR-1A-003.....	93	K220.....	119

K230	119	LLDPE02-3031C.....	132	MTLL007	13
K240	119	LLDPE02-3031M	132	MTLL013	13
K410	119	LNS	21, 37	MTLL025	13
K420	119	LP4.....	18	MTLL035	13
K430	119	LP34.....	18	MTLL055	13
K445	119	LT787	17	MTLL210	13
KFTL.....	20	LT878	17	MTLL220	13
KFTLL.....	20	Luers	12	MTLL230	13
KL210.....	120	M6P.....	124	MTLL240	13
KL220.....	120	M5210.....	124	MTLL250	13
KL230.....	120	M5220.....	124	MTLL360	13
KL410.....	120	M5230.....	124	MTLL410	13
KL420.....	120	M5240.....	124	MTLL420	13
KL430.....	120	M6210.....	124	MTLL430	13
KN.....	122	M6220.....	124	MTLL445	13
KP.....	122	M6230.....	124	MTLL450	13
KT210	121	M6240.....	124	MTLLP.....	18
KT220	121	M6250.....	124	MTLP.....	18
KT230	121	M6410.....	124	MTLS210.....	15
KT410	121	M6420.....	124	MTLS220.....	15
KT420	121	M6430.....	124	MTLS230.....	15
KT430	121	M6445.....	124	MTLS240.....	15
L10.....	27	M6450.....	124	MTLSP.....	18
L20.....	27	Manifold Connectors	34	Multiport Connectors	34
L30.....	27	MLB035.....	136	N004	23
L40.....	27	MLB055.....	136	N007	23
L50.....	27	MLB065.....	136	N007/004	24
L60.....	27	MLBC7	136	N013	23
L070.....	28	MLBC8	136	N013/004	24
L070/055.....	29	MLBP	137	N013/007	25
L080	28	MLRL004.....	14	N025	23
L080/070.....	29	MLRL007.....	14	N025/004	25
L090.....	28	MLRL013.....	14	N025/007	25
L100.....	28	MLRL025.....	14	N025/013	25
L100/080.....	29	MLRL035.....	14	N035	23
L110.....	28	MLRL055.....	14	N035/007	25
L110/100.....	29	MLRLB210.....	21	N035/013	25
L210.....	27	MLRLB220.....	21	N035/025	25
L220.....	27	MLRLB230.....	21	N055	23
L220/210.....	28	MLRLC	17	N055/013	25
L230.....	27	MLRLP.....	18	N055/025	25
L230/210.....	28	MLSL004.....	14	N055/035	26
L230/220.....	28	MLSL007.....	14	N065	23
L410.....	27	MLSL013.....	14	N065/025	26
L420.....	27	MLSL025.....	14	N065/035	26
L420/410.....	28	MLSL035.....	14	N065/055	26
L430.....	27	MLSL055.....	14	N070	24
L430/220.....	28	MLSLP.....	18	N070/013	26
L430/410.....	28	MQCF230.....	78	N070/055	26
L430/420.....	28	MQCF240.....	78	N080	24
L440.....	27	MQCF250.....	78	N080/055.....	26
L970.....	28	MQCMF230.....	78	N080/070	26
L975.....	28	MQCMF240.....	78	N090	24
L6110.....	28	MQCMF250.....	78	N090/055	26
LC34	17	MQCRM230	78	N090/070	26
LC78	17	MQCRM240	78	N100	24
LE87	17	MQCRM250	78	N100/080	26
LE7220.....	17	MSL004.....	15	N110	24
LE7230.....	18	MSL007.....	15	N110/100	26
LLDPE02-1019C.....	132	MSL013.....	15	N210	23
LLDPE02-1019M	132	MSL025.....	15	N220	23
LLDPE02-1031C.....	132	MSL035.....	15	N220/210	24
LLDPE02-1031M	132	MSL055.....	15	N230	23
LLDPE02-2031C.....	132	MTLCS	17	N230/210	24
LLDPE02-2031M	132	MTLL004	13	N230/220	24

N410	23	PIP50.....	36	S220	122
N420	23	PIP60.....	36	S230	122
N420/410.....	24	PIP070.....	36	S240	123
N430	23	PIP080.....	36	S250	123
N430/410.....	24	PIP090.....	36	S410	122
N430/420.....	25	PIP210.....	36	S420	122
N440	23	PIP220.....	36	S430	123
N670	24	PIP230.....	36	S445	123
N670/655.....	26	PIP240.....	36	S450	123
N680	24	PIP670.....	36	Sanitary Fittings	108
N680/655.....	26	PIP680.....	36	SBLF013	80
N680/670.....	26	PIP6100.....	36	SBLF035	80
N690	24	PIP6110.....	36	SBLM013	80
N690/655.....	26	Plug Connectors	36	SBLM035	80
N690/670.....	26	PM2N	115	SBLMFGB013.....	80
N970	24	PM3N	115	SBLMFGB035.....	80
N975	24	PM4N	82	SBLMFGB013.....	80
N6100	24	PMK210	37	SBLMFGB035.....	80
N6100/670.....	26	PMK220	37	SCV67FL/220.....	86
N6100/680.....	26	PMK220-210.....	37	SCV67FL/230.....	86
N6110	24	PMKN	21, 37	SCV06018.....	84
N6110/680.....	26	PMS230	37	SCV06246.....	84
N6110/6100.....	26	PMS230-210.....	37	SCV06248.....	84
Needlefree Swabable Valves	88	PMS230-220.....	37	SCV06249.....	85
NEFCV-50.....	41	PMSN	21, 37, 115	SCV06252.....	85
NEFCV-CM064	41	PVCA1190035.....	133	SCV06264.....	85
NEME21-50	41	PVCA1190049.....	133	SCV06265.....	85
NEME21-57	41	PVCA1190112.....	133	SCV06272.....	85
NEME31-50	41	PVCA1190126.....	133	SCV06308.....	85
NEMP-50	41	PVCA1190189.....	133	SCV06329.....	85
NEMP-57	41	PVCA1190203.....	133	SCV06397.....	85
NEMP-CM064	41	PVCA1190266.....	133	SCV06398.....	85
NEMPV-50	41	PVCA1190280.....	133	SCV08053.....	85
Neuraxial Connectors	40	PVCA1190343.....	133	SCV10090.....	85
NVF13VPC.....	88	PVCA1190357.....	133	SCV17090.....	85
NVF20PC.....	89	Quick Connect Couplings	42	SCV21053.....	85
NVF20VPC.....	88	Quick Connect Performance Ratings ...	141	SCV22050.....	85
NVF25PC.....	89	RB01	130	SCV23050.....	86
NVF28PC.....	89	Reducer Connectors	24	SCV67220.....	86
NVF37PC.....	89	Resin Sterilization Stability	148	SCV67220/230.....	86
NVF41PC.....	89	RMLLR.....	130	SCV67220/ML.....	86
NVFBGPC.....	89	RQCF655	79	SCV67230.....	86
NVFBV60PC	89	RQCF670	79	SCV67230/220.....	86
NVFBV66PC	89	RQCFC.....	79	SCV67230/ML.....	86
NVFMILL-001	89	RQCFL655	79	SDAFCNPS.....	106
NVFMILLCP	89	RQCFL670	79	SDAM655NPS	105
NVFMILLPC	89	RQCFLC.....	79	SDAM670NPS	106
NVFMILLS-001	89	RQCM655	79	SDAM680NPS	106
NVFM LPC.....	89	RQCM670	79	SDAMPNPS	106
NVFM LPCB.....	89	RQCMP.....	79	SDAVF655NPS	105
NVFSCMLLPC	89	RQXF670.....	79	SDAVF670NPS	106
NVFTPC.....	89	RQXF680.....	79	SDAVF680NPS	106
NV FY36PC.....	89	RQXFC	79	SDAVM655NPS	105
NV FY42S-001.....	89	RQXFL670.....	79	SDAVM670NPS	106
NV FY401CP.....	89	RQXFL680.....	79	SDAVM680NPS	106
NV FY401PC.....	89	RQXFLC	79	SFBP630	95
NVM31S-001.....	89	RQXM670.....	79	SFBP635/635.....	98
NVM41S-001.....	89	RQXM680.....	79	SFBP655	95
Opposing Barb Ports	98	RQXMP.....	79	SFBP655/655.....	98
Panel Mount Instructions.....	145	S10	122	SFBP670	95
PharmaLok™ Clamps	101	S20	122	SFBP670/670.....	98
PIP10.....	36	S30	122	SFBP680	95
PIP20.....	36	S40	123	SFBP680/680.....	98
PIP30.....	36	S50	123	SFBP6100	95
PIP40.....	36	S210	122	SFBP6100/6100.....	98

SFBP6110	95	T220/210	30	VBV1000NW	100
SFBPF1	96	T220/230	30	VLP4	18
SFBPF3	96	T220R210	30	VNG075WHT	99
SFBPF4	96	T230	29	VNG150WHT	99
SFM3C	97	T230/210	31	VNG200WHT	99
SFM4C	97	T230/220	31	VNG250WHT	99
SFM100/SFM075	110	T230R210	31	VNG300WHT	99
SFM150/SFM075	110	T230R220	31	VNG400WHT	99
SFM150/SFM100	110	T250	29	VNG600WH	99
SFM200/SFM075	110	T360/220	31	VNG800WH	99
SFM200/SFM150	110	T360/230	31	VPB0850400N	137
SFMN630	108	T360/250	31	VPM0480201N	137
SFMN655	108	T410	29	VPM1060200N	137
SFMN670	108	T410/420	30	VPM1080303N	137
SFMN680	108	T410/430	30	VPM2930601N	137
SFMN6100	108	T420	29	VPS3901074N	137
SFMNC	108	T420/410	30	VPS5069200N	137
SFMX655	109	T420/430	31	VPS5401001N	86
SFMX670	109	T420R410	30	VPS5401020N	86
SFMX680	109	T430	29	VPS5401036SN	86
SFMX6100	109	T430/410	31	VPS5401068N	86
SFMX6110	109	T430/420	31	X10	119
SFMXC	109	T430R410	31	X210	119
SFMXT680	111	T430R420	31	X220	119
SFMXT6110	111	T670	30	X230	119
SFTL	20	T670/655	32	X240	119
SFTLL	20	T680	30	X410	119
SILC7001098	133	T680/655	32	X420	119
SILC7001280	133	T680/670	32	X430	119
SILC7001422	133	T690	30	X445	119
SILC7001424	133	T970	30	XFTL	20
SILC7001657	133	T975	30	XFTLL	20
SILC7001839	133	T6100	30	XL210	120
SILC7002203	133	T6100/670	32	XL220	120
SILC7002593	133	T6100/680	32	XL230	120
SILC7002775	133	T6110	30	XL410	120
SILC7003139	133	T6110/6100	32	XL420	120
SILPG7021041	133	Tapered Seal	134	XL430	120
SILPG7021042	133	TECS116	123	XMLRL	19
SILPG7021043	133	TECS118	123	XMTLL	19
SILPG7021044	133	Tee Connectors	29	XN	122
SMLRL	19	Threaded Connectors	118	XP	122
SMTLL	19	Trademarks	149	XQCBM755	81
SP	123	TSC2F007	135	XQCBM770	81
Spaulding™ Aseptic Disconnects	105	TSC2F013	135	XQCLM755	81
T10	29	TSC2F035	135	XQCLM770	81
T013	29	TSC2F210	135	XQCM755	81
T20	29	TSC2FB230	135	XQCM770	81
T30	29	TSC2M013	135	XQCVBM755	82
T40	29	TSC2M035	135	XQCVBM770	82
T055	30	TSC2M210	135	XQCVLM755	82
T60	29	TSC2MP	135	XQCVM755	82
T070	30	Tube Fittings	22	XQCVM765	82
T070/055	32	Tubing	132	XQCVM770	82
T080	30	Typical Resin Properties	149	XQSBF755	81
T080/055	32	U10	37	XQSBF770	81
T080/070	32	U20	37	XQSF755	81
T100	30	U30	37	XQSF770	81
T100/080	32	UR00-1031C	133	XQSVBF755	82
T110	30	UR00-1031M	133	XQSVBF770	82
T110/100	32	UR00-2047C	133	XQSVF755	82
T210	29	UR00-2047M	133	XQSVF765	82
T210/220	30	UR00-3062C	133	XQSVF770	82
T210/230	30	UR00-3062M	133	XT210	121
T220	29	VBV0500NW	100	XT220	121

XT230.....	121
XT410.....	121
XT420.....	121
XT430.....	121
Y013	32
Y055	33
Y070	33
Y080	33
Y090	33
Y100	33
Y110	33
Y210	32
Y220	32
Y230	32
Y240	32
Y250	33
Y260	33
Y410	32
Y420	32
Y430	32
Y450	33
Y460	33
Y630/655	34
Y655/670	34
Y670	33
Y670/680	34
Y680	33
Y690	33
Y6100	33
Y6100/6110	34
Y6110	33
Y Connectors.....	32

About Nordson MEDICAL

Nordson MEDICAL is a global integrated solutions partner with more than 30 years of experience in design, engineering, and manufacturing of complex medical devices and components. Nordson MEDICAL specializes in components and services used in interventional and minimally invasive surgical products including catheters, balloons, extrusions, polyimide and composite tubing, heat shrink tubing, reinforced shafts, cleanroom injection molding, and finished device assembly and packaging. Visit Nordson MEDICAL at nordsonmedical.com.

©Nordson MEDICAL. All rights reserved. **PCV029 Rev. A** 2020-04



Value Plastics, Inc.
dba Nordson MEDICAL
805 West 71st Street
Loveland, CO 80538, USA

Phone: (970) 267-5200
Toll Free: (888) 404-5837

FMF-sales@nordsonmedical.com
www.nordsonmedical.com